

RQ-2019-07-22-67

Log-in Checklist

Login Station ID: #2

Shipping Method: FedEx UPS | HD | Greyhound

Date/Time of Receipt: 8-6-19/9:50

Cooler Temperature*	Samples Frozen YES	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
		HNO ₃	Formalin		Present?		* Intact?		
					Yes	No	Yes	No	
2.1				5		✓			

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C) without login confirmation from customer for samples requiring ice preservation. An NCR still is required.

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

****Chain of Custody/ Field Sheet(s) Included?** Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☐ No ☐

CONTAINER CHECK

****Evidence Tape on Bottles?** Yes ☒ No ☐ If yes, is it intact? Yes ☒ No ☐
****Caps on tight?** Yes ☒ No ☐ ****Containers intact?** Yes ☒ No ☐
****Sufficient Sample Volume:** Yes ☒ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ☐ No ☐ Init:

Chlorine detected? (One container checked per Field ID)

Yes ☐ No ☐ Init:

If "YES", list affected Field IDs and Test IDs below.

Field ID	List Test ID(s) where Chlorine was Detected, an NCR is required.

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

****Acid Preserved Samples:** pH ≤ 2.0? Yes ☐ No ☐ NA ☒ Init: SR
****W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?** Yes ☐ No ☐ NA ☒ Init: SR
Coolers Unpacked/Checked by: SR **Date:** 8-6-19 **NCRs (Y/N)?** N

Event ID: Algal Bloom Response - as needed
WQAP-2019-08-06-04 (BLOOM-RESP)

NCR #

Log-in Checklist

Samples with Incorrect pH Preservation

Field ID	Bottle Test ID(s)	pH	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap	Damaged Container	Evidence Tape Not Intact
		CK	BR	CK	BR

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_02_2018 EXP:4-29-2020
pH Test Strips 0 – 6 Lot #	213316BVEXP:05-15-21
pH Paper 0 – 3 Lot #	220416A EXP:07-30-20
pH Paper 0 – 13 Lot #	222516 EXP:08-15-2019
Chlorine Test Strips Lot #	9080 EXP:01-2022

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ☒ _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Florida Department of Environmental Protection
Central Laboratory Submittal Form

Algal Bloom Response - as needed

last revised Jan 25, 2018

Customer: WQAP

Requester: Cheryl A. Swanson

Field Report Prepared By: K. Panzner

Project ID: BLOOM-RESP

Collected By: S. Meyer K. Panzner

Sampling Agency: FDEP SWROC

Site Name: South Pasadena Bloom Response - Near Sunset Dr		☐ Comp ☒ Grab		Collection (comp begin or grab)		☒ Eastern ☐ Central		Composite end		Bottle Group(s)	
Date: 08/05/2019		Date 8/5/19		Time 1105		Date		Time		(See pg 2)	
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☒ mg/L ☐ % sat		Total Res. Chlorine (mg/L)				ABC	
Temp (C) 30.00		pH 7.44		Sample Depth 0.3		☐ ft ☒ m		Sp. Conductance (umho/cm) 46989			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					

Location		☐ Comp ☐ Grab		Collection (comp begin or grab)		☐ Eastern ☐ Central		Composite end		Bottle Group(s)	
Field ID		Date		Time		Date		Time			
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)					
Temp (C)		pH		Sample Depth		☐ ft ☐ m		Sp. Conductance (umho/cm)			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					

Location		☐ Comp ☐ Grab		Collection (comp begin or grab)		☐ Eastern ☐ Central		Composite end		Bottle Group(s)	
Field ID		Date		Time		Date		Time			
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)					
Temp (C)		pH		Sample Depth		☐ ft ☐ m		Sp. Conductance (umho/cm)			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					

Location		☐ Comp ☐ Grab		Collection (comp begin or grab)		☐ Eastern ☐ Central		Composite end		Bottle Group(s)	
Field ID		Date		Time		Date		Time			
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)					
Temp (C)		pH		Sample Depth		☐ ft ☐ m		Sp. Conductance (umho/cm)			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					

Location		☐ Comp ☐ Grab		Collection (comp begin or grab)		☐ Eastern ☐ Central		Composite end		Bottle Group(s)	
Field ID		Date		Time		Date		Time			
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)					
Temp (C)		pH		Sample Depth		☐ ft ☐ m		Sp. Conductance (umho/cm)			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					

Relinquished By: <i>Ryan</i>	Date/Time: 8/5/19 1700	Shipping Method: Fedex	Number of Coolers Shipped: 1	Received By: <i>SR</i>	Date/Time: 8-6-19/9:50
------------------------------	------------------------	------------------------	------------------------------	------------------------	------------------------

Group: A	Bottle Type:	PT-50ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
BLOOM-ID							
Group: B	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W							
Group: B	Bottle Type:	BG-250ML-WDMT H	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
W-MCYST-AA							
Group: B	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3							
W-NO2NO3							
W-S-A-TP							
W-TKN							
Group: B	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F							
Group: C	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	Lugols	Enter Number of Bottles Sent to Lab	1
DTY-QN-C							
PKD-QN-C							

Date of Request: 16-JUL-2019
 Created By: SWANSON_C On: 16-JUL-2019
 Modified By: WATTS_J On: 17-JUL-2019
 Customer: WQAP
 Project: BLOOM-RESP
 Division: Division of Environmental Assessment and Restoration
 District:
 Sampling Event: Algal Bloom Response - as needed

Send Coolers To:

Phone: 813-470-5951
 DEP SW ROC
 13051 North Telecom Parkway
 Temple Terrace, FL 33637-0926
 Attn: Matthew Bearden

Send Final Report To:

2600 Blair Stone Rd
 MS #3560
 Tallahassee, FL 32399-2400
 Attn: Julie Espy

Program Module Number:

Priority: 1
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 17-JUL-2019
 Biology Request Reviewed By: JASROTIA_P On: 16-JUL-2019
 Sampling Kit Required: YES Ship on: 19-JUL-2019
 Sampling Kit Shipped: YES On: 18-JUL-2019
 Sampling Kit Packed By: SHAIK_A
 Preservatives Needed: NO
 Date to Receive Samples: 22-JUL-2019
 Received By:

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Unless peak algal bloom season, may analyze toxin after consultation with Biology's Taxonomy Manager.

Bottle Groups A/B - separate event, priority 1
 Bottle Group C - separate event, priority 12

Bottle Group A - Priority 1
 For algal dominant taxon.

Bottle Group B - Priority 1
 For toxin testing, nutrients, and chlorophyll

Bottle Group C - Priority 12
 For algal full ID

There are multiple bottles per Group in case multiple sample targets or sample locations are needed. Samplers use best professional judgement to decide which samples are appropriate to collect for a particular response. All sets of bottles do not have to be filled at one time.

Bottle Group A (Biological) With 2 Samples:

Bottle Type: PT-50ML Number of Bottles: 2 Preserved With: ICE
 BLOOM-ID Template: DEFAULT (SOP-AB05) Assessment of dominant algal taxa in bloom or mat sample

Bottle Group B (Water) With 2 Samples:

Bottle Type: BP-1L Number of Bottles: 2 Preserved With: ICE
 CHLSUITE-W Template: DEFAULT (SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Type: BG-250ML-WDMTH Number of Bottles: 2 Preserved With: ICE
 W-MCYST-AA Template: DEFAULT (EPA 8321B) Microcystins in water matrices by HPLC/MS/MS

Bottle Type: P-500ML Number of Bottles: 2 Preserved With: ICE And H2SO4
 W-NH3 Template: DEFAULT (EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
 W-NO2NO3 Template: DEFAULT (EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
 W-S-A-TP Template: DEFAULT (EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
 W-TKN Template: DEFAULT (EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices

Bottle Type: P-125ML Number of Bottles: 2 Preserved With: FILTER-ICE
 W-PO4-F Template: DEFAULT (EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group C (Water) With 2 Samples:

Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: Lugols
DTY-QN-C	Template: DEFAULT	(SOP-AB05_1) No. of diatom taxa of quantitative phytoplankton sample.
PKD-QN-C	Template: DEFAULT	(SOP-AB05) No. of wet taxa of quantitative phytoplankton sample.

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.

Do Not Lift Using This Tag

ORIGIN ID: TLHA (813) 470-5772
BEN PASHWATER
FL DEP SOUTHWEST DISTRICT
13051 N TELECOM PKWY

TEMPLE TERRACE, FL 336370926
UNITED STATES US

SHIP DATE: 22JUL19
ACTWGT: 9.00 LB MAN
CAD: 0446970/CAFE3211

TO **JOHN WATTS**
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD


TALLAHASSEE FL 32399

(850) 245-8077

REF:

DEPT:

RMA: ||| ||| |||

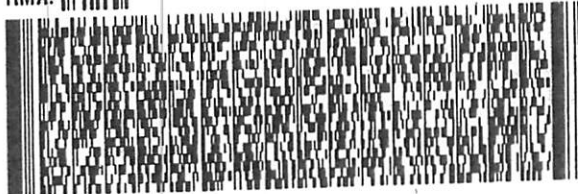

FedEx
Express



J18111806000700

FedEx

TRK# 4751 8782 0020
0221

RETURNING MON ED
TUE - 06 AUG 10:30A
PRIORITY OVERNIGHT

36 TLHA

32399
FL-US
TLH



FTD 5211412 05AUG19 MCFA 568C2/F551/0C8A

Cooler Packing Worksheet for Request: RQ-2019-07-22-67

Algal Bloom Response - as needed

Ship Cooler On: 7/19/2019

Customer/Project: WQAP/BLOOM-RESP

Assigned To: SHAIK_A

Requester: Cheryl A. Swanson

Priority: 1

813-470-5951
DEP SW ROC
13051 North Telecom Parkway
Temple Terrace, FL 33637-0926

Attn: Matthew Bearden

Comments:

No preservatives needed.

Unless peak algal bloom season, may analyze toxin after consultation with Biology's Taxonomy Manager.

Bottle Groups A/B - separate event, priority 1

Bottle Group C - separate event, priority 12

Bottle Group A - Priority 1

For algal dominant taxon.

Bottle Group B - Priority 1

For toxin testing, nutrients, and chlorophyll

Bottle Group C - Priority 12

For algal full ID

There are multiple bottles per Group in case multiple sample targets or sample locations are needed. Samplers use best professional judgement to decide which samples are appropriate to collect for a particular response. All sets of bottles do not have to be filled at one time.

Requested Analyses:

Bottle Group: A

of Sites: 2

Container ID: PT-50ML

Qty: 2 Preservation: ICE

Lot # 03775003

Description: Plastic Tube 50 mL

Analysis

BLOOM-ID

Description

Assessment of dominant algal taxa in bloom or mat sample

Bottle Group: B

of Sites: 2

Container ID: BP-1L

Qty: 2 Preservation: ICE

Lot # 1248230

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: BG-250ML-WDMTH

Qty: 2 Preservation: ICE

Lot # 88071670

Description: Wide Mouth Amber Glass Bottle 250 ml

Analysis

W-MCYST-AA

Description

Microcystins in water matrices by HPLC/MS/MS

Container ID: P-500ML

Qty: 2 Preservation: ICE And H2SO4

Lot # 88074433

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3
W-NO2NO3
W-S-A-TP
W-TKN

Description

Ammonia in aqueous matrices as mg N/L
Nitrite/Nitrate in aqueous matrices as mg N/L
Total Phosphorus in aqueous matrices as mg P/L
Total Kjeldahl Nitrogen in aqueous matrices

Container ID: P-125ML

Qty: 2 Preservation: FILTER-ICE

Lot # 88074433

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: C

of Sites: 2

Container ID: BP-1L

Qty: 2 Preservation: Lugols

Lot # 1248230

Description: Brown Plastic Bottle 1L

Analysis

DTY-QN-C
PKD-QN-C

Description

No. of diatom taxa of quantitative phytoplankton sample.
No. of wet taxa of quantitative phytoplankton sample.

Cooler Packed By: ASB

Number of Coolers: 1

Date: 7-18-19

Kit must also include:

- ☒ Field Sheets
- ☒ Temperature Control Bottle
- ☒ FedEx Bills, if applicable
- ☒ Plastic Bags/Zip Ties
- ☒ Acid Labels (sets of 24)

If Preservation Included:

ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____

PLEASE RETURN ALL COOLERS FOR RQ-2019-07-22-67