

Chemical Analysis Report

WQAP-2022-09-14-03

Florida DEP Laboratory
2600 Blair Stone Road
Tallahassee, FL 32399-2400
DOH Accreditation E31780

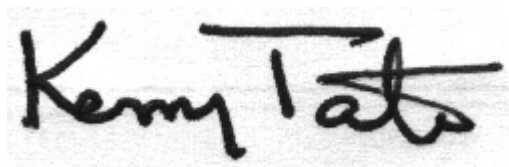
Event Description: **Algae Bloom Response NEROC-as needed**
Request ID: **RQ-2022-09-12-54**
Customer: **WQAP**
Project ID: **BLOOM-RESP**

Send Reports to:
FL Dept. of Environmental Protection
8800 Baymeadows Way W
Suite 100
Tallahassee, FL 32399-2400
Attn: Brian G Davis

For additional information please contact
Colin Wright, Ph.D.
Liang-Tsair Lin, Ph.D.
Kerry Tate, Ph.D.
Dr. rer. nat. Bettina Steinbock
Thekkekalathil Chandrasekhar, Ph.D, QA Officer
Phone (850) 245-8085

Certified by: Kerry Tate, Ph.D., Environmental Administrator

Date Certified: 28-SEP-2022 09:58

A handwritten signature in black ink that reads "Kerry Tate". The signature is written in a cursive, flowing style with a horizontal line extending from the end of the name.

Case Narrative

Unless otherwise noted, all samples included in this report were received in accordance with protocols referenced in Chapter 62-160, Florida Administrative Code (F.A.C.). Results published in this report pertain only to the samples as submitted to, and received by the laboratory. All times in this report are adjusted to the applicable Eastern Time Zone (EST or EDT).

Results for the following analytical groups are included in this report: Nutrients and Pesticides.

Scientific notation may be used in reporting very large or small values. Values reported using scientific notation will take the form of the following example: 1.3E+03, which is equivalent to 1.3×10^3 or 1300.

Unless otherwise noted, analytical values for soil and sediment samples are reported on a dry weight basis, and analytical values for waste and tissue samples are reported on a wet weight basis.

Results for TNI accredited tests met requirements established by The NELAC Institute. A double asterisk (**) is used to indicate an analyte/matrix/method for which the laboratory is not TNI accredited by the Florida Department of Health Environmental Laboratory Certification Program or where accreditation for that field of testing is not applicable.

Any significant anomalies or deviations from established protocols are documented in Non-Conformance Reports, which, where appropriate, are included within this analytical report. Additional comments related to specific analytical tests may be included as remarks following the analytical results for each sample. Such comments and remarks are for informational purposes only and are not intended to convey judgement about the usability of the reported data.

A quality control report on the performance of the test method for the submitted samples is included. Uncertainty associated with the analytical results contained in this report can be estimated from the reported quality assurance results and from published quality control acceptance limits for each analytical test. Matrix quality control results (matrix spike recoveries and matrix sample precision) pertain only to the matrix sample tested and do not necessarily reflect test method performance for other samples.

Typical matrix quality control (QC) measurements may include matrix spike recovery, matrix spike duplicate recovery, matrix spike precision and matrix sample precision. Not all matrix QC results may be available or reportable; where they are not an explanation is provided. Typical reasons for unavailable QC results include, but are not limited to, a) insufficient matrix sample to perform some or all QC measurements; b) analyte concentration in the sample replicated was too low for a meaningful measurement of precision and c) analyte concentration in the matrix sample spiked was too high (relative to the amount of analyte spiked) for a meaningful measurement of recovery. Where matrix QC results are unavailable, other method performance metrics (e.g., LCS recovery, LCS precision, surrogate recovery) may be used to assess performance of the method. Comments explaining any missing QC measurements are not intended to convey any adverse conclusions about the quality of the reported data.

Precision is reported as relative percent difference unless otherwise noted.

Quality Control codes as defined below may be used in this report to indicate results that are associated with one or more quality control elements which did not fall within established test method criteria. Such results may be qualified as estimates using a J qualifier as required by 62-160 F.A.C. Explanations are included in the report for any results that were reported as estimates for other reasons.

QC Codes used in this report may include:

- LCS – Recovery for the batch Laboratory Control Sample (LCS) was outside existing control limits;
- MS – Recovery for the batch matrix spike (MS) was outside existing control limits;
- CCV – Recovery for a continuing calibration verification (CCV) standard was outside existing control limits;
- SUR – Recovery of a surrogate (SUR) for associated analytes was outside existing control limits;
- RPD – The precision, measured as relative percent difference (RPD), of batch replicate measurements was outside existing control limits;
- RSD – The precision, measured as relative standard deviation (RSD), of batch replicate measurements was outside existing control limits;
- SMP – Sample - used precision derived from replicate analyses of a sample;

The following data qualifiers are used, where applicable, in this report as specified in 62-160 F.A.C.

- A - Value reported is the mean of two or more determinations.
- B - Results based on colony counts outside the acceptable range.
- I - The reported value is between the laboratory method detection limit and the laboratory practical quantitation limit.
- J - Estimated value and/or the analysis did not meet established quality control criteria.
- K - Actual value is known to be less than value given.
- L - Actual value is known to be greater than value given.
- N - Presumptive evidence of presence of material.
- O - Sampled, but analysis lost or not performed.
- Q - Sample held beyond normal holding time.
- T - Value reported is less than the criterion of detection.
- U - Material was analyzed for but not detected. The reported value is the method detection limit for the sample analyzed.
- V - Analyte was detected in both sample and method blank.
- X - Too few individuals to calculate SCI value.
- Y - The laboratory analysis was from an unpreserved or improperly preserved sample. The data may not be accurate.
- Z - Colonies were too numerous to count (TNTC).

Quality control information from overflow laboratories may not be included in this report. Please refer to the associated report from the overflow laboratory for additional information.

Sample Location: Lochloosa Lake @ Park

Collection Date/Time: 09/13/2022 11:00

Field ID: NEHAB0082

Matrix: W-SURF-FRH

Sample ID	Ref. Method	Component	Result	Code	Units	Batch ID	QC Codes
2356209	EPA 8321B	Anatoxin-a	0.29	I	ug/L	P419476	
		Cylindrospermopsin	0.10	U	ug/L	P419474	
		Desmethyl microcystin LR	0.25	U	ug/L	P419474	
		Microcystin HILR	0.25	U	ug/L	P419474	
		Microcystin HtyR	0.25	U	ug/L	P419474	
		Microcystin LA	0.10	U	ug/L	P419474	
		Microcystin LF	0.10	U	ug/L	P419474	
		Microcystin LR	0.25	U	ug/L	P419474	
		Microcystin LW	0.25	U	ug/L	P419474	
		Microcystin LY	0.10	U	ug/L	P419474	
		Microcystin RR	0.10	U	ug/L	P419474	
		Microcystin WR	0.50	U	ug/L	P419474	
		Neosaxitoxin	0.50	U	ug/L	P419476	
		Microcystin YR	0.25	U	ug/L	P419474	
		Saxitoxin	0.50	U	ug/L	P419476	
		Nodularin-R	0.10	U	ug/L	P419474	
2356215	EPA 350.1 Rev. 2.0	Ammonia-N	0.007		mg N/L	P419789	
	EPA 351.2 Rev. 2.0	Kjeldahl Nitrogen	1.8		mg N/L	P419806	
	EPA 353.2 Rev. 2.0	NO2NO3-N	0.004	U	mg N/L	P419990	
	EPA 365.1 Rev. 2.0	Total-P	0.075		mg P/L	P419690	
2356217	EPA 365.1 Rev. 2.0 dissolved	O-Phosphate-P	0.004	U	mg P/L	P419562	

Sample Location: Violet Lake @ 130 Violet Cir

Collection Date/Time: 09/13/2022 12:05

Field ID: NEHAB0087

Matrix: W-SURF-FRH

Sample ID	Ref. Method	Component	Result	Code	Units	Batch ID	QC Codes
2356210	EPA 8321B	Anatoxin-a	0.25	U	ug/L	P419580	
		Cylindrospermopsin	0.34	I	ug/L	P419585	
		Desmethyl microcystin LR	0.25	U	ug/L	P419585	
		Microcystin HILR	0.25	U	ug/L	P419585	
		Microcystin HtyR	0.25	U	ug/L	P419585	
		Microcystin LA	0.10	UJ	ug/L	P419585	MS
		Microcystin LF	0.10	UJ	ug/L	P419585	MS
		Microcystin LR	0.25	U	ug/L	P419585	
		Microcystin LW	0.25	U	ug/L	P419585	
		Microcystin LY	0.10	U	ug/L	P419585	
		Microcystin RR	0.10	U	ug/L	P419585	
		Microcystin WR	0.50	U	ug/L	P419585	
		Neosaxitoxin	0.50	U	ug/L	P419580	
		Microcystin YR	0.25	U	ug/L	P419585	
		Saxitoxin	0.50	U	ug/L	P419580	
		Nodularin-R	0.10	U	ug/L	P419585	
2356216	EPA 350.1 Rev. 2.0	Ammonia-N	0.002	U	mg N/L	P419642	
	EPA 351.2 Rev. 2.0	Kjeldahl Nitrogen	0.40		mg N/L	P419808	
	EPA 353.2 Rev. 2.0	NO2NO3-N	0.022		mg N/L	P419990	
	EPA 365.1 Rev. 2.0	Total-P	0.012		mg P/L	P419690	
2356218	EPA 365.1 Rev. 2.0 dissolved	O-Phosphate-P	0.004	U	mg P/L	P419562	

Ref. Method and Comment:

EPA 8321B: Refer to the Lab Analysis Report for an explanation of QC Codes.

Quality Assurance Report Method Blank Results

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419642

Component	Result	Code	Units
Ammonia-N	0.002	U	mg N/L

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419789

Component	Result	Code	Units
Ammonia-N	0.002	U	mg N/L

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419806

Component	Result	Code	Units
Kjeldahl Nitrogen	0.080	U	mg N/L

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419808

Component	Result	Code	Units
Kjeldahl Nitrogen	0.080	U	mg N/L

Reference Method: EPA 353.2 Rev. 2.0
Batch ID: P419990

Component	Result	Code	Units
NO2NO3-N	0.004	U	mg N/L

Reference Method: EPA 365.1 Rev. 2.0
Batch ID: P419690

Component	Result	Code	Units
Total-P	0.002	U	mg P/L

Reference Method: EPA 365.1 Rev. 2.0 dissolved
Batch ID: P419562

Component	Result	Code	Units
O-Phosphate-P	0.004	U	mg P/L

Reference Method: EPA 8321B
Batch ID: P419474

Component	Result	Code	Units
Cylindrospermopsin	0.10	U	ug/L
Desmethyl microcystin LR	0.25	U	ug/L
Microcystin HIR	0.25	U	ug/L
Microcystin HtyR	0.25	U	ug/L
Microcystin LA	0.10	U	ug/L
Microcystin LF	0.10	U	ug/L
Microcystin LR	0.25	U	ug/L
Microcystin LW	0.25	U	ug/L
Microcystin LY	0.10	U	ug/L
Microcystin RR	0.10	U	ug/L

Quality Assurance Report

Method Blank Results

Reference Method: EPA 8321B
Batch ID: P419474

Component	Result	Code	Units
Microcystin WR	0.50	U	ug/L
Microcystin YR	0.25	U	ug/L
Nodularin-R	0.10	U	ug/L

Reference Method: EPA 8321B
Batch ID: P419476

Component	Result	Code	Units
Anatoxin-a	0.25	U	ug/L
Neosaxitoxin	0.50	U	ug/L
Saxitoxin	0.50	U	ug/L

Reference Method: EPA 8321B
Batch ID: P419580

Component	Result	Code	Units
Anatoxin-a	0.25	U	ug/L
Neosaxitoxin	0.50	U	ug/L
Saxitoxin	0.50	U	ug/L

Reference Method: EPA 8321B
Batch ID: P419585

Component	Result	Code	Units
Cylindrospermopsin	0.10	U	ug/L
Desmethyl microcystin LR	0.25	U	ug/L
Microcystin HllR	0.25	U	ug/L
Microcystin HtyR	0.25	U	ug/L
Microcystin LA	0.10	U	ug/L
Microcystin LF	0.10	U	ug/L
Microcystin LR	0.25	U	ug/L
Microcystin LW	0.25	U	ug/L
Microcystin LY	0.10	U	ug/L
Microcystin RR	0.10	U	ug/L
Microcystin WR	0.50	U	ug/L
Microcystin YR	0.25	U	ug/L
Nodularin-R	0.10	U	ug/L

Quality Assurance Report Laboratory Control Sample Accuracy

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419642

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Ammonia-N	94.4		P	90 - 110

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419789

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Ammonia-N	98.8		P	90 - 110

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419806

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Kjeldahl Nitrogen	96.2		P	90 - 110

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419808

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Kjeldahl Nitrogen	104		P	90 - 110

Reference Method: EPA 353.2 Rev. 2.0
Batch ID: P419990

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
NO2NO3-N	97.5		P	90 - 110

Reference Method: EPA 365.1 Rev. 2.0
Batch ID: P419690

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Total-P	104		P	90 - 110

Reference Method: EPA 365.1 Rev. 2.0 dissolved
Batch ID: P419562

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
O-Phosphate-P	98.6		P	90 - 110

Reference Method: EPA 8321B
Batch ID: P419474

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Cylindrospermopsin	83.4		P	40 - 150
Desmethyl microcystin LR	96.0		P	40 - 150
Microcystin HIR	103		P	40 - 150
Microcystin HtyR	100		P	40 - 150
Microcystin LA	112		P	40 - 150
Microcystin LF	117		P	40 - 150
Microcystin LR	104		P	40 - 150
Microcystin LW	109		P	40 - 150
Microcystin LY	114		P	40 - 150
Microcystin RR	101		P	40 - 150

Quality Assurance Report Laboratory Control Sample Accuracy

Reference Method: EPA 8321B
Batch ID: P419474

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Microcystin WR	99.7		P	40 - 150
Microcystin YR	97.8		P	40 - 150
Nodularin-R	99.1		P	40 - 150

Reference Method: EPA 8321B
Batch ID: P419476

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Anatoxin-a	80.2		P	40 - 150
Neosaxitoxin	62.0		P	40 - 150
Saxitoxin	66.9		P	40 - 150

Reference Method: EPA 8321B
Batch ID: P419580

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Anatoxin-a	87.8		P	40 - 150
Neosaxitoxin	89.2		P	40 - 150
Saxitoxin	83.0		P	40 - 150

Reference Method: EPA 8321B
Batch ID: P419585

Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
Cylindrospermopsin	99.6		P	40 - 150
Desmethyl microcystin LR	117		P	40 - 150
Microcystin HllR	120		P	40 - 150
Microcystin HtyR	120		P	40 - 150
Microcystin LA	140		P	40 - 150
Microcystin LF	135		P	40 - 150
Microcystin LR	117		P	40 - 150
Microcystin LW	132		P	40 - 150
Microcystin LY	136		P	40 - 150
Microcystin RR	119		P	40 - 150
Microcystin WR	129		P	40 - 150
Microcystin YR	104		P	40 - 150
Nodularin-R	115		P	40 - 150

Quality Assurance Report Matrix Spike Accuracy

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419642

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356136	Ammonia-N	97.4	96.2	P/P	90 - 110

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419789

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356068	Ammonia-N	95.2	94.0	P/P	90 - 110

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419806

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356297	Kjeldahl Nitrogen	103	107	P/P	90 - 110

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419808

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356161	Kjeldahl Nitrogen	102	106	P/P	90 - 110

Reference Method: EPA 353.2 Rev. 2.0
Batch ID: P419990

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356215	NO2NO3-N	95.5	96.5	P/P	90 - 110

Reference Method: EPA 365.1 Rev. 2.0
Batch ID: P419690

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356216	Total-P	106	104	P/P	90 - 110

Reference Method: EPA 365.1 Rev. 2.0 dissolved
Batch ID: P419562

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356172	O-Phosphate-P	93.0	90.9	P/P	90 - 110

Reference Method: EPA 8321B
Batch ID: P419474

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2355681	Cylindrospermopsin	93.1	106	P/P	40 - 150
2355681	Desmethyl microcystin LR	126	121	P/P	40 - 150
2355681	Microcystin H1LR	121	121	P/P	40 - 150
2355681	Microcystin HtyR	126	122	P/P	40 - 150
2355681	Microcystin LA	123	117	P/P	40 - 150
2355681	Microcystin LF	132	137	P/P	40 - 150
2355681	Microcystin LR	124	122	P/P	40 - 150
2355681	Microcystin LW	124	124	P/P	40 - 150
2355681	Microcystin LY	126	124	P/P	40 - 150
2355681	Microcystin RR	127	121	P/P	40 - 150

Quality Assurance Report Matrix Spike Accuracy

Reference Method: EPA 8321B
Batch ID: P419474

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2355681	Microcystin WR	119	121	P/P	40 - 150
2355681	Microcystin YR	119	116	P/P	40 - 150
2355681	Nodularin-R	131	127	P/P	40 - 150

Reference Method: EPA 8321B
Batch ID: P419476

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2355681	Anatoxin-a	79.4	82.2	P/P	40 - 150
2355681	Neosaxitoxin	42.5	42.2	P/P	40 - 150
2355681	Saxitoxin	53.3	54.1	P/P	40 - 150

Reference Method: EPA 8321B
Batch ID: P419580

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356210	Anatoxin-a	105	99.1	P/P	40 - 150
2356210	Neosaxitoxin	73.9	72.3	P/P	40 - 150
2356210	Saxitoxin	63.4	64.0	P/P	40 - 150

Reference Method: EPA 8321B
Batch ID: P419585

Spiked Sample	Component	% Rec.1	% Rec.2	Pass/Fail	Control Limits
2356210	Cylindrospermopsin	99.0	97.6	P/P	40 - 150
2356210	Desmethyl microcystin LR	111	125	P/P	40 - 150
2356210	Microcystin HilR	114	122	P/P	40 - 150
2356210	Microcystin HtyR	110	125	P/P	40 - 150
2356210	Microcystin LA	130	151	P/F	40 - 150
2356210	Microcystin LF	135	151	P/F	40 - 150
2356210	Microcystin LR	110	114	P/P	40 - 150
2356210	Microcystin LW	117	126	P/P	40 - 150
2356210	Microcystin LY	130	140	P/P	40 - 150
2356210	Microcystin RR	112	128	P/P	40 - 150
2356210	Microcystin WR	109	112	P/P	40 - 150
2356210	Microcystin YR	103	109	P/P	40 - 150
2356210	Nodularin-R	105	115	P/P	40 - 150

Quality Assurance Report Precision

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419642

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356136	Ammonia-N	1.18	Spike	P	0 - 20

Reference Method: EPA 350.1 Rev. 2.0
Batch ID: P419789

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356068	Ammonia-N	0.880	Spike	P	0 - 20

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419806

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356297	Kjeldahl Nitrogen	2.42	Spike	P	0 - 20

Reference Method: EPA 351.2 Rev. 2.0
Batch ID: P419808

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356161	Kjeldahl Nitrogen	1.81	Spike	P	0 - 20

Reference Method: EPA 353.2 Rev. 2.0
Batch ID: P419990

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356215	NO2NO3-N	1.04	Spike	P	0 - 20

Reference Method: EPA 365.1 Rev. 2.0
Batch ID: P419690

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356216	Total-P	1.65	Spike	P	0 - 20

Reference Method: EPA 365.1 Rev. 2.0 dissolved
Batch ID: P419562

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356172	O-Phosphate-P	2.07	Spike	P	0 - 20

Reference Method: EPA 8321B
Batch ID: P419474

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2355681	Cylindrospermopsin	12.3	Spike	P	0 - 30
2355681	Desmethyl microcystin LR	3.72	Spike	P	0 - 30

Quality Assurance Report Precision

Reference Method: EPA 8321B
Batch ID: P419474

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2355681	Microcystin HILR	0.0497	Spike	P	0 - 30
2355681	Microcystin HtyR	2.79	Spike	P	0 - 30
2355681	Microcystin LA	5.48	Spike	P	0 - 30
2355681	Microcystin LF	3.95	Spike	P	0 - 30
2355681	Microcystin LR	1.19	Spike	P	0 - 30
2355681	Microcystin LW	0.0547	Spike	P	0 - 30
2355681	Microcystin LY	1.81	Spike	P	0 - 30
2355681	Microcystin RR	4.32	Spike	P	0 - 30
2355681	Microcystin WR	2.15	Spike	P	0 - 30
2355681	Microcystin YR	3.17	Spike	P	0 - 30
2355681	Nodularin-R	3.18	Spike	P	0 - 30

Reference Method: EPA 8321B
Batch ID: P419476

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2355681	Anatoxin-a	3.42	Spike	P	0 - 30
2355681	Neosaxitoxin	0.766	Spike	P	0 - 30
2355681	Saxitoxin	1.49	Spike	P	0 - 30

Reference Method: EPA 8321B
Batch ID: P419580

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356210	Anatoxin-a	5.52	Spike	P	0 - 30
2356210	Neosaxitoxin	2.22	Spike	P	0 - 30
2356210	Saxitoxin	1.07	Spike	P	0 - 30

Reference Method: EPA 8321B
Batch ID: P419585

Replicated Lab Sample	Component	% RSD/RPD	Sample/Spike/LCS*	Pass/Fail	Control Limits
2356210	Cylindrospermopsin	1.35	Spike	P	0 - 30
2356210	Desmethyl microcystin LR	12.1	Spike	P	0 - 30
2356210	Microcystin HILR	6.87	Spike	P	0 - 30
2356210	Microcystin HtyR	13.2	Spike	P	0 - 30
2356210	Microcystin LA	15.0	Spike	P	0 - 30
2356210	Microcystin LF	11.4	Spike	P	0 - 30
2356210	Microcystin LR	3.70	Spike	P	0 - 30
2356210	Microcystin LW	7.30	Spike	P	0 - 30
2356210	Microcystin LY	7.26	Spike	P	0 - 30
2356210	Microcystin RR	12.6	Spike	P	0 - 30
2356210	Microcystin WR	3.58	Spike	P	0 - 30
2356210	Microcystin YR	5.70	Spike	P	0 - 30
2356210	Nodularin-R	8.35	Spike	P	0 - 30

* Sample, spike and/or laboratory control sample precision (LCS) is reported.

Quality Assurance Report Surrogates

Lab Sample ID: 2356209
Field Sample ID: NEHAB0082

Reference Method	Surrogate	% Rec.	Pass/Fail	Control Limits
EPA 8321B	Microcystin LR ethyl-d5	104	P	30 - 160
EPA 8321B	Neosaxitoxin-15N7	55.7	P	30 - 160

Lab Sample ID: 2356210
Field Sample ID: NEHAB0087

Reference Method	Surrogate	% Rec.	Pass/Fail	Control Limits
EPA 8321B	Microcystin LR ethyl-d5	107	P	30 - 160
EPA 8321B	Neosaxitoxin-15N7	91.8	P	30 - 160

Quality Assurance Report Calibration Verification

Reference Method: EPA 365.1 Rev. 2.0 dissolved

Run ID: A114721

Included Lab Sample IDs: 2356217, 2356218

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
O-Phosphate-P	95.8	97.7	P/P	90 - 110
O-Phosphate-P	98.6	95.8	P/P	90 - 110

Reference Method: EPA 8321B

Run ID: A114726

Included Lab Sample IDs: 2356209

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Anatoxin-a	72.7	75.3	P/P	50 - 160
Neosaxitoxin	51.6	52.7	P/P	50 - 160
Saxitoxin	58.4	58.5	P/P	50 - 160

Reference Method: EPA 8321B

Run ID: A114729

Included Lab Sample IDs: 2356209

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Cylindrospermopsin	90.3	88.0	P/P	60 - 160
Desmethyl microcystin LR	95.8	95.2	P/P	60 - 160
Microcystin HIR	97.8	100	P/P	60 - 160
Microcystin HtyR	95.6	101	P/P	60 - 160
Microcystin LA	112	116	P/P	50 - 160
Microcystin LF	104	108	P/P	60 - 160
Microcystin LR	93.5	96.3	P/P	60 - 160
Microcystin LW	85.0	86.9	P/P	60 - 160
Microcystin LY	99.4	102	P/P	60 - 160
Microcystin RR	101	102	P/P	60 - 160
Microcystin WR	93.8	95.9	P/P	60 - 160
Microcystin YR	85.2	92.1	P/P	60 - 160
Nodularin-R	100	92.4	P/P	60 - 160

Reference Method: EPA 350.1 Rev. 2.0

Run ID: A114747

Included Lab Sample IDs: 2356216

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Ammonia-N	98.4	97.2	P/P	90 - 110

Reference Method: EPA 8321B

Run ID: A114756

Included Lab Sample IDs: 2356210

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Anatoxin-a	93.3	87.5	P/P	50 - 160
Neosaxitoxin	90.3	86.5	P/P	50 - 160
Saxitoxin	85.1	83.2	P/P	50 - 160

Quality Assurance Report Calibration Verification

Reference Method: EPA 8321B
Run ID: A114759
Included Lab Sample IDs: 2356210

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Cylindrospermopsin	96.4	97.7	P/P	60 - 160
Desmethyl microcystin LR	106	120	P/P	60 - 160
Microcystin HILR	113	118	P/P	60 - 160
Microcystin HtyR	110	125	P/P	60 - 160
Microcystin LA	131	148	P/P	50 - 160
Microcystin LF	139	153	P/P	60 - 160
Microcystin LR	102	119	P/P	60 - 160
Microcystin LW	122	138	P/P	60 - 160
Microcystin LY	122	141	P/P	60 - 160
Microcystin RR	108	122	P/P	60 - 160
Microcystin WR	110	122	P/P	60 - 160
Microcystin YR	109	115	P/P	60 - 160
Nodularin-R	102	113	P/P	60 - 160

Reference Method: EPA 350.1 Rev. 2.0
Run ID: A114804
Included Lab Sample IDs: 2356215

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Ammonia-N	100	99.7	P/P	90 - 110

Reference Method: EPA 365.1 Rev. 2.0
Run ID: A114815
Included Lab Sample IDs: 2356215

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Total-P	99.7	99.2	P/P	90 - 110

Reference Method: EPA 365.1 Rev. 2.0
Run ID: A114829
Included Lab Sample IDs: 2356216

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Total-P	99.5	100	P/P	90 - 110

Reference Method: EPA 351.2 Rev. 2.0
Run ID: A114856
Included Lab Sample IDs: 2356215, 2356216

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
Kjeldahl Nitrogen	95.3	99.9	P/P	90 - 110
Kjeldahl Nitrogen	99.8	99.6	P/P	90 - 110

Reference Method: EPA 353.2 Rev. 2.0
Run ID: A114889
Included Lab Sample IDs: 2356215, 2356216

Component	% Rec.1	% Rec.2	Pass/Fail*	Control Limits
NO2NO3-N	98.5	96.5	P/P	90 - 110
NO2NO3-N	99.0	98.5	P/P	90 - 110

Quality Assurance Report

Calibration Verification

* Pass/Fail determinations are made for each bracketing calibration verification check.

Control limits for initial calibration checks may be different from those for continuing checks, depending on method requirements.

Where they are different, both control limits are provided.

Quality Assurance Report Summary

Ref. Method	Analyte	LCS % Recovery	MS % Recovery		Precision LCS	SMP	MS
EPA 350.1 Rev. 2.0	Ammonia-N	94.4	97.4	96.2			1.18
	Ammonia-N	98.8	95.2	94.0			0.880
EPA 351.2 Rev. 2.0	Kjeldahl Nitrogen	96.2	103	107			2.42
	Kjeldahl Nitrogen	104	102	106			1.81
EPA 353.2 Rev. 2.0	NO ₂ NO ₃ -N	97.5	95.5	96.5			1.04
EPA 365.1 Rev. 2.0	Total-P	104	106	104			1.65
EPA 365.1 Rev. 2.0 dissolved	O-Phosphate-P	98.6	93.0	90.9			2.07
EPA 8321B	Anatoxin-a	80.2	79.4	82.2			3.42
	Anatoxin-a	87.8	105	99.1			5.52
	Cylindrospermopsin	83.4	93.1	106			12.3
	Cylindrospermopsin	99.6	99.0	97.6			1.35
	Desmethyl microcystin LR	96.0	126	121			3.72
	Desmethyl microcystin LR	117	111	125			12.1
	Microcystin HiIR	103	121	121			0.0497
	Microcystin HiIR	120	114	122			6.87
	Microcystin HtyR	100	126	122			2.79
	Microcystin HtyR	120	110	125			13.2
	Microcystin LA	112	123	117			5.48
	Microcystin LA	140	130	151			15.0
	Microcystin LF	117	132	137			3.95
	Microcystin LF	135	135	151			11.4
	Microcystin LR	104	124	122			1.19
	Microcystin LR	117	110	114			3.70
	Microcystin LW	109	124	124			0.0547
	Microcystin LW	132	117	126			7.30
	Microcystin LY	114	126	124			1.81
	Microcystin LY	136	130	140			7.26
	Microcystin RR	101	127	121			4.32
	Microcystin RR	119	112	128			12.6
	Microcystin WR	99.7	119	121			2.15
	Microcystin WR	129	109	112			3.58
	Microcystin YR	97.8	119	116			3.17
	Microcystin YR	104	103	109			5.70
	Neosaxitoxin	62.0	42.5	42.2			0.766
	Neosaxitoxin	89.2	73.9	72.3			2.22
	Nodularin-R	99.1	131	127			3.18
	Nodularin-R	115	105	115			8.35
	Saxitoxin	66.9	53.3	54.1			1.49
	Saxitoxin	83.0	63.4	64.0			1.07

Reference Method Descriptions

Method	Description	Associated Samples
EPA 350.1 Rev. 2.0	Ammonia in aqueous matrices as mg N/L	2356215, 2356216
EPA 351.2 Rev. 2.0	Total Kjeldahl Nitrogen in aqueous matrices	2356215, 2356216
EPA 353.2 Rev. 2.0	Nitrite/Nitrate in aqueous matrices as mg N/L	2356215, 2356216
EPA 365.1 Rev. 2.0	Total Phosphorus in aqueous matrices as mg P/L	2356215, 2356216
EPA 365.1 Rev. 2.0 dissolved	Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L	2356217, 2356218
EPA 8321B	Microcystins in water matrices by HPLC/MS/MS	2356209, 2356210
EPA 8321B	Saxitoxin in water matrices by HPLC/MS/MS	2356209, 2356210

Preparation and Analysis Log

Ref. Method	Received Date	Prep Date/Time	Prepared By	Analysis Date/Time	Analyzed By	Associated Samples
EPA 350.1 Rev. 2.0	09/14/2022			09/15/2022 12:09	Ping Hua	2356216
	09/14/2022			09/19/2022 12:49	Ping Hua	2356215
EPA 351.2 Rev. 2.0	09/14/2022	09/20/2022 11:51	Alexandra J Mattheus	09/21/2022 10:55	Alexandra J Mattheus	2356215
	09/14/2022	09/20/2022 11:52	Alexandra J Mattheus	09/21/2022 11:34	Alexandra J Mattheus	2356216
EPA 353.2 Rev. 2.0	09/14/2022			09/22/2022 09:48	Beverly Y. Sanders	2356215
	09/14/2022			09/22/2022 09:59	Beverly Y. Sanders	2356216
EPA 365.1 Rev. 2.0	09/14/2022	09/16/2022 15:34	Adam P Silver	09/19/2022 14:36	James L Waggerby	2356215
	09/14/2022	09/16/2022 15:34	Adam P Silver	09/20/2022 14:42	Simonne.V. Calhoun	2356216
EPA 365.1 Rev. 2.0 dissolved	09/14/2022			09/14/2022 15:20	Chandler E Cox	2356217
	09/14/2022			09/14/2022 15:26	Chandler E Cox	2356218
EPA 8321B	09/14/2022	09/14/2022 12:00	Manjita Shrestha	09/14/2022 17:31	Umesh Chiluwai	2356209
	09/14/2022	09/14/2022 13:00	Manjita Shrestha	09/14/2022 19:10	Manjita Shrestha	2356209
	09/14/2022	09/15/2022 12:20	Manjita Shrestha	09/15/2022 13:19	Manjita Shrestha	2356210
	09/14/2022	09/15/2022 12:20	Umesh Chiluwai	09/15/2022 15:23	Umesh Chiluwai	2356210