

Log-in Checklist

RQ- 2018-07-30-20

Login Station ID: #4

Shipping Method: FedEx | UPS | HD | Greyhound | _____

Date/Time of Receipt: 8/2/18 9:55

Cooler ID	Cooler Temperature*	No. Sample Containers in Cooler	Evidence Tape				Tracking # (fill out only if waybill copy is damaged)
			Present?		*Intact?		
			Yes	No	Yes	No	
RED	23°C	13		✓			
BLUE	1.6°C	15		✓			

*If the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX), or if Evidence Seal is damaged, identify the containers in the affected cooler(s) on back of this form.

**Note: If "No" is checked in any of the fields below, fill out the back of this form with additional details (Field ID, Test IDs, etc.) and generate an NCR.

**Chain of Custody/ Field Sheet(s) Included? Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☒ No ☐ Date _____

CONTAINER CHECK

**Evidence Tape on Bottles? Yes ☐ No ☒ If yes, is it intact? Yes ☐ No ☐

**Caps on tight? Yes ☒ No ☐ **Containers intact? Yes ☒ No ☐

**Sufficient Sample Volume: Yes ☐ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes ☒ No ☐ Init: JA

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.					
Analysis	Yes	No	Analysis	Yes	No
W-1,4 DIOX		☒	W-PC-AA-SF		
W-E8321-DI, -MS, -21		☒	W-NP-AA-SF		
W-BNA (-625, -R)			W-PDQUAT		
W-CARB-AA (-CAB-AA-R)			W-PCL-SQ-R		
W-U2060-MS			W-TR-TQ		
W-PCB-R			W-CT-CI-R		
W-PAH (-R)			W-PSNP-TQ		

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ☒ No ☐ NA ☐ Init: JA

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes ☐ No ☒ NA ☐ Init: JA

Coolers Unpacked/Checked by: JA Date: 8/2/18 NCRs (Y/N)? _____

Event ID: UE-DIST-2018-08-02-01 NCR # _____

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_04_2018
pH Paper 0.0 – 3.0 Lot #	230515 EXP:10/30/2018
pH Paper 0 – 13 Lot #	222516 EXP:08/15/2019
Chlorine Test Strips Lot #	7226 EXP:07/2020

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ☒

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

LSJR Downtown Fresh

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: NE-DIST

Requester: Jeremy M. Parrish

Field Report Prepared By:

Project ID: SMP

Collected By:

Sampling Agency:

WIN /Field ID:



20030752

TROUT RIVER AT RAMALLAH RD

Latitude

☐ dd

Longitude

☐ dms

WIN /Field ID:



20030751

TROUT RIVER AT GARDEN ST

Latitude

☐ dd

Longitude

☐ dms

Location

WIN /Field ID:



20030729

Field ID

WIN/ST

HOGANS CR AT BROAD ST.

Latitude

☐ dd

Longitude

☐ dms

Location

WIN /Field ID:



20030746

Field ID

WIN/ST

BLOCKHOUSE CR @ DUVAL RD

Latitude

☐ dd

Longitude

☐ dms

Relinquished By:

Date/Time

08-01-2018

Shipping Method

FedEx

Number of Coolers Shipped:

2

Received By:

Date/Time

8/2/18 9:55

☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08-01-2018 Time 0930	Eastern Central	Composite end Date Time	Bottle Group(s) (See pg 2)
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen 5.5	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
Temp (C) 24.8	pH 5.19	Sample Depth 0.15	☐ ft ☒ m	Sp. Conductance (umho/cm) 72
Salinity (ppt) 0.03	Comments			
☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08-01-2018 Time 1000	Eastern Central	Composite end Date Time	Bottle Group(s)
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen 6.16	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
Temp (C) 25.0	pH 5.87	Sample Depth 0.1	☐ ft ☒ m	Sp. Conductance (umho/cm) 82
Salinity (ppt) 0.04	Comments			
☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08-01-2018 Time 1255	Eastern Central	Composite end Date Time	Bottle Group(s)
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen 3.95	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
Temp (C) 28.5	pH 7.32	Sample Depth 0.3	☐ ft ☒ m	Sp. Conductance (umho/cm) 487
Salinity (ppt) 0.23	Comments			
☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08-01-2018 Time 1055	Eastern Central	Composite end Date Time	Bottle Group(s)
Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen 5.5	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
Temp (C) 27.0	pH 6.76	Sample Depth 0.10	☐ ft ☒ m	Sp. Conductance (umho/cm) 151
Salinity (ppt) 0.07	Comments			

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	CHLSUITE-W						
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	NE-ALS-ECQ						→ Sent to contract lab
Group: A	Bottle Type:	QPCR-P-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	PCR-FILTER						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	W-E8321-DI						
	W-E8321-MS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	2
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	2
	W-PO4-F						
Group: B	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						

LSJR Downtown Fresh

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: NE-DIST

Requester: Jeremy M. Parrish

Field Report Prepared By: A. Kalmbacher

Project ID: SMP

Collected By: P. O'Connor

Sampling Agency: FDEP/NE ROC

Field ID:  BLANK 08012018		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08-01-2018 Time 1125	☒ Eastern ☐ Central	Composite end Date Time	Bottle Group(s) (See pg 2)
Field Blank Location: NE ROC DI		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)
WIN /Field ID:  G2NE1159		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08-01-2018 Time 1135	☒ Eastern ☐ Central	Composite end Date Time	Bottle Group(s)
Blockhouse 100m abv Irma rd.		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Temp (C)	pH	Sample Depth	☐ ft ☒ m	Sp. Conductance (umho/cm)
WIN/STORET #		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Salinity (ppt)	Comments			
Location		☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	☐ Eastern ☐ Central	Composite end Date Time	Bottle Group(s)
Field ID		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
WIN/STORET #		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Salinity (ppt)	Comments			
Location		☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	☐ Eastern ☐ Central	Composite end Date Time	Bottle Group(s)
Field ID		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
WIN/STORET #		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Salinity (ppt)	Comments			

Relinquished By:

Date/Time

Shipping Method

Number of Coolers Shipped:

Received By:

Date/Time

A. Kalmbacher 08-01-2018

Fed Ex

2

A. Kalmbacher

8/2/18 9:55

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	NE-ALS-ECQ					
Group: A	Bottle Type:	QPCR-P-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	PCR-FILTER					
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-E8321-DI W-E8321-MS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC					

Group: B	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	W-TDS						
	W-TSS						
Group: B	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	CHLSUITE-W						
Group: B	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	NE-ALS-ECQ						→ Sent to contract lab
Group: B	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	2
	W-ICP						
	W-ICPMS						
Group: B	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	2
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: B	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	2
	W-PO4-F						
Group: C	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	NE-ALS-ECQ						→ Sent to contract lab
Group: C	Bottle Type:	QPCR-P-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	PCR-FILTER						
Group: C	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	W-E8321-DI						
	W-E8321-MS						
Group: D	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						

Group: D	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	W-TDS						
	W-TSS						
Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-ICP						
	W-ICPMS						
Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: D	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
	W-PO4-F						



SR#

CAS Contract

Distribution: White - Return to Originator, Yellow - Retained by Client

Date of Request: 06-JUL-2018
 Created By: PARRISH_J On: 06-JUL-2018
 Modified By: SWANSON_C On: 09-JUL-2018
 Customer: NE-DIST
 Project: SMP
 Division:
 District: Northeast District
 Sampling Event: LSJR Downtown Fresh

Send Coolers To:

Phone: 904-256-1700
 8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256-7577
 Attn: Jeremy Parrish

Send Final Report To:

8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256-7577
 Attn: Jeremy Parrish

Program Module Number:
 Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 06-JUL-2018
 Biology Request Reviewed By: SWANSON_C On: 09-JUL-2018
 Sampling Kit Required: YES Ship on: 20-JUL-2018
 Sampling Kit Shipped: YES On: 16-JUL-2018
 Sampling Kit Packed By: SHAIK_A
 Preservatives Needed: NO
 Date to Receive Samples: 30-JUL-2018
 Received By:

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Group A: Freshwater Kit, No Metals
 Group B: Freshwater Kit With Metals

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 2	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville
Bottle Type: QPCR-P-500ML	Number of Bottles: 2	Preserved With: ICE
PCR-FILTER	Template: DEFAULT	(SOP-PCR-4.0) Preparation of Samples by Filtration and held for qPCR Analysis
Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: ICE
W-E8321-DI	Template:	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 2 Samples:

Bottle Type: P-250ML-WI-THI	Number of Bottles: 2	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville

Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Calcium
Chromium
Iron
Magnesium
Nickel
Potassium
Sodium
Zinc

W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
---------	-----------	---

Arsenic
Cadmium
Copper
Lead
Silver

Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville

Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-FILTER	Template: DEFAULT	(SOP-PCR-4.0) Preparation of Samples by Filtration and held for qPCR Analysis

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template:	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L		Number of Bottles: 1	Preserved With: ICE
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices	
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm	
Color (true)			
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)	
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices	
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices	
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample	
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects	
Calcium			
Chromium			
Iron			
Magnesium			
Nickel			
Potassium			
Sodium			
Zinc			
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects	
Arsenic			
Cadmium			
Copper			
Lead			
Silver			
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L	
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L	
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L	
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices	
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices	
Bottle Type: P-125ML		Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L	

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.

Do Not Lift Using This Tag

ORIGIN: FLORIDA DEPT

ORIGIN ID: TLHA (904) 807-3300
JEREMY PARRISH
FL DEP NORTHEAST DISTRICT
8800 BAYMEADOWS WAY WEST
SUITE 100
JACKSONVILLE, FL 32256
UNITED STATES US

SHIP DATE: 16JUL18
ACTWGT: 30.00 LB MAN
CAD: 0446970/CAFE3210

TO **JOHN WATTS**
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD

551C2/8532/104C

TALLAHASSEE FL 32399

(850) 245-8077


REF:

DEPT:

INU:

PO:

RMA: ||| ||| |||

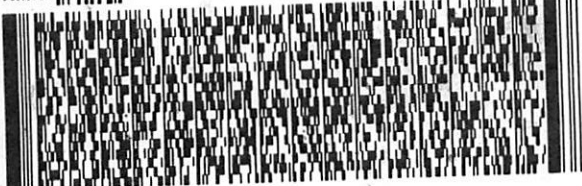

FedEx
Express



J1811180420010V

RETURNS MON-FRI
PRIORITY OVERNIGHT

TRK# 4385 5031 4431
0221

32399

FL - IIS

THU - 02 AUG 10:30A
PRIORITY OVERNIGHT

32399

FL-US

TLH

FedEx
TRK# 4385 5031 4431
0221

36 TLHA



FAD 936361 01AUG18 NRBA 546C1/3309/0C8A

Do Not Lift Using This Tag

ORIGIN: PRIORITY OVERNIGHT

ORIGIN ID: TLHA (904) 807-3300
JEREMY PARRISH
FL DEP NORTHEAST DISTRICT
8800 BAYMEADOWS WAY WEST
SUITE 100
JACKSONVILLE, FL 32256
UNITED STATES US

SHIP DATE: 16 JUL 18
ACTWGT: 30.00 LB MAN
CAD: 0446970/CAFE3210

TO **JOHN WATTS**
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD

551C2/8532/104C


TALLAHASSEE FL 32399

(850) 245-8077

REF:

DEPT:

RMA: ||| ||| |||

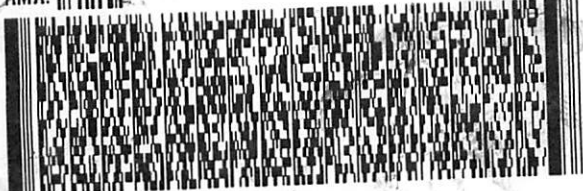

FedEx
Express



J1811180420014V

TRK# 4385 5031 4410
0221

RETURNS MON-FRI
PRIORITY OVERNIGHT

THU - 02 AUG 10:30A
PRIORITY OVERNIGHT

32399

JS

32399
FL-US
TLH



TLHA



061 01AUG18 NRBA 546C1/3309/0CBA

Cooler Packing Worksheet for Request: RQ-2018-07-30-20

LSJR Downtown Fresh

Ship Cooler On: 7/20/2018

Customer/Project: NE-DIST/SMP

Assigned To: GREENWOO_J

Requester: Jeremy M. Parrish

Priority: 3

904-256-1700
8800 Baymeadows Way West
Suite100
Jacksonville, FL 32256-7577
Attn: Jeremy Parrish

Please return
packing list with
sample submittal
forms.

Comments:

No preservatives needed.

Group A: Freshwater Kit, No Metals

Group B: Freshwater Kit With Metals

Requested Analyses:

Bottle Group: A

of Sites: 2

Container ID: P-1L

Qty: 2 Preservation: ICE

Lot # 00071006

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

W-CL-IC

W-COLOR

W-F

W-SO4-IC

W-TDS

W-TSS

Description

Alkalinity in aqueous matrices as mg CaCO₃/L

Turbidity in aqueous matrices

Chloride in aqueous matrices

True Color measured at 450 nm

Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Sulfate in aqueous matrices

Total Dissolved Solids in aqueous matrices

Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 2 Preservation: ICE

Lot # 1223624

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 2 Preservation: ICE

Lot # L18005

Description:

Analysis

NE-ALS-ECQ

Description

Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville

Container ID: QPCR-P-500ML

Qty: 2 Preservation: ICE

Lot # 55070587

Description: 500ml Sterile Nalgene

Affix PCR Instruction Label to container

Analysis

PCR-FILTER

Description

Preparation of Samples by Filtration and held for qPCR Analysis

Container ID: BG-500ML

Qty: 2 Preservation: ICE

Lot # 55071432

Description: Brown Glass Bottle 500ML

Affix BLUE dot to container.

Analysis

W-E8321-DI

W-E8321-MS

Description

Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Container ID: P-500ML

Qty: 2 Preservation: ICE And H2SO4

Lot # 55070684

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

W-NO2NO3

W-S-A-TP

W-TKN

W-TOC

Description

Ammonia in aqueous matrices as mg N/L

Nitrite/Nitrate in aqueous matrices as mg N/L

Total Phosphorus in aqueous matrices as mg P/L

Total Kjeldahl Nitrogen in aqueous matrices

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 2 Preservation: FILTER-ICE

Lot # 55070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: B

of Sites: 2

Container ID: P-1L

Qty: 2 Preservation: ICE

Lot # 0071006

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

W-CL-IC

W-COLOR

W-F

W-SO4-IC

W-TDS

W-TSS

Description

Alkalinity in aqueous matrices as mg CaCO3/L

Turbidity in aqueous matrices

Chloride in aqueous matrices

True Color measured at 450 nm

Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Sulfate in aqueous matrices

Total Dissolved Solids in aqueous matrices

Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 2 Preservation: ICE

Lot # 023623

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 2 Preservation: ICE

Lot # 418005

Description:

Analysis

NE-ALS-ECQ

Description

Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville

Container ID: P-500ML

Qty: 2 Preservation: HNO3

Lot # 0070684

Description: Plastic Bottle 500 mL

Affix RED nitric acid dot to container.

Analysis

W-ICP

W-ICPMS

Description

Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Container ID: P-500ML

Qty: 2 Preservation: ICE And H2SO4

Lot # 0070684

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

W-NO2NO3

Description

Ammonia in aqueous matrices as mg N/L

Nitrite/Nitrate in aqueous matrices as mg N/L

W-S-A-TP

Total Phosphorus in aqueous matrices as mg P/L

W-TKN

Total Kjeldahl Nitrogen in aqueous matrices

W-TOC

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 2 Preservation: FILTER-ICE

Lot # 00070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: C

of Sites: 1

Container ID: P-250ML-WI-THI

Qty: 1 Preservation: ICE

Lot # L16005

Description:

Analysis

NE-ALS-ECQ

Description

Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville

Container ID: QPCR-P-500ML

Qty: 1 Preservation: ICE

Lot # 00070587

Description: 500ml Sterile Nalgene

Affix PCR Instruction Label to container

Analysis

PCR-FILTER

Description

Preparation of Samples by Filtration and held for qPCR Analysis

Container ID: BG-500ML

Qty: 1 Preservation: ICE

Lot # 00071432

Description: Brown Glass Bottle 500ML

Affix BLUE dot to container.

Analysis

W-E8321-DI

Description

Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

W-E8321-MS

Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Group: D

of Sites: 1

Container ID: P-1L

Qty: 1 Preservation: ICE

Lot # 00071006

Description: Plastic Bottle 1L

Analysis

ALKALINITY

DescriptionAlkalinity in aqueous matrices as mg CaCO₃/L

TURBIDITY	Turbidity in aqueous matrices
W-CL-IC	Chloride in aqueous matrices
W-COLOR	True Color measured at 450 nm
W-F	Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Sulfate in aqueous matrices
W-TDS	Total Dissolved Solids in aqueous matrices
W-TSS	Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: P-500ML

Qty: 1 Preservation: HNO3

Lot # 00070684

Description: Plastic Bottle 500 mL

Affix RED nitric acid dot to container.

Analysis

W-ICP

W-ICPMS

Description

Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Container ID: P-500ML

Qty: 1 Preservation: ICE And H2SO4

Lot # 00070684

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

W-NO2NO3

W-S-A-TP

W-TKN

W-TOC

Description

Ammonia in aqueous matrices as mg N/L

Nitrite/Nitrate in aqueous matrices as mg N/L

Total Phosphorus in aqueous matrices as mg P/L

Total Kjeldahl Nitrogen in aqueous matrices

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 1 Preservation: FILTER-ICE

Lot # 00070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Cooler Packed By:

ARJ

Number of Coolers:

3

Date:

7-16-18

Kit must also include:

- ☒ Field Sheets
- ☒ Temperature Control Bottle
- ☒ FedEx Bills, if applicable
- ☒ Plastic Bags/Zip Ties
- ☒ Acid Labels (sets of 24)

If Preservation Included:

ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____

PLEASE RETURN ALL COOLERS FOR RQ-2018-07-30-20