

Log-in Checklist

RQ- 2018-09-10-16

Login Station ID: #4

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 9/27/18 9:30

Cooler Temperature*	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # (fill out only if waybill copy is damaged)
	HNO ₃	Formalin		Present?		*Intact?		
				Yes	No	Yes	No	
<u>0.6°C</u>			<u>5</u>		☒			

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form.

****Note:** If "No" is checked in any of the fields below, fill out the back of this form with additional details (Field ID, Test IDs, etc.) and generate an NCR.

**Chain of Custody/ Field Sheet(s) Included? Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☐ No ☒ Date

CONTAINER CHECK

****Evidence Tape on Bottles?** Yes ☐ No ☒ If yes, is it intact? Yes ☐ No ☐

****Caps on tight?** Yes ☒ No ☐ ****Containers intact?** Yes ☒ No ☐

****Sufficient Sample Volume:** Yes ☒ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes ☐ No ☒ Init: J.O.

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.					
Analysis	Yes	No	Analysis	Yes	No
W-1,4 DIOX			W-PC-AA-SF		
W-E8321-DI, -MS, -21			W-NP-AA-SF		
W-BNA (-625, -R)			W-PDQUAT		
W-CARB-AA (-CAB-AA-R)			W-PCL-SQ-R		
W-U2060-MS			W-TR-TQ		
W-PCB-R			W-CT-CI-R		
W-PAH (-R)			W-PSNP-TQ		

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

****Acid Preserved Samples:** pH ≤ 2.0? Yes ☒ No ☐ NA ☐ Init: J.O.

****W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?** Yes ☐ No ☒ NA ☐ Init: J.O.

Coolers Unpacked/Checked by: J.O. Date: 9/27/18 NCRs (Y/N)?

Event ID: 50-015T-2018-09-27-02 NCR #

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_04_2018
pH Paper 0.0 – 3.0 Lot #	230515 EXP: 10/30/2018
pH Paper 0 – 13 Lot #	222516 EXP: 08/15/2019
Chlorine Test Strips Lot #	7226 EXP: 07/2020

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ✓

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Tellywog

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: SO-DIST

Requester: Nicole Reistad

Field Report Prepared By: Gary Shord

Project ID: SMP

Collected By: Nicole ReistadSampling Agency: 8052 SKAC

Location	Field ID/WIN ID		☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>9/26/18</u> Time <u>1240</u>	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s) (See pg 2)	
Field ID	Location	<u>G3SD0100</u>	Matrix (type, e.g., Salt, Fresh, etc) <u>fresh</u>	Diss. Oxygen <u>3.43</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L) <u>—</u>	<u>B</u>	
WIN/STORE	Telegraph Swamp South		Temp (C) <u>27.34</u>	pH <u>6.84</u>	Sample Depth <u>0.3</u> ☐ ft ☒ m	Sp. Conductance (umho/cm) <u>376</u>		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.18</u>	Comments <u>coli sent to Sanders</u>			
Location	Field ID		☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>9/26/18</u> Time <u>1245</u>	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID	Location	<u>G3SD0100_DUP</u>	Matrix (type, e.g., Salt, Fresh, etc) <u>fresh</u>	Diss. Oxygen <u>3.43</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L) <u>—</u>	<u>BT C</u>	
WIN/STORE	WIN ID		Temp (C) <u>27.34</u>	pH <u>6.84</u>	Sample Depth <u>0.3</u> ☐ ft ☒ m	Sp. Conductance (umho/cm) <u>376</u>		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.18</u>	Comments <u>coli sent to Sanders</u>			
Location			☐ Comp ☐ Grab	Collection (comp begin or grab) Date _____ Time _____	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID			Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		
WIN/STORE #			Temp (C)	pH	Sample Depth ☐ ft ☐ m	Sp. Conductance (umho/cm)		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location			☐ Comp ☐ Grab	Collection (comp begin or grab) Date _____ Time _____	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID			Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		
WIN/STORE #			Temp (C)	pH	Sample Depth ☐ ft ☐ m	Sp. Conductance (umho/cm)		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			

Relinquished By: <u>Gary Shord</u>	Date/Time: <u>9/26/18 1430</u>	Shipping Method: <u>fedx</u>	Number of Coolers Shipped: <u>1</u>	Received By: <u>JS</u>	Date/Time: <u>9/27/18 9:30</u>
------------------------------------	--------------------------------	------------------------------	-------------------------------------	------------------------	--------------------------------

Group: A	Bottle Type: P-1L	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: A	Bottle Type: BP-1L	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
CHLSUITE-W					
Group: A	Bottle Type: P-250ML-WI-THI	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
SD-SND-ECQ					
Group: A	Bottle Type: P-500ML	# of Bottles: 2	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	0
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type: P-125ML	# of Bottles: 2	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	0
W-PO4-F					
Group: B	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: B	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0/1 1
CHLSUITE-W					
Group: B	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0/1 1

Group: B	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
SD-SND-ECQ <i>Ecoli sent to Sanders</i>					
Group: B	Bottle Type: P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab	1
W-ICP W-ICPMS <i>Lot# : NA8080070</i> <i>Exp. 03/21/19</i>					
Group: B	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC H2SO4 LOT# SA8101020 EXP: 04/11/2019					
Group: B	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F					
Group: C	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
SD-SND-ECQ <i>Ecoli sent to Sanders</i>					

Date of Request: 20-JUL-2018
 Created By: REISTAD_N On: 20-JUL-2018
 Modified By: WATTS_J On: 23-AUG-2018
 Customer: SO-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: South District
 Sampling Event: Tellywog

Send Coolers To:

Phone: 239-344-5718
 FLDEP South District Office
 2295 Victoria Ave. Suite 364
 Fort Myers, FL 33901-3381
 Attn: Nicole Reistad

Program Module Number:

Priority: 3

Event Status: S

Criminal Investigation: NO

Chemistry Request Reviewed By: WATTS_J On: 23-AUG-2018

Biology Request Reviewed By: SWANSON_C On: 21-AUG-2018

Sampling Kit Required: YES Ship on: 31-AUG-2018

Sampling Kit Shipped: YES On: 28-AUG-2018

Sampling Kit Packed By: GREENWOO_J

Preservatives Needed: NO

Date to Receive Samples: 10-SEP-2018

Received By:

Send Final Report To:

FLDEP South District Office
 2295 Victoria Ave. Suite 364
 Fort Myers, FL 33901-3381
 Attn: Nicole Reistad

Report Type: Final Only

FTP Data: NO

QC Report: YES

Date Log: NO

Authorization Log: NO

Comments: Group A: Freshwater Kit, No Metals (Kit A); E.coli x2
 (2 WBIDs)
 Group B: Freshwater Kit With Metals (Kit C); E.coli x1
 Group C: E.coli DUPLICATE x1

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 2	Preserved With: ICE
SD-SND-ECQ	Template: DEFAULT	(Colilert-18) Escherichia coli by Colilert 18 Quanti-Tray method by Sanders Laboratories in Fort Myers
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO ₄ -F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SD-SND-ECQ	Template: DEFAULT	(Colilert-18) Escherichia coli by Colilert 18 Quanti-Tray method by Sanders Laboratories in Fort Myers
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Chromium		
Iron		
Magnesium		
Nickel		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Arsenic		
Cadmium		
Copper		
Lead		
Silver		

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SD-SND-ECQ	Template: DEFAULT	(Colilert-18) Escherichia coli by Colilert 18 Quanti-Tray method by Sanders Laboratories in Fort Myers

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.

Do Not Lift Using This Tag

ORIGIN ID: TLHA (239) 344-5718
NICOLE REISTAD
FL DEP SOUTH DISTRICT OFFICE
2295 VICTORIA AVE SUITE 364

SHIP DATE: 28AUG18
ACTWGT: 15.00 LB MAN
CAD: 0446970/CAFE3210

FORT MYERS, FL 33901
UNITED STATES US

JOHN WATTS
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD

TALLAHASSEE FL 32399

(860) 246-8077


REF:

DEPT:

INU:

PO:

RMA: ||| ||| |||

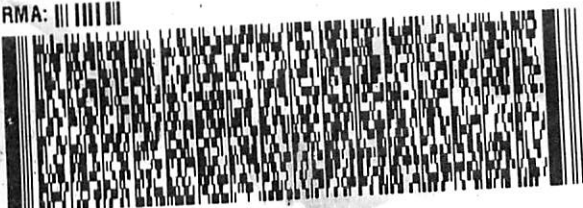

FedEx
Express



J161118042001BY

FedEx

TRK# 4385 5032 9651
0221

XH TLHA


THU - 27 SEP 10:30A
PRIORITY OVERNIGHT

32399

FL-US

TLH

EXP 04/19

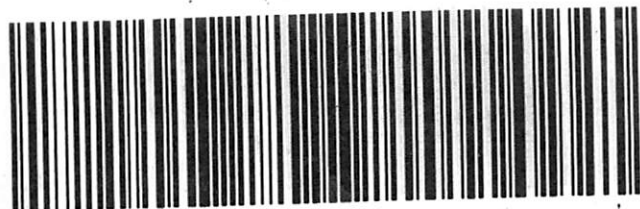

*286779 09/26 552J1/F78C/DCA5

Cooler Packing Worksheet for Request: RQ-2018-09-10-16

Tellywog

Ship Cooler On: 8/31/2018

Customer/Project: SO-DIST/SMP

Assigned To: GREENWOO_J

Requester: Nicole Reistad

Priority: 3

239-344-5718
FLDEP South District Office
2295 Victoria Ave. Suite 364
Fort Myers, FL 33901-3381
Attn: Nicole Reistad

Comments:

No preservatives needed.

~~Group A: Freshwater Kit, No Metals (Kit A); E.coli x2
(2 WBIDs)~~

Group B: Freshwater Kit With Metals (Kit C); E.coli x1

Group C: E.coli DUPLICATE x1

TELEGRAPH & E. coli Duplicate
ONLY NOW

Requested Analyses:

Bottle Group: A

of Sites: 2

Container ID: P-1L

Qty: 2 Preservation: ICE

Lot # 00571991

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

W-CL-IC

W-COLOR

W-F

W-SO4-IC

W-TDS

W-TSS

Description

Alkalinity in aqueous matrices as mg CaCO3/L

Turbidity in aqueous matrices

Chloride in aqueous matrices

True Color measured at 450 nm

Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Sulfate in aqueous matrices

Total Dissolved Solids in aqueous matrices

Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 2 Preservation: ICE

Lot # 1232362

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 2 Preservation: ICE

Lot # 418005

Description: 250 ml Sterilized with tablet

Analysis

SD-SND-ECQ

Description

Escherichia coli by Colilert 18 Quanti-Tray method by Sanders Laboratories in Fort Myers

Container ID: P-500ML

Qty: 2 Preservation: ICE And H2SO4

Lot # 00072044

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

W-NO2NO3

W-S-A-TP

W-TKN

W-TOC

Description

Ammonia in aqueous matrices as mg N/L

Nitrite/Nitrate in aqueous matrices as mg N/L

Total Phosphorus in aqueous matrices as mg P/L

Total Kjeldahl Nitrogen in aqueous matrices

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 2 Preservation: FILTER-ICE

Lot # 00070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: B

of Sites: 1

Container ID: P-1L

Qty: 1 Preservation: ICE

Lot # 00071991

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

W-CL-IC

W-COLOR

W-F

W-SO4-IC

W-TDS

W-TSS

Description

Alkalinity in aqueous matrices as mg CaCO3/L

Turbidity in aqueous matrices

Chloride in aqueous matrices

True Color measured at 450 nm

Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Sulfate in aqueous matrices

Total Dissolved Solids in aqueous matrices

Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 1 Preservation: ICE

Lot # 1232362

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 1 Preservation: ICE

Lot # 418005

Description: 250 ml Sterilized with tablet

Analysis

SD-SND-ECQ

Description

Escherichia coli by Colilert 18 Quanti-Tray method by Sanders Laboratories in Fort Myers

Container ID: P-500ML

Qty: 1 Preservation: HNO3

Lot # 00072044

Description: Plastic Bottle 500 mL

Affix RED nitric acid dot to container.

Analysis

W-ICP

W-ICPMS

Description

Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Container ID: P-500ML

Qty: 1 Preservation: ICE And H2SO4

Lot # 00072044

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

W-NO2NO3

W-S-A-TP

W-TKN

W-TOC

Description

Ammonia in aqueous matrices as mg N/L

Nitrite/Nitrate in aqueous matrices as mg N/L

Total Phosphorus in aqueous matrices as mg P/L

Total Kjeldahl Nitrogen in aqueous matrices

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 1 Preservation: FILTER-ICE

Lot # 00070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: C

of Sites: 1

Container ID: P-250ML-WI-THI

Qty: 1 Preservation: ICE

Lot # L18405

Description: 250 ml Sterilized with tablet

Analysis

SD-SND-ECQ

Description

Escherichia coli by Colilert 18 Quanti-Tray method by Sanders Laboratories in Fort Myers

Cooler Packed By: JE

Number of Coolers: 1

Date: 8/28/18

Kit must also include:

- ☒ Field Sheets
- ☒ Temperature Control Bottle
- ☒ FedEx Bills, if applicable
- ☒ Plastic Bags/Zip Ties
- ☐ Acid Labels (sets of 24)

If Preservation Included:

ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____

PLEASE RETURN ALL COOLERS FOR RQ-2018-09-10-16