

Log-in Checklist

RQ- 2018-03-19-15

Shipping Method: FedEx | UPS | HD | Greyhound | _____

Date/Time of Receipt: 3/23/18 9:11

Cooler ID	Cooler Temperature*	No. Sample Containers in Cooler	Evidence Tape				Tracking # (fill out only if waybill copy is damaged)
			Present?		*Intact?		
			Yes	No	Yes	No	
RED	3.2°C	14		✓			
RED	3.9°C	10		✓			

*If the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX), or if Evidence Seal is damaged, identify the containers in the affected cooler(s) on back of this form.

****Note:** If "No" is checked in any of the fields below, fill out the back of this form with additional details (Field ID, Test IDs, etc.) and generate an NCR.

**Chain of Custody/ Field Sheet(s) Included? Yes ☒ No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes _____ No ☒

CONTAINER CHECK

**Evidence Tape on Bottles? Yes _____ No ☒ If yes, is it intact? Yes _____ No _____

**Caps on tight? Yes ☒ No _____ **Containers intact? Yes ☒ No _____

**Sufficient Sample Volume: Yes _____ No _____

CHLORINE CHECK (Blue dot on container)

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.					
Analysis	Yes	No	Analysis	Yes	No
W-1,4 DIOX		☒	W-PC-AA-SF		
W-E8321-DI, -MS, -21		☒	W-NP-AA-SF		
W-BNA (-625, -R)		☒	W-PDQUAT		
W-CARB-AA (-CAB-AA-R)		☒	W-PCL-SQ-R		☒
W-U2060-MS			W-TR-TQ		☒
W-PCB-R			W-CT-CI-R		☒
W-PAH (-R)			W-PSNP-TQ		☒

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ☒ No _____ NA _____ Init: JD

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes _____ No ☒ NA ☒ Init: JD

Coolers Unpacked/Checked by: JD Date: 3/23/18 NCRs (Y/N)? N

Event ID: AE-NW-DIST-201803-23-01 NCR # _____

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ✓

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF; W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Additional Comments:

Central Laboratory Submittal Form

Customer: NW-DIST

Requester: Michael King

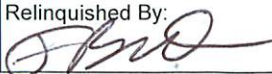
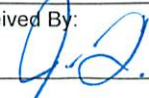
Field Report Prepared By:

Project ID: SMP

Collected By:

Sampling Agency:

Location	Field/WIN ID		G3NW0148	Location:	Black Creek Upstream of Hwy 20	☐ Comp ☒ Grab Collection (comp begin or grab) Date 3/22/2018 Time 9:54 Eastern Central	Composite end Date Time	Bottle Group(s) (See pg 2)			
Field ID	WIN/STORE			Matrix (type, e.g., Salt, Fresh, etc)	FRESH	Diss. Oxygen	9.43	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		
Latitude	dd	dms	dd	dms	Temp (C)	13.39	pH	4.18	Sample Depth 0.3 ☐ ft ☒ m	Sp. Conductance (umho/cm)	24
				Salinity (ppt)	0.01	Comments					
Location	Field/WIN ID		G3NW0133	Location:	Otter Creek upstream of Otter Creek Bridge Rd	☐ Comp ☒ Grab Collection (comp begin or grab) Date 3/22/18 Time 12:35 Eastern Central	Composite end Date Time	Bottle Group(s)			
Field ID	WIN/STORE			Matrix (type, e.g., Salt, Fresh, etc)	FRESH	Diss. Oxygen	8.32	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		
Latitude	dd	dms	dd	dms	Temp (C)	17.78	pH	6.14	Sample Depth 0.3 ☐ ft ☒ m	Sp. Conductance (umho/cm)	136
				Salinity (ppt)	0.06	Comments			NO SCI. Bact-Tis		
Location	Field ID	WIN/STORE #	Latitude	dd	dms	Longitude	dd	dms	Salinity (ppt)	Comments	
Location	Field ID	WIN/STORE #	Latitude	dd	dms	Longitude	dd	dms	Salinity (ppt)	Comments	
Location	Field ID	WIN/STORE #	Latitude	dd	dms	Longitude	dd	dms	Salinity (ppt)	Comments	

Relinquished By:	Date/Time	Shipping Method	Number of Coolers Shipped:	Received By:	Date/Time
	4/3/2018	Home & Deliver	2		3/23/18 9:11

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	CHLSUITE-W						
Group: A	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ECOLI-18QT						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	CH2CICOOH Pre-Preserved	Enter Number of Bottles Sent to Lab	1
	W-CARB-AA						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab	4
	W-CT-CI-R						
	W-PCL-SQ-R						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	W-E8321-DI						
	W-E8321-MS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-ICP						
	W-ICPMS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
	W-PO4-F						

Group: A	Bottle Type:	BG-1L	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab	4
W-PSNP-TQ							
Group: B	Bottle Type:	PT-50ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
ALGAL_ID							
Group: B	Bottle Type:	PJ-2L	# of Bottles: 2	Preservative:	Formalin	Enter Number of Bottles Sent to Lab	2
MI-FW-QLDC							
Group: C	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS							
Group: C	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W							
Group: C	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
ECOLI-18QT							
Group: C	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
W-ICP W-ICPMS							
Group: C	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC							
Group: C	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F							
Group: D	Bottle Type:	PT-50ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
ALGAL_ID							

Group: D	Bottle Type: MI-FW-QLDC	PJ-2L	# of Bottles: 2	Preservative: Formalin	Enter Number of Bottles Sent to Lab	0
Group: E	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
Group: E	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
Group: E	Bottle Type: NW-UWF-ECQ	TEST HAS NO BOTTLE	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
Group: E	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	0
Group: E	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	0
Group: F	Bottle Type: ALGAL_ID	PT-50ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
Group: F	Bottle Type: MI-FW-QLDC	PJ-2L	# of Bottles: 2	Preservative: Formalin	Enter Number of Bottles Sent to Lab	0
Group: G	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0

Date of Request: 15-FEB-2018
 Created By: KING_M On: 15-FEB-2018
 Modified By: WATTS_J On: 16-MAR-2018
 Customer: NW-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Northwest District
 Sampling Event: BIOASSESSMENT

Send Coolers To:

Phone: 850-595-0605
 160 Government Center
 Suite 208
 Pensacola, FL 32502-5794
 Attn: Michael King

Program Module Number: 1306
 Priority: 3
 Event Status: A
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 16-MAR-2018
 Biology Request Reviewed By: SWANSON_C On: 16-MAR-2018
 Sampling Kit Required: NO Ship on:
 Sampling Kit Shipped: YES On: 12-MAR-2018
 Sampling Kit Packed By: REYNOLDS_T
 Preservatives Needed: NO
 Date to Receive Samples: 19-MAR-2018
 Received By:

Send Final Report To:

160 Government Center
 Suite 208
 Pensacola, FL 32502-5794
 Attn: Cheryl Bunch

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Group A: (WBID 679) Freshwater Kit With Metals, Pesticides, Contract Lab Bacteria
 Group B: Biology
 Group C: (WBID 950) Freshwater Kit With Metals, Contract Lab Bacteria
 Group D: Biology
 Group E: (WBID 49A) Freshwater Kit, Contract Lab Bacteria
 Group F: Biology
 Group G: (WBID 951) Contract Lab Bacteria
 Group H: Biology

** Please ship 2 extra fedex express return labels and handle tags

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices
Bottle Type: BG-1L	Number of Bottles: 4	Preserved With: ICE
W-PSNP-TQ	Template: DEFAULT	(EPA 8270D) Organonitrogen and phosphorus pesticides (ultra trace level) in water matrices by GC/MS/MS
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE

Bottle Group A (Water) With 1 Samples:

Bottle Type: BP-1L CHLSUITE-W	Number of Bottles: 1 Template: DEFAULT	Preserved With: ICE (SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: ECOLI-18QT	Number of Bottles: 1 Template: DEFAULT	Preserved With: ICE (SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: BG-500ML W-CARB-AA	Number of Bottles: 1 Template: DEFAULT	Preserved With: CH2CICOOH Pre-Preserved (EPA 8321B) Carbamates in water matrices by HPLC/MS/MS
Bottle Type: BG-500ML W-CT-CI-R	Number of Bottles: 4 Template: DEFAULT	Preserved With: ICE (EPA 8276) Organochlorine pesticides in water matrices by GC/MS with Negative Chemical Ionization
W-PCL-SQ-R	Template: DEFAULT	(EPA 8270D) Organochlorine pesticides in water matrices by GC/MS-SIM
Bottle Type: BG-500ML W-E8321-DI	Number of Bottles: 1 Template: DEFAULT	Preserved With: ICE (EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: ADD-PYRACL	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML W-ICP	Number of Bottles: 1 Template:	Preserved With: HNO3 (EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium Chromium Iron Magnesium Nickel Potassium Sodium Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Arsenic Cadmium Copper Lead Silver		
Bottle Type: P-500ML W-NH3	Number of Bottles: 1 Template: DEFAULT	Preserved With: ICE And H2SO4 (EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML W-PO4-F	Number of Bottles: 1 Template: DEFAULT	Preserved With: FILTER-ICE (EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Biological) With 1 Samples:

Bottle Type: PT-50ML ALGAL_ID	Number of Bottles: 1 Template: DEFAULT	Preserved With: ICE (SOP-AB05) Qualitative assessment of algal taxa present without counts
Bottle Type: PJ-2L MI-FW-QLDC	Number of Bottles: 2 Template: DEFAULT	Preserved With: Formalin (SOP-IZ06) No. Taxa of FW macroinverts, qual samp, 20 dipnets.

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L ALKALINITY	Number of Bottles: 1 Template: DEFAULT	Preserved With: ICE (SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR Color (true)	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices
Bottle Type: BP-1L CHLSUITE-W	Number of Bottles: 1 Template: DEFAULT	Preserved With: ICE (SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Group C (Water) With 1 Samples:

Bottle Type:	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Chromium		
Iron		
Magnesium		
Nickel		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Arsenic		
Cadmium		
Copper		
Lead		
Silver		
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group D (Biological) With 1 Samples:

Bottle Type: PT-50ML	Number of Bottles: 1	Preserved With: ICE
ALGAL_ID	Template: DEFAULT	(SOP-AB05) Qualitative assessment of algal taxa present without counts
Bottle Type: PJ-2L	Number of Bottles: 2	Preserved With: Formalin
MI-FW-QLDC	Template: DEFAULT	(SOP-IZ06) No. Taxa of FW macroinverts, qual samp, 20 dipnets.

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: TEST HAS NO	Number of Bottles: 1	Preserved With: ICE
NW-UWF-ECQ	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method by UWF Overflow Laboratory for NWD
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group F (Biological) With 1 Samples:

Bottle Type: PT-50ML	Number of Bottles: 1	Preserved With: ICE
ALGAL_ID	Template: DEFAULT	(SOP-AB05) Qualitative assessment of algal taxa present without counts
Bottle Type: PJ-2L	Number of Bottles: 2	Preserved With: Formalin
MI-FW-QLDC	Template: DEFAULT	(SOP-IZ06) No. Taxa of FW macroinverts, qual samp, 20 dipnets.

Bottle Group G (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: TEST HAS NO	Number of Bottles: 1	Preserved With: ICE
NW-UWF-ECQ	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Collert 18 Quanti-Tray method by UWF Overflow Laboratory for NWD
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group H (Biological) With 1 Samples:

Bottle Type: PT-50ML	Number of Bottles: 1	Preserved With: ICE
ALGAL_ID	Template: DEFAULT	(SOP-AB05) Qualitative assessment of algal taxa present without counts
Bottle Type: PJ-2L	Number of Bottles: 2	Preserved With: Formalin
MI-FW-QLDC	Template: DEFAULT	(SOP-IZ06) No. Taxa of FW macroinverts, qual samp, 20 dipnets.

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.