

RQ 2018-09-03-53

Log-in Checklist

Login Station ID: #2

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 9-6-18/9:40

Cooler Temperature*	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # (fill out only if waybill copy is damaged)
	HNO ₃	Formalin		Present?		*Intact?		
				Yes	No	Yes	No	
0.4			8		✓			
0.8			10		✓			
0.6			14		✓			

* HNO₃ and Formalin preserved samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form.

****Note:** If "No" is checked in any of the fields below, fill out the back of this form with additional details (Field ID, Test IDs, etc.) and generate an NCR.

****Chain of Custody/ Field Sheet(s) Included?** Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☐ No ☒ Date

CONTAINER CHECK

****Evidence Tape on Bottles?** Yes ☐ No ☒ If yes, is it intact? Yes ☒ No ☐

****Caps on tight?** Yes ☒ No ☐ ****Containers intact?** Yes ☒ No ☐

****Sufficient Sample Volume:** Yes ☒ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes ☒ No ☐ Init: SR

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.						
Analysis	Yes	No	Analysis	Yes	No	
W-1,4 DIOX			W-PC-AA-SF			
W-E8321-DI, -MS, -21			W-NP-AA-SF			
W-BNA (-625, -R)			W-PDQUAT			
W-CARB-AA (-CAB-AA-R)			W-PCL-SQ-R			
W-U2060-MS			W-TR-TQ			
W-PCB-R			W-CT-CI-R			
W-PAH (-R)			W-PSNP-TQ			

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

****Acid Preserved Samples: pH ≤ 2.0?** Yes ☒ No ☐ NA ☐ Init: SR

****W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?** Yes ☐ No ☒ NA ☒ Init: SR

Coolers Unpacked/Checked by: SR **Date:** 9-6-18

NCRs (Y/N)? N

Event ID: SEDIST-2018-09-06-01

NCR #

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_02_2018
pH Paper 0.0 – 3.0 Lot #	230515 EXP:10-30-18
pH Paper 0 – 13 Lot #	222516 EXP:08-15-19
Chlorine Test Strips Lot #	7226 EXP: 07/2020

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ✓

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Customer: SE-DIST

Requester: Matt Sewell

Field Report Prepared By: Drew Liddick

Project IN. SMP

Collected By: Drew Liddick, Kristene Bonilla

Sampling Agency: 21FLWPB

Location	Location: Okeehetee Park Ramp		☐ Comp ☒ Grab		Collection (comp begin or grab)	Eastern Central	Composite end	Bottle Group(s)		
Field ID	WIN ID/Field ID:  G3SE0015		Date 09/05/2018		Time 1140		Date	Time	(See pg 2)	
WIN/STOI			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		A	
		Temp (C) 29.55		pH 7.77	Sample Depth 0.3	☐ ft ☒ m	Sp. Conductance (umho/cm) 170			
Latitude	26.655531	☒ dd ☐ dms	Longitude	-80.166361	☒ dd ☐ dms	Salinity (ppt) 0.06				Comments
Location	Location: Lake Clarke @ Pine Tree		☐ Comp ☒ Grab		Collection (comp begin or grab)	Eastern Central	Composite end	Bottle Group(s)		
Field ID	WIN ID/Field ID:  G3SE0001		Date 09/05/2018		Time 1240		Date	Time	B	
WIN/STORET #			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)			
		Temp (C) 28.65		pH 7.52	Sample Depth 0.3	☐ ft ☒ m	Sp. Conductance (umho/cm) 536		Comments	
Latitude	26.6540583	☒ dd ☐ dms	Longitude	-80.072781	☒ dd ☐ dms	Salinity (ppt) 0.26				
Location	Pine Lake		☐ Comp ☒ Grab		Collection (comp begin or grab)	Eastern Central	Composite end	Bottle Group(s)		
Field ID	G3SE0003		Date 09/05/18		Time 1315		Date	Time	E	
WIN/STORET #			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)			
		Temp (C) 29.50		pH 7.88	Sample Depth 0.3	☐ ft ☒ m	Sp. Conductance (umho/cm) 378		Comments	
Latitude	26.67683	☐ dd ☐ dms	Longitude	-80.07211	☒ dd ☐ dms	Salinity (ppt) 0.18				
Location	Location: Clear Lake near Australian		☐ Comp ☒ Grab		Collection (comp begin or grab)	Eastern Central	Composite end	Bottle Group(s)		
Field ID	Field ID:  CLEAR LAKE		Date 09/05/18		Time 1425		Date	Time	C	
WIN/STORET #			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)			
		Temp (C) 29.55		pH 8.17	Sample Depth 0.3	☐ ft ☒ m	Sp. Conductance (umho/cm) 365		Comments	
Latitude	26.717218	☒ dd ☐ dms	Longitude	-80.067152	☒ dd ☐ dms	Salinity (ppt) 0.17				

Relinquished By: FDEP	Date/Time 09/05/18 16:30	Shipping Method Fedex	Number of Coolers Shipped: 3	Received By: 	Date/Time 9-6-18/9140
--------------------------	-----------------------------	--------------------------	---------------------------------	---	--------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	SE-FLR-ECQ					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP					
	W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3					
	W-NO2NO3					
	W-S-A-TP					
	W-TKN					
	W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					

SMP- South Lakes

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: SE-DIST

Requester: Matt Sewell

Field Report Prepared By:

Drew Liddick

Project ID: SMP

Collected By:

Drew Liddick, Kristene Bonilla

Sampling Agency:

21FLWPB

Location	Location: Lake Mangonia North		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 09/05/18 Time 1510	Eastern Central	Composite end Date Time	Bottle Group(s) (See pg 2)
Field ID	WIN ID/Field ID: 		Matrix (type, e.g., Salt, Fresh, etc) Fresh		Diss. Oxygen 11.72	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
WIN/STORET #	G3SE0030		Temp (C) 32.65	pH 9.01	Sample Depth 0.3	☐ ft ☒ m	Sp. Conductance (umho/cm) 308
Latitude 26.746780	☒ dd ☐ dms	Longitude -80.083728	☒ dd ☐ dms	Salinity (ppt) 0.14	Comments		
Location			☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
WIN/STORET #			Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments		
Location			☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
WIN/STORET #			Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments		
Location			☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
WIN/STORET #			Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments		

Relinquished By: FDEP	Date/Time 09/05/18 16:30	Shipping Method Fedex	Number of Coolers Shipped: 23	Received By: SR	Date/Time 9-6-18/9:40
--------------------------	-----------------------------	--------------------------	----------------------------------	--------------------	--------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	SE-FLR-ECQ					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					

Group: B	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: SE-FLR-ECQ	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-CARB-AA	BG-500ML	# of Bottles: 1	Preservative:	CH2CICOOH Pre-Preserved	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-CT-CI-R W-PCL-SQ-R	BG-500ML	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>4</u>
Group: B	Bottle Type: W-E8321-DI W-E8321-MS	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-PSNP-TQ	BG-1L	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>4</u>
Group: C	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>

Group: C	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
W-TDS					
W-TSS					
Group: C	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W					
Group: C	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
SE-FLR-ECQ					
Group: C	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3					
W-NO2NO3					
W-S-A-TP					
W-TKN					
W-TOC					
Group: C	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F					
Group: D	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY					
TURBIDITY					
W-CL-IC					
W-COLOR					
W-F					
W-SO4-IC					
W-TDS					
W-TSS					
Group: D	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W					
Group: D	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
SE-FLR-ECQ					
Group: D	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3					
W-NO2NO3					
W-S-A-TP					
W-TKN					

Group: B	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: SE-FLR-ECQ	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-CARB-AA	BG-500ML	# of Bottles: 1	Preservative:	CH ₂ ClCOOH Pre-Preserved	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-CT-CI-R W-PCL-SQ-R	BG-500ML	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>4</u>
Group: B	Bottle Type: W-E8321-DI W-E8321-MS	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative:	HNO ₃	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-NH ₃ W-NO ₂ NO ₃ W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative:	ICE And H ₂ SO ₄	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-PO ₄ -F	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: B	Bottle Type: W-PSNP-TQ	BG-1L	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>4</u>
Group: C	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO ₄ -IC	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u>1</u>

Group: C	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
W-TDS					
W-TSS					
Group: C	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W					
Group: C	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
SE-FLR-ECQ					
Group: C	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3					
W-NO2NO3					
W-S-A-TP					
W-TKN					
W-TOC					
Group: C	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F					
Group: D	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY					
TURBIDITY					
W-CL-IC					
W-COLOR					
W-F					
W-SO4-IC					
W-TDS					
W-TSS					
Group: D	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W					
Group: D	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
SE-FLR-ECQ					
Group: D	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3					
W-NO2NO3					
W-S-A-TP					
W-TKN					

Group: D Bottle Type: P-500ML # of Bottles: 1 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 1

W-TOC

Group: D Bottle Type: P-125ML # of Bottles: 1 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 1

W-PO4-F

Group: E Bottle Type: P-1L # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1

ALKALINITY
TURBIDITY
W-CL-IC
W-COLOR
W-F
W-SO4-IC
W-TDS
W-TSS

Group: E Bottle Type: BP-1L # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1

CHLSUITE-W

Group: E Bottle Type: P-250ML-WI-THI # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1

SE-FLR-ECQ

Group: E Bottle Type: P-500ML # of Bottles: 1 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 1

W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC

Group: E Bottle Type: P-125ML # of Bottles: 1 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 1

W-PO4-F

Group: D Bottle Type: P-500ML # of Bottles: 1 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 1
W-TOC

Group: D Bottle Type: P-125ML # of Bottles: 1 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 1
W-PO4-F

Group: E Bottle Type: P-1L # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1
ALKALINITY
TURBIDITY
W-CL-IC
W-COLOR
W-F
W-SO4-IC
W-TDS
W-TSS

Group: E Bottle Type: BP-1L # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1
CHLSUITE-W

Group: E Bottle Type: P-250ML-WI-THI # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1
SE-FLR-ECQ

Group: E Bottle Type: P-500ML # of Bottles: 1 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 1
W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC

Group: E Bottle Type: P-125ML # of Bottles: 1 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 1
W-PO4-F

Date of Request: 15-AUG-2018
 Created By: SEWELL_M On: 15-AUG-2018
 Modified By: SWANSON_C On: 17-AUG-2018
 Customer: SE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Southeast District
 Sampling Event: SMP- South Lakes

Send Coolers To:

Phone: (772) 448-5883
 FL Dept. of Environmental Protection
 2199 South Rock Road
 Fort Pierce, FL 34945
 Attn: Matthew Sewell

Send Final Report To:

FL Dept. of Environmental Protection
 2199 South Rock Road
 Fort Pierce, FL 34945
 Attn: Matthew Sewell

Program Module Number: 1306
 Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 15-AUG-2018
 Biology Request Reviewed By: SWANSON_C On: 17-AUG-2018
 Sampling Kit Required: YES Ship on: 24-AUG-2018
 Sampling Kit Shipped: YES On: 20-AUG-2018
 Sampling Kit Packed By: GREENWOO_J
 Preservatives Needed: NO
 Date to Receive Samples: 03-SEP-2018
 Received By:

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Sampler Comments:
 Bottle Group A = 3245A
 Bottle Group B = 3245B
 Bottle Group C = 3245C1
 Bottle Group D = 3245C2
 Bottle Group E = 3254C4

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SE-FLR-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Chromium		
Iron		
Magnesium		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Nickel		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Arsenic		
Cadmium		
Copper		
Lead		
Silver		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: CH2CICOOH Pre-Preserved
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS
Bottle Type: BG-500ML	Number of Bottles: 4	Preserved With: ICE
W-CT-Cl-R	Template: DEFAULT	(EPA 8276) Organochlorine pesticides in water matrices by GC/MS with Negative Chemical Ionization
W-PCL-SQ-R	Template: DEFAULT	(EPA 8270D) Organochlorine pesticides in water matrices by GC/MS-SIM
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Chromium		
Iron		
Magnesium		
Nickel		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Arsenic		
Cadmium		
Copper		
Lead		
Silver		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: BG-1L	Number of Bottles: 4	Preserved With: ICE
W-PSNP-TQ	Template: DEFAULT	(EPA 8270D) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SE-FLR-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SE-FLR-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO ₄ -F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SE-FLR-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SE-FLR-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.



Do Not Lift Using This Tag

ORIGIN ID: TLHA (772) 467-5574
MATTHEW SEWELL
FLORIDA DEP-SOUTHEAST DISTRICT
2195 SOUTH ROCK ROAD

SHIP DATE: 18AUG18
ACTWGT: 15.00 LB MAN
CAD: 0446970/CAFE3210

FORT PIERCE, FL 34945
UNITED STATES US

TO **JOHN WATTS**
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD

TALLAHASSEE FL 32399

(850) 245-8077

REF:

DEPT:

FedEx

TRK# 4385 5032 5510
0221

THU - 06 SEP 10:30A
PRIORITY OVERNIGHT

36 TLHA

32399

FL-US TL



#164016 09/05 552J1/F78C/DCAS

3047/3309/104C

61/10 09/07/18

FedEx

FedEx

ORIGIN ID: TLHA (772) 467-5574
MATTHEW SEWELL
FLORIDA DEP-SOUTHEAST DISTRICT
2199 SOUTH ROCK ROAD
FORT PIERCE, FL 34945
UNITED STATES US

SHIP DATE: 20AUG18
ACTWGT: 15.00 LB MAN
CAD: 0446970/CAFE3210

TO **JOHN WATTS**
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD

TALLAHASSEE FL 32399

(860) 245-8077
INU:
PO:

REF:

DEPT:

RMA: ||| ||| |||



FedEx
Express



FedEx
TRK# 4385 5032 5873
0221

THU - 06 SEP 10:30A
PRIORITY OVERNIGHT

36 TLHA

RI
IT
32399
TLH
FL-US
19



#164016 09/05 552J1/F78C/DCAS

(Using This Tag

ORIGIN ID: TLHA (772) 467-5574
MATTHEW SEWELL
FLORIDA DEP-SOUTHEAST DISTRICT
2199 SOUTH ROCK ROAD

SHIP DATE: 20AUG18
ACTWGT: 15.00 LB MAN
CAD: 0446970/CAFE3210

FORT PIERCE, FL 34945
UNITED STATES US

TO **JOHN WATTS**
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD

TALLAHASSEE FL 32399

(850) 245-8077


REF:

DEPT:

RMA: III IIII III



FedEx
Express

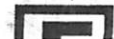

FedEx

TRK#
0221

4385 5032 5862

THU - 06 SEP 10:30A
PRIORITY OVERNIGHT

36 TLHA

32399

FL-US

TLH



#164016 09/05 552J1/F78C/DCAS

Cooler Packing Worksheet for Request: RQ-2018-09-03-53

SMP- South Lakes

Ship Cooler On: 8/24/2018

Customer/Project: SE-DIST/SMP

Assigned To: GREENWOO_J

Requester: Matt Sewell

Priority: 3

(772) 448-5883
FL Dept. of Environmental Protection
2199 South Rock Road
Fort Pierce, FL 34945

Attn: Matthew Sewell

Comments:

No preservatives needed.

Sampler Comments:

Bottle Group A = 3245A

Bottle Group B = 3245B

Bottle Group C = 3245C1

Bottle Group D = 3245C2

Bottle Group E = 3254C4

Requested Analyses:

Bottle Group: A

of Sites: 1

Container ID: P-1L

Qty: 1

Preservation: ICE

Lot # 00071878

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

W-CL-IC

W-COLOR

W-F

W-SO4-IC

W-TDS

W-TSS

Description

Alkalinity in aqueous matrices as mg CaCO3/L

Turbidity in aqueous matrices

Chloride in aqueous matrices

True Color measured at 450 nm

Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Sulfate in aqueous matrices

Total Dissolved Solids in aqueous matrices

Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 1

Preservation: ICE

Lot # 1225081

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI Qty: 1 Preservation: ICE Lot # LI8005
Description: 250 ml Sterilized with tablet

Analysis
SE-FLR-ECQ

Description
Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie

Container ID: P-500ML Qty: 1 Preservation: HNO3 Lot # 06071230
Description: Plastic Bottle 500 mL Affix RED nitric acid dot to container.

Analysis
W-ICP
W-ICPMS

Description
Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Container ID: P-500ML Qty: 1 Preservation: ICE And H2SO4 Lot # 00071230
Description: Plastic Bottle 500 mL Affix YELLOW sulfuric acid dot to container.

Analysis
W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC

Description
Ammonia in aqueous matrices as mg N/L
Nitrite/Nitrate in aqueous matrices as mg N/L
Total Phosphorus in aqueous matrices as mg P/L
Total Kjeldahl Nitrogen in aqueous matrices
Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML Qty: 1 Preservation: FILTER-ICE Lot # 00070886
Description: Plastic Bottle 125 mL

Analysis
W-PO4-F

Description
Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: B

of Sites: 1

Container ID: P-1L Qty: 1 Preservation: ICE Lot # 00071878
Description: Plastic Bottle 1L

Analysis
ALKALINITY
TURBIDITY

Description
Alkalinity in aqueous matrices as mg CaCO3/L
Turbidity in aqueous matrices

W-CL-IC	Chloride in aqueous matrices
W-COLOR	True Color measured at 450 nm
W-F	Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Sulfate in aqueous matrices
W-TDS	Total Dissolved Solids in aqueous matrices
W-TSS	Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 1 Preservation: ICE

Lot # 1215081

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 1 Preservation: ICE

Lot # 4805

Description: 250 ml Sterilized with tablet

Analysis

SE-FLR-ECQ

Description

Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie

Container ID: BG-500ML

Qty: 1 Preservation: CH2ClCOOH Pre-Preserved

Lot # 00071185

Description: Brown Glass Bottle 500ML

Affix BLUE dot to container.

Analysis

W-CARB-AA

Description

Carbamates in water matrices by HPLC/MS/MS

Container ID: BG-500ML

Qty: 4 Preservation: ICE

Lot # 00070815

Description: Brown Glass Bottle 500ML

Affix BLUE dot to container.

Analysis

W-CT-Cl-R

W-PCL-SQ-R

Description

Organochlorine pesticides in water matrices by GC/MS with Negative Chemical Ionization

Organochlorine pesticides in water matrices by GC/MS-SIM

Container ID: BG-500ML

Qty: 1 Preservation: ICE

Lot # 00070815

Description: Brown Glass Bottle 500ML

Affix BLUE dot to container.

Analysis

W-E8321-DI

W-E8321-MS

Description

Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Container ID: P-500ML

Qty: 1

Preservation: HNO3

Lot # 00071230

Description: Plastic Bottle 500 mL

Affix RED nitric acid dot to container.

Analysis

W-ICP

W-ICPMS

Description

Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Container ID: P-500ML

Qty: 1

Preservation: ICE And H2SO4

Lot # 00071230

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

W-NO2NO3

W-S-A-TP

W-TKN

W-TOC

Description

Ammonia in aqueous matrices as mg N/L

Nitrite/Nitrate in aqueous matrices as mg N/L

Total Phosphorus in aqueous matrices as mg P/L

Total Kjeldahl Nitrogen in aqueous matrices

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 1

Preservation: FILTER-ICE

Lot # 00070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Container ID: BG-1L

Qty: 4

Preservation: ICE

Lot # 00071178

Description: Brown Glass Bottle 1L

Affix BLUE dot to container.

Analysis

W-PSNP-TQ

Description

Organonitrogen and phosphorus pesticides (ultra trace level) in water matrices by GC/MS/MS

Bottle Group: C

of Sites: 1

Container ID: P-1L

Qty: 1

Preservation: ICE

Lot # 00071878

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

Description

Alkalinity in aqueous matrices as mg CaCO3/L

Turbidity in aqueous matrices

W-CL-IC	Chloride in aqueous matrices
W-COLOR	True Color measured at 450 nm
W-F	Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Sulfate in aqueous matrices
W-TDS	Total Dissolved Solids in aqueous matrices
W-TSS	Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 1 Preservation: ICE

Lot # 1215081

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 1 Preservation: ICE

Lot # 118605

Description: 250 ml Sterilized with tablet

Analysis

SE-FLR-ECQ

Description

Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie

Container ID: P-500ML

Qty: 1 Preservation: ICE And H2SO4

Lot # 00071230

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC

Description

Ammonia in aqueous matrices as mg N/L
Nitrite/Nitrate in aqueous matrices as mg N/L
Total Phosphorus in aqueous matrices as mg P/L
Total Kjeldahl Nitrogen in aqueous matrices
Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 1 Preservation: FILTER-ICE

Lot # 00070866

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: D

of Sites: 1

Container ID: P-1L

Qty: 1 Preservation: ICE

Lot # 00071078

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

W-CL-IC

W-COLOR

W-F

W-SO4-IC

W-TDS

W-TSS

Description

Alkalinity in aqueous matrices as mg CaCO3/L

Turbidity in aqueous matrices

Chloride in aqueous matrices

True Color measured at 450 nm

Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Sulfate in aqueous matrices

Total Dissolved Solids in aqueous matrices

Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 1 Preservation: ICE

Lot # 1215081

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 1 Preservation: ICE

Lot # L18005

Description: 250 ml Sterilized with tablet

Analysis

SE-FLR-ECQ

Description

Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie

Container ID: P-500ML

Qty: 1 Preservation: ICE And H2SO4

Lot # 00071230

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

W-NO2NO3

W-S-A-TP

W-TKN

W-TOC

Description

Ammonia in aqueous matrices as mg N/L

Nitrite/Nitrate in aqueous matrices as mg N/L

Total Phosphorus in aqueous matrices as mg P/L

Total Kjeldahl Nitrogen in aqueous matrices

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 1 Preservation: FILTER-ICE

Lot # 06070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group: E

of Sites: 1

Container ID: P-1L

Qty: 1 Preservation: ICE

Lot # 60071878

Description: Plastic Bottle 1L

Analysis

ALKALINITY

TURBIDITY

W-CL-IC

W-COLOR

W-F

W-SO4-IC

W-TDS

W-TSS

Description

Alkalinity in aqueous matrices as mg CaCO3/L

Turbidity in aqueous matrices

Chloride in aqueous matrices

True Color measured at 450 nm

Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Sulfate in aqueous matrices

Total Dissolved Solids in aqueous matrices

Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Container ID: BP-1L

Qty: 1 Preservation: ICE

Lot # 1215081

Description: Brown Plastic Bottle 1L

Analysis

CHLSUITE-W

Description

Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Container ID: P-250ML-WI-THI

Qty: 1 Preservation: ICE

Lot # 618005

Description: 250 ml Sterilized with tablet

Analysis

SE-FLR-ECQ

Description

Escherichia coli by Colilert 18 Quanti-Tray method by Flowers Chemical Laboratory in Port St. Lucie

Container ID: P-500ML

Qty: 1 Preservation: ICE And H2SO4

Lot # 00071230

Description: Plastic Bottle 500 mL

Affix YELLOW sulfuric acid dot to container.

Analysis

W-NH3

Description

Ammonia in aqueous matrices as mg N/L

W-NO2NO3

W-S-A-TP

W-TKN

W-TOC

Nitrite/Nitrate in aqueous matrices as mg N/L

Total Phosphorus in aqueous matrices as mg P/L

Total Kjeldahl Nitrogen in aqueous matrices

Nonpurgeable Organic Carbon in aqueous matrices

Container ID: P-125ML

Qty: 1 Preservation: FILTER-ICE

Lot # 00070886

Description: Plastic Bottle 125 mL

Analysis

W-PO4-F

Description

Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Cooler Packed By: Q.12

Number of Coolers: 3

Date: 8/20/18

Kit must also include:

- ☒ Field Sheets
- ☒ Temperature Control Bottle
- ☒ FedEx Bills, if applicable
- ☒ Plastic Bags/Zip Ties
- ☐ Acid Labels (sets of 24)

If Preservation Included:

ID <u>MCAA</u>	Exp <u>3/26/19</u>	Lot # <u>91013</u>
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____
ID _____	Exp _____	Lot # _____

PLEASE RETURN ALL COOLERS FOR RQ-2018-09-03-53