

Log-in Checklist

RQ- 2019-03-11-14

Login Station ID: #5

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 3/15 /19 @ 9:55

Cooler Temperature*	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # (fill out only if waybill copy is damaged)
	HNO ₃	Formalin		Present?		*Intact?		
				Yes	No	Yes	No	
<u>0.9</u>			<u>12</u>		☒			<u>4751 8777 1820</u>
<u>1.7</u>			<u>14</u>		☒			<u>4751 8777 1830</u>

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form.

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

**Chain of Custody/ Field Sheet(s) Included? Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☒ No ☐ Date 3-15-19

CONTAINER CHECK

**Evidence Tape on Bottles? Yes ☐ No ☒ If yes, is it intact? Yes ☒ No ☐
 **Caps on tight? Yes ☒ No ☐ **Containers intact? Yes ☒ No ☐
 **Sufficient Sample Volume: Yes ☒ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes ☐ No ☒ Init: KK

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.					
Analysis	Yes	No	Analysis	Yes	No
W-1,4 DIOX			W-PC-AA-SF		
W-E8321-DI, -MS, -21			W-NP-AA-SF		
W-BNA (-625, -R)			W-PDQUAT		
W-CARB-AA (-CAB-AA-R)			W-PCL-SQ-R		
W-U2060-MS			W-TR-TQ		
W-PCB-R			W-CT-CI-R		
W-PAH (-R)			W-PSNP-TQ		
W-PFAS-MS					

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ☒ No ☐ NA ☐ Init: KK

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes ☐ No ☒ NA ☐ Init: KK

Coolers Unpacked/Checked by: KK Date: 3/15 /19 NCRs (Y/N)? N

Event ID: UEC 2019
NE-DIST-2019-03-15-01 (SMP)
 Date Received: 15-MAR-2019 09:55 NCR #

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_05-2018
pH Paper 0.0 – 3.0 Lot #	
pH Paper 0 – 13 Lot #	222516EXP:9/15/19
Chlorine Test Strips Lot #	091417V EXP:9/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No X

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

UEC 2019

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: NE-DIST

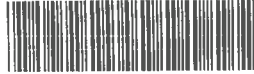
Requester: Jeremy M. Parrish

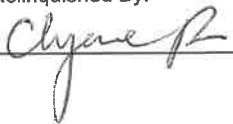
Field Report Prepared By: C. Ruben

Project ID: SMP

Collected By: C. Ruben & A. Kalmbacher

Sampling Agency: FDEP NE ROC

WIN ID/Field ID:  G5NE0026		☐ Comp ☒ Grab		Collection (comp begin or grab) Date 3/14/19 Time 1045		☒ Eastern ☐ Central	Composite end Date — Time —		Bottle Group(s) (See pg 2)	
Palm Coast Canal 0.6mi down from road		Matrix (type, e.g., Salt, Fresh, etc) Fresh		Diss. Oxygen 3.13		☒ mg/L ☐ % sat		Total Res. Chlorine (mg/L) —		C
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Temp (C) 20.4		pH 7.05		Sample Depth 0.05		☐ ft ☒ m		
Field ID:  G5NE0026 DUP		Salinity (ppt) 0.18		Comments						
Location: Palm Coast Canal 0.6 mi down from road		☐ Comp ☒ Grab		Collection (comp begin or grab) Date 3/14/19 Time 1050		☒ Eastern ☐ Central	Composite end Date — Time —		Bottle Group(s)	
WIN ID:  G5NE0026		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)		C
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Temp (C)		pH		Sample Depth		☐ ft ☐ m		
		Salinity (ppt)		Comments						
WIN ID/Field ID:  27010060		☐ Comp ☒ Grab		Collection (comp begin or grab) Date 3/14/19 Time 0950		☒ Eastern ☐ Central	Composite end Date — Time —		Bottle Group(s)	
CRACKER BR. @ CR 204		Matrix (type, e.g., Salt, Fresh, etc) Fresh		Diss. Oxygen 4.1		☒ mg/L ☐ % sat		Total Res. Chlorine (mg/L) —		C
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Temp (C) 19.4		pH 6.65		Sample Depth 0.3		☐ ft ☒ m		
		Salinity (ppt) 0.14		Comments						
WIN ID/Field ID:  27010050		☐ Comp ☒ Grab		Collection (comp begin or grab) Date 3/14/19 Time 0920		☒ Eastern ☐ Central	Composite end Date — Time —		Bottle Group(s)	
MOSES CR AT US 1		Matrix (type, e.g., Salt, Fresh, etc) Fresh		Diss. Oxygen 5.4		☒ mg/L ☐ % sat		Total Res. Chlorine (mg/L) —		A
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms		Temp (C) 18.7		pH 7.10		Sample Depth 0.05		☐ ft ☒ m		
		Salinity (ppt) 0.17		Comments						

Relinquished By: 	Date/Time: 3/14/19 1500	Shipping Method: Fedex	Number of Coolers Shipped: 2	Received By: 	Date/Time: 3/15/19/9:55
--	-------------------------	------------------------	------------------------------	--	-------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	CHLSUITE-W						
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1 to ALS
	NE-ALS-ECQ						
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-ICP						
	W-ICPMS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
	W-PO4-F						
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						

Group: B	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS	P-1L	# of Bottles: 3	Preservative: ICE	Enter Number of Bottles Sent to Lab	3
Group: C	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 3	Preservative: ICE	Enter Number of Bottles Sent to Lab	3
Group: C	Bottle Type: NE-ALS-ECQ	P-250ML-WI-THI	# of Bottles: 3	Preservative: ICE	Enter Number of Bottles Sent to Lab	3 to ACS
Group: C	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 3	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	3
Group: C	Bottle Type:	P-125ML	# of Bottles: 3	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	3

Group: C	Bottle Type:	P-125ML	# of Bottles: 3	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab
----------	--------------	---------	-----------------	--------------------------	-------------------------------------

W-PO4-F

3

Group: D	Bottle Type:	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab
----------	--------------	------	-----------------	-------------------	-------------------------------------

ALKALINITY
TURBIDITY
W-CL-IC
W-COLOR
W-F
W-SO4-IC
W-TDS
W-TSS

1

Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab
----------	--------------	---------	-----------------	--------------------	-------------------------------------

W-ICP
W-ICPMS

1

Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab
----------	--------------	---------	-----------------	-----------------------------	-------------------------------------

W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC

1

Group: D	Bottle Type:	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab
----------	--------------	---------	-----------------	--------------------------	-------------------------------------

W-PO4-F

1

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: NE-DIST

Project ID: SMP

Requester: Jeremy M. Parrish

Field Report Prepared By: Cheyenne RibenCollected By: C. Riben & A. KalmbacherSampling Agency: FDEP

Location: LAKE VEDRA AT SOLANA RD.		☐ Comp ☒ Grab		Collection (comp begin or grab) Date 3/14/19 Time 1350		☒ Eastern ☐ Central		Composite end Date Time		Bottle Group(s) (See pg 2)	
FIELD ID:  FB 03142019		Matrix (type, e.g., Salt, Fresh, etc) D7 water		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)		D	
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Temp (C)		pH		Sample Depth ☐ ft ☐ m			
WIN ID/Field ID:  27010180		Matrix (type, e.g., Salt, Fresh, etc) Fresh		Diss. Oxygen		☒ mg/L ☐ % sat		Total Res. Chlorine (mg/L)		B	
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Temp (C) 21.0		pH 8.93		Sample Depth 0.3 ☒ ft ☐ m			
Location		☐ Comp ☐ Grab		Collection (comp begin or grab) Date Time		☒ Eastern ☐ Central		Composite end Date Time		Bottle Group(s)	
Field ID		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)			
WIN/STORET #		Temp (C)		pH		Sample Depth ☐ ft ☐ m		Sp. Conductance (umho/cm)			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					
Location		☐ Comp ☐ Grab		Collection (comp begin or grab) Date Time		☒ Eastern ☐ Central		Composite end Date Time		Bottle Group(s)	
Field ID		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)			
WIN/STORET #		Temp (C)		pH		Sample Depth ☐ ft ☐ m		Sp. Conductance (umho/cm)			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					
Location		☐ Comp ☐ Grab		Collection (comp begin or grab) Date Time		☒ Eastern ☐ Central		Composite end Date Time		Bottle Group(s)	
Field ID		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)			
WIN/STORET #		Temp (C)		pH		Sample Depth ☐ ft ☐ m		Sp. Conductance (umho/cm)			
Latitude ☐ dd ☐ dms		Longitude ☐ dd ☐ dms		Salinity (ppt)		Comments					

Relinquished By: <u>Cheyenne Riben</u>	Date/Time: 3/14/19 1300	Shipping Method: Fedex	Number of Coolers Shipped: 2	Received By: <u>KH</u>	Date/Time: 3.15.19 9:55
--	-------------------------	------------------------	------------------------------	------------------------	-------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	NE-ALS-ECQ					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					

Date of Request: 29-JAN-2019
 Created By: PARRISH_J On: 29-JAN-2019
 Modified By: JASROTIA_P On: 18-FEB-2019
 Customer: NE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Northeast District
 Sampling Event: UEC 2019

Send Coolers To:

Phone: 904-807-3300
 8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Program Module Number:
 Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 18-FEB-2019
 Biology Request Reviewed By: JASROTIA_P On: 18-FEB-2019
 Sampling Kit Required: YES Ship on: 01-MAR-2019
 Sampling Kit Shipped: YES On: 25-FEB-2019
 Sampling Kit Packed By: KLUG_K
 Preservatives Needed: NO
 Date to Receive Samples: 11-MAR-2019
 Received By:

Send Final Report To:

8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Group A: Nutrients, Metals, Bacteria
 Group B: Nutrients, Metals,
 Group C: Nutrients, Bacteria
 Group D: Nutrients,metals

NE-ALS-ECQ

Packing Notes:

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Group C (Water) With 3 Samples:

Bottle Type: P-1L	Number of Bottles: 3	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 3	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 3	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville
Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 3	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-500ML

W-ICP

Template:

Number of Bottles: 1

Preserved With: HNO3

(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Barium

Calcium

Iron

Magnesium

Potassium

Sodium

Zinc

W-ICPMS

Template:

(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

*Boron

Aluminum

Antimony

Arsenic

Cadmium

Chromium

Copper

Lead

Nickel

Selenium

Silver

Thallium

Bottle Type: P-500ML

Number of Bottles: 1

Preserved With: ICE And H2SO4

W-NH3

Template: DEFAULT

(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L

W-NO2NO3

Template: DEFAULT

(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L

W-S-A-TP

Template: DEFAULT

(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L

W-TKN

Template: DEFAULT

(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices

W-TOC

Template: DEFAULT

(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML

Number of Bottles: 1

Preserved With: FILTER-ICE

W-PO4-F

Template: DEFAULT

(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.