

RQ- 2019-04-08-59

Log-in Checklist

Login Station ID: #2Shipping Method: FedEx | UPS | HD | GreyhoundDate/Time of Receipt: 4-10-19/9:40

Cooler Temperature*	Samples Frozen YES	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
		HNO ₃	Formalin		Present?		*Intact?		
					Yes	No	Yes	No	
<u>2.4</u>				<u>5</u>		☒			

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. DOH coolers can exceed 2.0 °C above standard temperature without customer notification/confirmation.

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

**Chain of Custody/ Field Sheet(s) Included? Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☐ No ☒

CONTAINER CHECK

**Evidence Tape on Bottles? Yes ☒ No ☒ If yes, is it intact? Yes ☐ No ☐

**Caps on tight? Yes ☒ No ☐ **Containers intact? Yes ☒ No ☐

**Sufficient Sample Volume: Yes ☒ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes ☐ No ☐ Init:

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.					
Analysis	Yes	No	Analysis	Yes	No
W-1,4 DIOX			W-PC-AA-SF		
W-E8321-DI, -MS, -21			W-NP-AA-SF		
W-BNA (-625, -R)			W-PDQUAT		
W-CARB-AA (-CAB-AA-R)			W-PCL-SQ-R		
W-U2060-MS			W-TR-TQ		
W-PCB-R			W-CT-CI-R		
W-PAH (-R)			W-PSNP-TQ		
W-PFAS-MS					

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ☒ No ☐ NA ☐ Init: SR

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes ☐ No ☒ NA ☒ Init: SR

Coolers Unpacked/Checked by: SR Date: 4-10-19 NCRs (Y/N)? N

Event ID: Continuous Monitoring Event
NE-DIST-2019-04-10-01 (SMP)

NCR #

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_02_2018EXP:4/27/19
pH Paper 0.0 – 3.0 Lot #	220416A EXP:7-30-20
pH Paper 0 – 13 Lot #	222516 EXP:08-15-19
Chlorine Test Strips Lot #	8212 EXP:05-2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ✓

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF; W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Continuous Monitoring Event

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: NE-DIST

Requester: Brian Davis

Field Report Prepared By: B. Davis

Project ID: SMP

Collected By: B. Davis, P. O'ConnorSampling Agency: DEAR NEV

Location <u>Hopkins Creek at Kings Rd</u>	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>4/9/19</u> Time <u>1430</u>	Eastern Central	Composite end Date Time	Bottle Group(s) (See pg 2)
Field ID	Matrix (type, e.g., Salt, Fresh, etc) <u>S-F</u>	Diss. Oxygen <u>6.89</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	<u>C</u>
WIN/STORET # <u>G2NE0766</u>	Temp (C) <u>24.2</u>	pH <u>7.5</u>	Sample Depth <u>0.3</u>	Sp. Conductance (umho/cm) <u>586</u>	
Latitude ☐ dd ☐ dms	Longitude ☐ dd ☐ dms	Salinity (ppt) <u>21.74</u>	Comments		

Location	☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORET #	Temp (C)	pH	Sample Depth	Sp. Conductance (umho/cm)	
Latitude ☐ dd ☐ dms	Longitude ☐ dd ☐ dms	Salinity (ppt)	Comments		

Location	☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORET #	Temp (C)	pH	Sample Depth	Sp. Conductance (umho/cm)	
Latitude ☐ dd ☐ dms	Longitude ☐ dd ☐ dms	Salinity (ppt)	Comments		

Location	☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORET #	Temp (C)	pH	Sample Depth	Sp. Conductance (umho/cm)	
Latitude ☐ dd ☐ dms	Longitude ☐ dd ☐ dms	Salinity (ppt)	Comments		

Relinquished By: <u>B. Davis</u>	Date/Time <u>4/9/19 1630</u>	Shipping Method <u>FedEx</u>	Number of Coolers Shipped: <u>1</u>	Received By: <u>SR</u>	Date/Time <u>4-10-19 9:40</u>
----------------------------------	---------------------------------	---------------------------------	-------------------------------------	------------------------	----------------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-BR-IC					
	W-F					
	W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	NE-ALS-ENQ					
Group: A	Bottle Type:	QPCR-P-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	PCR-HF183					
Group: A	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-E8321-DI					
	W-E8321-MS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP					
	W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3					
	W-NO2NO3					
	W-S-A-TP					
	W-TKN					
	W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-BR-IC					
	W-F					
	W-TSS					

Group: B	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
CHLSUITE-W					
Group: B	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
NE-ALS-ENQ					
Group: B	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	0
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: B	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	0
W-PO4-F					
Group: C	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY TURBIDITY W-BR-IC W-F W-TSS					
Group: C	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W					
Group: C	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
NE-ALS-ENQ					
Group: C	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-F-F W-PO4-F					
Group: C	Bottle Type: P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab	1
W-ICP W-ICPMS					
Group: C	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3 W-NO2NO3 W-S-A-TP					

Group: C

Bottle Type:

P-500ML

of Bottles: 1

Preservative: ICE And H₂SO₄

Enter Number of Bottles Sent to Lab

1

W-TKN

W-TOC

Date of Request: 03-APR-2019
 Created By: DAVIS_B On: 03-APR-2019
 Modified By: JASROTIA_P On: 04-APR-2019
 Customer: NE-DIST
 Project: SMP
 Division:
 District: Northeast District
 Sampling Event: Continuous Monitoring Event

Send Coolers To:

Phone: 904-256-1568
 8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Brian G Davis

Program Module Number:
 Priority: 3
 Event Status: A
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 04-APR-2019
 Biology Request Reviewed By: JASROTIA_P On: 04-APR-2019
 Sampling Kit Required: NO Ship on:
 Sampling Kit Shipped:
 Sampling Kit Packed By:
 Preservatives Needed: NO
 Date to Receive Samples: 08-APR-2019
 Received By:

Send Final Report To:

8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Brian G Davis

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Continuous Monitoring SMP
 All Groups MARINE, do not send coolers.
 Group A: Nutrients, Metals, Tracers, Markers, Bacteria
 Group B: Nutrients, Bacteria
 Group C: Nutrients, Metals, Bacteria
 Each group will be sent in through separate sampling events over the next 3 months.

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
NE-ALS-ENQ	Template: DEFAULT	(ASTM D650399) Enterococci by Enterolert Quanti-Tray Method by ALS Environmental Laboratory in Jacksonville
Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Manganese		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
NE-ALS-ENQ	Template: DEFAULT	(ASTM D650399) Enterococci by Enterolert Quanti-Tray Method by ALS Environmental Laboratory in Jacksonville
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L W-TSS	Template: DEFAULT	Number of Bottles: 1 (SM 2540 D-97)	Preserved With: ICE Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L CHLSUITE-W	Template: DEFAULT	Number of Bottles: 1 (SM 10200 H (mod.))	Preserved With: ICE Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI NE-ALS-ENQ	Template: DEFAULT	Number of Bottles: 1 (ASTM D650399)	Preserved With: ICE Enterococci by Enterolert Quanti-Tray Method by ALS Environmental Laboratory in Jacksonville
Bottle Type: P-125ML W-F-F	Template: DEFAULT	Number of Bottles: 1 (SM 4500 F-C-97)	Preserved With: FILTER-ICE Fluoride, dissolved, in filtered aqueous matrices, without distillation (not valid for NPDES monitoring)
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-500ML W-ICP	Template:	Number of Bottles: 1 (EPA 200.7 Rev. 4.4)	Preserved With: HNO3 Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium			
Calcium			
Iron			
Manganese			
Potassium			
Sodium			
Zinc			
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4)	Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum			
Antimony			
Arsenic			
Cadmium			
Chromium			
Lead			
Nickel			
Selenium			
Silver			
Thallium			
Bottle Type: P-500ML W-NH3	Template: DEFAULT	Number of Bottles: 1 (EPA 350.1 Rev. 2.0)	Preserved With: ICE And H2SO4 Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0)	Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0)	Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00)	Nonpurgeable Organic Carbon in aqueous matrices

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.

Do Not Lift Using This Tag

ORIGIN ID: TLHA (904) 807-3300
JEREMY PARRISH
FL DEP NORTHEAST DISTRICT
8800 BAYMEADOWS WAY WEST
SUITE 100
JACKSONVILLE, FL 32256
UNITED STATES US

SHIP DATE: 07FEB19
ACTWGT: 10.00 LB MAN
CHD: 0448970/CAFE3211


TO **JOHN WATTS**
FLORIDA DEP/CHEMISTRY SECTION
2600 BLAIRSTONE RD

TALLAHASSEE FL 32399

(850) 245-8077
REF: 1001
DOI:

DEPT:

RMA: |||||

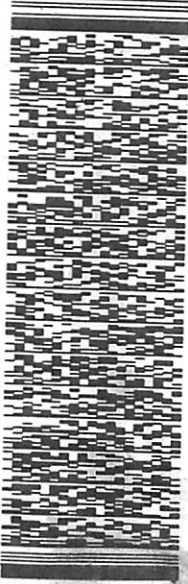

FedEx
Express



TRK# 4751 8776 4910
0221

FedEx

4751 8776 4910

RETURNS MON - FRI
PRIORITY OVERNIGHT
WED - 10 APR 10:30A
PRIORITY OVERNIGHT

36 TLHA

32399
FL-US
TLH



10 3690753 09APR19 NRBA 553C1/D7E5/0C8A

R1 **856** 1 **B**
FZ 10:30 4910 04.10