

Log-in Checklist

RQ- 2019-04-15-17

Login Station ID: #5

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 4/18 /19 @ 10:00

Cooler Temperature*	Samples Frozen	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
0.7				16		✓			4751 8778 4613
2.3				14		✓			4751 8778 4624

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. **DOH coolers can exceed 2.0 °C above standard temperature without customer notification/confirmation.**

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

**Chain of Custody/ Field Sheet(s) Included? Yes ✓ No

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes No ✓

CONTAINER CHECK

**Evidence Tape on Bottles? Yes No ✓ If yes, is it intact? Yes No

**Caps on tight? Yes ✓ No **Containers intact? Yes ✓ No

**Sufficient Sample Volume: Yes ✓ No

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes No ✓ Init: KK

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.					
Analysis	Yes	No	Analysis	Yes	No
W-1,4 DIOX			W-PC-AA-SF		
W-E8321-DI, -MS, -21			W-NP-AA-SF		
W-BNA (-625, -R)			W-PDQUAT		
W-CARB-AA (-CAB-AA-R)			W-PCL-SQ-R		
W-U2060-MS			W-TR-TQ		
W-PCB-R			W-CT-CI-R		
W-PAH (-R)			W-PSNP-TQ		
W-PFAS-MS					

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ✓ No NA Init: KK

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes No NA ✓ Init: KK

Coolers Unpacked/Checked by: KK Date: 4/18 /19 NCRs (Y/N)? N

Event ID: LSJR North Boat
NE-DIST-2019-04-18-01 (SMP)

Page 1 of 2, 03/25/2019 Date Received: 18-APR-2019 10:00

NCR #

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	
pH Paper 0.0 – 3.0 Lot #	220416A EXP:7/30/20
pH Paper 0 – 13 Lot #	222516 EXP:9/15/2019
Chlorine Test Strips Lot #	091417V EXP:9/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No X

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: NE-DIST

Project ID: SMP

Requester: Jeremy M. Parrish

Field Report Prepared By:

Collected By:

Pat O'Connor, Nicole Reistad

Sampling Agency:

Nicole Reistad
FDEP - NE ROC

Location	WIN ID/Field ID:		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 04/17/19 Time 1125	Eastern Central	Composite end Date	Time	Bottle Group(s) (See pg 2)
Field ID		20030129	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		A
WIN/STORE	RIBAUT R MONCRIEF RD BR		Temp (C)	pH	Sample Depth	ft m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location	WIN ID/Field ID:		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 04/17/19 Time 1300	Eastern Central	Composite end Date	Time	Bottle Group(s)
Field ID		20030115	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		C
WIN/STORE	MONCRIEF CREEK NEAR MOUTH		Temp (C)	pH	Sample Depth	ft m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location	WIN ID/Field ID:		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 04/17/19 Time 1410	Eastern Central	Composite end Date	Time	Bottle Group(s)
Field ID		20030114	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		C
WIN/STORE #	MONCRIEF CR LEM TURNER RD BR/ NORWOOD AVE		Temp (C)	pH	Sample Depth	ft m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location			☐ Comp ☐ Grab	Collection (comp begin or grab) Date	Eastern Central	Composite end Date	Time	Bottle Group(s)
Field ID			Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		
WIN/STORE #			Temp (C)	pH	Sample Depth	ft m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			

Relinquished By:	Date/Time	Shipping Method	Number of Coolers Shipped:	Received By:	Date/Time
N. Reistad	04/17/19	FED Ex	2 (total)	NR	4.18.19/10:00

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab
	ALKALINITY TURBIDITY W-BR-IC W-F W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab
	NE-ALS-ENQ					
Group: A	Bottle Type:	QPCR-P-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab
	PCR-FILTER					
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab
	W-E8321-DI W-E8321-MS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab
	ALKALINITY TURBIDITY W-BR-IC W-F W-TSS					
Group: B	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab
	CHLSUITE-W					
Group: B	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab

Group: B	Bottle Type: P-250ML-WI-THI NE-ALS-ENQ	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
Group: B	Bottle Type: P-500ML W-ICP W-ICPMS	# of Bottles: 2	Preservative: HNO3	Enter Number of Bottles Sent to Lab	2
Group: B	Bottle Type: P-500ML W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	# of Bottles: 2	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	2
Group: B	Bottle Type: P-125ML W-PO4-F	# of Bottles: 2	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	2
Group: C	Bottle Type: P-1L ALKALINITY TURBIDITY W-BR-IC W-F W-TSS	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	2
Group: C	Bottle Type: BP-1L CHLSUITE-W	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	2
Group: C	Bottle Type: P-250ML-WI-THI NE-ALS-ENQ	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	0
Group: C	Bottle Type: P-500ML W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	# of Bottles: 2	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	2
Group: C	Bottle Type: P-125ML W-PO4-F	# of Bottles: 2	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	2
Group: D	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1

Group: D	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-BR-IC						
	W-F						
	W-TSS						
Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-ICP						
	W-ICPMS						
Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: D	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
	W-PO4-F						

Request Number: RQ-2019-04-15-17

Florida Department of Environmental Protection Central Laboratory Submittal Form

Page 1 of 4

LSJR North Boat

last revised Jan 25, 2018

Customer: NE-DIST

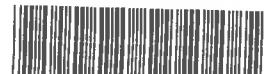
Requester: Jeremy M. Parrish

Field Report Prepared By: Nicole Reistad

Project ID: SMP

Collected By: Pat O'Connor, Nicole Reistad

Sampling Agency: FDEP-NE ROC

Location	Location: TROUT R @ I-95 BURT MAXWELL DOCK	☐ Comp ☒ Grab	Collection (comp begin or grab) Date 04/17/19 Time 0925	Eastern Central	Composite end Date Time	Bottle Group(s) (See pg 2)
Field ID	FIELD ID:  FB 04172019	Matrix (type, e.g., Salt, Fresh, etc) Field Blank		Diss. Oxygen ☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	D
WIN/STOR		Temp (C)	pH	Sample Depth ☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments	
Location	WIN ID/Field ID:  20030121	☐ Comp ☒ Grab	Collection (comp begin or grab) Date 04/17/19 Time 0900	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID	TROUT R @ I-95 BURT MAXWELL DOCK	Matrix (type, e.g., Salt, Fresh, etc) Salt		Diss. Oxygen 7.20	☒ mg/L ☐ % sat	AB PO
WIN/STORE		Temp (C) 23.5	pH 7.75	Sample Depth 0.3	Sp. Conductance (umho/cm) 16001	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) 9.36	Comments	
Location	Field ID:  20030121 DUP	☐ Comp ☒ Grab	Collection (comp begin or grab) Date 04/17/19 Time 0915	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID	Location: TROUT R @ I-95 BURT MAXWELL DOCK	Matrix (type, e.g., Salt, Fresh, etc) Field Replicate		Diss. Oxygen ☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	AB PO
WIN/STORET	WIN ID:  20030121	Temp (C)	pH	Sample Depth 0.3	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments	
Location	WIN ID/Field ID:  20030133	☐ Comp ☒ Grab	Collection (comp begin or grab) Date 04/17/19 Time 1050	Eastern Central	Composite end Date Time	Bottle Group(s)
Field ID	RIBAUT RIVER AT US-1	Matrix (type, e.g., Salt, Fresh, etc) Salt		Diss. Oxygen 7.33	☒ mg/L ☐ % sat	A
WIN/STORET		Temp (C) 23.5	pH 7.59	Sample Depth 0.3	Sp. Conductance (umho/cm) 5191	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) 2.79	Comments	

Relinquished By: N. Reistad	Date/Time 04/17/19	Shipping Method Fed Ex	Number of Coolers Shipped: 2 (total)	Received By: HK	Date/Time 4.18.19 / 10:00
--------------------------------	-----------------------	---------------------------	---	--------------------	------------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	ALKALINITY TURBIDITY W-BR-IC W-F W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	CHLSUITE-W						
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	NE-ALS-ENQ						
Group: A	Bottle Type:	QPCR-P-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	PCR-FILTER						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	W-E8321-DI W-E8321-MS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	2
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC						
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	2
	W-PO4-F						
Group: B	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	ALKALINITY TURBIDITY W-BR-IC W-F W-TSS						
Group: B	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	CHLSUITE-W						
Group: B	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0

Date of Request: 08-MAR-2019
 Created By: PARRISH_J On: 08-MAR-2019
 Modified By: JASROTIA_P On: 25-MAR-2019
 Customer: NE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Northeast District
 Sampling Event: LSJR North Boat

Send Coolers To:

Phone: 904-256-1700
 8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Program Module Number:

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 19-MAR-2019
 Biology Request Reviewed By: JASROTIA_P On: 25-MAR-2019
 Sampling Kit Required: YES Ship on: 05-APR-2019
 Sampling Kit Shipped: YES On: 03-APR-2019
 Sampling Kit Packed By: REYNOLDS_T
 Preservatives Needed: NO
 Date to Receive Samples: 15-APR-2019
 Received By:

Send Final Report To:

8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments:

Group A: Marinewater
 Nutrients, Tracers, Markers, Bacteria

Group B: Marinewater
 Nutrients, Metals, Bacteria

Group C: Marinewater
 Nutrients, Bacteria,

Group D: Marinewater
 Nutrients, Metals

NE-ALS-ENQ

Packing Notes:

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 2	Preserved With: ICE
NE-ALS-ENQ	Template: DEFAULT	(ASTM D650399) Enterococci by Enterolert Quanti-Tray Method by ALS Environmental Laboratory in Jacksonville
Bottle Type: QPCR-P-500ML	Number of Bottles: 2	Preserved With: ICE
PCR-FILTER	Template: DEFAULT	(SOP-PCR-4.0) Preparation of Samples by Filtration and held for qPCR Analysis
Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: ICE

Bottle Group A (Water) With 2 Samples:

Bottle Type: BG-500ML		Number of Bottles: 2	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B)	Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B)	Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML		Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0)	Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0)	Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0)	Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00)	Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML		Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 2 Samples:

Bottle Type: P-250ML-WI-THI		Number of Bottles: 2	Preserved With: ICE
NE-ALS-ENQ	Template: DEFAULT	(ASTM D650399)	Enterococci by Enterolert Quanti-Tray Method by ALS Environmental Laboratory in Jacksonville
Bottle Type: P-500ML		Number of Bottles: 2	Preserved With: HNO3
W-ICP	Template: S-METALS	(EPA 200.7 Rev. 4.4)	Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron			
Manganese			
Zinc			
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4)	Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum			
Antimony			
Arsenic			
Cadmium			
Chromium			
Copper			
Lead			
Nickel			
Selenium			
Silver			
Thallium			
Bottle Type: P-500ML		Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0)	Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0)	Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0)	Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00)	Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML		Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L		Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97)	Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0)	Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0)	Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97)	Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97)	Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L		Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.))	Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Group C (Water) With 2 Samples:

Bottle Type: P-1L		Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97)	Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0)	Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0)	Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97)	Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Bottle Group C (Water) With 2 Samples:

Bottle Type: P-1L W-TSS	Template: DEFAULT	Number of Bottles: 2 (SM 2540 D-97)	Preserved With: ICE Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L CHLSUITE-W	Template: DEFAULT	Number of Bottles: 2 (SM 10200 H (mod.))	Preserved With: ICE Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI NE-ALS-ENQ	Template: DEFAULT	Number of Bottles: 2 (ASTM D650399)	Preserved With: ICE Enterococci by Enterolert Quanti-Tray Method by ALS Environmental Laboratory in Jacksonville
Bottle Type: P-500ML W-NH3	Template: DEFAULT	Number of Bottles: 2 (EPA 350.1 Rev. 2.0)	Preserved With: ICE And H2SO4 Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0)	Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0)	Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00)	Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML W-PO4-F	Template: DEFAULT	Number of Bottles: 2 (EPA 365.1 Rev. 2.0)	Preserved With: FILTER-ICE Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L ALKALINITY	Template: DEFAULT	Number of Bottles: 1 (SM 2320 B-97)	Preserved With: ICE Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0)	Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0)	Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97)	Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97)	Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: P-500ML W-ICP	Template: S-METALS	Number of Bottles: 1 (EPA 200.7 Rev. 4.4)	Preserved With: HNO3 Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron			
Manganese			
Zinc			
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4)	Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum			
Antimony			
Arsenic			
Cadmium			
Chromium			
Copper			
Lead			
Nickel			
Selenium			
Silver			
Thallium			
Bottle Type: P-500ML W-NH3	Template: DEFAULT	Number of Bottles: 1 (EPA 350.1 Rev. 2.0)	Preserved With: ICE And H2SO4 Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0)	Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0)	Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00)	Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML W-PO4-F	Template: DEFAULT	Number of Bottles: 1 (EPA 365.1 Rev. 2.0)	Preserved With: FILTER-ICE Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.