

RQ-2019-04-22.33

Log-in Checklist

Login Station ID: #2

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 4-24-19/15:32

Cooler Temperature*	Samples Frozen YES	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
		HNO ₃	Formalin		Present?		*Intact?		
					Yes	No	Yes	No	
2.7				23		✓			

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. **DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C) without login confirmation from customer for samples requiring ice preservation. An NCR still is required.**

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

****Chain of Custody/ Field Sheet(s) Included?** Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☐ No ☐

CONTAINER CHECK

****Evidence Tape on Bottles?** Yes ☐ No ☒ **If yes, is it intact?** Yes ☐ No ☐

****Caps on tight?** Yes ☒ No ☐ ****Containers intact?** Yes ☒ No ☐

****Sufficient Sample Volume:** Yes ☒ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes ☐ No ☐ Init:

Chlorine detected? (One container checked per Field ID) Yes ☐ No ☐ Init:

If "YES", list affected Field IDs and Test IDs below.

Field ID	List Test ID(s) where Chlorine was Detected, an NCR is required.

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

****Acid Preserved Samples: pH ≤ 2.0?** Yes ☒ No ☐ NA ☐ Init: SR

****W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?** Yes ☐ No ☐ NA ☒ Init: SR

Coolers Unpacked/Checked by: SR **Date:** 4-24-19 **NCRs (Y/N)?** N

Event ID: Run 12
TLH-ROC-2019-04-24-01 (SMP) **NCR #**

Log-in Checklist

Samples with Incorrect pH Preservation

Field ID	Bottle Test ID(s)	pH	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_02_2018 EXP:4-27-19
pH Test Strips 0 – 6 Lot #	213316VEXP:05-15-2019
pH Paper 0 – 3 Lot #	220416A EXP:07-30-20
pH Paper 0 – 13 Lot #	222516 EXP:08-15-2019
Chlorine Test Strips Lot #	8212 EXP:05-2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ✓

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Run 12

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: TLH-ROC

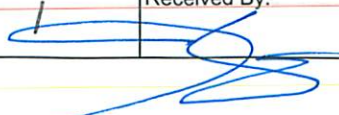
Requester: Michael Eckles

Field Report Prepared By: Jon Woods

Project ID: SMP

Collected By: Woods/ButchettSampling Agency: 8034

Location	Field/WIN ID → Location ↓	 G2TLHR0038	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>4/24/19</u> Time <u>1100</u>	<u>Eastern</u> Central	Composite end Date Time	Bottle Group(s) (See pg 2)	
Field ID				Matrix (type, e.g., Salt, Fresh, etc.)	Diss. Oxygen <u>5.10</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STOR		 CROOKED RIVER @ POST 46		Temp (C) <u>21.9</u>	pH <u>6.56</u>	Sample Depth <u>0.3</u>	Sp. Conductance (umho/cm) <u>448.6</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.22</u>	Comments			
Location	Field/WIN ID → Location ↓	 G2TLHR0025	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>4/24/19</u> Time <u>1220</u>	<u>Eastern</u> Central	Composite end Date Time	Bottle Group(s)	
Field ID				Matrix (type, e.g., Salt, Fresh, etc.)	Diss. Oxygen <u>7.05</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STOR		 Alligator Harbor near clam beds		Temp (C) <u>22.1</u>	pH <u>7.90</u>	Sample Depth <u>0.3</u>	Sp. Conductance (umho/cm) <u>474.5</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>30.90</u>	Comments			
Location	Field/WIN ID → Location ↓	 G2TLHR0029	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>4/24/19</u> Time <u>1350</u>	<u>Eastern</u> Central	Composite end Date Time	Bottle Group(s)	
Field ID				Matrix (type, e.g., Salt, Fresh, etc.)	Diss. Oxygen <u>7.36</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STOR		 DICKERSON BAY SOUTH		Temp (C) <u>22.6</u>	pH <u>8.10</u>	Sample Depth <u>0.3</u>	Sp. Conductance (umho/cm) <u>40624</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>26.01</u>	Comments			
Location	Field/WIN ID → Location ↓	 G2TLHR0028	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>4/24/19</u> Time <u>1400</u>	<u>Eastern</u> Central	Composite end Date Time	Bottle Group(s)	
Field ID				Matrix (type, e.g., Salt, Fresh, etc.)	Diss. Oxygen <u>7.13</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STOR		 DICKERSON BAY NORTH		Temp (C) <u>22.4</u>	pH <u>8.17</u>	Sample Depth <u>0.3</u>	Sp. Conductance (umho/cm) <u>41152</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>26.39</u>	Comments			

Relinquished By: <u>Jon Woods</u>	Date/Time <u>4/24/19 1529</u>	Shipping Method <u>HO</u>	Number of Coolers Shipped: <u>1</u>	Received By: 	Date/Time <u>4/24/19 1532</u>
--------------------------------------	----------------------------------	------------------------------	--	---	----------------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ECOLI-18QT					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP					
	W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3					
	W-NO2NO3					
	W-S-A-TP					
	W-TKN					
	W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab <u> / </u>
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					

Group: B Bottle Type: BP-1L # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1
CHLSUITE-W

Group: B Bottle Type: P-250ML-WO-THI # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1
ECOLI-18QT

Group: B Bottle Type: P-500ML # of Bottles: 1 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 1
W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC
8344670
exp. 12/10/19

Group: B Bottle Type: P-125ML # of Bottles: 1 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 1
W-PO4-F

Group: C Bottle Type: P-1L # of Bottles: 3 Preservative: ICE Enter Number of Bottles Sent to Lab 3
ALKALINITY
TURBIDITY
W-BR-IC
W-F
W-TSS

Group: C Bottle Type: BP-1L # of Bottles: 3 Preservative: ICE Enter Number of Bottles Sent to Lab 3
CHLSUITE-W

Group: C Bottle Type: P-250ML-WO-THI # of Bottles: 3 Preservative: ICE Enter Number of Bottles Sent to Lab 3
ENT-24-QT

Group: C Bottle Type: P-250ML-WO-THI # of Bottles: 3 Preservative: ICE Enter Number of Bottles Sent to Lab 3
FCOLI-MF

Group: C Bottle Type: P-500ML # of Bottles: 3 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 3
W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC
8344070
exp. 12/10/19

Group: C Bottle Type: P-125ML # of Bottles: 3 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 3
W-PO4-F

Printed: 4/24/2019 3:37:34 PM

Page 1 of 3

Date of Request: 04-APR-2019
 Created By: ECKLES_M On: 04-APR-2019
 Modified By: JASROTIA_P On: 08-APR-2019
 Customer: TLH-ROC
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Not Applicable
 Sampling Event: Run 12

Send Coolers To:

Phone: 850-245-7586

2600 Blair Stone Rd

Tallahassee, FL 32399

Attn: Michael Eckles

Send Final Report To:

2600 Blair Stone Rd

Tallahassee, FL 32399

Attn: Michael Eckles

Program Module Number: 0100
 Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 08-APR-2019
 Biology Request Reviewed By: JASROTIA_P On: 08-APR-2019
 Sampling Kit Required: YES Ship on: 18-APR-2019
 Sampling Kit Shipped: YES On: 15-APR-2019
 Sampling Kit Packed By: KLUG_K
 Preservatives Needed: NO
 Date to Receive Samples: 22-APR-2019
 Received By:

Report Type: Final Only

FTP Data: NO

QC Report: YES

Date Log: NO

Authorization Log: NO

Comments: Group A WBID 1242
 Group B WBID 1251
 Group C WBID 1256 & 1239(2 Sites)
 No fedex labels

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
, ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Group C (Water) With 3 Samples:

Bottle Type: P-1L	Number of Bottles: 3	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices

Bottle Group C (Water) With 3 Samples:

Bottle Type: P-1L	Number of Bottles: 3	Preserved With: ICE
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 3	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 3	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-250ML-WO-THI	Number of Bottles: 3	Preserved With: ICE
FCOLI-MF	Template: DEFAULT	(SM 9222 D) Fecal coliforms by membrane filter method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 3	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.