

Log-in Checklist

RQ- 2019-06-17-41

Login Station ID: #5

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 8/ 5 /19 @ 14:22

Cooler Temperature*	Samples Frozen	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
2.4				5		✓			

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. **DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C) without login confirmation from customer for samples requiring ice preservation. An NCR still is required.**

**Chain of Custody/ Field Sheet(s) Included? Yes ✓ No

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes No ✓

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

CONTAINER CHECK

**Evidence Tape on Bottles? Yes No ✓ If yes, is it intact? Yes No

**Caps on tight? Yes ✓ No **Containers intact? Yes ✓ No

**Sufficient Sample Volume: Yes ✓ No

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes No ✓ Init: KK

Chlorine detected? (One container checked per Field ID) Yes No Init:

If "YES", list affected Field IDs and Test IDs below.

Field ID	List TestID(s) where Chlorine was Detected, an NCR is required.

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ✓ No NA Init: KK

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes No NA ✓ Init: KK

Coolers Unpacked/Checked by: KK Date: 8/ 5 /19 NCRs (Y/N)? N

Event ID: Run 3
TLH-ROC-2019-08-05-01 (SMP)
Date Received: 05-AUG-2019 14:22
 NCR #

Log-in Checklist

Samples with Incorrect pH Preservation

Field ID	Bottle Test ID(s)	pH	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_07_2018 EXP:12/28/19
pH Test Strips 0 – 6 Lot #	
pH Paper 0 – 3 Lot #	220416A EXP:7/30/20
pH Paper 0 – 13 Lot #	222516 EXP:9/15/19
Chlorine Test Strips Lot #	091417V EXP:9/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes X No X

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF; W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Run 3

last revised Jan 25, 2018

Customer: TLH-ROC

Requester: Michael Eckles

Field Report Prepared By:

J. Burtchett

Project ID: SMP

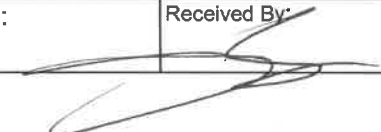
Collected By:

J. Burtchett, J. Woods

Sampling Agency:

8034

Location		☐ Comp ☒ Grab		Collection (comp begin or grab) Date <u>08/05/19</u> Time <u>1220</u>		☒ Eastern ☐ Central		Composite end Date <u>-</u> Time <u>-</u>		Bottle Group(s) (See pg 2)	
Field	Field/WIN ID →	 G1TLHR0081		Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>		Diss. Oxygen <u>5.83</u>		☒ mg/L ☐ % sat		Total Res. Chlorine (mg/L) <u>-</u>	
WIN/	Location ↓	<u>Bailey Mill Creek @ Wilson Rd</u>		Temp (C) <u>25.4</u>		pH <u>7.14</u>		Sample Depth <u>0.1</u>		☐ ft ☒ m	
Latitude		☐ dd ☐ dms		Salinity (ppt) <u>0.10</u>		Comments					
Location		☐ Comp ☐ Grab		Collection (comp begin or grab) Date Time		☐ Eastern ☐ Central		Composite end Date Time		Bottle Group(s)	
Field ID				Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)	
WIN/STORET #				Temp (C)		pH		Sample Depth		☐ ft ☐ m	
Latitude		☐ dd ☐ dms		Longitude		☐ dd ☐ dms		Salinity (ppt)		Comments	
Location		☐ Comp ☐ Grab		Collection (comp begin or grab) Date Time		☐ Eastern ☐ Central		Composite end Date Time		Bottle Group(s)	
Field ID				Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)	
WIN/STORET #				Temp (C)		pH		Sample Depth		☐ ft ☐ m	
Latitude		☐ dd ☐ dms		Longitude		☐ dd ☐ dms		Salinity (ppt)		Comments	
Location		☐ Comp ☐ Grab		Collection (comp begin or grab) Date Time		☐ Eastern ☐ Central		Composite end Date Time		Bottle Group(s)	
Field ID				Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L ☐ % sat		Total Res. Chlorine (mg/L)	
WIN/STORET #				Temp (C)		pH		Sample Depth		☐ ft ☐ m	
Latitude		☐ dd ☐ dms		Longitude		☐ dd ☐ dms		Salinity (ppt)		Comments	

Relinquished By: <u>JB, JW</u>	Date/Time <u>8/5/19 1427</u>	Shipping Method <u>HO</u>	Number of Coolers Shipped: <u>1</u>	Received By: 	Date/Time <u>8/5/19 1427</u>
-----------------------------------	---------------------------------	------------------------------	--	---	---------------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-BR-IC W-F W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WO-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ENT-24-QT					
Group: A	Bottle Type:	QPCR-P-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	PCR-FILTER					
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-E8321-DI W-E8321-MS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC					

Group: B	Bottle Type:	P-1L	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-TDS					
	W-TSS					
Group: B	Bottle Type:	BP-1L	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: B	Bottle Type:	P-250ML-WO-THI	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ECOLI-18QT					
Group: B	Bottle Type:	QPCR-P-500ML	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	PCR-FILTER					
Group: B	Bottle Type:	BG-500ML	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-E8321-DI					
	W-E8321-MS					
Group: B	Bottle Type:	P-500ML	# of Bottles: 3	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3					
	W-NO2NO3					
	W-S-A-TP					
	W-TKN					
	W-TOC					
Group: B	Bottle Type:	P-125ML	# of Bottles: 3	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: C	Bottle Type:	P-1L	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					
Group: C	Bottle Type:	BP-1L	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: C	Bottle Type:	P-250ML-WO-THI	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____

Group: C	Bottle Type: P-250ML-WO-THI ECOLI-18QT	# of Bottles: 3	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: P-500ML W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	# of Bottles: 3	Preservative: ICE And H2SO4 9010036 41/10/20	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: P-125ML W-PO4-F	# of Bottles: 3	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
Group: D	Bottle Type: P-1L ALKALINITY TURBIDITY W-BR-IC W-F W-TSS	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	
Group: D	Bottle Type: BP-1L CHLSUITE-W	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	
Group: D	Bottle Type: P-250ML-WO-THI ENT-24-QT	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	
Group: D	Bottle Type: P-500ML W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	
Group: D	Bottle Type: P-125ML W-PO4-F	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	
Group: E	Bottle Type: P-1L ALKALINITY TURBIDITY W-BR-IC W-F	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	

Group: E	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab
W-TSS				
Group: E	Bottle Type: P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab
ENT-24-QT				
Group: E	Bottle Type: QPCR-P-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab
PCR-FILTER				
Group: E	Bottle Type: P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab
W-ICP				
W-ICPMS				
Group: E	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab
W-NH3				
W-NO2NO3				
W-S-A-TP				
W-TKN				
W-TOC				
Group: E	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab
W-PO4-F				

Date of Request: 10-JUN-2019
 Created By: ECKLES_M On: 10-JUN-2019
 Modified By: JASROTIA_P On: 10-JUN-2019
 Customer: TLH-ROC
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Not Applicable
 Sampling Event: Run 3

Send Coolers To:

Phone: 850-245-7586
 2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles

Send Final Report To:

2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles

Program Module Number: 0100
 Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 10-JUN-2019
 Biology Request Reviewed By: JASROTIA_P On: 10-JUN-2019
 Sampling Kit Required: YES Ship on: 17-JUN-2019
 Sampling Kit Shipped: YES On: 17-JUN-2019
 Sampling Kit Packed By: SHAIK_A
 Preservatives Needed: NO
 Date to Receive Samples: 17-JUN-2019
 Received By: KLUG_K On: 02-JUL-2019

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Group A WBID 3310A & DUP

 Group B WBID 3512, 965, & 977

 Group C WBID 3518, 919, & 3419A

 Goup D WBID 8029

 Group E Blank

 No fedex labels

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 2	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: QPCR-P-500ML	Number of Bottles: 2	Preserved With: ICE
PCR-FILTER	Template: DEFAULT	(SOP-PCR-4.0) Preparation of Samples by Filtration and held for qPCR Analysis
Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template:	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: HNO3
W-ICP	Template: S-METALS	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: HNO3
W-ICP	Template: S-METALS	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron		
Manganese		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 3 Samples:

Bottle Type: BP-1L	Number of Bottles: 3	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 3	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: QPCR-P-500ML	Number of Bottles: 3	Preserved With: ICE
PCR-FILTER	Template: DEFAULT	(SOP-PCR-4.0) Preparation of Samples by Filtration and held for qPCR Analysis
Bottle Type: BG-500ML	Number of Bottles: 3	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template:	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 3	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 3	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Group C (Water) With 3 Samples:

Bottle Group C (Water) With 3 Samples:

Bottle Type: P-1L	Number of Bottles: 3	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 3	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 3	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 3	Preserved With: FILTER-ICE
W-PO ₄ -F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO ₄ -F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-FILTER	Template: DEFAULT	(SOP-PCR-4.0) Preparation of Samples by Filtration and held for qPCR Analysis

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects	
Barium			
Calcium			
Iron			
Manganese			
Potassium			
Sodium			
Zinc			
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects	
Aluminum			
Antimony			
Arsenic			
Cadmium			
Chromium			
Copper			
Lead			
Nickel			
Selenium			
Silver			
Thallium			
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L	
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L	
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L	
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices	
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices	
Bottle Type: P-125ML		Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L	

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.