

Log-in Checklist

RQ-2019-08-12-34

Login Station ID: #1

Shipping Method: FedEx | UPS | HD Greyhound

Date/Time of Receipt: 8/13/19 1430

Cooler Temperature*	Samples Frozen YES	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
		HNO ₃	Formalin		Present?		*Intact?		
					Yes	No	Yes	No	
2.4				13		✓			

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C) without login confirmation from customer for samples requiring ice preservation. An NCR still is required.

**Chain of Custody/ Field Sheet(s) Included? Yes ✓ No

Micro-Biology Overflow Chain of Custody/ FieldSheet(s) Included? Yes No ✓

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

CONTAINER CHECK

**Evidence Tape on Bottles? Yes No ✓ If yes, is it intact? Yes No

**Caps on tight? Yes ✓ No **Containers intact? Yes ✓ No

**Sufficient Sample Volume: Yes No

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes No ✓ Init:

Chlorine detected? (One container checked per Field ID) Yes No Init:

If "YES", list affected Field IDs and Test IDs below.

Field ID	List TestID(s) where Chlorine was Detected, an NCR is required.

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ✓ No NA Init:

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes No NA ✓ Init:

Coolers Unpacked/Checked by: Date: 8/13/19 NCRs (Y/N)?

Event ID: Run 2
TLH-ROC-2019-08-13-01 (SMP) NCR #

Log-in Checklist

Samples with Incorrect pH Preservation

Field ID	Bottle Test ID(s)	pH	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_01_2018 EXP: 04/29/2020
pH Test Strips 0 – 6 Lot #	213316AV EXP: 05/15/2020
pH Paper 0 – 3 Lot #	
pH Paper 0 – 13 Lot #	215517 EXP: 06/01/2020
Chlorine Test Strips Lot #	8261 EXP: 07/30/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ✓

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF; W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Run 2

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: TLH-ROC

Requester: Michael Eckles

Field Report Prepared By:

Blair, Harden

Project ID: SMP

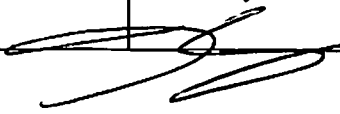
Collected By:

Blair, Harden

Sampling Agency:

8034

Field/WIN ID →		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08/13/2019 Time 1025	Eastern Central	Composite end Date	Time	Bottle Group(s) (See pg 2)	
Location ↓	Steinhatchee River @ S Canal Rd	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		C	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Temp (C)	pH	Sample Depth		Sp. Conductance (umho/cm)
				Salinity (ppt)	Comments			
Field/WIN ID →		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 08/13/2019 Time 1045	Eastern Central	Composite end Date	Time	Bottle Group(s)	
Location ↓	Steinhatchee River 40m N of Yearling Rd	Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		B	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Temp (C)	pH	Sample Depth		Sp. Conductance (umho/cm)
				Salinity (ppt)	Comments			
Location		☐ Comp ☐ Grab	Collection (comp begin or grab) Date	Time	Eastern Central	Composite end Date	Time	Bottle Group(s)
Field ID		Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)			
WIN/STORET #		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location		☐ Comp ☐ Grab	Collection (comp begin or grab) Date	Time	Eastern Central	Composite end Date	Time	Bottle Group(s)
Field ID		Matrix (type, e.g., Salt, Fresh, etc)	Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)			
WIN/STORET #		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			

Relinquished By:	Date/Time	Shipping Method	Number of Coolers Shipped:	Received By:	Date/Time
Harden/Blair	08/13/19	Hand Delivered			8/13/19 1430

Group: A	Bottle Type:	P-1L	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-BR-IC W-F W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WO-THI	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ENT-24-QT					
Group: A	Bottle Type:	P-250ML-WO-THI	# of Bottles: 3	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	FCOLI-MF					
Group: A	Bottle Type:	P-500ML	# of Bottles: 3	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 3	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 3	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: B	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____

Group: B	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab _____
Group: B	Bottle Type: ECOLI-18QT	P-250ML-WO-THI	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab _____
Group: B	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 2	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
Group: B	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 2	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab _____
Group: C	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab _____
Group: C	Bottle Type: ECOLI-18QT	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab _____
Group: C	Bottle Type: W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab _____
Group: C	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
Group: C	Bottle Type:	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab _____

Group: C	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab _____
Group: D	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab _____
Group: D	Bottle Type: ECOL-18QT	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab _____
Group: D	Bottle Type: W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab _____
Group: D	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
Group: D	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab _____

Date of Request: 25-JUL-2019
 Created By: ECKLES_M On: 25-JUL-2019
 Modified By: JASROTIA_P On: 29-JUL-2019
 Customer: TLH-ROC
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Not Applicable
 Sampling Event: Run 2

Send Coolers To:

Phone: 850-245-7586
 2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles

Send Final Report To:

2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles

Program Module Number: 0100
 Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 25-JUL-2019
 Biology Request Reviewed By: JASROTIA_P On: 29-JUL-2019
 Sampling Kit Required: YES Ship on: 08-AUG-2019
 Sampling Kit Shipped: YES On: 05-AUG-2019
 Sampling Kit Packed By: KLUG_K
 Preservatives Needed: NO
 Date to Receive Samples: 12-AUG-2019
 Received By: SHAIK_A On: 13-AUG-2019

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Group A WBID 3699A 2 sites & DUP
 Group B WBID 3556 & 3573B
 Group C WBID 3573A
 Group D Blank
 No fedex labels

Bottle Group A (Water) With 3 Samples:

Bottle Type: P-1L	Number of Bottles: 3	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 3	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 3	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-250ML-WO-THI	Number of Bottles: 3	Preserved With: ICE
FCOLI-MF	Template: DEFAULT	(SM 9222 D) Fecal coliforms by membrane filter method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: HNO ₃
W-ICP	Template: S-METALS	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron		
Manganese		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Bottle Group A (Water) With 3 Samples:

Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: HNO3
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Aluminum
Antimony
Arsenic
Cadmium
Chromium
Copper
Lead
Nickel
Selenium
Silver
Thallium

Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 3	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 2 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 2	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L W-TSS	Template: DEFAULT	Number of Bottles: 1 (SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample	Preserved With: ICE
Bottle Type: P-250ML-WO-THI ECOLI-18QT	Template: DEFAULT	Number of Bottles: 1 (SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems	Preserved With: ICE
Bottle Type: P-500ML W-ICP	Template:	Number of Bottles: 1 (EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects	Preserved With: HNO3
Barium Calcium Iron Magnesium Potassium Sodium Zinc			
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects	
*Boron Aluminum Antimony Arsenic Cadmium Chromium Copper Lead Nickel Selenium Silver Thallium			
Bottle Type: P-500ML W-NH3	Template: DEFAULT	Number of Bottles: 1 (EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L	Preserved With: ICE And H2SO4
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L	
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L	
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices	
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices	
Bottle Type: P-125ML W-PO4-F	Template: DEFAULT	Number of Bottles: 1 (EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L	Preserved With: FILTER-ICE

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L ALKALINITY	Template: DEFAULT	Number of Bottles: 1 (SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L	Preserved With: ICE
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices	
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices	
W-COLOR Color (true)	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm	
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)	
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices	
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices	
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample	
Bottle Type: P-250ML-WO-THI ECOLI-18QT	Template: DEFAULT	Number of Bottles: 1 (SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems	Preserved With: ICE
Bottle Type: P-500ML W-ICP	Template:	Number of Bottles: 1 (EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects	Preserved With: HNO3
Barium Calcium Iron Magnesium			

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects	
Potassium			
Sodium			
Zinc			
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects	
*Boron			
Aluminum			
Antimony			
Arsenic			
Cadmium			
Chromium			
Copper			
Lead			
Nickel			
Selenium			
Silver			
Thallium			
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L	
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L	
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L	
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices	
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices	
Bottle Type: P-125ML		Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L	

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.