

Log-in Checklist

RQ- 2019-08-19-24

Login Station ID: #5

Shipping Method: ~~Next~~ | UPS ~~(HD)~~ Greyhound

Date/Time of Receipt: 9/ 5 /19 @16:03

Cooler Temperature*	Samples Frozen	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
1.1				22		✓			
2.4				13		✓			

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. **DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C) without login confirmation from customer for samples requiring ice preservation. An NCR still is required.**

**Chain of Custody/ Field Sheet(s) Included? Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☐ No ☒

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

CONTAINER CHECK

**Evidence Tape on Bottles? Yes ☐ No ☒ If yes, is it intact? Yes ☐ No ☐

**Caps on tight? Yes ☒ No ☐ **Containers intact? Yes ☒ No ☐

**Sufficient Sample Volume: Yes ☒ No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ☐ No ☒ Init: KK

Chlorine detected? (One container checked per Field ID)

Yes ☐ No ☐ Init:

If "YES", list affected Field IDs and Test IDs below.

Field ID	List TestID(s) where Chlorine was Detected, an NCR is required.

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ☒ No ☐ NA ☐ Init: KK

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes ☐ No ☐ NA ☒ Init: KK

Coolers Unpacked/Checked by: KK Date: 9/ 5 /19 NCRs (Y/N)? N

Event ID: Run 12
TLH-ROC-2019-09-05-02 (SMP)
 Date Received: 05-SEP-2019 16:03 NCR #

Log-in Checklist

Samples with Incorrect pH Preservation

Field ID	Bottle TestID(s)	pH	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle TestID(s)

Bottles Not Intact

Field ID	Bottle TestID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle TestID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	REC_07_2018 EXP:12/28/19
pH Test Strips 0 – 6 Lot #	
pH Paper 0 – 3 Lot #	220416A EXP:7/30/20
pH Paper 0 – 13 Lot #	222516 EXP:9/15/19
Chlorine Test Strips Lot #	091417V EXP:9/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No X _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF; W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: TLH-ROC

Requester: Michael Eckles

Field Report Prepared By:

J. Woods

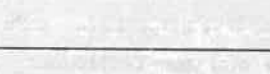
Project ID: SMP

Collected By:

Woods/Spencer

Sampling Agency:

8034

Location		☐ Comp	Collection (comp begin or grab)		☒ Eastern	Composite end		Bottle Group(s) (See pg 2)	
Field ID	Field/WIN ID →  G2TLHR0029	☒ Grab	Date 9/5/19	Time 1000	Central	Date	Time		
WIN/S	Location ↓  DICKSON BAY SOUTH	Matrix (type, e.g., Salt/Fresh, etc)		Diss. Oxygen 6.02	☒ mg/L	Total Res. Chlorine (mg/L)			
Latitude	☐ dms	Temp (C) 28.7	pH 7.97	Sample Depth 0.3	☐ ft	Sp. Conductance (umho/cm) 28604			
		☐ dd	Salinity (ppt) 17.55	Comments					C
Location		☐ Comp	Collection (comp begin or grab)		☒ Eastern	Composite end			
Field ID	Field/WIN ID →  G2TLHR0028	☒ Grab	Date 9/5/19	Time 1015	Central	Date	Time		
WIN/STOR	Location ↓  DICKSON BAY NORTH	Matrix (type, e.g., Salt/Fresh, etc)		Diss. Oxygen 5.90	☒ mg/L	Total Res. Chlorine (mg/L)			
Latitude	☐ dms	Temp (C) 28.7	pH 7.95	Sample Depth 0.3	☐ ft	Sp. Conductance (umho/cm) 32996		C	
		☐ dd	Salinity (ppt) 20.56	Comments					
Location		☐ Comp	Collection (comp begin or grab)		☒ Eastern	Composite end			
Field ID	Field/WIN ID →  G2TLHR0025	☒ Grab	Date 9/5/19	Time 1110	Central	Date	Time		
WIN/STOR	Location ↓  Alligator Harbor near clam beds	Matrix (type, e.g., Salt/Fresh, etc)		Diss. Oxygen 6.00	☒ mg/L	Total Res. Chlorine (mg/L)		C	
Latitude	☐ dms	Temp (C) 28.8	pH 7.97	Sample Depth 0.3	☐ ft	Sp. Conductance (umho/cm) 46864			
		☐ dd	Salinity (ppt) 30.36	Comments					
Location		☐ Comp	Collection (comp begin or grab)		☒ Eastern	Composite end			
Field ID	Field/WIN ID →  G2TLHR0025	☒ Grab	Date 9/5/19	Time 1225	Central	Date	Time	A	
WIN/STOR	Location ↓  Alligator Harbor near clam beds	Matrix (type, e.g., Salt/Fresh, etc)		Diss. Oxygen 4.97	☒ mg/L	Total Res. Chlorine (mg/L)			
Latitude	☐ dms	Temp (C) 26.0	pH 7.97	Sample Depth 0.15	☐ ft	Sp. Conductance (umho/cm) 42			
		☐ dd	Salinity (ppt) 0.02	Comments					

Relinquished By:	Date/Time	Shipping Method	Number of Coolers Shipped:	Received By:	Date/Time
JW	9/5/19	HD	2	[Signature]	9/5/19 1603

2. row, 3. row, 3 ENT
9/5/19 @ 16:54 (A)

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	<u>1</u>
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	<u>X</u>
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	<u>207</u> <u>2</u>
	ECOLI-18QT					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab	<u>X</u>
	W-ICP W-ICPMS 9129070					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	<u>X</u>
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC 9010030					
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	<u>X</u>
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	<u>X</u>
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					

Group: B	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab
CHLSUITE-W				
Group: B	Bottle Type: P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab
ECOLI-18QT				
Group: B	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC				
Group: B	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab
W-PO4-F				
Group: C	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>3</u>
ALKALINITY TURBIDITY W-BR-IC W-F W-TSS				
Group: C	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>3</u>
CHLSUITE-W				
Group: C	Bottle Type: P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>3 6</u>
ENT-24-QT				
Group: C	Bottle Type: P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>3</u>
FCOLI-MF				
Group: C	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab <u>3</u>
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC				
Group: C	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab <u>3</u>
W-PO4-F				

9010030
9/5/19

Run 12

Central Laboratory Submittal Form

last revised Jan 25, 2018

Customer: TLH-ROC

Requester: Michael Eckles

Field Report Prepared By: J. Woods

Project ID: SMP

Collected By: Woods / spences

Sampling Agency:

Location		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 9/5/19 Time 1225		Eastern Central	Composite end Date	Time	Bottle Group(s) (See pg 2)
Field ID	Field/WIN ID → G2TLHR0006	Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen 4.97	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		A
WIN/ST	Location ↓ Gully Branch 80m north of Gully Branch Rd	Temp (C) 26.0	pH 4.50	Sample Depth 0.15	☐ ft ☒ m	Sp. Conductance (umho/cm) 42		
Latitude	☐ dms	☐ dd ☐ dms Salinity (ppt) 0.02		Comments				
Location		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 9/5/19 Time 1315		Eastern Central	Composite end Date	Time	Bottle Group(s)
Field	Field/WIN ID → G2TLHR0038	Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen 3.37	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		B
WIN/	Location ↓ CROOKED RIVER @ POST 46	Temp (C) 29.9	pH 5.64	Sample Depth 0.6	☐ ft ☒ m	Sp. Conductance (umho/cm) 113		
Latitude	☐ dms	☐ dd ☐ dms Salinity (ppt) 0.05		Comments				
Location		☐ Comp ☐ Grab	Collection (comp begin or grab) Date		Eastern Central	Composite end Date	Time	Bottle Group(s)
Field ID		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		
WIN/STORET #		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location		☐ Comp ☐ Grab	Collection (comp begin or grab) Date		Eastern Central	Composite end Date	Time	Bottle Group(s)
Field ID		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)		
WIN/STORET #		Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)		
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			

Relinquished By: JW	Date/Time 9/5/19	Shipping Method HD	Number of Coolers Shipped: 1	Received By: [Signature]	Date/Time 9/5/19 1603
------------------------	---------------------	-----------------------	---------------------------------	-----------------------------	--------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	CHLSUITE-W						
Group: A	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	ECOLI-18QT						
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-ICP W-ICPMS 9.29070 exp. 4/10/20						
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC 9010030 exp. 5/9/20						
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
	W-PO4-F						
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS						

Group: B	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: ECOLI-18QT	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	2
Group: B	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4 9010030 exp. 5/9/20	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: ALKALINITY TURBIDITY W-BR-IC W-F W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	
Group: C	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	
Group: C	Bottle Type: ENT-24-QT	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	
Group: C	Bottle Type: FCOLI-MF	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	
Group: C	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	
Group: C	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	

Date of Request: 25-JUL-2019
 Created By: ECKLES_M On: 25-JUL-2019
 Modified By: JASROTIA_P On: 31-JUL-2019
 Customer: TLH-ROC
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Not Applicable
 Sampling Event: Run 12

Send Coolers To:

Phone: 850-245-7586
 2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles

Send Final Report To:

2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles

Program Module Number: 0100
 Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 30-JUL-2019
 Biology Request Reviewed By: JASROTIA_P On: 31-JUL-2019
 Sampling Kit Required: YES Ship on: 15-AUG-2019
 Sampling Kit Shipped: YES On: 15-AUG-2019
 Sampling Kit Packed By: REYNOLDS_T
 Preservatives Needed: NO
 Date to Receive Samples: 19-AUG-2019
 Received By:

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Group A WBID 1242
 Group B WBID 1251
 Group C WBID 1256

No fedex labels

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L			Number of Bottles: 1	Preserved With: ICE
W-BR-IC	Template: DEFAULT	(EPA 300.0) Bromide in aqueous matrices		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)		
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample		
Bottle Type: BP-1L			Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry		
Bottle Type: P-250ML-WO-THI			Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems		
Bottle Type: P-250ML-WO-THI			Number of Bottles: 1	Preserved With: ICE
FCOLI-MF	Template: DEFAULT	(SM 9222 D) Fecal coliforms by membrane filter method - non-chlorinated water systems		
Bottle Type: P-500ML			Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L		
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L		
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L		
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices		
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices		
Bottle Type: P-125ML			Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L		

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.