

RQ- 2020-03-23-08

Login Checklist

Login Station ID: #2

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 3/26/2020 @ 9:45

Cooler Temperature*	Samples Frozen	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
		HNO ₃	Formalin		Present?		*Intact?		
	YES				Yes	No	Yes	No	
3.8				14		✓			127840127909

*All samples received on ice unless otherwise noted. If the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), received without ice, or if Evidence Seal is damaged; identify the containers in the affected cooler(s) on back of this form. HNO₃ and Formalin preserved samples do **NOT** have a temperature requirement. **DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C) without login confirmation from customer for samples requiring ice preservation; An NCR still is required.**

**Chain of Custody/ Field Sheet(s) Included?

Yes ✓

No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ✓

No _____

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

CONTAINER CHECK

**Evidence Tape on Bottles? Yes ✓

No _____

If yes, is it intact? Yes ✓

No _____

**Caps on tight? Yes ✓

No _____

**Containers intact? Yes ✓

No _____

**Sufficient Sample Volume: Yes ✓

No _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ✓ No _____

Init: IM

Chlorine detected? (One container checked per Field ID)

Yes _____ No ✓

Init: IM

If "YES", list affected Field IDs and Test IDs below.

Field ID	List TestID(s) where Chlorine was Detected, an NCR is required.

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0?

Yes ✓

No _____

NA _____

Init: IM

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes _____

No _____

NA ✓

Init: IM

Coolers Unpacked/Checked by: IM

Date: 3/26/2020

NCRs (Y/N)?

N

Event ID: LSJR West 2020
NE-DIST-2020-03-26-01 (SMP)

NCR #

Date Received: 26-MAR-2020 09:45

Login Checklist

Samples with Incorrect pH Preservation

Field ID	Bottle Test ID(s)	pH	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	IR TEMP GUN #2 exp:11/5/20
pH Test Strips 0 – 6 Lot #	213316AV Exp:5/15/20
pH Paper 0 – 3 Lot #	
pH Paper 0 – 13 Lot #	208717 EXP: 3/31/20
Chlorine Test Strips Lot #	8261 EXP:7/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes_____ No_____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes_____ If No_____

no, were samples shipped on dry ice? Yes_____ No_____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes_____ No_____

Date and time samples preserved: _____

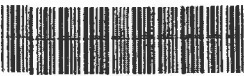
Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: NE-DIST

Project ID: SMP

Requester: Jeremy M. ParrishCollected By: K. McDonaldField Report Prepared By: A. KalmbacherSampling Agency: FDEP/NE ROC

WIN ID/Field ID:  20030549	☐ Comp ☒ Grab Collection (comp begin or grab) Date <u>03-25-2020</u> Time <u>0920</u>	☒ Eastern ☐ Central Composite end Date _____ Time _____	Bottle Group(s) (See pg 2) C
YELLOW WATER CR AT SR 228	Matrix (type, e.g., Salt, <u>Fresh</u> , etc)	Diss. Oxygen <u>1.49</u> ☒ mg/L ☐ % sat Total Res. Chlorine (mg/L)	
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms	Temp (C) <u>20.5</u> pH <u>6.34</u> Sample Depth <u>0.3</u> ☐ ft ☒ m Sp. Conductance (umho/cm) <u>149.6</u>		
Salinity (ppt) <u>0.07</u>	Comments		
WIN ID/Field ID:  20030751	☐ Comp ☒ Grab Collection (comp begin or grab) Date <u>03-25-2020</u> Time <u>1010</u>	☒ Eastern ☐ Central Composite end Date _____ Time _____	Bottle Group(s) B
TROUT RIVER AT GARDEN ST	Matrix (type, e.g., Salt, <u>Fresh</u> , etc)	Diss. Oxygen <u>6.5</u> ☒ mg/L ☐ % sat Total Res. Chlorine (mg/L)	
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms	Temp (C) <u>20.6</u> pH <u>6.36</u> Sample Depth <u>0.05</u> ☐ ft ☒ m Sp. Conductance (umho/cm) <u>136</u>		
Salinity (ppt) <u>0.06</u>	Comments		
WIN ID/Field ID:  G2NE1183	☐ Comp ☒ Grab Collection (comp begin or grab) Date <u>03-25-2020</u> Time <u>1100</u>	☒ Eastern ☐ Central Composite end Date _____ Time _____	Bottle Group(s) D
WILLS BR at Watgate rd.	Matrix (type, e.g., Salt, <u>Fresh</u> , etc)	Diss. Oxygen <u>1.83</u> ☒ mg/L ☐ % sat Total Res. Chlorine (mg/L)	
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms	Temp (C) <u>24.2</u> pH <u>6.89</u> Sample Depth <u>0.3</u> ☐ ft ☒ m Sp. Conductance (umho/cm) <u>3898</u>		
Salinity (ppt) <u>2.07</u>	Comments		
Location: WILLS BR at Watgate rd.	☐ Comp ☒ Grab Collection (comp begin or grab) Date <u>03-25-2020</u> Time <u>1105</u>	☒ Eastern ☐ Central Composite end Date _____ Time _____	Bottle Group(s) E
FIELD ID:  FB 03252020	Matrix (type, e.g., Salt, <u>Fresh</u> , etc)	Diss. Oxygen ☐ mg/L ☐ % sat Total Res. Chlorine (mg/L)	
Latitude ☐ dd ☐ dms Longitude ☐ dd ☐ dms	Temp (C) _____ pH _____ Sample Depth ☐ ft ☐ m Sp. Conductance (umho/cm) _____		
Salinity (ppt) _____	Comments: Field Blank		

Relinquished By: A. Kalmbacher

Date/Time

03-25-2020

Shipping Method

FedEx

Number of Coolers Shipped:

1Received By: IM

Date/Time

3/26/20 9:45

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	CHLSUITE-W						
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	NE-ALS-ECQ						
Group: A	Bottle Type:	QPCR-P-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	PCR-BACR						
	PCR-DG3						
	PCR-HF183						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	W-E8321-DI						
	W-E8321-MS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	0
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	0
	W-PO4-F						
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						

Group: B	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
W-F W-SO4-IC W-TDS W-TSS					
Group: B	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1 + 0
CHLSUITE-W					
Group: B	Bottle Type: P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1 Sent to ALS lab
NE-ALS-ECQ					
Group: B	Bottle Type: BG-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
W-E8321-DI W-E8321-MS					
Group: B	Bottle Type: P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab	1
W-HARD W-ICP W-ICPMS					
Group: B	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: B	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	0
W-PO4-F					
Group: C	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					

Group: C	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>AK-10</u>
Group: C	Bottle Type: W-HARD W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab <u>1</u>
Group: C	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab <u>1</u>
Group: C	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab <u>0</u>
Group: D	Bottle Type: ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: D	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>0</u>
Group: D	Bottle Type: NE-ALS-ECQ	P-250ML-WI-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>1</u> <i>→ Sent to ALS lab</i>
Group: D	Bottle Type: W-E8321-DI W-E8321-MS	BG-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab <u>1</u>
Group: D	Bottle Type: W-NH3 W-NO2NO3	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab <u>1</u>

Group: D Bottle Type: P-500ML # of Bottles: 1 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 1

W-S-A-TP
W-TKN
W-TOC

Group: D Bottle Type: P-125ML # of Bottles: 1 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 6

W-PO4-F

Group: E Bottle Type: P-1L # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1

ALKALINITY
TURBIDITY
W-CL-IC
W-COLOR
W-F
W-SO4-IC
W-TDS
W-TSS

Group: E Bottle Type: BG-500ML # of Bottles: 1 Preservative: ICE Enter Number of Bottles Sent to Lab 1

W-E8321-DI
W-E8321-MS

Group: E Bottle Type: P-500ML # of Bottles: 1 Preservative: HNO3 Enter Number of Bottles Sent to Lab 1

W-HARD
W-ICP
W-ICPMS

Group: E Bottle Type: P-500ML # of Bottles: 1 Preservative: ICE And H2SO4 Enter Number of Bottles Sent to Lab 1

W-NH3
W-NO2NO3
W-S-A-TP
W-TKN
W-TOC

Group: E Bottle Type: P-125ML # of Bottles: 1 Preservative: FILTER-ICE Enter Number of Bottles Sent to Lab 0

W-PO4-F

Date of Request: 17-FEB-2020
 Created By: PARRISH_J On: 17-FEB-2020
 Modified By: JASROTIA_P On: 05-MAR-2020
 Customer: NE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Northeast District
 Sampling Event: LSJR West 2020

Send Coolers To:

Phone: 904-256-1693
 8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Send Final Report To:

8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Program Module Number:

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 05-MAR-2020
 Biology Request Reviewed By: JASROTIA_P On: 05-MAR-2020
 Sampling Kit Required: YES Ship on: 13-MAR-2020
 Sampling Kit Shipped: YES On: 12-MAR-2020
 Sampling Kit Packed By: MCCOY_I
 Preservatives Needed: YES
 Date to Receive Samples: 23-MAR-2020
 Received By:

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: Group A:(Surf.-Fresh) Nutrients, Bacteria, Tracers, Markers
 Group B:(Surf.-Fresh) Nutrients, Metals, Bacteria, tracers
 Group C: (Surf.-Fresh)Nutrients, Metals
 Group D:(Surf.-Fresh)Nutrients, Bacteria, Tracers
 Group E: (Blank) Nutrients, Metals, Tracers

NE-ALS-ECQ

Packing Notes:

1 - case of sulfuric acid/1 set of labels

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 2	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) E. coli by Colilert/QT by ALS in Jacksonville
Bottle Type: QPCR-P-500ML	Number of Bottles: 2	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction

Bottle Group A (Water) With 2 Samples:

Bottle Type: QPCR-P-500ML	Number of Bottles: 2	Preserved With: ICE
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction
Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) E. coli by Colilert/QT by ALS in Jacksonville
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) E. coli by Colilert/QT by ALS in Jacksonville

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO ₃
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Barium
Calcium
Iron
Magnesium
Potassium
Sodium
Zinc

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

*Boron
 Aluminum
 Antimony
 Arsenic
 Cadmium
 Chromium
 Copper
 Lead
 Nickel
 Selenium
 Silver
 Thallium

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.