

Login Checklist

RQ- ⁵⁶2020 - 07-02-71

Login Station ID: # 3

Shipping Method: ^{APR 8-11-20} FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 07-24-2020 / 10:45

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with:		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
		YES	HNO ₃		Formalin	Present?		*Intact?	
	Yes			No		Yes	No		
1.3C				APR-54		X			1278 4011 1108

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes X No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No X

CONTAINER CHECK

Evidence Tape on Bottles? Yes _____ No _____ NA X Evidence Tape Intact? Yes _____ No _____

If Criminal, photographs taken? Yes _____ No _____ Containers Intact? Yes X No _____

Caps on tight? Yes X No _____

Sufficient Sample Volume: Yes X No _____

If NO is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes _____ No X Init: ABS

Chlorine detected? (One container checked per Field ID)

Yes _____ No _____ Init: ABS

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0? Yes X No _____ NA _____ Init: ABS

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes _____ No _____ NA X Init: ABS

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: ABS

Date: 07-24-2020

Event contains NCRs (Y/N)? N

Event ID: CEN-DIST - 2020 - 07 - 24 - 0 2

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	ITG #3 EXP:11-05-2020
pH Test Strips 0 – 6	213316 BV EXP:05/15/2021
pH Paper 0 – 3	XXXXXXXXXX
pH Paper 0 – 13	231019 EXP:NOV-01-2022
Chlorine Test Strips	9080 / 2022 / 01

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: ABS

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Request Number: RQ-2020-07-06-71

Bearhead/Terrace/Theresa/Hourglass

Florida Department of Environmental Protection Central Laboratory Submittal Form

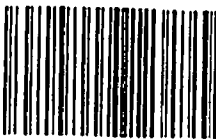
Page 1 of 2

last revised Jan 25, 2018

Customer: CEN-DIST

Project ID: SMP

Requester: Natalie AyalaField Report Prepared By: Katie MarchCollected By: Leif BermanSampling Agency: 212CEN

Location	WIN ID / Field ID → Hourglass Lake @ Center		 G4CE0178		☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>7/23/20</u> Time <u>1306</u>	☒ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s) (See pg 2) A	
Field ID					Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen <u>6.98</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L) _____	
WIN/STORE					Temp (C) <u>30.3</u>	pH <u>7.04</u>	Sample Depth <u>0.300</u>	☐ ft ☒ m	Sp. Conductance (umho/cm) <u>184.7</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.09</u>	Comments					
Location					☐ Comp ☐ Grab	Collection (comp begin or grab) Date _____ Time _____	☐ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID					Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORET #					Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments					
Location					☐ Comp ☐ Grab	Collection (comp begin or grab) Date _____ Time _____	☐ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID					Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORET #					Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments					
Location					☐ Comp ☐ Grab	Collection (comp begin or grab) Date _____ Time _____	☐ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID					Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORET #					Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments					

Relinquished By: <u>K March</u>	Date/Time: <u>7/23/20 1600</u>	Shipping Method: <u>FedEx</u>	Number of Coolers Shipped: <u>1</u>	Received By: <u>[Signature]</u>	Date/Time: <u>07-24-20 1145</u>
---------------------------------	--------------------------------	-------------------------------	-------------------------------------	---------------------------------	---------------------------------

Group: A

Bottle Type:

of Bottles: 0

Preservative:

Enter Number of Bottles Sent to Lab

1

ALKALINITY
CHLSUITE-W
TURBIDITY
W-CL-IC
W-COLOR
W-F
W-NH3
W-NO2NO3
W-PO4-F
W-S-A-TP
W-SO4-IC
W-TDS
W-TKN
W-TOC
W-TSS

Date of Request: 02-JUL-2020
 Created By: AYALA_N On: 02-JUL-2020
 Modified By: JASROTIA_P On: 06-JUL-2020
 Customer: CEN-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Central District
 Event Description: Bearhead/Terrace/Theresa/Hourglass

Send Coolers To:

Phone: 407-897-4195
 FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Natalia Ayala
 Cooler Ship Email: Natalia.Ayala@Floridadep.gov

Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 06-JUL-2020
 Biology Request Reviewed By: JASROTIA_P On: 06-JUL-2020
 Sampling Kit Required: NO Ship on:
 Sampling Kit Shipped:
 Sampling Kit Packed By:
 Preservatives Needed: NO
 Date to Receive Samples: 06-JUL-2020
 Logged In By: SHAIK_A On: 24-JUL-2020

Send Final Report To:

FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Natalia Ayala
 Report Notification Email: Natalia.Ayala@Floridadep.gov

Comments: All Bottle Group A

Bottle Group A (Water) With 4 Samples:

Bottle Type: P-1L	Number of Bottles: 4	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 4	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-500ML	Number of Bottles: 4	Preserved With: ICE And H₂SO₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 385.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 4	Preserved With: FILTER-ICE
W-PO ₄ -F	Template: DEFAULT	(EPA 385.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L