

RQ-2020-08-03-60

Login Checklist

Login Station ID: #4

Shipping Method: FedEx UPS | HD | Greyhound

Date/Time of Receipt: 08/06/2020 @ 10:05

*Cooler Temperature	Samples Frozen?	* All Samples In Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
2.2°C				71 49 9		x			1278 4015 7018

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No ☒

CONTAINER CHECK

Evidence Tape on Bottles?

Yes _____ No _____ NA ☒

Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No _____

Containers intact? Yes ☒ No _____

Caps on tight?

Yes ☒ No _____

Sufficient Sample Volume:

Yes ☒ No _____If NO is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes _____ No ☒ Init: Tyler

Chlorine detected? (One container checked per Field ID)

Yes _____ No _____ Init: _____

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes ☒ No _____ NA _____ Init: Tyler**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes _____ No _____ NA ☒ Init: Tyler

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: Tyler Reynolds

Date:

08/06/2020

Event contains NCRs (Y/N)? N

Event ID:

Ellen/Searcy

CEN-DIST-2020-08-06-01 (SMP)

Page 1 of 2, 0

Date Received: 06-AUG-2020 10:05

REF-4-6-3

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	#4 Temp Gun EXP:11/05/2020
pH Test Strips 0 – 6	
pH Paper 0 – 3	
pH Paper 0 – 13	231019 EXP:11/01/2022
Chlorine Test Strips	8261 EXP:07/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 8/5/2020 10:40
Customer: CEN-DIST
RQ: RQ-2020-08-03-60
Collected By: Leif Boman, Victoria Schwartz
Sampling Agency: FDEP
Project ID: SMP


PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 1

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
Lake Searcy @ Center	 G2CE0137
Lake Ellen @ Center	 G2CE0060

RELINQUISHED BY(SIGNATURE):

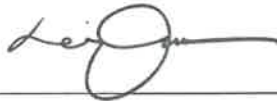

DATE/TIME: 08/05/2020 16:00

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By: Tyler Reynolds

Rec'd/Inspected Date: 08/06/2020

10:05



RQ-2020-2020-08-03-60

Customer: WQAP

Project ID: SMP

Collected By: Leif Boman, Victoria Schwartz

Field Report Prepared By: Leif Boman

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID:  G2CE0137										
Location: Lake Searcy @ Center										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPT)
8/5/2020	10:40 EST	5.53	74.7	31.1	7.02	0.3	0.9 ☐ VOB (S)	2.4	150.8	0.07
Central Lab please use top row for login purposes										
	10:46 EST	2.78	37.3	30.2	6.63	1.909			151.4	0.07
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☒ H2SO4*		☒ <2	1	B
		Lot #		Exp:						
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO3*		☐ <2		
		Lot #		Exp:						
Filtered Nutrients (P125ML)		☒ W-PO4-F ☒ Field Filtered 0.45 µm Filter				☒ Ice			1	B
		Filter Lot #								
		Syringe Lot #								
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	B
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CH-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	B
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice				
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☐ E Coli ☐ F Coli ☐ Enterococci ☐ Contract Lab				☐ Ice				
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☐ Ice				
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice				
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ W-E8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R ☐ W-PCL-SQ-R ☐ W-CARB-AA (5 ml of MCAA** Buffer Solution)				☐ Ice				
						☐ Ice				
Name: Leif Boman Victoria Schwartz					Signature:					Date:

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution



RQ-2020-2020-08-03-60

Customer: WQAP

Project ID: SMP

Collected By: Leif Boman, Victoria Schwartz

Field Report Prepared By: Leif Boman

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID:  G2CE0060										
Location: Lake Ellen @ Center										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPTH)
8/5/2020	11:25 EST	5.88	80.2	32	6.21	0.305	1.3 ☐ VOB (S)	6.8	98.1	0.04
Central Lab please use top row for login purposes										
	11:28 EST	0.1	1.1	19.4	6.03	6.311			121.3	0.06
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☐ H2SO4*		☒ <2	1	A
		Lot #		SA9354050		Exp:		12/20/2020		
Metals (P-500ML)		☒ W-ICPMS ☒ W-ICP ☒ W-HARD				☒ Ice ☒ HNO3*		☒ <2	1	A
		Lot #		NA9301040		Exp:		10/28/2020		
Filtered Nutrients (P125ML)		☒ W-PO4-F ☒ Field Filtered 0.45 µm Filter				☒ Ice			1	A
		Filter Lot #		16930280						
		Syringe Lot #		9149785						
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	A
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CH-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	A
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice				
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☐ E Coli ☐ F Coli ☐ Enterococci ☐ Contract Lab				☐ Ice				
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☐ Ice				
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice				
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ W-E8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R ☐ W-PCL-SQ-R ☐ W-CARB-AA (5 ml of MCAA** Buffer Solution)				☐ Ice				
						☐ Ice				
Name: Leif Boman Victoria Schwartz				Signature:				Date:		

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution

Date of Request: 31-JUL-2020
 Created By: BOMAN_L On: 31-JUL-2020
 Modified By: JASROTIA_P On: 03-AUG-2020
 Customer: CEN-DIST
 Project: SMP
 Division:
 District: Central District
 Event Description: Ellen/Searcy

Send Coolers To:

Phone: 407-897-4178
 FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Natalie Ayala
 Cooler Ship EMail: natalie.ayala@dep.state.fl.us

Priority: 3
 Event Status: A
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 03-AUG-2020
 Biology Request Reviewed By: JASROTIA_P On: 03-AUG-2020
 Sampling Kit Required: NO Ship on:
 Sampling Kit Shipped:
 Sampling Kit Packed By:
 Preservatives Needed: NO
 Date to Receive Samples: 03-AUG-2020
 Logged In By:

Send Final Report To:

FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Kalina Warren
 Report Notification EMail: natalie.ayala@dep.state.fl.us

Comments: Bottle Group A: Ellen
 Bottle Group B: Searcy

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
Total Alkalinity		
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
Turbidity		
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
Chloride		
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
Fluoride		
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
Sulfate		
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
TDS		
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
TSS		
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Chlorophyll-a, Corrected		
Chlorophyll-a, Uncorrected		
Phaeophytin-a		
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO ₃
Hardness		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

*Boron
Aluminum
Antimony
Arsenic
Cadmium
Chromium
Copper
Lead
Nickel
Selenium
Silver
Thallium

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
Ammonia-N		
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
NO2NO3-N		
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
Total-P		
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
Kjeldahl Nitrogen		
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Organic Carbon		

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
O-Phosphate-P		

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
Total Alkalinity		
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
Turbidity		
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
Chloride		
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
Fluoride		
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
Sulfate		
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
TDS		
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
TSS		
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Group B (Water) With 1 Samples:

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Chlorophyll-a, Corrected		
Chlorophyll-a, Uncorrected		
Phaeophytin-a		
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
Ammonia-N		
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
NO2NO3-N		
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
Total-P		
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
Kjeldahl Nitrogen		
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Organic Carbon		
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
O-Phosphate-P		

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.

