

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 11/5/2020 9:25

*Cooler Temperature	Samples Frozen?	* All Samples In Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
2.9				9		✓			929380068760

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO3 and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes ☐ No ☒

CONTAINER CHECK

Evidence Tape on Bottles?

Yes ☐ No ☐ NA ☒

Evidence Tape intact? Yes ☐ No ☐

If Criminal, photographs taken?

Yes ☐ No ☐

Containers intact? Yes ☒ No ☐

Caps on tight?

Yes ☒ No ☐

Sufficient Sample Volume:

Yes ☒ No ☐

If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ☐ No ☒ Init: IM

Chlorine detected? (One container checked per Field ID)

Yes ☐ No ☐ Init: _____

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes ☒ No ☐ NA ☐ Init: IM

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes ☐ No ☐ NA ☒ Init: IM

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: IM

Date: 11/5/2020

Event contains NCRs (Y/N)? ☒

Event ID: John's Lake Outlet/Minnie
CEN-DIST-2020-11-05-01 (SMP)

Page 1 of 2, 06, 11/5/2020, 09:25 Date Received: 05-NOV-2020 09:25

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	IR TEMP GUN #2 exp:10/28/21
pH Test Strips 0 – 6	
pH Paper 0 – 3	
pH Paper 0 – 13	231019 exp: 11/1/22
Chlorine Test Strips	8261 exp: 07/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 11/4/2020 **Collected By:** Terry Riordan, Leif Boman
Customer: CEN-DIST **Sampling Agency:** FDEP
RQ: RQ-2020-10-26-14 **Project ID:** SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 1

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
Blue Lake @ Center	 G2CE0232
Johns Lake @ Conservation Area Dock	 G1CE0111

RELINQUISHED BY(SIGNATURE):



DATE/TIME:

11/04/2020

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By:

IM

Rec'd/Inspected Date:

11/5/20



RQ-2020-10-26-14

Customer: WQAP

Project ID: SMP

Collected By: Terry Riordan, Leif Boman

Field Report Prepared By: Terry Riordan

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID:  G1CE0111										
Location: Johns Lake @ Conservation Area Dock										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPT)
11/4/2020	13:47 EST	6.22	72.9	23.3	6.91	0.304	0.6 ☐ VOB (S)	2.5	193.6	0.09
Central Lab please use top row for login purposes										
	13:50 EST	4.5	52.1	22.5	6.74	2.032			193.8	0.09
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☒ H2SO4*		☒ <2	1	A
		Lot #		SA9354050		Exp:		12/20/2020		
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO3*		☐ <2		
		Lot #				Exp:				
Filtered Nutrients (P125ML)		☒ W-P04-F ☒ Field Filtered 0.45 µm Filter				☒ Ice			1	A
		Filter Lot #		16930280						
		Syringe Lot #		9149785						
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	A
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CHC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	A
Physical/Aggregate (Marine) (P-1L)		☒ Alkalinity / W-CHC / W-TSS / W-TDS / W-COLOR				☒ Ice			1	A
Macroinvertebrates (PJ-2L)		☐ MI-FW-QUDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☐ E Coli ☐ F Coli ☐ Enterococci ☐ Contract Lab				☐ Ice				
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☐ Ice				
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice				
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ W-E8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R ☐ W-PCL-SQ-R ☐ W-CAR8-AA (5 ml of MCAA** Buffer Solution)				☐ Ice				
						☐ Ice				
Name: Terry Riordan Leif Boman					Signature: 					Date: 11/04/2020

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution



RQ-2020-10-26-14

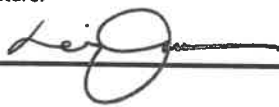
Customer: WQAP

Project ID: SMP

Collected By: Leif Boman, Terry Riordan

Field Report Prepared By: Leif Boman

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID: 										
G2CE0232										
Location: Blue Lake @ Center										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPT)
11/4/2020	11:46 EST	5.57	64.2	22.4	6.41	0.3	0.6 ☐ VOB (S)	2.9	127.7	0.06
↑ Central Lab please use top row for login purposes ↑										
	11:51 EST	5.33	61.5	22.3	6.4	2.479			127.7	0.06
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☒ H2SO4*		☒ <2	1	B
		Lot #		A9354050		Exp:		12/20/2020		
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO3*		☐ <2		
		Lot #				Exp:				
Filtered Nutrients (P125ML)		☒ W-PC4-F ☒ Field Filtered 0.45 µm Filter				☒ Ice			1	B
		Filter Lot #		16930280						
		Syringe Lot #		9149785						
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☐ Ice			1	B
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CH-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	B
Physical/Aggregate (Marine) (P-1L)		☒ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☒ Ice			1	B
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☐ E Coli ☐ F Coli ☐ Enterococci				☐ Ice				
		☐ Contract Lab								
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3				☐ Ice				
		☐ PCR-GULL2 ☐ PCR-GFD								
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice				
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ W-E8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R				☐ Ice				
		☐ W-PCL-SQ-R ☐ W-CARB-AA (5 ml of MCAA** Buffer Solution)				☐ Ice				
						☐ Ice				
Name: Leif Boman Terry Riordan				Signature: 				Date: 11/04/2020		

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution

Date of Request: 22-SEP-2020
 Created By: BOMAN_L On: 22-SEP-2020
 Modified By: JASROTIA_P On: 06-OCT-2020
 Customer: CEN-DIST
 Project: SMP
 Division:
 District: Central District
 Event Description: John's Lake Outlet/Minnie

Send Coolers To:

Phone: 407-897-4178
 FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Natalie Ayala
 Cooler Ship EMail: natalie.ayala@dep.state.fl.us

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 06-OCT-2020
 Biology Request Reviewed By: JASROTIA_P On: 06-OCT-2020
 Sampling Kit Required: YES Ship on: 16-OCT-2020
 Sampling Kit Shipped: YES On: 12-OCT-2020
 Sampling Kit Packed By: REYNOLDS_T
 Preservatives Needed: NO
 Date to Receive Samples: 26-OCT-2020
 Logged In By:

Send Final Report To:

FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Kalina Warren
 Report Notification EMail: natalie.ayala@dep.state.fl.us

Comments: Bottle Group A: Johns
 Bottle Group B: Minnie

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-1L		Number of Bottles: 1	Preserved With: ICE
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices	
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices	
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)	
Color (true)			
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)	
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices	
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices	
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample	
Bottle Type: BP-1L		Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry	
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L	
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L	
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L	
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices	
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices	
Bottle Type: P-125ML		Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L	