

Login Checklist

RQ- 2020-11-16-19

Login Station ID: #5

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 11/ 17 /20 @ 9:45

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
0.3				15		✓			9293 8008 1974

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No ☒

CONTAINER CHECK

Evidence Tape on Bottles?

Yes _____ No _____ NA ☒

Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No _____

Containers intact? Yes ☒ No _____

Caps on tight?

Yes ☒ No _____

Sufficient Sample Volume:

Yes ☒ No _____

If NO is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes _____ No ☒ Init: hkh

Chlorine detected? (One container checked per Field ID)

Yes _____ No _____ Init: _____

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes ☒ No _____ NA _____ Init: hkh

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes _____ No _____ NA ☒ Init: hkh

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: hkh

Date: 11/ 17 /20

Event contains NCRs (Y/N)? N

Event ID: SJR ABV Puzzle/Three Forks/Little Saw

CEN-DIST-2020-11-17-02 (SMP)

—Date Received: 17-NOV-2020 09:45

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	Temp Gun #5 EXP:10/28/2021
pH Test Strips 0 – 6	5/15/21 LOT#213316BV
pH Paper 0 – 3	
pH Paper 0 – 13	231019 EXP:11/01/2022
Chlorine Test Strips	8261 EXP:07/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

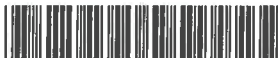
FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 11/16/2020 **Collected By:** Natalie Ayala, Victoria Schwartz
Customer: CEN-DIST **Sampling Agency:** FDEP
RQ: RQ-2020-11-16-19 **Project ID:** SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 1

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
Sjr At Inlet To Little Lake Sawgrass	 G3CE0017
Sjr Upstream Of Lake Hellen Blazes	 G3CE0018

RELINQUISHED BY(SIGNATURE):

Victoria Schwartz

DATE/TIME: 11/16/20 16:00

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By: *KR*

Rec'd/Inspected Date: *11-17-20/0945*



RQ-2020-11-16-19

Customer: WQAP

Project ID: SMP

Collected By: Natalie Ayala, Victoria Schwartz

Field Report Prepared By: Natalie Ayala

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID: 										
G3CE0017										
Location: Sjr At Inlet To Little Lake Sawgrass										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPT)
11/16/2020	13:36 EST	0.67	8.5	25.7	6.8	0.3	0.4 ☐ VOB (S)	1.6 (J)	418	0.2
Central Lab please use top row for login purposes										
	13:40 EST	0.51	6.3	25.6	6.54	1.107			407.8	0.19
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☒ H2SO4*		☒ <2	1	C
		Lot # 112544		Exp: 12/20/20						
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO3*		☐ <2		
		Lot #		Exp:						
Filtered Nutrients (P125ML)		☒ W-PO4-F ☒ Field Filtered 0.45 µm Filter				☒ Ice			1	C
		Filter Lot # 16930280								
		Syringe Lot # 9149785								
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	C
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CH-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	C
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice				
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☐ E Coli ☐ F Coli ☐ Enterococci ☐ Contract Lab				☐ Ice				
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☐ Ice				
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice				
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ WE8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R ☐ W-PCL-SQ-R ☐ W-CARB-AA (5 ml of MCAA** Buffer Solution)				☐ Ice				
						☐ Ice				
Name: Natalie Ayala Victoria Schwartz					Signature: <i>Victoria Schwartz</i>					Date: 11/16/20

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution



RQ-2020-11-16-19

Customer: WQAP

Project ID: SMP

Collected By: Natalie Ayala, Victoria Schwartz

Field Report Prepared By: Natalie Ayala

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID: 										
G3CE0018										
Location: Sjr Upstream Of Lake Hellen Blazes										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPTH)
11/16/2020	13:04 EST	0.27	3.3	25.7	6.52	0.302	0.5 ☐ VOB (S)	1.6 (J)	321.5	0.15
Central Lab please use top row for login purposes										
	13:15 EST	0.15	1.8	25.7	6.27	1.113			324.8	0.15
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☒ H2SO4*		☒ <2	1	B
		Lot #		112544		Exp:		12/20/20		
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO3*		☐ <2		
		Lot #				Exp:				
Filtered Nutrients (P125ML)		☒ W-PO4-F ☒ Field Filtered 0.45 µm Filter				☒ Ice			1	B
		Filter Lot #		16930280						
		Syringe Lot #		9149785						
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	B
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CH-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	B
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice				
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☒ E Coli ☐ F Coli ☐ Enterococci ☐ Contract Lab				☒ Ice			1	B
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☐ Ice				
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice				
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ W-E8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R ☐ W-PCL-SQ-R ☐ W-CARB-AA (5 ml of MCAA** Buffer Solution)				☐ Ice				
						☐ Ice				
Name: Natalie Ayala Victoria Schwartz				Signature: <i>Victoria Schwartz</i>				Date: 11/16/20		

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution

Date of Request: 23-OCT-2020
 Created By: BOMAN_L On: 23-OCT-2020
 Modified By: JASROTIA_P On: 26-OCT-2020
 Customer: CEN-DIST
 Project: SMP
 Division:
 District: Central District
 Event Description: SJR ABV Puzzle/Three Forks/Little Sawgrass

Send Coolers To:

Phone: 407-897-4178
 FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Natalie Ayala
 Cooler Ship EMail: natalie.ayala@dep.state.fl.us

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 26-OCT-2020
 Biology Request Reviewed By: JASROTIA_P On: 26-OCT-2020
 Sampling Kit Required: YES Ship on: 06-NOV-2020
 Sampling Kit Shipped: YES On: 02-NOV-2020
 Sampling Kit Packed By: SHAIK_A
 Preservatives Needed: NO
 Date to Receive Samples: 16-NOV-2020
 Logged In By:

Send Final Report To:

FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Kalina Warren
 Report Notification EMail: natalie.ayala@dep.state.fl.us

Comments: Bottle Group A: SJR
 Bottle Group B: Three Forks
 Bottle Group C: Little Sawgrass

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO ₃
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Beryllium		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L			Number of Bottles: 1	Preserved With: ICE
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)		
Color (true)				
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)		
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices		
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices		
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample		
Bottle Type: BP-1L			Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry		
Bottle Type: P-500ML			Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L		
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L		
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L		
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices		
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices		
Bottle Type: P-125ML			Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L		

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.