

RQ- 2020-07-06-52

Login Checklist

Login Station ID: #2

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 7/10/2020 10:50

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
1.7				11		✓			1278 4014 9075
3.6				16		✓			1278 4014 9053
6.6				14		✓			1278 4014 9031
3.9				6					1278 4014 9064

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# 8345
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes ✓ No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No _____

CONTAINER CHECK

Evidence Tape on Bottles?

Yes _____ No _____ NA ✓

Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No _____

Containers intact? Yes ✓ No _____

Caps on tight?

Yes ✓ No _____

Sufficient Sample Volume:

Yes ✓ No _____If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ✓ No _____ Init: IM

Chlorine detected? (One container checked per Field ID)

Yes _____ No ✓ Init: IM

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes ✓ No _____ NA _____ Init: IM

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes _____ No _____ NA ✓ Init: IM

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: IM

Date: 7/10/2020

Event contains NCRs (Y/N)? Y

Event ID:

2020 SMP : Myakka

SO-DIST-2020-07-10-01 (SMP)

Date Received: 10-JUL-2020 10:50

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Infrared Thermometer ID:	IR TEMP GUN #2 exp:11/5/20
pH Test Strips 0 – 6	
pH Paper 0 – 3	
pH Paper 0 – 13	231019 exp: 11/1/22
Chlorine Test Strips	8261 exp: 07/2021

Additional Comments:

Cooler 3 arrived with melted ice - IM 7/10/20

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: SO-DIST

Requester: Miranda CarrollField Report Prepared By: Miranda Carroll

Project ID: SMP

Collected By: Gary SnorekSampling Agency: FPEP SROC

Location	Field ID/WIN ID		☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>7/9/20</u> Time <u>0945</u>	☒ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s) (See pg 2)	
Field ID	Location	<u>G3SD0157</u>		Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>	Diss. Oxygen <u>6.55</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STOR				Temp (C) <u>27.8</u>	pH <u>7.29</u>	Sample Depth <u>0.1</u>	Sp. Conductance (umho/cm) <u>710</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.34</u>	Comments			<u>A</u>
Location	Field ID/WIN ID		☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>7/9/20</u> Time <u>1035</u>	☒ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID	Location	<u>G3SD0083</u>		Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>	Diss. Oxygen <u>6.95</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STC				Temp (C) <u>26.7</u>	pH <u>7.42</u>	Sample Depth <u>0.05</u>	Sp. Conductance (umho/cm) <u>1197</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.59</u>	Comments			<u>B</u>
Location	Field ID/WIN ID		☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>7/9/20</u> Time <u>1155</u>	☒ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID	Location	<u>G3SD0133</u>		Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>	Diss. Oxygen <u>5.50</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STC				Temp (C) <u>31.7</u>	pH <u>7.54</u>	Sample Depth <u>0.3</u>	Sp. Conductance (umho/cm) <u>386.4</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.18</u>	Comments			<u>A</u>
Location	Field ID/WIN ID		☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>7/9/20</u> Time <u>1310</u>	☒ Eastern ☐ Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID	Location	<u>G3SD0107</u>		Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>	Diss. Oxygen <u>5.13</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STOR				Temp (C) <u>32.3</u>	pH <u>7.28</u>	Sample Depth <u>0.3</u>	Sp. Conductance (umho/cm) <u>362.1</u>	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.17</u>	Comments			<u>C</u>

Relinquished By: [Signature]

Date/Time

7/9/20 1500

Shipping Method

FedEx

Number of Coolers Shipped:

4Received By: IM

Date/Time

7/10/20 10:50

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	CHLSUITE-W						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	CH2CICOOH Pre-Preserved	Enter Number of Bottles Sent to Lab	2
	W-CARB-AA						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 6	Preservative:	ICE	Enter Number of Bottles Sent to Lab	6
	W-CT-CI-R						
	W-PCL-SQ-R						
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	W-E8321-DI						
	W-E8321-MS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	2
	W-HARD						
	W-ICP						
	W-ICPMS						
					Lot # NA9354060		
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	2
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
					Lot # SA9301030		
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	2
	W-PO4-F						
Group: A	Bottle Type:	BG-1L	# of Bottles: 6	Preservative:	ICE	Enter Number of Bottles Sent to Lab	6
	W-PSNP-TQ						

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: SO-DIST

Requester: Miranda Carroll

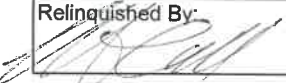
Field Report Prepared By:

Project ID: SMP

Collected By:

Sampling Agency:

Location	Field ID/WIN ID		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 7/9/20 Time 1405	☒ Eastern ☐ Central	Composite end Date Time	Bottle Group(s) (See pg 2)	
Field ID	Location	G3SD0158	Matrix (type, e.g., Salt, Fresh, etc) Fresh		Diss. Oxygen 7.11	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORE	Howard Creek NE Branch		Temp (C) 27.6	pH 7.27	Sample Depth 0.05	☐ ft ☒ m	Sp. Conductance (umho/cm) 333.7	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) 0.16	Comments			D
Location	Field ID	WIN/STORET #	☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	☐ Eastern ☐ Central	Composite end Date Time	Bottle Group(s)	
			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
			Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location	Field ID	WIN/STORET #	☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	☐ Eastern ☐ Central	Composite end Date Time	Bottle Group(s)	
			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
			Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			
Location	Field ID	WIN/STORET #	☐ Comp ☐ Grab	Collection (comp begin or grab) Date Time	☐ Eastern ☐ Central	Composite end Date Time	Bottle Group(s)	
			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
			Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			

Relinquished By:	Date/Time	Shipping Method	Number of Coolers Shipped:	Received By:	Date/Time
	7/9/20 1500	FedEx	4	IM	7/10/20 10:50

Group: A	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	CH2CICOOH Pre-Preserved	Enter Number of Bottles Sent to Lab _____
	W-CARB-AA					
Group: A	Bottle Type:	BG-500ML	# of Bottles: 6	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-CT-CI-R W-PCL-SQ-R					
Group: A	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-E8321-DI W-E8321-MS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-HARD W-ICP W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 2	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: A	Bottle Type:	BG-1L	# of Bottles: 6	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	W-PSNP-TQ					

Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS						
Group: B	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	CHLSUITE-W						
Group: B	Bottle Type:	QPCR-P-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	PCR-BACR PCR-DG3 PCR-HF183						
Group: B	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	SD-BES-ECQ						
Group: B	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	CH2CICOOH Pre-Preserved	Enter Number of Bottles Sent to Lab	1
	W-CARB-AA						
Group: B	Bottle Type:	BG-500ML	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
	W-CT-CI-R W-PCL-SQ-R						
Group: B	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	
	W-E8321-DI W-E8321-MS						
Group: B	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-HARD W-ICP W-ICPMS						
	Lot # NA 9354060						
Group: B	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3 W-NO2NO3 W-S-A-TP						
	Lot # SA 9301030						

Group: B	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-TKN W-TOC <i>See previous page</i>					
Group: B	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F					
Group: B	Bottle Type: BG-1L	# of Bottles: 2	Preservative: ICE	Enter Number of Bottles Sent to Lab	2
W-PSNP-TQ					
Group: C	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS					
Group: C	Bottle Type: BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W					
Group: C	Bottle Type: P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC <i>Lot #SA9 301030</i>					
Group: C	Bottle Type: P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F					
Group: D	Bottle Type: P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC					

Group: D	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	W-TDS						
	W-TSS						
Group: D	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	CHLSUITE-W						
Group: D	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	SD-BES-ECQ						
Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-HARD						
	W-ICP						
	W-ICPMS						
Group: D	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: D	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
	W-PO4-F						

Lot # NA 9354060

Lot # SA9301030

Date of Request: 23-JUN-2020
 Created By: CARROLL_M On: 23-JUN-2020
 Modified By: JASROTIA_P On: 24-JUN-2020
 Customer: SO-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: South District
 Event Description: 2020 SMP : Myakka

Send Coolers To:

Phone: 239-344-5718
 FLDEP South District Office
 2295 Victoria Ave. Suite 364
 Fort Myers, FL 33901-3381
 Attn: Miranda Carroll
 Cooler Ship EMail: Miranda.Carroll@FloridaDEP.gov

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 24-JUN-2020
 Biology Request Reviewed By: JASROTIA_P On: 24-JUN-2020
 Sampling Kit Required: YES Ship on: 29-JUN-2020
 Sampling Kit Shipped: YES On: 29-JUN-2020
 Sampling Kit Packed By: MCCOY_I
 Preservatives Needed: NO
 Date to Receive Samples: 06-JUL-2020
 Logged In By:

Send Final Report To:

FLDEP South District Office
 2295 Victoria Ave. Suite 364
 Fort Myers, FL 33901-3381
 Attn: Miranda Carroll
 Report Notification EMail:
 Nicole.Reistad@FloridaDEP.gov; Kirby.Wolfe@dep.state.fl.us; Gary.Snorek@dep.state.fl.us; Miranda.Carroll@FloridaDEP.gov; Chris.Hormuth@FloridaDEP.gov

Comments: Group A: (Horse, Upper Myakka) Freshwater Metals, Pesticides x2
 Group B: (Brandy) Freshwater Metals, E.coli, Markers, Pesticides
 Group C: (Lower Myakka) Freshwater Nutrients
 Group D: (Howard) Freshwater Metals, E. coli

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: CH ₂ ClCOOH Pre-Preserved
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS
Bottle Type: BG-500ML	Number of Bottles: 6	Preserved With: ICE
W-CT-CI-R	Template: DEFAULT	(EPA 8276) Organochlorine pesticides in water matrices by GC/MS with Negative Chemical Ionization
W-PCL-SQ-R	Template:	(EPA 8270E) Organochlorine pesticides in water matrices by GC/MS-SIM
Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: HNO ₃
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO ₃
Hardness		

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: BG-1L	Number of Bottles: 6	Preserved With: ICE
W-PSNP-TQ	Template: MASTERLIST	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SD-BES-ECQ	Template: DEFAULT	(SM 9223 Quanti-Tray) E. coli by Colilert/QT by BEA South in North Port
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: CH2CICOOH Pre-Preserved
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS

Bottle Group B (Water) With 1 Samples:

Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: ICE
W-CT-CI-R	Template: DEFAULT	(EPA 8276) Organochlorine pesticides in water matrices by GC/MS with Negative Chemical Ionization
W-PCL-SQ-R	Template:	(EPA 8270E) Organochlorine pesticides in water matrices by GC/MS-SIM
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: BG-1L	Number of Bottles: 2	Preserved With: ICE
W-PSNP-TQ	Template: MASTERLIST	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
SD-BES-ECQ	Template: DEFAULT	(SM 9223 Quanti-Tray) E. coli by Colilert/QT by BEA South in North Port
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L