

McCoy, Ian

From: King, Michael
Sent: Wednesday, September 2, 2020 4:09 PM
To: McCoy, Ian
Cc: Bunch, Cheryl; Ayres, Joshua
Subject: Re: Late Arriving Cooler

Ian,

Thank you for the notice. Please run samples with the NCR and I will make a new RQ for the qualified samples.

Thank you,

Mike K.

From: "McCoy, Ian" <Ian.McCoy@FloridaDEP.gov>
Date: Wednesday, September 2, 2020 at 3:05 PM
To: "King, Michael" <Michael.King@dep.state.fl.us>, "Bunch, Cheryl" <Cheryl.Bunch@dep.state.fl.us>
Cc: "Ayres, Joshua" <Joshua.Ayres@dep.state.fl.us>
Subject: Late Arriving Cooler

Good afternoon, we received a cooler today that contained an RQ that was sampled on 8/31/20. The RQ is "RQ-2020-08-24-19", the samples inside were at temp but outside of hold time for the short hold samples. I'm going to create an NCR for it after login. Data will be qualified as appropriate, would you like to do any resampling of this event?

If you have any more questions, feel free to ask. Have a great day!

Ian McCoy
Laboratory Technician IV
850-245-8062
Florida Department of Environmental Protection

RQ-2020-08-24-19

Login Checklist

Login Station ID: #2

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 9/2/2020 14:50

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
1.6				5		✓			9117 6381 0855

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 °C (above 10.0 °C for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# 8408

Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No ☐

CONTAINER CHECK

Evidence Tape on Bottles?

Yes ☐ No ☐ NA ☒Evidence Tape intact? Yes ☐ No ☐

If Criminal, photographs taken?

Yes ☐ No ☐Containers intact? Yes ☒ No ☐

Caps on tight?

Yes ☒ No ☐

Sufficient Sample Volume:

Yes ☒ No ☐If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ☐ No ☒ Init: IM

Chlorine detected? (One container checked per Field ID)

Yes ☐ No ☐ Init: _____

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes ☒ No ☐ NA ☐ Init: IM

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes ☐ No ☐ NA ☒ Init: IM

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: IM

Date: 9/2/2020

Event contains NCRs (Y/N)? ☒

Event ID: Holmes Run
 — NW-DIST-2020-09-02-02 (SMP)
 Page 1 of 2, Date Received: 02-SEP-2020 14:50

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	IR TEMP GUN #2 exp:11/5/20
pH Test Strips 0 – 6	
pH Paper 0 – 3	
pH Paper 0 – 13	231019 exp: 11/1/22
Chlorine Test Strips	8261 exp: 07/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 8/31/2020 **Collected By:** Frank Butera, Myranda Clark
Customer: NW-DIST **Sampling Agency:** FDEP
RQ: RQ-2020-08-24-19 **Project ID:** SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 1

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
Baggett Creek Upstream Of Hwy 90	 G4NW0336

RELINQUISHED BY(SIGNATURE): *Myranda Clark*

DATE/TIME: 8/31/2020

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By:

IM

Rec'd/Inspected Date:


9/2/20

CALIBRATION AND VERIFICATION LOG (FDEP SOP FT 1000-FT 1500, FD 1000-FD 4000)

Meter ID: 154179

RQ: 2020-09-24-19

Project: SMP

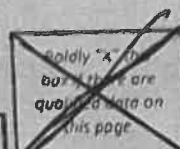

Notes: (1) Always wait for meter to stabilize before recording any readings.

(2) Report all digits displayed. Do not round before reporting measurements. (See special instructions for depth).

(3) For Calibrations, record calibrated meter reading. Do not record initial meter reading before calibration.

*DO Failed
J qualified*

Temperature (Quarterly) FT 1400

Date of Last Temperature Verification 6/4/20

DO DEP SOP FT 1500	Name	Date	Time CT-ET	Temp °C	Baro- meter mmHg	D.O. Chart mg/L	Meter D.O. mg/L	% DO	Probe Charge	Probe Gain	Pass / Fail	Lab / Field
Calibr.	<u>M. Clark</u>	<u>6/3/20</u>	<u>7:54</u>	<u>27.7</u>	<u>759.5</u>	<u>7.9</u>	<u>7.90</u>	<u>100.2</u>			<u>P/F</u>	<u>L/F</u>
ICV			<u>8:00</u>	<u>27.1</u>	<u>759.6</u>	<u>8.00</u>	<u>8.00</u>	<u>100.7</u>			<u>P/F</u>	<u>L/F</u>
CCV			<u>16:20</u>	<u>26.0</u>							<u>P/F</u>	<u>L/F</u>
CCV				<u>23.4</u>	<u>758.4</u>	<u>8.5</u>	<u>9.09</u>	<u>106.5</u>			<u>P/F</u>	<u>L/F</u>

DO Acceptance criteria from Table ± 0.3 mg/L

Rapid-Pulse Sensors: DO Gain Range 0.7 to 1.4, DO Charge Range 25-75.

Optical: DO gain range 0.85 to 1.15 (Pro DSS 0.75 to 1.50); DO charge N/A. Steady-state & Galvanic Sensors: DO Gain & Charge N/A.

Spec. Cond. FT 1200	Name	Date	Time CT-ET	Lot #	Expir. Date	Standard μ mhos/cm	Meter Reading μ mhos/cm	Pass / Fail	Lab / Field
Calibr.	<u>M. Clark</u>	<u>6/3/20</u>	<u>8:03</u>	<u>200124B</u>	<u>2/21</u>	<u>1000</u>	<u>1001</u>	<u>P/F</u>	<u>L/F</u>
ICV		<u>6/3/20</u>	<u>8:09</u>	<u>200124A</u>	<u>2/21</u>	<u>100</u>	<u>103.7</u>	<u>P/F</u>	<u>L/F</u>
CCV			<u>15:10</u>	<u>200124A</u>	<u>2/21</u>	<u>100</u>	<u>102.8</u>	<u>P/F</u>	<u>L/F</u>
CCV								<u>P/F</u>	<u>L/F</u>

Conductivity Acceptance criteria $\pm 5\%$

pH DEP SOP FT 1100	Name	Date	Time CT-ET	Lot #	Expir. Date	pH Buffer SU	Temp °C	Meter reading SU	mV	Pass / Fail	Lab / Field
Calibr.	<u>M. Clark</u>	<u>6/3/20</u>	<u>8:11</u>	<u>2904B02</u>	<u>3/21</u>	<u>7.0</u>	<u>27.4</u>	<u>7.00</u>	<u>-10.4</u>	<u>P/F</u>	<u>L/F</u>
Calibr.			<u>8:14</u>	<u>203620</u>	<u>6/22</u>	<u>4.0</u>	<u>27.1</u>	<u>4.01</u>	<u>-28.50</u>	<u>P/F</u>	<u>L/F</u>
Calibr.			<u>8:18</u>	<u>4904E09</u>	<u>9/20</u>	<u>10.0</u>	<u>26.7</u>	<u>9.99</u>	<u>-33.82</u>	<u>P/F</u>	<u>L/F</u>
ICV			<u>8:19</u>	<u>203620</u>	<u>6/22</u>	<u>4.0</u>	<u>26.7</u>	<u>4.08</u>	<u>-33.4</u>	<u>P/F</u>	<u>L/F</u>
CCV			<u>15:16</u>	<u>2901B02</u>	<u>3/21</u>	<u>7.0</u>	<u>24.7</u>	<u>6.98</u>	<u>-183.0</u>	<u>P/F</u>	<u>L/F</u>
CCV										<u>P/F</u>	<u>L/F</u>

pH Acceptance criteria ± 0.2 SU, mV pH 7 Range 0 ± 50 , mV pH 4 Range $+180 \pm 50$, mV pH 10 Range -180 ± 50 .

If mV are recorded: slope from 7 to 10 _____, slope from 4 to 7 _____ (both must be between 165 and 180 mV)

Does meter have a depth sensor that will be used to measure total depth & sample depth? YES / NO / NA (not surf. water project)

If YES, complete daily Calibr. & ICV below and list date of last quarterly depth verification: _____

If NO, what will be used? (Secchi Disk Line / Sonar) Unique ID: _____; Date of last verification: _____

Depth Sensor (Daily Calibration & ICV)	Name	Date	Time CT-ET	Calibrated Value (0.00 or Offset), meters	ICV Value, meters	Pass / Fail	Lab / Field
Pressure mode in air	<u>M. Clark</u>	<u>6/3/20</u>	<u>8:20</u>	<u>0.00</u>	<u>0.00</u>	<u>P/F</u>	<u>L/F</u>

Report two decimal places. Round numbers ≤ 4 down, ≥ 5 up. ICV acceptance criteria $\pm 5\%$ or ± 0.05 m, whichever is greater.

COMMENTS:



RQ-2020-08-24-19

Customer: WQAP

Project ID: SMP

Collected By: Frank Butera, Myranda Clark

Field Report Prepared By: Frank Butera

Sampling Agency: FDEP

Place Label Here						Comments: Microbiology sample(s) submitted to contract laboratory				
WIN ID/Field ID:										
		G4NW0336								
Location:		Baggett Creek Upstream Of Hwy 90								
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPTH)
8/31/2020	13:00 CEN	7.57	93.3	25.9	6.71	0.2	0.4 ☒ VOB (S)	0.4	34.1	0.01
↑ Central Lab please use top row for login purposes ↑										
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☒ H2SO4*		☒ <2	1	B
		Lot #	SA9354050		Exp:	12/20				
Metals (P-500ML)		☒ W-ICPMS ☒ W-ICP ☒ W-HARD				☒ Ice ☒ HNO3*		☒ <2	1	B
		Lot #	NA9301040		Exp:	10/20				
Filtered Nutrients (P125ML)		☒ W-PO4-F ☒ Field Filtered 0.45 µm Filter				☒ Ice			1	B
		Filter Lot #	16953247							
		Syringe Lot #	0799299							
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	B
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CH-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	B
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice				
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☒ E Coli ☐ F Coli ☐ Enterococci ☒ Contract Lab				☒ Ice			1	B
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☐ Ice				
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice				
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ W-E8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R ☐ W-PCL-SQ-R ☐ W-CARB-AA (5 ml of MCAA** Buffer Solution)				☐ Ice				
						☐ Ice				
Name: Frank Butera Myranda Clark				Signature: <i>Myranda Clark</i>				Date: 08/31/2020		

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution



FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION STRATEGIC MONITORING PROGRAM FIELD SHEET

January 2020

Meter Passed End
of Day CCV
Yes / No

WIN Station Name: Baggett Creek Upstream Of Hwy 90

Date: 8/31/2020

WIN/Field ID: G4NW0336

WBID: 275

ENR:

Latitude: 30.72808

County: OKALOOSA

Org ID: 21FLPNS

Longitude: -86.6603

Receiving Body of Water: Pensacola

RQ- 2020-08-24-19

Sampling Team Members		WQ Measurements	Sample Collection	Documentation	Preservation	Preservation Check
☒	Frank Butera		☒		☒	☒
☒	Myranda Clark	☒		☒		
☐						
☐						

Sampling Equipment	
☒ Sample Container	☐ Carboy
☐ Dipper	☐ Van Dorn
☒ Syringe	☐ Pole
Depth measured with: Meter	
Secchi Equipment ID:	
Collection Method	
☒ Wading	☐ Canoe/Kayak
☐ From Shore	☐ Electric Motor
☐ Other	☐ Gasoline Motor

Water Level: ☐ Normal ☐ Low ☒ High ☐ Dry ☐ Pooled ☐ Flooded		Meter ID# 154179	
Photos: No	Flow: FLOW	Tide:	Tide Times: High: Low:
Biological Samples Collected: ☒ None ☐ HA ☐ SCI ☐ RPS ☐ Algae ID ☐ LVS ☐ LVI Other:			
SBIO ID Numbers:			
Notes: MYRANDA CLARK WAS ON SAMPLING TEAM AND COMPLETED WQ MEASUREMENTS AND DOCUMENTATION			

Duplicate Sample

Time Zone: CEN	Time: N/A	Field ID: N/A	(WIN ID_DUP)
WIN DMT Activity ID:			

Blank

Time Zone:	Time: N/A	Field ID: N/A	(Blank Type_DDMMYYYY)
DI Source: N/A	Location: N/A		
DI Type: N/A	Equipment ID: N/A		

LIMS EVENT ID: WIN Import ID:

Enter Field Parameters, Verify Station ID In LIMS DMT: QC Station ID, Field Parameters In LIMS DMT:

Scan/Save Field Sheets Migrate Data to WIN: Verify Data In WIN:

Document Reference # G510

Project FDEP

Laboratory Wetlands Research Laboratory

Sampler's initials

Method of transport. Drop off on ice by automobile

Samples delivered by HB of FDEP

Sample received by: Cholde of WRL

Page 1 of 1

Date/Time 8/3/20, 7:15 PM

Date/Time: 08/31/2020 1:1412

Instructions: Fill in all fields except "Lab ID" on sample transport form completely. Each sample container must use a separate line. Place this form in a waterproof bag and place it with the samples for transport. If a field log or other chain of custody is attached, this information may be referenced on this form and the log or COC attached.

[illegible]

Comments sample container ID 20024#...022

RQ-2020-08-24-19

Contact organization FDEP

Contact name Michael King

Phone: (850) 595-8300

Address 160 W Governmental Street suite 208 Pensacola, FL 32502

Sampler's name Frank Butera

¹Matrix: SW = surface water, PW = pore water, GW = ground water, S = soil/sediment, A =

²Type codes: G = grab, C = composite, X =

Recorded using IRT-2

Document Reference # 05.1.0

FDEP

Laboratory: Wetlands Research Laboratory

Sampler's initials :

Drop off

of

of

Page _1_ of _1_

Date/Time: 8/5/20 FLD 1912

Date/Time: 08/31/2020 1412

If a field log or other chain of custody is attached, that information may be referenced on this form and the log or COC attached.

[illegible]

Comments sample container lot# 620024#..022

RQ-2020-08-24-19

Contact organization FDEP

Contact name Michael King

Phone: (850) 595-8300

Address 160 W Governmental Street, suite 208, Pensacola, FL 32502

Sampler's name: Frank Butera

¹Matrix: SW = surface water, PW = pore water, GW = ground water, S = soil/sediment, X =

²Type codes: G = grab, C = composite, X =

³ Recorded using IRT-2

1 CHL
9/2/2020

Date of Request: 30-JUL-2020
 Created By: KING_M On: 30-JUL-2020
 Modified By: JASROTIA_P On: 04-AUG-2020
 Customer: NW-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Northwest District
 Event Description: Holmes Run

Send Coolers To:

Phone: 850-595-0605
 160 W Government Street
 Suite 208
 Pensacola, FL 32502-5794
 Attn: Michael King
 Cooler Ship EMail: michael.king@dep.state.fl.us

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 04-AUG-2020
 Biology Request Reviewed By: JASROTIA_P On: 04-AUG-2020
 Sampling Kit Required: YES Ship on: 14-AUG-2020
 Sampling Kit Shipped: YES On: 07-AUG-2020
 Sampling Kit Packed By: MCCOY_I
 Preservatives Needed: NO
 Date to Receive Samples: 24-AUG-2020
 Logged In By:

Send Final Report To:

160 W Government Street
 Suite 208
 Pensacola, FL 32502-5794
 Attn: Cheryl Bunch
 Report Notification EMail:
 michael.king@dep.state.fl.us;cheryl.bunch@dep.state.fl.us

Comments: Group A: (WBID 142) Freshwater Kit With Metals, Pesticides, Contract Lab Bacteria
 Group B: (WBID 239) Freshwater Kit, Metals, Contract Lab Bacteria
 Group C: (WBID 54) , Tracers Markers

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BG-1L	Number of Bottles: 4	Preserved With: ICE
W-PSNP-TQ	Template: DEFAULT	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: TEST HAS NO	Number of Bottles: 1	Preserved With: ICE
NW-UWF-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) E. coli by Collert-18/QT by UWF in Pensacola
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: CH ₂ CICOOH Pre-Preserved
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS
Bottle Type: BG-500ML	Number of Bottles: 4	Preserved With: ICE
W-CT-CI-R	Template: DEFAULT	(EPA 8276) Organochlorine pesticides in water matrices by GC/MS with Negative Chemical Ionization
W-PCL-SQ-R	Template:	(EPA 8270E) Organochlorine pesticides in water matrices by GC/MS-SIM

Bottle Group A (Water) With 1 Samples:

Bottle Type: BG-500ML		Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B)	Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B)	Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B)	Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness			
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects	
Barium			
Calcium			
Iron			
Magnesium			
Potassium			
Sodium			
Zinc			
W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects	
*Boron			
Aluminum			
Antimony			
Arsenic			
Barium			
Beryllium			
Cadmium			
Chromium			
Cobalt			
Copper			
Lead			
Manganese			
Molybdenum			
Nickel			
Selenium			
Silver			
Thallium			
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0)	Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0)	Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0)	Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011)	Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML		Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0)	Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-1L		Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011)	Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0)	Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1)	Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011)	True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)			
W-F	Template: DEFAULT	(SM 4500 F-C-2011)	Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1)	Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011)	Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011)	Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L		Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.))	Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: TEST HAS NO		Number of Bottles: 1	Preserved With: ICE
NW-UWF-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray)	E. coli by Colilert-18/QT by UWF in Pensacola
Bottle Type: P-500ML		Number of Bottles: 1	Preserved With: HNO3

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group C (Water) With 1 Samples:

Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.