

RQ-2020-02-20-35

Login Checklist

Login Station ID: #1

Shipping Method: FedEx | UPS ☒ HD | Greyhound

Date/Time of Receipt: 9/3/2020 0805

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
					Yes	No	Yes	No	
4.0		HNO ₃	Formalin	5		✓			

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No ☒

CONTAINER CHECK

Evidence Tape on Bottles?

Yes _____ No _____ NA ☒

Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No _____

Containers intact? Yes ☒ No _____

Caps on tight?

Yes ☒ No _____

Sufficient Sample Volume:

Yes ☒ No _____If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ☒ No _____

Init: _____

Chlorine detected? (One container checked per Field ID)

Yes _____ No ☒

Init: _____

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes ☒ No _____

NA _____

Init: _____

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes _____ No _____

NA ☒

Init: _____

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: _____

Date: 9/3/2020

Event contains NCRs (Y/N)? ☒ N

Event ID: _____

Run 1

TLH-ROC-2020-09-03-01 (SMP)

Page 1 of 2, Date Received: 03-SEP-2020 08:05

DEC-4-4-2

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	#4 Temp Gun EXP:11/05/2020
pH Test Strips 0 – 6	N/A
pH Paper 0 – 3	N/A
pH Paper 0 – 13	231019 EXP:11/01/2022
Chlorine Test Strips	8261 EXP:07/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 9/2/2020 **Collected By:** Mike Eckles, Terry "Chris" Green
Customer: TLH-DIST **Sampling Agency:** FDEP
RQ: RQ-2020-07-20-35 **Project ID:** SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 1

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
Bailey Mill Creek @ Wilson Rd	 G1TLHR0081

RELINQUISHED BY(SIGNATURE):




DATE/TIME:

9/2/2020

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By:

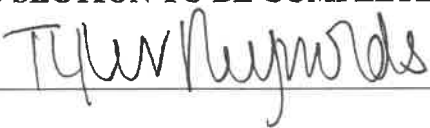

Rec'd/Inspected Date:

09/03/2020



FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
STRATEGIC MONITORING PROGRAM FIELD SHEET

January 2020

Meter Passed End
of Day CCV
Yes / No

WIN Station Name: Bailey Mill Creek @ Wilson Rd Date: 9/2/2020

WIN/Field ID: G1TLHR0081 WBID: 3419A ENR: Latitude: 30.2749

County: JEFFERSON Org ID: 21FLTLHR Longitude: -83.96667

Receiving Body of Water: Suwannee RQ- 2020-07-20-35

Sampling Team Members		WQ Measurements	Sample Collection	Documentation	Preservation	Preservation Check
☒	Mike Eckles	☒	☒	☒	☒	☒
☒	Terry "Chris" Green	☒	☒	☒	☒	☒
☐		☐	☐	☐	☐	☐
☐		☐	☐	☐	☐	☐

Sampling Equipment	
☒ Sample Container	☐ Carboy
☐ Dipper	☐ Van Dorn
☐ Syringe	☐ Pole
Depth measured with: Secchi	
Secchi Equipment ID: #1	
Collection Method	
☒ Wading	☐ Canoe/Kayak
☐ From Shore	☐ Electric Motor
☐ Other	☐ Gasoline Motor

Water Level: ☒ Normal ☐ Low ☐ High ☐ Dry ☐ Pooled ☐ Flooded		Meter ID# 153484	
Photos: NO	Flow: FLOW	Tide:	Tide Times: High: Low:
Biological Samples Collected: ☒ None ☐ HA ☐ SCI ☐ RPS ☐ Algae ID ☐ LVS ☐ LVI Other:			
SBIO ID Numbers:			
Notes:			

Duplicate Sample			
Time Zone: EST	Time: N/A	Field ID: N/A	(WIN ID_DUP)
WIN DMT Activity ID:			

Blank			
Time Zone:	Time: N/A	Field ID: N/A	(Blank Type_DDMMYYYY)
DI Source: N/A	Location: N/A		
DI Type: N/A	Equipment ID: N/A		

LIMS EVENT ID: WIN Import ID:

Enter Field Parameters, Verify Station ID In LIMS DMT: (Initials/Date) QC Station ID, Field Parameters In LIMS DMT: (Initials/Date)

Scan/Save Field Sheets (Initials/Date) Migrate Data to WIN: (Initials/Date) Verify Data In WIN: (Initials/Date)



RQ-2020-07-20-35

Customer: WQAP

Project ID: SMP

Collected By: Mike Eckles, Terry "Chris" Green

Field Report Prepared By: Mike Eckles

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID: 										
G1TLHR0081										
Location: Bailey Mill Creek @ Wilson Rd										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPTH)
9/2/2020	10:50 EST	5.09	62	25.4	6.39	0.3	0.3 ☐ VOB (S)	0.7	108.5	0.05
Central Lab please use top row for login purposes										
EST										
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☒ H2SO4*		☒ <2	1	E
		Lot #		Sa038110		Exp:		2/7/21		
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO3*		☐ <2		
		Lot #				Exp:				
Filtered Nutrients (P125ML)		☒ W-PO4-F ☐ Field Filtered 0.45 µm Filter				☒ Ice			1	E
		Filter Lot #		1696882						
		Syringe Lot #		9240133						
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	E
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CL-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	E
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice				
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☐ E Coli ☐ F Coli ☐ Enterococci				☐ Ice				
		☐ Contract Lab								
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3				☐ Ice				
		☐ PCR-GULL2 ☐ PCR-GFD ☐								
Tracers (BG-500ML)		☒ W-E8321-DI ☒ W-E8321-MS				☒ Ice			1	E
Pesticides (BG-500ML) (BG-1L)		☐ W-E8321-DI ☐ W-E8321-MS ☐ W-PSNP-TQ ☐ W-CT-CL-R				☐ Ice				
		☐ W-PCL-SQ-R ☐ W-CARB-AA (5 ml of MCAA** Buffer Solution)								
						☐ Ice				
Name: Mike Eckles Terry "Chris" Green						Signature: 			Date: 9/2/2020	

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution

CALIBRATION AND VERIFICATION LOG (for DEP SCF FT 1000-FT 1500, FD 1000-FD 4000)

Depth "N" data
Note if there are
qualified data on
this page.

Water ID: 152417

NO: 780-83-70-85

Project: SMP

Notes: (1) Always wait for meter to stabilize before recording any readings.

(2) Report all digits displayed. Do not round before reporting measurements. (See special instructions for depth).

(3) For Calibrations, record calibrated meter reading. Do not record initial meter reading before calibration.

Temperature (Quarterly) FT 1000

Date of Last Temperature Verification

ID	Name	Date	Time	Temp	Baro-	D.O.	Meter	% DO	Probe	Probe	Pass	Lab
DEP SCF			CT-ET	°C	meter	Chart	D.O.		Charge	Gain	/	/
FT 1000					mmHg	mg/L	mg/L				Fail	Field
Calib.	A. J. Eklun	9/1/20	8:55	21.7	752.7	8.2	8.29	100.0	-	101	P/F	Q/F
ICV	A. J. Eklun	9/2/20	8:27	21.3	752.7	8.1	8.2	100.7	-	-	Q/F	Q/F
CCV	A. J. Eklun	9/2/20	12:30	20.4	757.9	7.4	7.57	101.0	-	-	P/F	L/F
CCV												

DO Acceptance criteria from Table 1 0.3 mg/L

Rapid-Pulse Sensors: DO Gain Range 0.7 to 1.4; DO Charge Range 25-75.

Optical DO gain range 0.85 to 1.15 (Pre DO 0.75 to 1.50); DO charge N/A. Steady-state & Galvanic Sensors: DO Gain & Charge N/A.

Spec.	Name	Date	Time	Lot #	Expir.	Standard	Meter	Pass	Lab
Cont.			CT-ET		Date	umhos/cm	Reading	/	/
FT 1000							umhos/cm	Fail	Field
Calib.	A. J. Eklun	9/1/20	5:00	1911200	12/1/20	1000	998	Q/F	L/F
ICV	A. J. Eklun	9/2/20	9:05	2911083	11/2/20	100	101.0	Q/F	Q/F
CCV	A. J. Eklun	9/2/20	13:35	1911200	12/1/20	1000	975	Q/F	Q/F
CCV								P/F	L/F

Conductivity Acceptance criteria $\pm 5\%$

pH	Name	Date	Time	Lot #	Expir.	pH	Temp	Meter	mV	Pass	Lab
DEP SCF			CT-ET		Date	Buffer	°C	reading		/	/
FT 1000						SU		SU		Fail	Field
Calib.	A. J. Eklun	9/2/20	9:10	2005078	10/21	7.01	22.1	7.01	21.8	P/F	Q/F
Calib.	A. J. Eklun	9/2/20	9:12	81120A	6/21	4.00	22.4	4.00	19.4	Q/F	N/F
Calib.	A. J. Eklun	9/2/20	9:15	191120C	6/21	10.03	22.6	10.03	20.7	Q/F	Q/F
ICV	A. J. Eklun	9/2/20	9:20	191120C	6/21	7.0	22.5	7.0	20.3	Q/F	Q/F
CCV	A. J. Eklun	9/2/20	13:10	191120A	6/21	4.0	23.5	3.9	19.2	Q/F	Q/F
CCV										P/F	L/F

pH Acceptance criteria ± 0.2 SU; mV pH 7 Range 0 ± 50 ; mV pH 4 Range $\pm 180 \pm 50$; mV pH 10 Range -180 ± 50 ;

If mV are recorded: slope from 7 to 10; slope from 4 to 7; [both must be between 165 and 180 mV]

Does meter have a depth sensor that will be used to measure total depth & sample depth? YES / NO / NA (not surf. water project)

IF YES, complete daily Calib. & ICV below and list date of last quarterly depth verification:

IF NO, what will be used? Secchi Disk Line / Sonar Unique ID: _____ Date of last verification: _____

Depth Sensor	Name	Date	Time	Calibrated	ICV Value	Pass /	Lab /
(Daily Calibration & ICV)			CT-ET	Value (0.00 or Offset), meters	meters	Fail	Field
Pressure mode in air						P/F	L/F

Report two decimal places. Round numbers ≤ 4 down, ≥ 5 up. ICV acceptance criteria $\pm 5\%$ or ± 0.05 m, whichever is greater.

COMMENTS:

Date of Request: 30-JUN-2020
 Created By: ECKLES_M On: 30-JUN-2020
 Modified By: JASROTIA_P On: 02-JUL-2020
 Customer: TLH-ROC
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Not Applicable
 Event Description: Run 1

Send Coolers To:

Phone: 850-245-7586
 2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles
 Cooler Ship EMail: michael.eckles@dep.state.fl.us

Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 02-JUL-2020
 Biology Request Reviewed By: JASROTIA_P On: 02-JUL-2020
 Sampling Kit Required: YES Ship on: 23-JUL-2020
 Sampling Kit Shipped: YES On: 22-JUL-2020
 Sampling Kit Packed By: SHAIK_A
 Preservatives Needed: NO
 Date to Receive Samples: 20-JUL-2020
 Logged In By: KLUG_K On: 02-SEP-2020

Send Final Report To:

2600 Blair Stone Rd

 Tallahassee, FL 32399
 Attn: Michael Eckles
 Report Notification EMail: michael.eckles@dep.state.fl.us

Comments: Group A WBID 3314 & 3473B
 Group B WBID 3473C, 3518, Dup
 Group C WBID 3518 & 3318 BacT only
 Group D WBID 3318
 Group E 3322A
 Group F Blank
 No fedex labels

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 2	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: HNO ₃

Bottle Group A (Water) With 2 Samples:

Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 3 Samples:

Bottle Type: P-500ML	Number of Bottles: 3	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 3	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 3	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 3	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 3	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group F (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

