

Login Checklist

RQ- 2021-08-02-23

Login Station ID: #5

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 11/ 2 /21 @ 9:33

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
0.5				7		✓			5225 5541 1588

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO3 and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No ☒

CONTAINER CHECK

Evidence Tape on Bottles?

Yes _____ No _____ NA ☒

Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No _____

Containers intact? Yes ☒ No _____

Caps on tight?

Yes ☒ No _____

Sufficient Sample Volume:

Yes ☒ No _____

If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes _____ No ☒ Init: kh

Chlorine detected? (One container checked per Field ID)

Yes _____ No _____ Init: _____

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH< 2.0?

Yes _____ No _____ NA ☒ Init: kh

**W-1-4-DIOX (Green dot on container) preserved to pH<4.0?

Yes _____ No _____ NA ☒ Init: kh

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: kh

Date: 11/ 2 /21

Event contains NCRs (Y/N)? N

Event ID: Tomoka Basin/IRL (2963A1)
CEN-DIST-2021-11-02-01 (SMP)

Date Received: 02-NOV-2021 09:33

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	IR TEMP GUN #5 EXP:10-19-22
pH Test Strips 0 – 6	Lot#210620 EXP:4-15-23
pH Paper 0 – 3	
pH Paper 0 – 13	Lot#210620 EXP:4-15-23
Chlorine Test Strips	Lot#0070 EXP:2-23

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 11/1/2021
Customer: CEN-DIST
RQ: RQ-2021-08-02-23

Collected By: Terry Riordan, Elizabeth Bates
Sampling Agency: FDEP
Project ID: SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB


Number of Coolers Shipped: 1

Delivery Method: FedEx

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
IRL @ Castaway Point Park	 G5CE0155

RELINQUISHED BY(SIGNATURE):

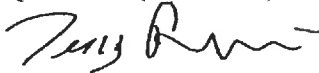

DATE/TIME: 11/1/2021
1600 EST

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By: *HK*

Rec'd/Inspected

Date: *11-2-21/0933*



RQ-2021-08-02-23

Collected By: Terry Riordan, Elizabeth Bates

Customer: FDEP

Field Report Prepared by: Terry Riordan, Elizabeth Bates

Project ID: SMP

Sampling Agency: FDEP

 WIN ID/Field ID: 
 WIN ID/Field ID: G5CE0155

Comments:

Location:

Location: IRL @ Castaway Point Park

Matrix: Marine

(Check One)

☒ Grab☐ Composite

Date and Time {EST}	DO (mg/L)	DO (%SAT)	Temp. ©	pH	Sample Depth (m)	Secchi Depth	Total Depth (m)	Sp COND (µmhos /cm)	Salinity (PPT)
11/1/2021 1112	5.73	75.8	23.5	7.65	0.3	1.3	1.3	32033	19.98

Parameter Suite	Parameter Methods	Preservation	pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)	W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP	Ice H2SO4*	<2		
	Lot # Exp:				
Metals (P-500ML)	W-ICPMS W-ICP W-HARD	Ice HNO3*	<2		
	Lot # Exp:				
Filtered Nutrients (P125ML)	W-PO4-F Field Filtered 0.45 um Filter	Ice			
	Filter Lot # Syringe Lot #				
Chlorophyll (BP-1L)	CHLSUITE-W	Ice			
Physical/Aggregate (Fresh) (P-1L)	Alkalinity / W-F / W-TSS / TURBIDITY / W-CL-IC / W-COLOR / W-SO4-IC / W- TDS	Ice			
Physical/Aggregate (Marine) (P-1L)	Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY	Ice			
Macroinvertebrates (PJ-2L)	MI-FW-QLDC Formalin Lot #	Formalin			
Periphyton (PT-50ML)	ALGAE-ID	Ice			
Microbiology (Contract Lab Bottle)	E Coli F Coli X Enterococci Contract Lab	x Ice		2	B
Molecular (QPCR-P-500ML)	PCR-HF183 PCR-BacR PCR-DG3 PCR-GULL2 PCR-GFD	Ice			
Tracers (BG-500ML)	W-E8321-DI W-E8321-MS	Ice			
Pesticides (BG-500ML) (BG-1L)	W-E8321-DI W-E8321-MS W-PSNP-TQ W-PCL-TQ-R W-CARB-AA (5 ml of MCAA** Buffer Solution) W-CT-CL-R	Ice			
Phytoplankton	DTY-QN-C PKD-QN-BV	Ice			

Date of Request: 02-JUL-2021
 Created By: BOMAN_L On: 02-JUL-2021
 Modified By: JASROTIA_P On: 08-JUL-2021
 Customer: CEN-DIST
 Project: SMP
 Division:
 District: Central District
 Event Description: Tomoka Basin/IRL (2963A1)

Send Coolers To:

Phone: 407-897-4181
 FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Leif Boman
 Cooler Ship EMail: Leif.Boman@FloridaDEP.gov

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 08-JUL-2021
 Biology Request Reviewed By: JASROTIA_P On: 08-JUL-2021
 Sampling Kit Required: YES Ship on: 23-JUL-2021
 Sampling Kit Shipped: YES On: 21-JUL-2021
 Sampling Kit Packed By: MCCOY_I
 Preservatives Needed: NO
 Date to Receive Samples: 02-AUG-2021
 Logged In By:

Send Final Report To:

3319 Maguire Blvd., Suite 232

 Orlando, FL 32803-3767
 Attn: Matt Bearden
 Report Notification EMail: Matthew.Bearden@FloridaDEP.gov

Comments: Bottle Group A: Tomoka
 Bottle Group B: IRL

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-ICP	Template: S-METALS	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron		
Manganese		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Silver

Thallium

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems