

Login Checklist

RQ-2022-04-04-39

Login Station ID: #4

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 04/21/22 @ 9:21

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
1.9				5		X			5225 5545 6491

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO3 and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No ☐

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes ☒ No ☐

CONTAINER CHECK

Evidence Tape on Bottles?

Yes ☐ No ☐ NA ☒

Evidence Tape intact? Yes ☐ No ☐

If Criminal, photographs taken?

Yes ☐ No ☐

Containers intact? Yes ☒ No ☐

Caps on tight?

Yes ☒ No ☐

Sufficient Sample Volume:

Yes ☒ No ☐

If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ☒ No ☐ Init: RR

Chlorine detected? (One container checked per Field ID)

Yes ☐ No ☒ Init: RR

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes ☒ No ☐ NA ☐ Init: RR

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes ☐ No ☐ NA ☒ Init: RR

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: RR

Date: 04/21/22 Event contains NCRs (Y/N)? N

Event ID: Pellicer Creek Tribs
NE-DIST-2022-04-21-01 (SMP)

- Date Received: 21-APR-2022 09:21

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	#4 Temp Gun EXP:10/19/2022
pH Test Strips 0 – 6	223819AV EXP:08/30/2023
pH Paper 0 – 3	N/A
pH Paper 0 – 13	207621 EXP:03/15/2024
Chlorine Test Strips	0070 EXP:02/2023

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 4/20/2022
Customer: NE-DIST
RQ: RQ-2022-04-04-38
Collected By: Jeremy Parrish, Kimberly Appleby
Sampling Agency: FDEP
Project ID: SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 1

Delivery Method: FedEx

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
Stevens Branch East Branch at CR204	 G5NE0600
STEVENS BR @ CR 204	 27010062

RELINQUISHED BY(SIGNATURE):



DATE/TIME: 4/20/2022
1600EST

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By: 

Rec'd/Inspected Date: 4-21-2022

CALIBRATION AND VERIFICATION LOG (FDEP SOP FT 1000-FT 1500, FD 1000-FD 4000)

Boldly "X" this box if there are qualified data on this page.

Meter ID: Brule RQ: 2022-04-04-38 Project: SMP Pellicer

Tris

- Notes: (1) Always wait for meter to stabilize before recording any readings.
 (2) Report all digits displayed. Do not round before reporting measurements. (See special instructions for depth).
 (3) For Calibrations, record calibrated meter reading. Do not record initial meter reading before calibration.

Temperature (Quarterly) FT 1400

Date of Last Temperature Verification 2/3/22

DO DEP SOP FT 1500	Name	Date	Time CT-ET	Temp °C	Baro-meter mmHg	D.O. Chart mg/L	Meter D.O. mg/L	% DO	Probe Charge	Probe Gain	Pass / Fail	Lab / Field
Calibr.	Kim Appleby	4/20/22	836	20.3	768.0	9.1	-	99.9	-	1.031	☒ F	☒ F
ICV	↓	↓	837	20.4	768.1	9.1	9.13	101.1	-	1.031	☒ F	☒ F
CCV	J.P. Amsh	4/20/22	1330	21.4	768	8.9	8.99	101.5	-	1.031	☒ F	☒ F
CCV											P / F	L / F

DO Acceptance criteria from Table ± 0.3 mg/L.

Rapid-Pulse Sensors: DO Gain Range 0.7 to 1.4; DO Charge Range 25-75.

Optical: DO gain range 0.85 to 1.15 (Pro DSS 0.75 to 1.50); DO charge N/A. Steady-state & Galvanic Sensors: DO Gain & Charge N/A.

Spec. Cond. FT 1200	Name	Date	Time CT-ET	Lot #	Expir. Date	Standard μ hos/cm	Meter Reading μ hos/cm	Pass / Fail	Lab / Field
Calibr.	Kim Appleby	4/20/22	841	22022SA	3/23	1000	1000	☒ F	☒ F
ICV	↓	↓	843	210917A	9/22	100	100.5	☒ F	☒ F
CCV	J.P. Amsh	4/20/22	1338	211122B	12/22	1413	1419	☒ F	☒ F
CCV								P / F	L / F

Conductivity Acceptance criteria $\pm 5\%$

pH DEP SOP FT 1100	Name	Date	Time CT-ET	Lot #	Expir. Date	pH Buffer SU	Temp °C	Meter reading SU	mV	Pass / Fail	Lab / Field
Calibr.	Kim Appleby	4/20/22	846	2112B19	12/23	7.01	21.9	7.00	-13.7	☒ F	☒ F
Calibr.	↓	↓	848	2201D88	1/24	4.00	21.7	4.00	163.5	☒ F	☒ F
Calibr.	↓	↓	851	2202A63	7/23	10.02	21.7	10.03	-185.0	☒ F	☒ F
ICV	↓	↓	953	2112B19	12/23	7.01	21.8	7.07	-14.9	☒ F	☒ F
CCV	J.P. Amsh	4/20/22	1333	2201D88	1/24	4.00	21.9	4.00	159.8	☒ F	☒ F
CCV										P / F	L / F

pH Acceptance criteria ± 0.2 SU; mV pH 7 Range 0 ± 50 ;

mV pH 4 Range $+180 \pm 50$;

mV pH 10 Range -180 ± 50 ;

If mV are recorded: slope from 7 to 10 171.3, slope from 4 to 7 178.3 (both must be between 165 and 180 mV)

Does meter have a depth sensor that will be used to measure total depth & sample depth? YES / NO / NA (not surf. water project)

If YES, complete daily Calibr. & ICV below and list date of last quarterly depth verification: 2/7/22

If NO, what will be used? (circle one) Secchi Disk Line / Sonar Unique ID: _____; Date of last verification: _____

Depth Sensor (Daily Calibration & ICV)	Name	Date	Time CT-ET	Calibrated Value (0.00 or Offset), meters	ICV Value, meters	Pass / Fail	Lab / Field
Pressure mode in air	Kim Appleby	4/20/22	854	0.00	0.00	☒ F	☒ F

Report two decimal places. Round numbers ≤ 4 down, ≥ 5 up. ICV acceptance criteria $\pm 5\%$ or ± 0.05 m, whichever is greater.

COMMENTS:



RQ-2022-04-04-38

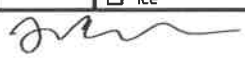
Customer: NE-DIST

Project ID: SCI BENCHMARK

Collected By: Jeremy Parrish, Kimberly Appleby

Field Report Prepared By: Jeremy Parrish

Sampling Agency: FDEP

Place Label Here						Comments: Microbiology sample(s) submitted to contract laboratory.					
WIN ID/Field ID: 											
G5NE0600											
Location: Stevens Branch East Branch at CR204											
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE					
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPTH)	
4/20/2022	1050 EST	7.28	75.3	17.1	6.59	0.3	0.4 ☐ VOB (S)	0.47	150	0.07	
	EST						Central Lab please use top row for login purposes				
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group	
Preserved Nutrients (P-500ML)		☐ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☐ Ice ☐ 250 ⁺		☐ <2			
		Lot #		Exp:							
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO ₃ ⁺		☐ <2			
		Lot #		Exp:							
Filtered Nutrients (P125ML)		☐ W-PO4-F ☐ Field Filtered 0.45 µm Filter				☐ Ice					
		Filter Lot #				Syringe Lot #					
BOD (BP-1L)		☐ OV-BODUNIN				☐ Ice					
Chlorophyll (BP-1L)		☐ CHLSUITE-W				☐ Ice					
Physical/Aggregate (Fresh) (P-1L)		☐ Alkalinity / W-F / W-TSS / TURBIDITY / W-CHIC / W-COLOR / W-SO4-IC / W-TDS				☐ Ice					
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice					
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin					
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice					
Microbiology (Contract Lab Bottle)		☒ Coli ☐ Coli ☐ Enterococci ☒ Contract Lab				☒ Ice			1	E	
Molecular (QPCR-P-500ML)		☒ PCR-HF183 ☒ PCR-BacR ☒ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☒ Ice			1	E	
Tracers (BG-500ML)		☒ W-E8321-DI ☒ W-E8321-MS				☒ Ice		MCAA Exp.	1	E	
Pesticides (BG-500ML) (BG-1L)		W-E8321-DI WE8321-MS W-PSNP-TQ W-CT-CL-R W-PCL-TQ-R W-CARB-AA (5 ml of MCAA** Buffer Solution) W-TR-SQ-R W-PCB-TQ-R				☐ Ice		MCAA Lot#			
Phytoplankton/Periphyton		DTY-QN-C PKD-QN-C				☐ Ice					
Name: Jeremy Parrish Kimberly Appleby					Signature: 					Date: 04/20/2022	

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution

**FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION****SCI BENCHMARK FIELD SHEET**

January 2022

**Meter Passed End
of Day CCV**
Yes / No**WIN Station Name:** Stevens Branch East Branch at CR204**Date:** 4/20/2022**WIN/Field ID:** G5NE0600**WBID:** 2551A**ENR:** -**Latitude:** 29.66**County:** ST. JOHNS**Org ID:** 21FLA**Longitude:** -81.33**Receiving Body of Water:** Pringle Branch**RQ-** 2022-04-04-38

Sampling Team Members		WQ Measurements	Sample Collection	Documentation	Preservation	Preservation Check
☒	Jeremy Parrish	x		x		
☒	Kimberly Appleby		x		x	x
☐						
☐						

Sampling Equipment	
☒ Sample Container	☐ Carboy
☐ Dipper	☐ Van Dorn
☐ Syringe	☒ Pole
Depth measured with: Meter	
Secchi Equipment ID:	
Collection Method	
☐ Wading	☐ Canoe/Kayak
☒ From Shore	☐ Electric Motor
☐ Other _____	☐ Gasoline Motor

Water Level: ☐ Normal ☐ Low ☒ High ☐ Dry ☐ Pooled ☐ Flooded		Meter ID# BRUCE
Photos: NO	Flow: FLOW	Tide:
Biological Samples Collected: ☒ None ☐ HA ☐ SCI ☐ RPS ☐ Algae ID ☐ LVS ☐ LVI Other:		Tide Times: High: Low:
SBIO ID Numbers:		
Notes: HIGH RECENT RAIN		

Duplicate Sample

Time Zone: EST	Time:	Field ID: N/A	(WIN ID_DUP)
WIN DMT Activity ID:			

Blank

Time Zone:	Time:	Field ID:	(Blank Type_DDMMYYYY)
DI Source:	Location:		
DI Type:	Equipment ID:		
DI Container (Name/ID):		DI Container Type:	

LIMS EVENT ID: WIN Import ID:

Enter Field Parameters, Verify Station ID In LIMS DMT: (Initials/Date) QC Station ID, Field Parameters In LIMS DMT: (Initials/Date)

Scan/Save Field Sheets (Initials/Date) Migrate Data to WIN: (Initials/Date) Verify Data In WIN: (Initials/Date)



RQ-2022-04-04-38

Customer: NE-DIST

Project ID: STRATEGIC MONITORING PLAN

Collected By: Jeremy Parrish, Kimberly Appleby

Field Report Prepared By: Jeremy Parrish

Sampling Agency: FDEP

Place Label Here						Comments:				
WIN ID/Field ID: 										
27010062										
Location: STEVENS BR @ CR 204										
Matrix: Fresh						(Check one) ☒ GRAB ☐ COMPOSITE				
Date	Time	DO (mg/L)	DO (%SAT)	Temp. (C°)	pH	Sample Depth (m)	Secchi Disk (m)	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPT)
4/20/2022	1105 EST	5.3	56	17.9	5.4	0.1	0.2 ☒ VOB (S)	0.2	118	0.06
	EST						Central Lab please use top row for login purposes			
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group
Preserved Nutrients (P-500ML)		☒ W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				☒ Ice ☐ 2SO4*		☒ <2	1	D
		Lot #		Sa1103080		Exp:		12/01/22		
Metals (P-500ML)		☐ W-ICPMS ☐ W-ICP ☐ W-HARD				☐ Ice ☐ HNO3*		☐ <2		
		Lot #				Exp:				
Filtered Nutrients (P125ML)		☐ W-PO4-F ☐ Field Filtered 0.45 µm Filter				☐ Ice				
		Filter Lot #				Syringe Lot #				
BOD (BP-1L)		☐ OV-BODUNIN				☐ Ice				
Chlorophyll (BP-1L)		☒ CHLSUITE-W				☒ Ice			1	D
Physical/Aggregate (Fresh) (P-1L)		☒ Alkalinity / W-F / W-TSS / TURBIDITY / W-CL-IC / W-COLOR / W-SO4-IC / W-TDS				☒ Ice			1	D
Physical/Aggregate (Marine) (P-1L)		☐ Alkalinity / W-BR-IC / W-F / W-TSS / TURBIDITY				☐ Ice				
Macroinvertebrates (PJ-2L)		☐ MI-FW-QLDC				☐ Formalin				
Periphyton (PT-50ML)		☐ ALGAE-ID				☐ Ice				
Microbiology (Contract Lab Bottle)		☒ Coli ☐ Coli ☐ Enterococci ☒ Contract Lab				☒ Ice			1	D
Molecular (QPCR-P-500ML)		☐ PCR-HF183 ☐ PCR-BacR ☐ PCR-DG3 ☐ PCR-GULL2 ☐ PCR-GFD ☐				☐ Ice				
Tracers (BG-500ML)		☐ W-E8321-DI ☐ W-E8321-MS				☐ Ice		MCAA Exp.		
Pesticides (BG-500ML) (BG-1L)		W-E8321-DI W-E8321-MS W-PSNP-TQ W-CT-CL-R W-PCL-TQ-R W-CARB-AA (5 ml of MCAA** Buffer Solution) W-TR-SQ-R W-PCB-TQ-R				☐ Ice		MCAA Lot#		
Phytoplankton/Periphyton		DTY-QN-C PKD-QN-C				☐ Ice				
Name: Jeremy Parrish Kimberly Appleby						Signature: 			Date: 04/20/2022	

*2.0 ml of Preservative used unless noted otherwise.

**MCAA: Monochloroacetic acid buffer solution

**FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION****STRATEGIC MONITORING PLAN FIELD SHEET**

January 2022

Meter Passed End

of Day CCV

Yes / No

WIN Station Name: STEVENS BR @ CR 204

Date: 4/20/2022

WIN/Field ID: 27010062

WBID: 2551A

ENR: -

Latitude: 29.66

County: ST. JOHNS

Org ID: 21FLA

Longitude: -81.34

Receiving Body of Water: Pringle Branch

RQ- 2022-04-04-38

Sampling Team Members		WQ Measurements	Sample Collection	Documentation	Preservation	Preservation Check
☒	Jeremy Parrish	x		x		
☒	Kimberly Appleby		x		x	x
☐						
☐						

Sampling Equipment	
☒ Sample Container	☐ Carboy
☐ Dipper	☐ Van Dorn
☐ Syringe	☒ Pole
Depth measured with: Meter	
Secchi Equipment ID:	
Collection Method	
☐ Wading	☐ Canoe/Kayak
☒ From Shore	☐ Electric Motor
☐ Other	☐ Gasoline Motor

Water Level: ☐ Normal ☐ Low ☒ High ☐ Dry ☐ Pooled ☐ Flooded		Meter ID# BRUCE	
Photos: NO	Flow: FLOW	Tide:	Tide Times: High: Low:
Biological Samples Collected: ☒ None ☐ HA ☐ SCI ☐ RPS ☐ Algae ID ☐ LVS ☐ LVI Other:			
SBIO ID Numbers:			
Notes: HIGH RECENT RAIN			

Duplicate Sample

Time Zone: EST	Time:	Field ID: N/A	(WIN ID_DUP)
WIN DMT Activity ID:			

Blank

Time Zone:	Time:	Field ID:	(Blank Type_DDMMYYYY)
DI Source:	Location:		
DI Type:	Equipment ID:		
DI Container (Name/ID):		DI Container Type:	

LIMS EVENT ID: WIN Import ID:

Enter Field Parameters, Verify Station ID In LIMS DMT: (Initials/Date) QC Station ID, Field Parameters In LIMS DMT: (Initials/Date)

Scan/Save Field Sheets (Initials/Date) Migrate Data to WIN: (Initials/Date) Verify Data In WIN: (Initials/Date)



DIVERSIFIED
ENVIRONMENTAL LABORATORIES, INC.

CHAIN OF CUSTODY RECORD

Project Number:

(Lab Use Only)

3653 Regent Boulevard, Suite # 509

Jacksonville, FL 32224

P. (904) 807-9625 E. info@delilab.com

F. (904) 807-9627

Page 1 of 1

Company Name <u>Florida Department of Environmental Protection</u>										Container Information (Qty, Size, & Type)										Sample Matrix Codes:																																							
Address <u>8800 Baymeadows Way W</u>										NPW NPW NPW										A. Air P. Petroleum B. Bulk S. Soil/Sludge W. Wipe O. Other PW. Potable Water NPW. Non-Potable Water																																							
City <u>Jacksonville</u> State <u>FL</u> Zip <u>32256</u>										Container Preservation Codes										Sample Kit or Container Lot #																																							
Attn.										6 6 6																																																	
Tel <u>(904) 256-1700</u> E-mail or Fax <u>Dep-overflaminc@dep.state.fl.us</u>										Parameters										Preservation Codes:																																							
Project Name <u>Polluter Tribs</u>										E. Coli - Quantitray Enterococci - Quantitray Coliform, Fecal - MF 3 Dilutions										0. None 5. Na2S2O3 1. H2SO4 6. Crushed Ice (4C) 2. HNO3 7. C6H8O6 3. HCl 8. NH4Cl 4. NaOH 9. Other: (List)																																							
Project Number <u>RA-2022-04-04-38</u>																				Thermal Preservation:																																							
Project Location																				In Crushed Ice? ☒ Yes / ☐ No																																							
Sampler Name/Signature <u>Kim Appleby</u>																				Cooler Temp: <u>0.2</u> °C																																							
																				Sample Cont. Temp: <u>2.4</u> °C																																							
																				Turnaround Requested:																																							
																				☒ Standard (10 Business Days)																																							
																				☐ Semi-Rush (3-5 Business Days)																																							
																				☐ Rush (1-2 Business Days)																																							
																				Turnaround times are estimated & subject to change. Surcharges apply for all Rush samples.																																							
																				Billing Information:																																							
																				☐ Invoice with Report																																							
																				☐ Invoice to Accounting																																							
																				Req. #: _____																																							
																				P.O. #: _____																																							
Special Instructions:										Transfer Numbers:										Relinquished By										Date / Time										Received By										Date / Time									
										1										Kim Appleby										4/20/22 1245																				04/20/22 1245									



Certification #E821059

Lab General Form 059, Rev. #1, Rev. Date: 02/07/22, Effective Date: 02/07/22.

Visit us online & order kits at: www.delilab.com

Date of Request: 11-MAR-2022
 Created By: PARRISH_J On: 11-MAR-2022
 Modified By: AYRES_J On: 13-APR-2022
 Customer: NE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Northeast District
 Event Description: Pellicer Creek Tribs

Send Coolers To:

Phone: 904-256-1693
 8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish
 Cooler Ship EMail: Jeremy.M.Parrish@dep.state.fl.us

Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: AYRES_J On: 13-APR-2022
 Biology Request Reviewed By: AYRES_J On: 13-APR-2022
 Sampling Kit Required: NO Ship on:
 Sampling Kit Shipped: YES On: 28-MAR-2022
 Sampling Kit Packed By: ROBERTS_Z
 Preservatives Needed: NO
 Date to Receive Samples: 04-APR-2022
 Logged In By: SHAIK_A On: 14-APR-2022

Send Final Report To:

8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish
 Report Notification EMail: Jeremy.M.Parrish@floridadep.gov

Comments: Group A:(Surf.Fresh) Nutrients, Metals,Bacteria
 Group B:(Surf.Fresh) Nutrient, Metal, Bacteria, Tracer, Marker
 Group C:(Surf.Fresh)Nutrient, Bacteria, Tracer, PCR
 Group D:(Surf.Fresh)Nutrient, Bacteria
 Group E:(Surf.Fresh)Bacteria, Tracer, Marker
 Group F:(Surf.Fresh) Metals, Bacteria

Packing Notes:**Bottle Group A (Water) With 1 Samples:**

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
*Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-HARD	Template: DEFAULT	(SM 2340 B-2011) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO ₃

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340 B-2011) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template: SMP-FRESH	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Strontium		
Tin		
Titanium		
Vanadium		
Zinc		
W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Barium		
Beryllium		
Cadmium		
Chromium		
Cobalt		
Copper		
Lead		
Manganese		
Molybdenum		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Group B (Water) With 1 Samples:

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340 B-2011) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template: SMP-FRESH	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Strontium		
Tin		
Titanium		
Vanadium		
Zinc		
W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

*Boron
Aluminum
Antimony
Arsenic
Barium
Beryllium
Cadmium
Chromium
Cobalt
Copper
Lead
Manganese
Molybdenum
Nickel
Selenium
Silver
Thallium

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
*Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
*Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
*Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Group E (Water) With 2 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 2	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: QPCR-P-500ML	Number of Bottles: 2	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction

Bottle Group E (Water) With 2 Samples:

Bottle Type: QPCR-P-500ML	Number of Bottles: 2	Preserved With: ICE
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction

Bottle Type: BG-500ML	Number of Bottles: 2	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Group F (Water) With 1 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340 B-2011) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template: SMP-FRESH	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Calcium
Iron
Magnesium
Potassium
Sodium
Strontium
Tin
Titanium
Vanadium
Zinc

W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
---------	-------------------	---

*Boron
Aluminum
Antimony
Arsenic
Barium
Beryllium
Cadmium
Chromium
Cobalt
Copper
Lead
Manganese
Molybdenum
Nickel
Selenium
Silver
Thallium

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.

