

Login Checklist

RQ- 2022-10-03-27

Login Station ID: 5

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 11 30 / 22 @ 10:04

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
0.4				1		X			5859 47248342

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO3 and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included? Yes X No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes _____ No X

CONTAINER CHECK

Evidence Tape on Bottles? Yes _____ No _____ NA X Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken? Yes _____ No X Containers intact? Yes X No _____

Caps on tight? Yes X No _____

Sufficient Sample Volume: Yes X No _____

If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes _____ No X Init: C-P

Chlorine detected? (One container checked per Field ID) Yes _____ No _____ Init: _____

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH< 2.0? Yes _____ No _____ NA X Init: C-P

**W-1-4-DIOX (Green dot on container) preserved to pH< 4.0? Yes _____ No _____ NA X Init: C-P

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: C-P Date: 11/ 30 /22 Event contains NCRs (Y/N)? N

Event ID: CEN-DIST-2022-11-30-02 (SMP)
Date Received: 30-NOV-2022 10:04

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	221426058 10/17/23
pH Test Strips 0 – 6	223819AV 08/30/23
pH Paper 0 – 3	NA
pH Paper 0 – 13	210620 4/15/2023
Chlorine Test Strips	0070 02/2023

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 11/29/2022 **Collected By:** Leif Boman, Katie March
Customer: CEN-DIST **Sampling Agency:** FDEP
RQ: RQ-2022-10-03-27 **Project ID:** SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 1

Delivery Method: FedEx

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
SECOND CREEK AT S.R. 320	 G3CE0066
FB	 FB_11292022
Florence Lake @ Center	 G3CE0012

LB 11/29/22

LB 11/29/22


RELINQUISHED BY(SIGNATURE):



DATE/TIME: 11/29/2022
1600 EST

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By:

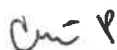

Rec'd/Inspected

Date: 11-30-22



RQ-2022-10-03-27	Collected By: Leif Boman, Katie March
Customer: FDEP	Field Report Prepared by: Leif Boman, Katie March
Project ID: SMP	Sampling Agency: FDEP
WIN ID/Field ID: WIN ID/Field ID: G3CE0012	Comments: Please log in PKD-QN-BV & DTY-QN-C under a separate event with priority 12, Bottle Group - B.

Location: Florence Lake @ Center

Matrix: Fresh

(Check One)

☒ Grab

☐ Composite

Date and Time {EST}	DO (mg/L)	DO (%SAT)	Temp. °C	pH	Sample Depth (m)	Secchi Depth	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPT)					
11/29/2022 1300	1.35	15.8	23	6.75	0.3	0.5	1.4	917	0.45					
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group				
Preserved Nutrients (P-500ML)	X	W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				X	Ice	X	H2SO4*					
		Lot #	SA2229010	Exp:	08/30/2023			X	<2	1	A			
Metals (P-500ML)	X	W-ICPMS				X	Ice	X	HNO3*					
		Lot #	NA1105100	Exp:	06/07/2023			X	<2	1	A			
Filtered Nutrients (P125ML)	X	W-PO4-F				X	Ice			1	A			
		Filter Lot #	17413308	Syringe Lot #	2090604									
BOD (BP-1L)		OV-BODUNIN					Ice							
Chlorophyll (BP-1L)	X	CHLSUITE-W				X	Ice			1	A			
Physical/Aggregate (Fresh) (P-1L)	X	Alkalinity / W-F / W-TSS / TURBIDITY / W-Cl-IC / W-COLOR / W-SO4-IC / W-TDS				X	Ice			1	A			
Physical/Aggregate (Marine) (P-1L)		Alkalinity / W-F / W-TSS / TURBIDITY					Ice							
Macroinvertebrates (PJ-2L)		MI-FW-QLDC	Formalin Lot #				Formalin							
Periphyton (PT-50ML)		ALGAE-ID					Ice							
Microbiology (Contract Lab Bottle)		E Coli					Ice							
		Contract Lab												
Molecular (QPCR-P-500ML)		PCR-HF183					Ice							
		PCR-GULL2												
Tracers (BG-500ML)		W-E8321-DI					Ice							
Pesticides (BG-500ML) (BG-1L)		W-E8321-DI	W-E8321-MS		W-PSNP-TQ		Ice	MCAA Exp.	MCAA Lot#					
		W-PCL-TQ-R	W-CARB-AA (5 ml of MCAA Buffer Solution)		W-CT-CL-R									
		W-TR-SQ-R	W-PCB-TQ-R											
Phytoplankton	X	DTY-QN-C		X	PKD-QN-BV		X	Ice		1	B			
DOC		W-DOC	Field Filtered 0.45 um Filter				Ice		H2SO 4*					
			Acid Preservative Expiration Date							<2				
			Acid Preservative Lot Number											

Date of Request: 29-AUG-2022
 Created By: BEARDEN_M On: 29-AUG-2022
 Modified By: AYRES_J On: 16-SEP-2022
 Customer: CEN-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Central District
 Event Description: Florence Lake and Faulk Canal

Send Coolers To:

Phone: 407-897-4123
 FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Matthew Bearden
 Cooler Ship Email: matthew.bearden@floridadep.gov

Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: AYRES_J On: 16-SEP-2022
 Biology Request Reviewed By: JASROTIA_P On: 30-AUG-2022
 Sampling Kit Required: YES Ship on: 23-SEP-2022
 Sampling Kit Shipped: YES On: 21-SEP-2022
 Sampling Kit Packed By: KLUG_K
 Preservatives Needed: NO
 Date to Receive Samples: 07-OCT-2022
 Logged In By: SHAIK_A On: 07-OCT-2022

Send Final Report To:

FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Matthew Bearden
 Report Notification Email:
 matthew.bearden@floridadep.gov; victoria.schwartz@floridadep.gov

Comments: Please log in PKD-QN-BV & DTY-QN-C under a separate event with priority 12, Bottle Group - B

Florence: A and B
 Faulk: C

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-HARD	Template: DEFAULT	(SM 2340 B-2011) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO ₃
Hardness		
W-ICP	Template: SMP-FRESH	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Strontium		
Tin		

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template: SMP-FRESH	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Titanium		
Vanadium		
Zinc		
W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum		
Antimony		
Arsenic		
Barium		
Beryllium		
Boron		
Cadmium		
Chromium		
Cobalt		
Copper		
Lead		
Manganese		
Molybdenum		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: Lugols
DTY-QN-C	Template: DEFAULT	(SOP-AB05_1) No. of diatom taxa of quantitative phytoplankton sample.
PKD-QN-BV	Template: DEFAULT	(SOP-AB05) No. of wet taxa of quantitative phytoplankton sample w/speciation and cell counts.

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.