

Login Checklist

RQ- 2023-01-02-58

Login Station ID: 5

Shipping Method FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 1/11/23 @ 9:30

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
0.3				3		X			585947318142
1.8				12		X			585947317374

COOLER CHECK

*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes X No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No X

CONTAINER CHECK

Evidence Tape on Bottles?

Yes _____ No _____ NA X Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No X Containers intact? Yes X No _____

Caps on tight?

Yes X No _____

Sufficient Sample Volume:

Yes X No _____

If NO is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes X No _____ Init: C-P

Chlorine detected? (One container checked per Field ID)

Yes _____ No X Init: C-P

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH< 2.0?

Yes X No _____ NA _____ Init: C-P

**W-1-4-DIOX (Green dot on container) preserved to pH< 4.0?

Yes _____ No _____ NA X Init: C-P

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: C-P

Date: 1/11/23

Event contains NCRs (Y/N)? N

EventID:

CEN-DIST-2023-01-11-01 (SMP)

Date Received: 11-JAN-2023 09:30

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	221426058 10/17/23
pH Test Strips 0 – 6	223819AV 08/30/23
pH Paper 0 – 3	NA
pH Paper 0 – 13	210620 4/15/2023
Chlorine Test Strips	0070 02/2023

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 1/10/2023
Customer: CEN-DIST
RQ: RQ-2023-01-02-58
Collected By: Terry Riordan, Jessie Atchison
Sampling Agency: FDEP
Project ID: SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB

Number of Coolers Shipped: 2

Delivery Method: FedEx

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
Tiger Lake near Center	 G4CE0206

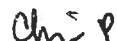
RELINQUISHED BY(SIGNATURE):



DATE/TIME: 1/10/2023

1600 EST

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By: 

Rec'd/Inspected

Date: 1-11-23



RQ-2023-01-02-58	Collected By: Terry Riordan, Jessie Atchison
Customer: FDEP	Field Report Prepared by: Terry Riordan, Jessie Atchison
Project ID: SMP	Sampling Agency: FDEP
WIN ID/Field ID:  WIN ID/Field ID: G4CE0206	Comments: Bottom meter readings were lost when survey was submitted

Location: Tiger Lake near Center

Matrix: Fresh

(Check One)

☒ Grab☐ Composite

Date and Time {EST}	DO (mg/L)	DO (%SAT)	Temp. °C	pH	Sample Depth (m)	Secchi Depth	Total Depth (m)	Sp COND (µmhos /cm)	Salinity (PPTH)	
1/10/2023 1300	11.37	125	20	8.54	0.3	0.5	2.5	104.1	0.05	
Parameter Suite	Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group	
Preserved Nutrients (P-500ML)	X	W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP	X	Ice	X	H2SO4*		<2	1	A
		Lot # SA2229010	Exp:	08/30/2023						
Metals (P-500ML)	X	W-ICPMS	X	W-ICP	X	W-HARD	X	Ice	X	HNO3*
		Lot # NA1105100	Exp:	06/07/2023				<2	1	A
Filtered Nutrients (P125ML)	X	W-PO4-F	X	Field Filtered 0.45 um Filter	X	Ice			1	A
		Filter Lot # 17413308	Syringe Lot # 2090604							
BOD (BP-1L)		OV-BODUNIN				Ice				
Chlorophyll (BP-1L)	X	CHLSUITE-W			X	Ice			1	A
Physical/Aggregate (Fresh) (P-1L)	X	Alkalinity / W-F / W-TSS / TURBIDITY / W-CL-IC / W-COLOR / W-SO4-IC / W-TDS			X	Ice			1	A
Physical/Aggregate (Marine) (P-1L)		Alkalinity / W-F / W-TSS / TURBIDITY				Ice				
Macroinvertebrates (PJ-2L)		MI-FW-QLDC	Formalin Lot #			Formalin				
Periphyton (PT-50ML)		ALGAE-ID				Ice				
Microbiology (Contract Lab Bottle)		E Coli	F Coli	Enterococci		Ice				
		Contract Lab								
Molecular (QPCR-P-500ML)		PCR-HF183	PCR-BacR	PCR-DG3		Ice				
		PCR-GULL2	PCR-GFD							
Tracers (BG-500ML)		W-E8321-DI	W-E8321-MS			Ice				
Pesticides (8G-500ML) (8G-1L)	X	W-E8321-DI	X	W-E8321-MS	X	W-PSNP-TQ		MCAA Exp.	MCAA Lot#	
		W-PCL-TQ-R	X	W-CARB-AA (5 ml of MCAA Buffer Solution)	X	W-CT-CL-R				10
	X	W-TR-SQ-R	X	W-PCB-TQ-R			02/04/2023	128343		A
Phytoplankton	X	DTY-QN-C	X	PKD-QN-BV	X	Ice			1	D
DOC		W-DOC		Field Filtered 0.45 um Filter		Ice		H2SO4*		
				Acid Preservative Expiration Date					<2	
				Acid Preservative Lot Number						

Date of Request: 29-DEC-2022
 Created By: BOMAN_L On: 29-DEC-2022
 Modified By: JASROTIA_P On: 03-JAN-2023
 Customer: CEN-DIST
 Project: SMP
 Division:
 District: Central District
 Event Description: Tiger/Henry/Smart

Send Coolers To:

Phone: 407-897-4181
 FL Dept. of Environmental Protection
 3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Leif Boman
 Cooler Ship EMail: Leif.Boman@FloridaDEP.gov

Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: AYRES_J On: 03-JAN-2023
 Biology Request Reviewed By: JASROTIA_P On: 03-JAN-2023
 Sampling Kit Required: NO Ship on:
 Sampling Kit Shipped:
 Sampling Kit Packed By:
 Preservatives Needed: NO
 Date to Receive Samples: 02-JAN-2023
 Logged In By: KREUCH_B On: 05-JAN-2023

Send Final Report To:

3319 Maguire Blvd., Suite 232
 Orlando, FL 32803-3767
 Attn: Jessie Atchison
 Report Notification EMail: Jessie.Atchison@FloridaDEP.gov

Comments: Please log in PKD-QN-BV & DTY-QN-C under a separate event with priority 12, Bottle Group - D

Bottle Group A: Tiger
 Bottle Group B: Henry
 Bottle Group C: Smart
 Bottle Group D: Tiger Phyto

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: CH ₂ ClCOOH and Ice
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: BG-500ML FILL TO	Number of Bottles: 4	Preserved With: ICE
W-PCB-TQ-R	Template: DEFAULT	(EPA 8270E) PCBs in water matrices by GC/MS/MS
W-PCL-TQ-R	Template: DEFAULT	(EPA 8270E) Organochlorine pesticides, pyrethroids and tracers in water matrices by GC/MS/MS
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: BG-1L	Number of Bottles: 4	Preserved With: ICE
W-PSNP-TQ	Template: DEFAULT	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS
W-TR-SQ-R	Template: DEFAULT	(EPA 8270E) Analysis of Sunscreen Compounds and Flame Retardants in Water by GC/MS-SIM

Bottle Group B (Water) With 1 Samples:

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Group C (Water) With 1 Samples:

Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS

Bottle Group D (Water) With 1 Samples:

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: Lugols
DTY-QN-C	Template: DEFAULT	(SOP-AB05_1) No. of diatom taxa of quantitative phytoplankton sample.
PKD-QN-BV	Template: DEFAULT	(SOP-AB05) No. of wet taxa of quantitative phytoplankton sample w/speciation and cell counts.

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.