

RQ-2023-12-04-46

Login Checklist

Login Station ID: #4

Shipping Method: ~~FedEx~~ | UPS | HD | Greyhound

Date/Time of Receipt: 12/ 1 /2023@09:22

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
1.3				12					648484606043

COOLER CHECK*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes X No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No X**CONTAINER CHECK**

Evidence Tape on Bottles?

Yes _____ No _____ NA X Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No _____ Containers intact? Yes X No _____

Caps on tight?

Yes X No _____

Sufficient Sample Volume:

Yes X No _____If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____**CHLORINE CHECK REQUIRED?** (Blue dot on container)Yes X No _____ Init: TL

Chlorine detected? (One container checked per Field ID)

Yes _____ No X Init: TL

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH < 2.0?

Yes X No _____ NA _____ Init: TL

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes _____ No _____ NA X Init: TL

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: TL

Date: 12/ 1 /2023

Event contains NCRs (Y/N)? N

Event ID: SW-DIST-2023-12-01-01

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	IR TEMP GUN #4 Exp:10/16/2024
pH Test Strips 0 – 6	LOT# 223819AV Exp: AUG 30, 2023
pH Paper 0 – 3	N/A
pH Paper 0 – 13	LOT# 207621 Exp: March 15, 2024
Chlorine Test Strips	LOT# 2313 Exp: July 2025

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 11/30/2023
Customer: SW-DIST
Collected By: Megan Maloney, Hannah Henning
Sampling Agency: FDEP
RQ: RQ-2023-12-04-46
Project ID: SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB


Number of Coolers Shipped: 2

Delivery Method: FedEx

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
LAKE ROUSSEAU AT CENTER	 G4SW0118

RELINQUISHED BY(SIGNATURE):

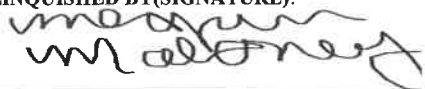

DATE/TIME: 11/30/2023
1700 EST

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By: TL

Rec'd/Inspected 12/1/23
Date:



RQ-2023-12-04-46

Collected By: Megan Maloney, Hannah Henning

Customer: FDEP

Field Report Prepared by: Megan Maloney,
Hannah Henning

Project ID: SMP

Sampling Agency: FDEP

Comments:

WIN ID/Field ID:



WIN ID/Field ID: G4SW0118

Location: LAKE ROUSSEAU AT CENTER

Matrix: Fresh

(Check One)

☒ Grab☐ Composite

Date and Time {EST}	DO (mg/L)	DO (%SAT)	Temp. °C	pH	Sample Depth (m)	Secchi Depth	Total Depth (m)	Sp COND (µmhos /cm)	Salinity (PPTH)
11/30/2023 1050	10	101.2	16	7.91	0.3	0.7	0.7	315.6	0.15

Parameter Suite	Parameter Methods				Preservation		pH Check		# Bottles sent to Lab	Bottle Group		
Preserved Nutrients (P-500ML)	X	W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP			X	Ice	X	H2SO4*	X	<2	1	A
		Lot #		SA3193040	Exp:	07/31/2024						
Metals (P-500ML)		W-ICPMS		W-ICP	W-HARD	Ice		HNO3*		<2		
		Lot #		Exp:								
Filtered Nutrients (P125ML)		W-PO4-F		Field Filtered 0.45 µm Filter		Ice						
BOD (BP-1L)		OV-BODUNIN				Ice						
Chlorophyll (BP-1L)	X	CHLSUITE-W				Ice						
Physical/Aggregate (Fresh) (P-1L)	X	Alkalinity / W-F / W-TSS / TURBIDITY / W-CL-IC / W-COLOR / W-SO4-IC / W-TDS				Ice						
Physical/Aggregate (Marine) (P-1L)		Alkalinity / W-F / W-TSS / TURBIDITY				Ice						
Macroinvertebrates (PJ-2L)		MI-FW-QLDC	Formalin Lot #			Formalin						
Periphyton (PT-50ML)		ALGAE-ID				Ice						
Microbiology (Contract Lab Bottle)		E Coli	F Coli	Enterococci		Ice						
		Contract Lab										
Molecular (QPCR-P-500ML)		PCR-HF183	PCR-BacR	PCR-DG3		Ice						
		PCR-GULL2	PCR-GFD									
Tracers (BG-500ML)		W-E8321-DI	W-E8321-MS			Ice						
Pesticides (BG-500ML) (BG-1L)	X	W-E8321-DI	X	W-E8321-MS	X	W-PSNP-TQ	x	Ice	MCAA Exp.	MCAA Lot#	9	A
	X	W-PCL-TQ-R	W-CARB-AA (5 ml of MCAA Buffer Solution)		X	W-CT-CL-R						
	X	W-TR-SQ-R	X	W-PCB-TQ-R								
Phytoplankton (BP-1L)		DTY-QN-C		PKD-QN-BV		Ice						
DOC (P-500ML)		W-DOC	Field Filtered 0.45 µm Filter		Ice		H2SO4*		<2			
	Acid Preservative Expiration Date											
	Acid Preservative Lot Number											

Date of Request: 09-NOV-2023
 Created By: PASWATER_B On: 09-NOV-2023
 Modified By: JASROTIA_P On: 16-NOV-2023
 Customer: SW-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Southwest District
 Event Description: Run 1

Send Coolers To:

Phone: 813-632-7600
 FL Dept. of Environmental Protection
 13051 N. Telecom Parkway
 Temple Terrace, FL 33637-0926
 Attn: Ben Paswater
 Cooler Ship EMail: benjamin.paswater@dep.state.fl.us

Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: AYRES_J On: 16-NOV-2023
 Biology Request Reviewed By: JASROTIA_P On: 16-NOV-2023
 Sampling Kit Required: YES Ship on: 29-NOV-2023
 Sampling Kit Shipped: YES On: 21-NOV-2023
 Sampling Kit Packed By: SHAIK_A
 Preservatives Needed: NO
 Date to Receive Samples: 04-DEC-2023
 Logged In By: LEBRUN_T On: 01-DEC-2023

Send Final Report To:

FL Dept. of Environmental Protection
 13051 N. Telecom Parkway
 Temple Terrace, FL 33637-0926
 Attn: Ben Paswater
 Report Notification EMail:
 benjamin.paswater@dep.state.fl.us;sarah.meyer@dep.state.fl.us
 ;Devin.Mathias@Floridadep.gov

Comments:**Bottle Group A (Water) With 1 Samples:**

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10150 B) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: BG-500ML FILL TO	Number of Bottles: 4	Preserved With: ICE
W-CT-CI-R	Template: DEFAULT	(EPA 8276) Chlordane and Toxaphene in water matrices by GC/MS with Negative Chemical Ionization
W-PCB-TQ-R	Template: DEFAULT	(EPA 8270E) PCBs in water matrices by GC/MS/MS
W-PCL-TQ-R	Template:	(EPA 8270E) Organochlorine pesticides, pyrethroids and tracers in water matrices by GC/MS/MS
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: BG-500ML	Number of Bottles: 4	Preserved With: ICE

Bottle Group A (Water) With 1 Samples:

Bottle Type: BG-500ML	Number of Bottles: 4	Preserved With: ICE
W-PSNP-TQ	Template: DEFAULT	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS
W-TR-SQ-R	Template:	(EPA 8270E) Analysis of Sunscreen Compounds and Flame Retardants in Water by GC/MS-SIM

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
FCOLI-MF	Template: DEFAULT	(SM 9222 D) Fecal coliforms by membrane filter method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10150 B) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.