

RQ-2023-09-11-27

Login Checklist

Login Station ID: #4

Shipping Method: ~~FedEx~~ | UPS | HD | Greyhound

Date/Time of Receipt: 10/ 6 /2023@09:04

*Cooler Temperature	Samples Frozen?	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Waybill Tracking #
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
2.8				8		X			585947320705
0.7				6		X			648484582281

COOLER CHECK*All samples received on ice unless otherwise noted. HNO₃ and Formalin preserved samples do NOT have a temperature requirement.

- If the temperature of a cooler exceeds 6.0 oC (above 10.0 oC for W-1-4-DIOX and microbiology) or received without ice, identify affected samples in an NCR. NCR# _____
- DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C for W-1-4-DIOX) without login confirmation from customer for samples requiring ice preservation; identify affected samples in an NCR. NCR# _____
- If cooler or samples are received with a damaged Evidence Seal; identify the affected samples in an NCR. NCR# _____
- If coolers or samples are received late; identify the affected samples in an NCR. NCR# _____

Chain of Custody/ Field Sheet(s) Included?

Yes X No _____

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included?

Yes _____ No X**CONTAINER CHECK**

Evidence Tape on Bottles?

Yes _____ No _____ NA X Evidence Tape intact? Yes _____ No _____

If Criminal, photographs taken?

Yes _____ No _____

Containers intact? Yes X No _____

Caps on tight?

Yes X No _____

Sufficient Sample Volume:

Yes X No _____If **NO** is checked above, generate an NCR listing affected sample IDs in its appropriate category. NCR# _____**CHLORINE CHECK REQUIRED?** (Blue dot on container)Yes X No _____ Init: TL

Chlorine detected? (One container checked per Field ID)

Yes _____ No X Init: TL

If chlorine is detected, generate an NCR listing affected sample IDs . NCR# _____

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH< 2.0?

Yes X No _____ NA _____ Init: TL

**W-1-4-DIOX (Green dot on container) preserved to pH< 4.0?

Yes _____ No _____ NA X Init: TL

If samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

Coolers Unpacked/Checked by: TLDate: 10/ 6 /2023Event contains NCRs (Y/N)? NEvent ID: SE-DIST-2023-10-06-02

Login Checklist

Login Preservation Equipment

Lot # & Exp: Date

Additional Comments:

Infrared Thermometer ID:	IR TEMP GUN #4 Exp:10/17/2023
pH Test Strips 0 – 6	LOT# 223819AV Exp: AUG 30, 2023
pH Paper 0 – 3	N/A
pH Paper 0 – 13	LOT# 207621 Exp: March 15, 2024
Chlorine Test Strips	LOT# 2313 Exp: July 2025

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ Init: _____

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ No _____

If no, were samples shipped on dry ice? Yes _____ No _____

If S-VOC-MS samples were not preserved correctly, generate an NCR listing affected sample IDs . NCR# _____

FLORIDA DEPARTMENT OF ENVIRONMENTAL PROTECTION
CHAIN OF CUSTODY FORM

Date: 10/5/2023 **Collected By:** Jordan Skaggs, Natalie Howard
Customer: SE-DIST **Sampling Agency:** FDEP
RQ: RQ-2023-09-11-27 **Project ID:** SMP

PLEASE RETURN THE ORIGINAL OF THIS FORM WITH SAMPLE SUBMITTAL
FORMS AND THE COOLER PACKING WORKSHEET TO THE LAB


Number of Coolers Shipped: 2

Delivery Method: FedEx

FIELD ID (Labels) submitted under this RQ for this collection date.

Site Location	Field ID/WIN ID
SEVENMILE SLOUGH AT PEAVINE TRAIL	 G4SE0153
OAK CREEK AT 68	 G4SE0080

RELINQUISHED BY(SIGNATURE):

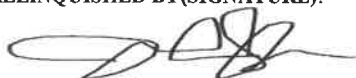

DATE/TIME: 10/5/2023
1600 EST

THIS SECTION TO BE COMPLETED BY LABORATORY

Rec'd/Inspected By: TL

Rec'd/Inspected Date: 10/6/23



RQ-2023-09-11-27

Customer: FDEP

Project ID: SMP

Collected By: Jordan Skaggs, Natalie Howard

Field Report Prepared by: Jordan Skaggs, Natalie Howard

Sampling Agency: FDEP

Comments:

WIN ID/Field ID:



WIN ID/Field ID: G4SE0153

Location: SEVENMILE SLOUGH AT PEAVINE TRAIL

Matrix: Fresh

(Check One)

☒ Grab☐ Composite

Date and Time {EST}	DO (mg/L)	DO (%SAT)	Temp. °C	pH	Sample Depth (m)	Secchi Depth	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPTH)
10/5/2023 1030	2.53	31.5	26.5	5.53	0.3	0.6	1.3	60.1	0.03

Parameter Suite	Parameter Methods				Preservation		pH Check		# Bottles sent to Lab	Bottle Group		
Preserved Nutrients (P-500ML)	X	W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP			X	Ice	X	H2SO4*	X	<2	1	C
		Lot #		SA3193140	Exp:	07/31/2024						
Metals (P-500ML)		W-ICPMS		W-ICP	W-HARD	Ice	HNO3*		<2			
		Lot #			Exp:							
Filtered Nutrients (P125ML)		W-PO4-F		Field Filtered 0.45 um Filter		Ice						
		Filter Lot #		Syringe Lot #								
BOD (BP-1L)		OV-BODUNIN				Ice						
Chlorophyll (BP-1L)	X	CHLSUITE-W				X	Ice			1	C	
Physical/Aggregate (Fresh) (P-1L)	X	Alkalinity / W-F / W-TSS / TURBIDITY / W-Cl-IC / W-COLOR / W-SO4-IC / W-TDS				X	Ice			1	C	
Physical/Aggregate (Marine) (P-1L)		Alkalinity / W-F / W-TSS / TURBIDITY					Ice					
Macroinvertebrates (PJ-2L)		MI-FW-QLDC	Formalin Lot #			Formalin						
Periphyton (PT-50ML)		ALGAE-ID				Ice						
Microbiology (Contract Lab Bottle)	X	E Coli	F Coli	Enterococci	X	Ice				1	C	
		Contract Lab										
Molecular (QPCR-P-500ML)		PCR-HF183	PCR-BacR	PCR-DG3		Ice						
		PCR-GULL2	PCR-GFD									
Tracers (BG-500ML)		W-E8321-DI	W-E8321-MS	W-PSNP-TQ		Ice						
Pesticides (BG-500ML) (BG-1L)		W-E8321-DI	W-E8321-MS	W-PSNP-TQ		Ice	MCAA Exp.	MCAA Lot#				
		W-PCL-TQ-R	W-CARB-AA (5 ml of MCAA Buffer Solution)	W-CT-CL-R								
		W-TR-SQ-R	W-PCB-TQ-R									
Phytoplankton (BP-1L)		DTY-QN-C	PKD-QN-BV			Ice						
DOC (P-500ML)		W-DOC	Field Filtered 0.45 um Filter			Ice	H2SO4*		<2			
			Acid Preservative Expiration Date									
			Acid Preservative Lot Number									



RQ-2023-09-11-27

Customer: FDEP

Project ID: SMP

Collected By: Jordan Skaggs, Natalie Howard

Field Report Prepared by: Jordan Skaggs, Natalie Howard

Sampling Agency: FDEP

Comments:

WIN ID/Field ID:



WIN ID/Field ID: G4SE0080

Location: OAK CREEK AT 68

Matrix: Fresh

(Check One)

☒ Grab☐ Composite

Date and Time {EST}	DO (mg/L)	DO (%SAT)	Temp. °C	pH	Sample Depth (m)	Secchi Depth	Total Depth (m)	Sp COND (µmhos/cm)	Salinity (PPT)			
10/5/2023 1115	3.06	39	27.3	6.01	0.3	0.5	1.8	49.6	0.02			
Parameter Suite		Parameter Methods				Preservation		pH Check	# Bottles sent to Lab	Bottle Group		
Preserved Nutrients (P-500ML)		W-NH3 / W-NO2NO3 / W-TKN / W-TOC / W-S-A-TP				Ice		H2SO4*	<2			
		Lot #		Exp:								
Metals (P-500ML)		W-ICPMS				W-HARD		HNO3*	<2			
		Lot #		Exp:								
Filtered Nutrients (P125ML)		W-PO4-F				Field Filtered 0.45 um Filter						
		Filter Lot #		Syringe Lot #		Ice						
BOD (BP-1L)		OV-BODUNIN				Ice						
Chlorophyll (BP-1L)		CHLSUITE-W				Ice						
Physical/Aggregate (Fresh) (P-1L)		Alkalinity / W-F / W-TSS / TURBIDITY / W-CI-IC / W-COLOR / W-SO4-IC / W-TDS				Ice						
Physical/Aggregate (Marine) (P-1L)		Alkalinity / W-F / W-TSS / TURBIDITY				Ice						
Macroinvertebrates (PJ-2L)		MI-FW-QLDC				Formalin Lot #						
Periphyton (PT-50ML)		ALGAE-ID				Ice						
Microbiology (Contract Lab Bottle)		X E Coli				F Coli		Enterococci	1	D		
		Contract Lab				x		Ice				
Molecular (QPCR-P-500ML)		X PCR-HF183				X PCR-BacR		X PCR-DG3	1	D		
		PCR-GULL2				PCR-GFD						
Tracers (BG-500ML)		W-E8321-DI				W-E8321-MS		Ice				
Pesticides (BG-500ML) (BG-1L)		X W-E8321-DI	X	W-E8321-MS		X	W-PSNP-TQ		07/13/2024	143355		
		X W-PCL-TQ-R	X	W-CAR8-AA (5 ml of MCAA Buffer Solution)		X	W-CT-CL-R					
		X W-TR-SQ-R	X	W-PCB-TQ-R								
Phytoplankton (BP-1L)		DTY-QN-C				PKD-QN-BV		Ice				
DOC (P-500ML)		W-DOC				Field Filtered 0.45 um Filter		Ice	<2			
						Acid Preservative Expiration Date						
						Acid Preservative Lot Number						

Date of Request: 25-JUL-2023
 Created By: SKAGGS_J On: 25-JUL-2023
 Modified By: JASROTIA_P On: 22-AUG-2023
 Customer: SE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Southeast District
 Event Description: SMP - Basinger

Send Coolers To:

Phone: (772) 448-5884
 FL Dept. of Environmental Protection
 2199 South Rock Road
 Fort Pierce, FL 34945
 Attn: Jordan Skaggs
 Cooler Ship EMail: jordan.skaggs@FloridaDEP.gov

Priority: 3
 Event Status: R
 Criminal Investigation: NO
 Chemistry Request Reviewed By: AYRES_J On: 22-AUG-2023
 Biology Request Reviewed By: JASROTIA_P On: 22-AUG-2023
 Sampling Kit Required: YES Ship on: 01-SEP-2023
 Sampling Kit Shipped: YES On: 23-AUG-2023
 Sampling Kit Packed By: KREUCH_B
 Preservatives Needed: YES
 Date to Receive Samples: 22-SEP-2023
 Logged In By: KREUCH_B On: 26-SEP-2023

Send Final Report To:

FL Dept. of Environmental Protection
 2199 South Rock Road
 Fort Pierce, FL 34945
 Attn: Andrew Luering
 Report Notification EMail:
 jordan.skaggs@floridadep.gov;andrew.luering@floridadep.gov

Comments: #####

Buffer Solution for W-CARB-AA Bottle
 1 vial per bottle

#####

Sampler Comments:

A
 G4SE0142
 - Nutrients, E. coli, Pest Kit

B
 G4SE0144
 - Nutrients, Metals, E. coli

C
 G4SE0153
 - Nutrients, E. coli

D
 G4SE0080
 - E. coli, PCR, Pest Kit

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10150 B) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: CH2CICOOH and Ice
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS
Bottle Type: BG-500ML FILL TO	Number of Bottles: 4	Preserved With: ICE
W-CT-CI-R	Template: DEFAULT	(EPA 8276) Chlordane and Toxaphene in water matrices by GC/MS with Negative Chemical Ionization
W-PCB-TQ-R	Template: DEFAULT	(EPA 8270E) PCBs in water matrices by GC/MS/MS
W-PCL-TQ-R	Template: DEFAULT	(EPA 8270E) Organochlorine pesticides, pyrethroids and tracers in water matrices by GC/MS/MS
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: BG-500ML	Number of Bottles: 4	Preserved With: ICE
W-PSNP-TQ	Template: DEFAULT	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS
W-TR-SQ-R	Template:	(EPA 8270E) Analysis of Sunscreen Compounds and Flame Retardants in Water by GC/MS-SIM

Bottle Group B (Water) With 1 Samples:

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10150 B) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340 B-2011) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template: SMP-FRESH	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Strontium		
Tin		
Titanium		
Vanadium		
Zinc		
W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum		
Antimony		
Arsenic		
Barium		
Beryllium		
Boron		
Cadmium		

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICPMS	Template: DEFAULT	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Chromium		
Cobalt		
Copper		
Lead		
Manganese		
Molybdenum		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(DEP SOP: NU-094-1) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500-F- C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2015) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2015) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10150 B) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems

Bottle Type: QPCR-P-500ML	Number of Bottles: 1	Preserved With: ICE
---------------------------	----------------------	---------------------

Bottle Group D (Water) With 1 Samples:

Bottle Type: QPCR-P-500ML		Number of Bottles: 1	Preserved With: ICE
PCR-BACR	Template: DEFAULT	(SOP-PCR-1.4) Detection and quantification of Ruminant specific Bacteroidetes BacR genetic marker with quantitative polymerase chain reaction	
PCR-DG3	Template: DEFAULT	(SOP-PCR-1.5) Detection and quantification of canine-specific DG3 Bacteroides genetic marker with quantitative polymerase chain reaction	
PCR-HF183	Template: DEFAULT	(SOP-PCR-1.0) Detection and quantification of human-specific HF183 Bacteroides genetic marker with quantitative polymerase chain reaction	
Bottle Type: BG-500ML		Number of Bottles: 1	Preserved With: CH ₂ ClCOOH and Ice
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS	
Bottle Type: BG-500ML FILL TO		Number of Bottles: 2	Preserved With: ICE
W-CT-CI-R	Template: DEFAULT	(EPA 8276) Chlordane and Toxaphene in water matrices by GC/MS with Negative Chemical Ionization	
W-PCB-TQ-R	Template: DEFAULT	(EPA 8270E) PCBs in water matrices by GC/MS/MS	
W-PCL-TQ-R	Template: DEFAULT	(EPA 8270E) Organochlorine pesticides, pyrethroids and tracers in water matrices by GC/MS/MS	
Bottle Type: BG-500ML		Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS	
W-E8321-MS	Template: SMP	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS	
Bottle Type: BG-500ML		Number of Bottles: 2	Preserved With: ICE
W-PSNP-TQ	Template: DEFAULT	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS	
W-TR-SQ-R	Template:	(EPA 8270E) Analysis of Sunscreen Compounds and Flame Retardants in Water by GC/MS-SIM	

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.