

Log-in Checklist

RQ- 2019-05-06-18

Login Station ID: #5

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 5/ 8 /19 @ 10:10

Cooler Temperature*	Samples Frozen	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
1.4 -0.1				14 10		✓ ✓			4751 8779 0608 4751 8779 0619

* HNO₃ and Formalin persevered samples do **NOT** have a temperature requirement. For all others, if the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), or if Evidence Seal is damaged. Identify the containers in the affected cooler(s) on back of this form. **DOH coolers can exceed 2.0 °C above standard temperature without customer notification/confirmation.**

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

**Chain of Custody/ Field Sheet(s) Included? Yes ✓ No

Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes No ✓

CONTAINER CHECK

**Evidence Tape on Bottles? Yes No ✓ If yes, is it intact? Yes No

**Caps on tight? Yes ✓ No **Containers intact? Yes ✓ No

**Sufficient Sample Volume: Yes ✓ No

CHLORINE CHECK REQUIRED? (Blue dot on container) Yes No ✓ Init: KK

Check the Box marked "Yes" or "No" for the Presence of Chlorine. If Cl detected, an NCR is required.					
Analysis	Yes	No	Analysis	Yes	No
W-1,4 DIOX			W-PC-AA-SF		
W-E8321-DI, -MS, -21			W-NP-AA-SF		
W-BNA (-625, -R)			W-PDQUAT		
W-CARB-AA (-CAB-AA-R)			W-PCL-SQ-R		
W-U2060-MS			W-TR-TQ		
W-PCB-R			W-CT-CI-R		
W-PAH (-R)			W-PSNP-TQ		
W-PFAS-MS					

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0? Yes ✓ No NA Init: KK

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0? Yes No NA ✓ Init: KK

Coolers Unpacked/Checked by: KK

Date: 5/ 8 /19

NCRs (Y/N)? N

Event ID: Nassau West
NE-DIST-2019-05-08-01 (SMP)

NCR #

Page 1 of 2, 03/25/2018 Date Received: 08-MAY-2019 10:10

Log-in Checklist

Samples with Incorrect pH Preservation or presence of Chlorine

Field ID	Bottle Test ID(s)	pH / Cl	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle Test ID(s)

Bottles Not Intact

Field ID	Bottle Test ID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle Test ID(s)	Action Taken

Additional Comments:

Infrared Thermometer ID:	
pH Paper 0.0 – 3.0 Lot #	220416A EXP:7/30/20
pH Paper 0 – 13 Lot #	222516 EXP:9/15/2019
Chlorine Test Strips Lot #	091417V EXP:9/2021

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes_____ No X

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes_____ If No_____

no, were samples shipped on dry ice? Yes_____ No_____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check (W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes_____ No_____

Date and time samples preserved: _____

Florida Department of Environmental Protection Central Laboratory Submittal Form

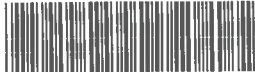
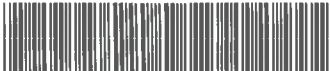
Customer: NE-DIST

Requester: Jeremy M. Parrish

Field Report Prepared By: Jeremy Parrish

Project ID: SMP

Collected By: Jeremy Parrish / Cheyenne RivenSampling Agency: FDEP NE ROC

Location	WIN ID/Field ID:		19010021	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>5/7/19</u> Time <u>950</u>	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s) (See pg 2)	
Field ID	ST MARYS S PRONG @ US 90			Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>	Diss. Oxygen <u>4.8</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	A	
WIN/STOR				Temp (C) <u>22.8</u>	pH <u>7.60</u>	Sample Depth <u>0.3</u>	☐ ft ☒ m		Sp. Conductance (umho/cm) <u>106.107</u>
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.05</u>	Comments				
Location	WIN ID/Field ID:		19010136	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>5/7/19</u> Time <u>1015</u>	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID	DAUGHTERY BR @ CREWS RD			Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>	Diss. Oxygen <u>6.3</u>	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	A	
WIN/STOR				Temp (C) <u>21.6</u>	pH <u>7.8</u>	Sample Depth <u>0.09</u>	☐ ft ☒ m		Sp. Conductance (umho/cm) <u>140</u>
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>0.07</u>	Comments				
Location	WIN ID		19010136	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>5/7/19</u> Time <u>1020</u>	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID	Field ID		19010136DUP	Matrix (type, e.g., Salt, Fresh, etc) <u>Fresh</u>	Diss. Oxygen <u>—</u>	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	A	
WIN/STOR	DAUGHTERY BR @ CREWS RD DUP			Temp (C) <u>—</u>	pH <u>—</u>	Sample Depth <u>0.09</u>	☐ ft ☒ m		Sp. Conductance (umho/cm)
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>—</u>	Comments				
Location	Field ID:		SMP050719FB	☐ Comp ☒ Grab	Collection (comp begin or grab) Date <u>5/7/19</u> Time <u>1025</u>	Eastern Central	Composite end Date _____ Time _____	Bottle Group(s)	
Field ID	Field Blank 050719			Matrix (type, e.g., Salt, Fresh, etc) <u>Field Blank</u>	Diss. Oxygen <u>—</u>	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	C	
WIN/STO				Temp (C) <u>—</u>	pH <u>—</u>	Sample Depth <u>—</u>	☐ ft ☐ m		Sp. Conductance (umho/cm)
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) <u>—</u>	Comments <u>Add w/PO4 please JP</u>				

Relinquished By:	Date/Time	Shipping Method	Number of Coolers Shipped:	Received By:	Date/Time
<u>Jeremy Parrish</u>	<u>5/7/19 1600</u>	<u>Fedex</u>	<u>2</u>	<u>KH</u>	<u>5.8.19 110:10</u>

Group: A	Bottle Type:	P-1L	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab	4
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						4
Group: A	Bottle Type:	BP-1L	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab	
	CHLSUITE-W						
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	NE-ALS-ECQ						
Group: A	Bottle Type:	P-500ML	# of Bottles: 4	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	4
	W-ICP						
	W-ICPMS						
Group: A	Bottle Type:	P-500ML	# of Bottles: 4	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	4
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: A	Bottle Type:	P-125ML	# of Bottles: 4	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	4
	W-PO4-F						
Group: B	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	NE-ALS-ECQ						
Group: C	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Nassau West

last revised Jan 25, 2018

Customer: NE-DIST

Requester: Jeremy M. Parrish

Field Report Prepared By:

Jeremy Parrish

Project ID: SMP

Collected By:

Jeremy Parrish / Cheyenne Riden

Sampling Agency:

FDEP MROC

Location	WIN ID/Field ID:		19010023	ST MARYS S PRONG @ CR 23C	☐ Comp ☒ Grab	Collection (comp begin or grab)	Eastern	Composite end	Bottle
Field ID					Date	5/7/19	Time	1045	Group(s)
WIN/STORET					Matrix (type, e.g., Salt, Fresh, etc)	Fresh	Diss. Oxygen	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	0.06	Temp (C)	22.6	pH	7.87
						Sample Depth	0.3	☐ ft ☒ m	Sp. Conductance (umho/cm)
									126
Comments									

Location	Field ID	WIN/STORET #	Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	☐ dd ☐ dms	Comments

Location	Field ID	WIN/STORET #	Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	☐ dd ☐ dms	Comments

Location	Field ID	WIN/STORET #	Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	☐ dd ☐ dms	Comments

Relinquished By:	Date/Time	Shipping Method	Number of Coolers Shipped:	Received By:	Date/Time
JParrish	5/7/19 1600	Fedex	2	KH	5-8-19 / 10:10

Group: A	Bottle Type:	P-1L	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-250ML-WI-THI	# of Bottles: 4	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	NE-ALS-ECQ					
Group: A	Bottle Type:	P-500ML	# of Bottles: 4	Preservative:	HNO3	Enter Number of Bottles Sent to Lab _____
	W-ICP					
	W-ICPMS					
Group: A	Bottle Type:	P-500ML	# of Bottles: 4	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3					
	W-NO2NO3					
	W-S-A-TP					
	W-TKN					
	W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 4	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-250ML-WI-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	NE-ALS-ECQ					
Group: C	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					

Group: C	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	W-TDS						
	W-TSS						

Group: C	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
	W-ICP						
	W-ICPMS						

Group: C	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						

Group C

Bottle type

W-P04-F

125ml # bottle\$ | preserve ICE Enter # bottles sent lab

~~add~~
Please
Add

Date of Request: 22-MAR-2019
 Created By: PARRISH_J On: 22-MAR-2019
 Modified By: JASROTIA_P On: 17-APR-2019
 Customer: NE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Northeast District
 Sampling Event: Nassau West

Send Coolers To:

Phone: 904-807-3300
 8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Send Final Report To:

8800 Baymeadows Way West
 Suite 100
 Jacksonville, FL 32256
 Attn: Jeremy Parrish

Program Module Number:

Priority: 3

Event Status: S

Criminal Investigation: NO

Chemistry Request Reviewed By: WATTS_J On: 17-APR-2019

Biology Request Reviewed By: JASROTIA_P On: 17-APR-2019

Sampling Kit Required: YES Ship on: 26-APR-2019

Sampling Kit Shipped: YES On: 22-APR-2019

Sampling Kit Packed By: KLUG_K

Preservatives Needed: NO

Date to Receive Samples: 06-MAY-2019

Received By:

Report Type: Final Only

FTP Data: NO

QC Report: YES

Date Log: NO

Authorization Log: NO

Comments: Group A: Nutrients, Metals, Bacteria
 Group B: Bacteria

Packing Notes:
 None

Bottle Group A (Water) With 4 Samples:

Bottle Type: P-1L	Number of Bottles: 4	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO ₄ -IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 4	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WI-THI	Number of Bottles: 4	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville
Bottle Type: P-500ML	Number of Bottles: 4	Preserved With: HNO ₃
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		

Bottle Group A (Water) With 4 Samples:

Bottle Type: P-500ML	Number of Bottles: 4	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 4	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 4	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-250ML-WI-THI	Number of Bottles: 1	Preserved With: ICE
NE-ALS-ECQ	Template: DEFAULT	(SM 9223 B Quanti-Tray) Escherichia coli by Colilert Quanti-Tray method by ALS Environmental Laboratory in Jacksonville

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-00) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-97) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B) True Color measured at 450 nm
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-97) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-97) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-97) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects

Barium
Calcium
Iron
Magnesium
Potassium
Sodium
Zinc

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-500ML

Number of Bottles: 1

Preserved With: HNO3

W-ICPMS

Template:

(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects

*Boron
Aluminum
Antimony
Arsenic
Cadmium
Chromium
Copper
Lead
Nickel
Selenium
Silver
Thallium

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.