

RQ- 2020-03-16-32

Login Checklist

Login Station ID: #2

Shipping Method: FedEx | UPS | HD | Greyhound

Date/Time of Receipt: 3/18/2020 @ 10:35

Cooler Temperature*	Samples Frozen	* All Samples in Cooler Preserved with		Number of Sample Containers in Cooler	Evidence Tape				Tracking # Damaged waybills only
					Present?		*Intact?		
	YES	HNO ₃	Formalin		Yes	No	Yes	No	
2.2				10		✓			1278 4012 4818
3.8				22		✓			1278 4012 4829
2.1				20 18 IM 3/18/20					1278 4012 4807

*All samples received on ice unless otherwise noted. If the temperature of a cooler is above 6.0 °C (or above 10.0 °C for W-1-4-DIOX and microbiology), received without ice, or if Evidence Seal is damaged; identify the containers in the affected cooler(s) on back of this form. HNO₃ and Formalin preserved samples do **NOT** have a temperature requirement. **DOH coolers can exceed 2.0 °C above the standard temperature (6.0 or 10.0 °C) without login confirmation from customer for samples requiring ice preservation; An NCR still is required.**

**Chain of Custody/ Field Sheet(s) Included?

Yes ☒No ☐Micro-Biology Overflow Chain of Custody/ Field Sheet(s) Included? Yes ☐No ☒

****Note:** If "No" is checked below, fill out the back of this form (Field ID, Test IDs, etc.) and generate an NCR.

CONTAINER CHECK

Evidence Tape on Bottles? Yes ☐No ☐If yes, is it intact? Yes ☐No ☐Caps on tight? Yes ☒No ☐**Containers intact? Yes ☒No ☐**Sufficient Sample Volume: Yes ☒No ☐

CHLORINE CHECK REQUIRED? (Blue dot on container)

Yes ☒No ☐

Init: IM

KML

Chlorine detected? (One container checked per Field ID)

Yes ☐No ☒

Init: IM

IM

If "YES", list affected Field IDs and Test IDs below.

Field ID	List TestID(s) where Chlorine was Detected, an NCR is required.

PRESERVATION CHECK (EXCLUDING SAMPLES SENT TO OVERFLOW LABS, "OV-")

**Acid Preserved Samples: pH ≤ 2.0?

Yes ☒No ☐NA ☐

Init: IM

KML

**W-1-4-DIOX (Green dot on container) preserved to pH < 4.0?

Yes ☐No ☐NA ☐

Init: IM

Coolers Unpacked/Checked by: IM

KML

Date: 3/18/2020

NCRs (Y/N)?

N

Event ID:

SMP- Martin 2/5

SE-DIST-2020-03-18-01 (SMP)

NCR #

Date Received: 18-MAR-2020 10:35

IM 3/18/20

Login Checklist

Samples with Incorrect pH Preservation

Field ID	Bottle TestID(s)	pH	Action Taken

Samples in Coolers with Incorrect Preservation or Damaged Evidence Tape

Cooler ID	Field ID	Bottle TestID(s)

Bottles Not Intact

Field ID	Bottle TestID(s)	Loose Cap	Damaged Cap		Damaged Container		Evidence Tape Not Intact
			CK	BR	CK	BR	

Samples with Insufficient Volume for Analyses

Field ID	Bottle TestID(s)	Action Taken

Infrared Thermometer ID:	IR TEMP GUN #2 exp:11/5/20
pH Test Strips 0 – 6 Lot #	213316AV Exp:5/15/20
pH Paper 0 – 3 Lot #	
pH Paper 0 – 13 Lot #	208717 EXP: 3/31/20
Chlorine Test Strips Lot #	8261 EXP:7/2021

Additional Comments:

ENCORE SAMPLES

S-VOC-MS samples in Encores must be frozen within 48 hours of collection.

Were Encore samples received? Yes _____ No ✓

Date and time placed in freezer _____

Were all samples placed in freezer within 48 hours of collection? Yes _____ If No _____

no, were samples shipped on dry ice? Yes _____ No _____

FOR CLEAN LAB USE ONLY – Total Mercury Preservation Check(W-HG-AF;W-HG-AF-F):

Samples preserved within 48 hours? Yes _____ No _____

Date and time samples preserved: _____

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: SE-DIST

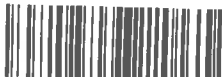
Requester: Jordan Skaggs

Field Report Prepared By: Kristene Taylor

Project ID: SMP

Collected By: Alex Miranda

Sampling Agency: 21FLWPB

Loca	Location: Roebuck Creek @ Locks		☐ Comp	Collection (comp begin or grab)	☒ Grab	Date 03/17/20	Time 10:15	Eastern Central	Composite end	Date	Time	Bottle Group(s)
Field	WIN ID/Field ID:		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☒ mg/L		Total Res. Chlorine		E	
WIN			Marine		3.28		☐ % sat		(mg/L)			
Latitude		☐ dd ☐ dms	Longitude		Temp (C) 21.5		pH 7.24		Sample Depth 0.15			Sp. Conductance (umho/cm) 649
Salinity (ppt)		☐ dd ☐ dms	0.32		Comments							
Loca	Location: South Fork SLR @ Hosford Park		☐ Comp	Collection (comp begin or grab)	☒ Grab	Date 03/17/20	Time 1040	Eastern Central	Composite end	Date	Time	Bottle Group(s)
Field	WIN ID/Field ID:		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☒ mg/L		Total Res. Chlorine		D	
WIN			marine		5.75		☐ % sat		(mg/L)			
Latitude		☐ dd ☐ dms	Longitude		Temp (C) 24.7		pH 7.52		Sample Depth 0.3			Sp. Conductance (umho/cm) 4508
Salinity (ppt)		☐ dd ☐ dms	2.40		Comments							
Loca	Location: Danforth Creek @ Martin Hwy		☐ Comp	Collection (comp begin or grab)	☒ Grab	Date 03/17/20	Time 11:10	Eastern Central	Composite end	Date	Time	Bottle Group(s)
Field	WIN ID/Field ID:		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☒ mg/L		Total Res. Chlorine		G	
WIN/S			fresh		4.77		☐ % sat		(mg/L)			
Latitude		☐ dd ☐ dms	Longitude		Temp (C) 22.9		pH 7.28		Sample Depth 0.10			Sp. Conductance (umho/cm) 640
Salinity (ppt)		☐ dd ☐ dms	0.31		Comments							
Location	Location: Danforth Creek @ Martin Hwy		☐ Comp	Collection (comp begin or grab)	☒ Grab	Date 03/17/20	Time 11:15	Eastern Central	Composite end	Date	Time	Bottle Group(s)
Field ID	WIN ID/Field ID:		Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen		☐ mg/L		Total Res. Chlorine		G	
WIN/STC							☐ % sat		(mg/L)			
Latitude		☐ dd ☐ dms	Longitude		Temp (C)		pH		Sample Depth			Sp. Conductance (umho/cm)
Salinity (ppt)		☐ dd ☐ dms			Comments Duplicate							

Relinquished By: FDEP Kristene Taylor	Date/Time 03/17/20 14:15	Shipping Method fedex	Number of Coolers Shipped: 3	Received By: IM	Date/Time 3/18/20 10:35
--	-----------------------------	--------------------------	---------------------------------	--------------------	----------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	0
	CHLSUITE-W						
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	0
	W-NH3						
	W-NO2NO3						
	W-S-A-TP						
	W-TKN						
	W-TOC						
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	0
	W-PO4-F						
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ALKALINITY						
	TURBIDITY						
	W-CL-IC						
	W-COLOR						
	W-F						
	W-SO4-IC						
	W-TDS						
	W-TSS						
Group: B	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	CHLSUITE-W						
Group: B	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
	ECOLI-18QT						
Group: B	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1

Florida Department of Environmental Protection

Central Laboratory Submittal Form

Customer: SE-DIST

Project ID: SMP

Requester: Jordan Skaggs

Field Report Prepared By: Kristene Taylor

Collected By: Alexander Miranda

Sampling Agency: 2IFLWAB

Location	Location: Bessey Creek @ Murphy		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 03/17/20 Time 11:45	☒ Eastern ☐ Central	Composite end Date / Time	Bottle Group(s) (See pg 2)	
Field ID	WIN ID/Field ID:  G2SE0017		Matrix (type, e.g., Salt, Fresh, etc) marine		Diss. Oxygen 3.08	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/S			Temp (C) 24.8	pH 7.20	Sample Depth 0.3	☐ ft ☒ m	Sp. Conductance (umho/cm) 15730	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) 9.19	Comments			F
Location	Location: Willoughby Creek		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 03/17/20 Time 12:20	☒ Eastern ☐ Central	Composite end Date / Time	Bottle Group(s)	
Field ID	WIN ID/Field ID:  G2SE0034		Matrix (type, e.g., Salt, Fresh, etc) marine		Diss. Oxygen 4.03	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/S			Temp (C) 25.9	pH 7.56	Sample Depth 0.25	☐ ft ☒ m	Sp. Conductance (umho/cm) 48793	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) 0.21	Comments			C
Location	Location: Warner Creek @24th		☐ Comp ☒ Grab	Collection (comp begin or grab) Date 03/17/20 Time 13:35	☒ Eastern ☐ Central	Composite end Date / Time	Bottle Group(s)	
Field ID	WIN ID/Field ID:  G2SE0038		Matrix (type, e.g., Salt, Fresh, etc) fresh		Diss. Oxygen 3.58	☒ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STC			Temp (C) 25.6	pH 7.60	Sample Depth 0.1	☐ ft ☒ m	Sp. Conductance (umho/cm) 410	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt) 0.20	Comments			B
Location			☐ Comp ☐ Grab	Collection (comp begin or grab) Date / Time	☐ Eastern ☐ Central	Composite end Date / Time	Bottle Group(s)	
Field ID			Matrix (type, e.g., Salt, Fresh, etc)		Diss. Oxygen	☐ mg/L ☐ % sat	Total Res. Chlorine (mg/L)	
WIN/STORET #			Temp (C)	pH	Sample Depth	☐ ft ☐ m	Sp. Conductance (umho/cm)	
Latitude	☐ dd ☐ dms	Longitude	☐ dd ☐ dms	Salinity (ppt)	Comments			

Relinquished By: FDEP Kristene Taylor	Date/Time 03/17/20 14:15	Shipping Method Fedex	Number of Coolers Shipped: 3	Received By: IM	Date/Time 3/18/20 10:35
--	-----------------------------	--------------------------	---------------------------------	--------------------	----------------------------

Group: A	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					
Group: A	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: A	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab _____
	W-NH3					
	W-NO2NO3					
	W-S-A-TP					
	W-TKN					
	W-TOC					
Group: A	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab _____
	W-PO4-F					
Group: B	Bottle Type:	P-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ALKALINITY					
	TURBIDITY					
	W-CL-IC					
	W-COLOR					
	W-F					
	W-SO4-IC					
	W-TDS					
	W-TSS					
Group: B	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	CHLSUITE-W					
Group: B	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____
	ECOLI-18QT					
Group: B	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab _____

Group: B	Bottle Type: W-E8321-DI W-E8321-MS	BG-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: W-HARD W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
Group: B	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: ALKALINITY TURBIDITY W-BR-IC W-F W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: ENT-24-QT	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: W-CARB-AA	BG-500ML	# of Bottles: 1	Preservative: CH2CICOOH Pre-Preserved	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: W-CT-CI-R W-PCL-SQ-R	BG-500ML	# of Bottles: 4	Preservative: ICE	Enter Number of Bottles Sent to Lab	4
Group: C	Bottle Type: W-E8321-DI	BG-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1

Group: C	Bottle Type: W-E8321-MS	BG-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: W-ICP W-ICPMS	P-500ML	# of Bottles: 1	Preservative: HNO3	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
Group: C	Bottle Type: W-PSNP-TQ	BG-1L	# of Bottles: 4	Preservative: ICE	Enter Number of Bottles Sent to Lab	4
Group: D	Bottle Type: ALKALINITY TURBIDITY W-BR-IC W-F W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: D	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: D	Bottle Type: ENT-24-QT	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: D	Bottle Type: W-E8321-DI W-E8321-MS	BG-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: D	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1

Group: D	Bottle Type: W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
Group: D	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
Group: E	Bottle Type: ALKALINITY TURBIDITY W-BR-IC W-F W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: E	Bottle Type: CHLSUITE-W	BP-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: E	Bottle Type: ENT-24-QT	P-250ML-WO-THI	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: E	Bottle Type: W-E8321-DI W-E8321-MS	BG-500ML	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1
Group: E	Bottle Type: W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC	P-500ML	# of Bottles: 1	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
Group: E	Bottle Type: W-PO4-F	P-125ML	# of Bottles: 1	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	1
Group: F	Bottle Type: ALKALINITY TURBIDITY W-BR-IC W-F W-TSS	P-1L	# of Bottles: 1	Preservative: ICE	Enter Number of Bottles Sent to Lab	1

Group: F	Bottle Type:	BP-1L	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
CHLSUITE-W							
Group: F	Bottle Type:	P-250ML-WO-THI	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
ENT-24-QT							
Group: F	Bottle Type:	BG-500ML	# of Bottles: 1	Preservative:	ICE	Enter Number of Bottles Sent to Lab	1
W-E8321-DI W-E8321-MS							
Group: F	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	HNO3	Enter Number of Bottles Sent to Lab	1
W-ICP W-ICPMS							
Group: F	Bottle Type:	P-500ML	# of Bottles: 1	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	1
W-NH3 W-NO2NO3 W-S-A-TP W-TKN W-TOC							
Group: F	Bottle Type:	P-125ML	# of Bottles: 1	Preservative:	FILTER-ICE	Enter Number of Bottles Sent to Lab	1
W-PO4-F							
Group: G	Bottle Type:	P-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
ALKALINITY TURBIDITY W-CL-IC W-COLOR W-F W-SO4-IC W-TDS W-TSS							
Group: G	Bottle Type:	BP-1L	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
CHLSUITE-W							
Group: G	Bottle Type:	P-250ML-WO-THI	# of Bottles: 2	Preservative:	ICE	Enter Number of Bottles Sent to Lab	2
ECOLI-18QT							
Group: G	Bottle Type:	P-500ML	# of Bottles: 2	Preservative:	ICE And H2SO4	Enter Number of Bottles Sent to Lab	2

Group: G	Bottle Type:	P-500ML	# of Bottles: 2	Preservative: ICE And H2SO4	Enter Number of Bottles Sent to Lab	<u>2</u>
	W-NH3					
	W-NO2NO3					
	W-S-A-TP					
	W-TKN					
	W-TOC					

Group: G	Bottle Type:	P-125ML	# of Bottles: 2	Preservative: FILTER-ICE	Enter Number of Bottles Sent to Lab	<u>2</u>
	W-PO4-F					

Date of Request: 19-FEB-2020
 Created By: SKAGGS_J On: 19-FEB-2020
 Modified By: SWANSON_C On: 27-FEB-2020
 Customer: SE-DIST
 Project: SMP
 Division: Division of Environmental Assessment and Restoration
 District: Southeast District
 Sampling Event: SMP- Martin 2/5

Send Coolers To:

Phone: (772) 448-5884
 FL Dept. of Environmental Protection
 2199 South Rock Road
 Fort Pierce, FL 34945
 Attn: Jordan Skaggs

Program Module Number:

Priority: 3
 Event Status: S
 Criminal Investigation: NO
 Chemistry Request Reviewed By: WATTS_J On: 26-FEB-2020
 Biology Request Reviewed By: SWANSON_C On: 27-FEB-2020
 Sampling Kit Required: YES Ship on: 06-MAR-2020
 Sampling Kit Shipped: YES On: 02-MAR-2020
 Sampling Kit Packed By: SHAIK_A
 Preservatives Needed: YES
 Date to Receive Samples: 18-MAR-2020
 Received By:

Send Final Report To:

FL Dept. of Environmental Protection
 2199 South Rock Road
 Fort Pierce, FL 34945
 Attn: Drew Liddick

Report Type: Final Only
 FTP Data: NO
 QC Report: YES
 Date Log: NO
 Authorization Log: NO

Comments: 1 - case of nitric acid/1 set of labels
 1 - case of sulfuric acid/1 set of labels
 5 - return shipping labels/handle tags

Sampler Comments:
 3194G-Group A-G2SE0014
 3194W-Group B-G2SE0038
 3208B-Group C-G2SE0034
 3210C-Group D-G2SE0043
 3210E-Group E-G2SE0029
 3211-Group F-G2SE0017
 3215-Group G-G2SE0018+G2SE0018_DUP

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices

Bottle Group A (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group B (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-HARD	Template: DEFAULT	(SM 2340B) Hardness by calculation, in aqueous matrices, reported as mg/L CaCO3
Hardness		
W-ICP	Template:	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Barium		
Calcium		
Iron		
Magnesium		
Potassium		
Sodium		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
*Boron		
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Group C (Water) With 1 Samples:

Bottle Group C (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: CH ₂ ClCOOH Pre-Preserved
W-CARB-AA	Template: DEFAULT	(EPA 8321B) Carbamates in water matrices by HPLC/MS/MS
Bottle Type: BG-500ML	Number of Bottles: 4	Preserved With: ICE
W-CT-CI-R	Template: DEFAULT	(EPA 8276) Organochlorine pesticides in water matrices by GC/MS with Negative Chemical Ionization
W-PCL-SQ-R	Template: DEFAULT	(EPA 8270E) Organochlorine pesticides in water matrices by GC/MS-SIM
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO ₃
W-ICP	Template: S-METALS	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron		
Manganese		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H ₂ SO ₄
W-NH ₃	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO ₂ NO ₃	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO ₄ -F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L
Bottle Type: BG-1L	Number of Bottles: 4	Preserved With: ICE
W-PSNP-TQ	Template: MASTERLIST	(EPA 8270E) Organonitrogen and phosphorus pesticides in water matrices by GC/MS/MS

Bottle Group D (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO ₃ /L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE

Bottle Group D (Water) With 1 Samples:

Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group E (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group F (Water) With 1 Samples:

Bottle Type: P-1L	Number of Bottles: 1	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-BR-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Bromide in aqueous matrices
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 1	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 1	Preserved With: ICE
ENT-24-QT	Template: DEFAULT	(Enterolert/QT) Enterococci by Enterolert Quanti-Tray method - non-chlorinated water systems
Bottle Type: BG-500ML	Number of Bottles: 1	Preserved With: ICE
W-E8321-DI	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
W-E8321-MS	Template: DEFAULT	(EPA 8321B) Pesticides, herbicides, and other chemicals in water matrices by HPLC/MS/MS
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3

Bottle Group F (Water) With 1 Samples:

Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: HNO3
W-ICP	Template: S-METALS	(EPA 200.7 Rev. 4.4) Total Recoverable Metals analysis using ICP emission spectroscopy for aqueous samples supporting Clean Water Act Projects
Iron		
Manganese		
Zinc		
W-ICPMS	Template:	(EPA 200.8 Rev. 5.4) Total Recoverable Metals analysis using ICP-MS for aqueous samples supporting Clean Water Act Projects
Aluminum		
Antimony		
Arsenic		
Cadmium		
Chromium		
Copper		
Lead		
Nickel		
Selenium		
Silver		
Thallium		
Bottle Type: P-500ML	Number of Bottles: 1	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 1	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

Bottle Group G (Water) With 2 Samples:

Bottle Type: P-1L	Number of Bottles: 2	Preserved With: ICE
ALKALINITY	Template: DEFAULT	(SM 2320 B-2011) Alkalinity in aqueous matrices as mg CaCO3/L
TURBIDITY	Template: DEFAULT	(EPA 180.1 Rev. 2.0) Turbidity in aqueous matrices
W-CL-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Chloride in aqueous matrices
W-COLOR	Template: DEFAULT	(SM 2120 B-2011) True Color measured at 450 nm (not valid for NPDES monitoring)
Color (true)		
W-F	Template: DEFAULT	(SM 4500 F-C-2011) Fluoride in aqueous matrices, without distillation (not valid for NPDES monitoring)
W-SO4-IC	Template: DEFAULT	(EPA 300.0 Rev. 2.1) Sulfate in aqueous matrices
W-TDS	Template: DEFAULT	(SM 2540 C-2011) Total Dissolved Solids in aqueous matrices
W-TSS	Template: DEFAULT	(SM 2540 D-2011) Total Suspended Solids in aqueous matrices using a maximum of 250 mL of sample
Bottle Type: BP-1L	Number of Bottles: 2	Preserved With: ICE
CHLSUITE-W	Template: DEFAULT	(SM 10200 H (mod.)) Phytoplankton chlorophyll-a corrected, uncorrected, and phaeophytin by spectrophotometry
Bottle Type: P-250ML-WO-THI	Number of Bottles: 2	Preserved With: ICE
ECOLI-18QT	Template: DEFAULT	(SM 9223 Quanti-Tray) Escherichia coli by Colilert 18 Quanti-Tray method - non-chlorinated water systems
Bottle Type: P-500ML	Number of Bottles: 2	Preserved With: ICE And H2SO4
W-NH3	Template: DEFAULT	(EPA 350.1 Rev. 2.0) Ammonia in aqueous matrices as mg N/L
W-NO2NO3	Template: DEFAULT	(EPA 353.2 Rev. 2.0) Nitrite/Nitrate in aqueous matrices as mg N/L
W-S-A-TP	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Total Phosphorus in aqueous matrices as mg P/L
W-TKN	Template: DEFAULT	(EPA 351.2 Rev. 2.0) Total Kjeldahl Nitrogen in aqueous matrices
W-TOC	Template: DEFAULT	(SM 5310 B-2011) Nonpurgeable Organic Carbon in aqueous matrices
Bottle Type: P-125ML	Number of Bottles: 2	Preserved With: FILTER-ICE
W-PO4-F	Template: DEFAULT	(EPA 365.1 Rev. 2.0) Ortho-phosphate, dissolved, in filtered, aqueous matrices as mg P/L

* - The laboratory is not NELAP certified for this analyte/method, or certification is not applicable.