Preparation of Preliminary Algal Assay Procedure (PAAP) Media for Freshwater Algae Cultures

1. SCOPE AND APPLICATION

    1.1. This method is based on Miller, W. E., Greene, J. C., and Shiroyama, T. 1978. "The *Selenastrum capricornutum* Printz Algal Assay Bottle Test". EPA600/9-78-018. U. S. EPA, Corvallis, Oregon. 126 p. It is used to prepare culture media for *Raphidocelis subcapitata* (formerly *Selenastrum capricornutum*) from stock solutions.

2. SUMMARY OF METHOD

    2.1. Aliquots of pre-made nutrient stock solutions are added to autoclaved deionized water. The solution is pH-adjusted as needed, filtered, and stored in the walk-in refrigerator. Batch information is recorded in the Laboratory Information Management System (LIMS).

    2.2. Interferences: Contaminated items including glassware, filters, deionized water and nutrient solutions can cause unwanted growth in the algal cultures.

3. EQUIPMENT AND SUPPLIES

    3.1. 2 L Volumetric flask, clean and autoclaved

    3.2. Fisherbrand Accumet AB315 pH/mV Meter, calibrated per SOP TA-06.43

    3.3. 30 mL disposable plastic vial

    3.4. 2 mL Disposable pipettes, 8

    3.5. 2 L Filter assembly: filter flask, filter holder, and filter funnel. Refer to SOP TA-08.01.

    3.6. 0.45 µm Membrane filters, washed per TA-08.01

    3.7. Parafilm®

    3.8. Computer with access to the Laboratory Information Management System (LIMS) per SOP TA-06.42

4. REAGENTS AND STANDARDS

    4.1. Deionized Water, autoclaved 30 min at 121 ºC and 15 psi

    4.2. PAAP stock solutions (prepared in advance per SOP TA-09.09):

        4.2.1. Sodium nitrate

        4.2.2. Magnesium chloride

        4.2.3. Calcium chloride

        4.2.4. Trace salts

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

4.2.5. Magnesium sulfate

4.2.6. Potassium phosphate

4.2.7. Sodium bicarbonate

4.2.8. EDTA

5. PREPARATION

5.1. Get PAAP stock solutions from the walk-in refrigerator.

5.2. Rinse the volumetric flask with autoclaved DI water. Fill to approximately half volume with autoclaved DI water.

5.3. Add 2 mL of each PAAP stock solution to the 2 L volumetric flask.

5.3.1. Use a separate pipette for each stock solution to avoid cross contamination of the stock solutions.

5.3.2. Aliquots of each stock solution can be poured into 30 mL plastic cups for ease of transfer or pipetted directly from the flask with a sterile glass pipette. Do not return unused solution to the flask when done.

5.3.3. Add 1 mL of each stock solution per liter of media that is being prepared.

5.4. Cover the flask with Parafilm®. Swirl and shake the flask. Remove the Parafilm® and finish filling to the 2 L mark with autoclaved DI water. Replace the Parafilm® and invert the flask several times.


5.5. Pour approximately 30 mL of the mixed medium into the disposable plastic vial and measure the pH.

5.5.1. Adjust final media pH to 7.5 ± 0.2 with NaOH or HCl as required. Refer to SOP TA-08.06.

5.6. Label with PAAP Media, preparation date/batch number, analyst initials, final pH measurement, LIMS ID, and expiration date 3 months from preparation date.

5.7. Filter the pH-adjusted media through a washed 0.45 µm membrane filter. Refer to SOP TA-08.01. See figure below.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

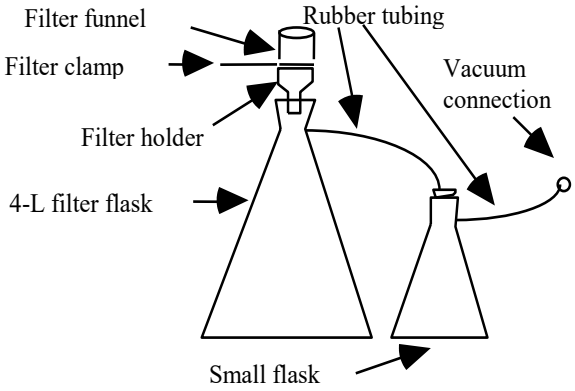

5.8. Record batch information in the LIMS Standard Prep Tracker.

5.9. Store flask in walk-in refrigerator.

6. QUALITY CONTROL

6.1. After three months, discard unused media and prepare new PAAP media.

7. SAFETY/HAZARDOUS WASTE MANAGEMENT

7.1. Proper laboratory protection, including lab hood or ventilation system, lab coat, gloves and safety glasses, is required when working with chemicals and unprocessed samples.

7.2. Refer to the FDEP Laboratory Quality Assurance Manual.

7.3. Refer to the FDEP Health and Safety Manual.

7.4. Read MSDS information for reagent handling and disposal.

8. REFERENCES

8.1. FDEP Health and Safety Manual

8.2. FDEP Quality Assurance Manual

8.3. Miller, W. E., Greene, J. C., and Shiroyama, T. 1978. "The *Selenastrum capricornutum* Printz Algal Assay Bottle Test". EPA600/9-78-018. U. S. EPA, Corvallis, Oregon. 126 p.

8.4. FDEP SOP TA-06.43, Calibration and Use of Fisherbrand™ accumet™ AB315 pH/mV Meter

8.5. FDEP SOP TA-08.01, Initial Algal Assay Sample handling

8.6. FDEP SOP TA-08.06, pH adjustment of Samples: Algal Assay and Toxicity bioassay

8.7. FDEP SOP TA-09.09, Preparation of PAAP-Media Stock Solutions.

8.8. FDEP SOP TA-06.42, LIMS Entry for Laboratory Parameters

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

8.9. MSDS information sheets.

APPENDIX OF SIGNIFICANT CHANGES

09/25/2008 Added Scope and Application, Summary of the Method, Interferences, References and Appendix of Changes.

02/09/2015 Formatted for 508 compliance.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.