PURIFICATION OF FILTERED SAMPLES USING THE QIAGEN DNEASY POWERSOIL® PRO KIT

1. SCOPE AND APPLICATION

1.1. This method is based on the Qiagen DNeasy PowerSoil® Pro kit user manual and describes a procedure for recovering genomic DNA (gDNA) from samples that contain PCR inhibitors. This method is used to purify samples treated with propidium monoazide (PMA) and to re-analyze samples when inhibition is detected after performing the crude lysate method (as described in SOP PCR-4.1).

2. SUMMARY OF METHOD

2.1. Environmental samples (water filters) are added to bead tubes and homogenized on a bead beater for rapid and thorough homogenization. Cell lysis occurs by mechanical shearing and chemical interactions. Total gDNA is bound on a silica membrane in a spin column. DNA is washed and eluted from the spin column to yield purified DNA.

2.2. Interferences: Purification removes PCR inhibitors but may result in a low yield of DNA. Ethanol in Solution C5 can interfere with qPCR. It is important to ensure that the Spin Filter is thoroughly centrifuged to remove any traces of Solution C5 and that the Spin Filter is removed carefully to avoid splashing any Solution C5 onto the Spin Filter.

2.3. Definitions:

2.3.1. qPCR: Quantitative Polymerase Chain Reaction is a procedure in which the amount of DNA amplified in each PCR cycle is measured using fluorescent dyes to determine the original amount of DNA in the sample.

2.3.2. Target sequence: The DNA sequence in the Gene of Interest (GOI) that is amplified by the primers/probe in each PCR cycle. The primers/probe series is designed to attach specifically to this sequence.

3. EQUIPMENT AND SUPPLIES

3.1. Equipment and Supplies –

Equipment/Supply Name	Notes
Gloves, Nitriles	N/A
Tube Racks	N/A
Forceps	N/A
Micropipettes	N/A
Sterile PCR Grade Filter Tips	N/A
CBS Optimizer PCR Workstation #2	S/N 007632
UV Sterilizer Unit	N/A
Fisher Mini Vortex	N/A
Centrifuge	N/A
Fisher Scientific Bead Mill 24	S/N 24G2336
Qiagen DNeasy PowerSoil® Pro Kit	Catalog Number: 47014

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

PowerBead Pro Tubes	N/A
MB Spin Columns	N/A
2 mL Collection Tubes	N/A
1.5 mL Elution Tubes	N/A
2 mL Microcentrifuge Tubes	N/A
PCR#5, Thermo Fridge/Freezer	S/N 0121174701121218. Maintained at 4 ± 2°C and -20 ± 10°C, Room B246
PowerVac® Manifold	Catalog Number 11991. Optional
MilliporeSigma™ Isopore™ Membrane Filters, 0.4 µm, 47 mm	Catalog Number HTTP04700
MilliporeSigma™ Glass Fiber Filters with Binder, 2 µm (100 filters)	Catalog Number AP2004700.

4. REAGENTS AND STANDARDS

4.1. Reagents and Standards –

Reagent/Standard Name	Notes
Qiagen DNeasy PowerLyzer® Pro Kit	Catalog Number 47014
	Solution CD1 – Store at room temperature (15–30°C)
	Solution CD2 – Store at 2-8°C upon arrival.
	Solution CD3 – Store at room temperature (15–30°C)
	Solution EA – Store at room temperature (15–30°C)
	Solution C5 – Store at room temperature (15–30°C)
	Solution C6 – Store at room temperature (15–30°C)
1E4 SKETA Standard (Salmon DNA)	Prepare per SOP PCR-3.1
E6 HFG2GD spike	Prepare per SOP PCR-3.4
E7 BRD3 spike	Prepare per SOP PCR-3.5
E7 C2H2 spike	Prepare per SOP PCR-3.7
E7 HoF spike	Prepare per SOP PCR-3.9
Sterile Deionized (DI) Water	N/A
Bleach solution, 10% (v/v)	N/A
80% Ethanol	Prepare per SOP IZ-09

5. CALIBRATION AND STANDARDIZATION

5.1. See DEP Laboratory Quality Manual, Section 5.5.

6. PREPARATION

6.1 Print two labels per sample and control ID using the Brady Workstation, and label Power Bead and final elution tubes. Use a permanent marker to write sample and control IDs on all other tubes. Refer to SOP PCR-5.0 for printing sample and control labels.


6.1 Clean work area before and after use. Wipe down all surfaces with 10% bleach and wait 15 minutes before rinsing with sterile water.

6.2 Samples and standards should not be removed from storage until extraction process begins. All empty extraction tubes can be labeled and set up the day before extraction. Seal the tubes in a rack with a sterile plastic cover.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

6.3 Count number of samples and controls to obtain number of extraction tubes needed from the Power Soil Pro kitIf tubes are set up at the same time as extraction, they can be placed directly into the Argosflexirack per the diagram below.

6.4 Table diagram of PowerSoil Pro kit extraction tubes

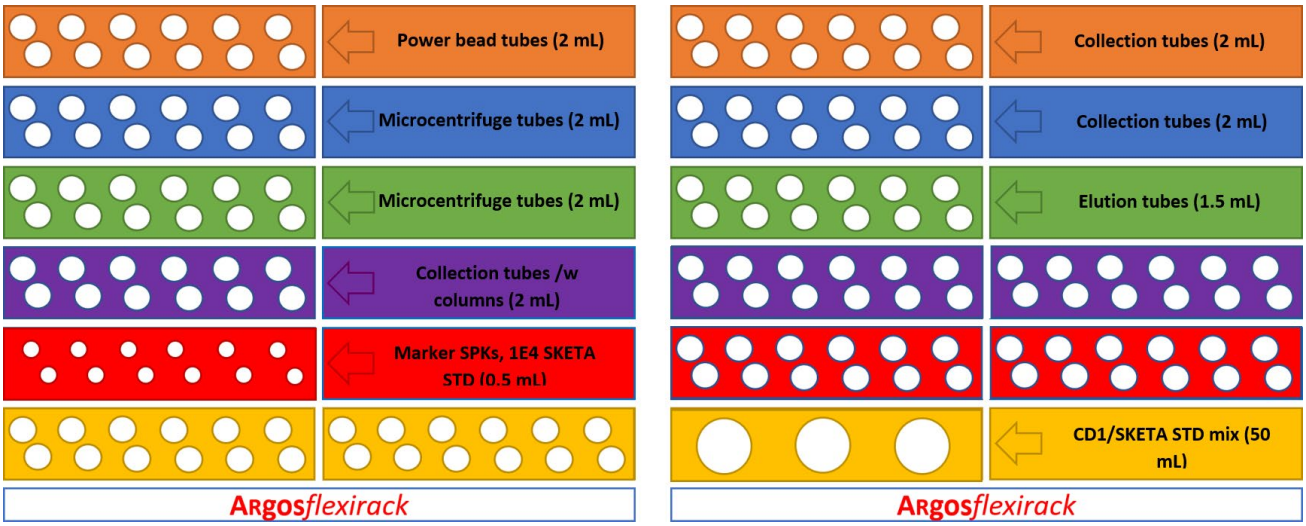

7. SAMPLE PURIFICATION

7.1. Place forceps in PCR Workstation #2 or UV sterilizer and turn on UV lights. Sterilize forceps for at least 2 minutes.

7.2. Retrieve samples. method blank, 1E4 SKETA, and extraction controls from -80 freezer. Use sterile forceps to transfer a filter to a PowerBead Pro Tube labeled for each sample or for controls (see Section 8).

NOTE: Keep forceps in a small beaker of ethanol during the purification process. Dip in the ethanol and shake off excess before grasping each filter. If forceps are accidentally dropped or placed on counter, re-sterilize in PCR Workstation #2. Make sure no ethanol residue from the forceps touches the filter.

7.3. Count the number of tubes. In a 50 mL conical tube, add 800 µL of Solution CD1 X total number of tubes plus one, and 2 µL of 1E4 SKETA standard X total number of tubes plus one. Invert several times to mix well.

Example: CD1 volume = (10 (tubes) samples and controls) +1 x 800 µL = 8,800 µL

1E4 SKETA standard volume = (10 (tubes) samples and controls) +1 x 2 µL of 1E4 SKETA standard = 22 µL

7.4 Fold a sterile filter and place into each control tube labeled EB and EC.

7.5. With a 1000 µL micropipette, dispense 802 µL of Solution CD1 to each PowerBead Pro Tube.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

7.6. Make sure tubes are sealed tightly. Vortex briefly on speed 7 to mix.

7.7.

Place the tubes in the bead beater. Follow arrows on bead beater for correct placement to balance. Select setting for 45 seconds at a speed setting of 4.

7.8. Centrifuge at 15000 x g for 1 min to pellet the debris.

7.9. With a 1000 µL micropipette, carefully transfer the supernatant to a new 2 mL microcentrifuge tube for each sample. Expect 500–600 µL. The supernatant may still contain some particles.

7.10. With a 200 µL micropipette, add 200 µL of Solution CD2 (stored in fridge) into each tube and vortex at speed 7 for 5 seconds.

7.11. Centrifuge at 15000 x g for 1 min to pellet the debris.

7.12. Carefully, without disturbing the pellet, transfer up to 700 µL of the supernatant to a new 2 mL microcentrifuge tube for each sample. Expect 500-600 µL.

7.13. With a 1000 µL micropipette, add 600 µL of Solution CD3 into each tube and vortex at speed 7 for 5 seconds.

NOTE: *If Solution CD3 has precipitated, heat at 60°C until precipitate dissolves.*

7.14. With a 1000 µL micropipette, transfer 650 µL of the lysate to an MB Spin Column for each sample and centrifuge at 15000 x g for 1 min.

7.15. Discard the flow through.

7.16. Repeat steps 6.13 and 6.14 until all the solution is centrifuged through the MB Spin Column.

7.17. Carefully transfer the MB Spin Column into a clean 2 mL Collection Tube. Avoid splashing any flow-through onto the MB Spin Column.

7.18. With a 1000 µL micropipette, add 500 µL of Solution EA into each MB Spin Column.

7.19. Centrifuge at 15000 x g for 1 min.

7.20. Discard the flow-through and place the MB Spin Column back into the same 2 mL Collection Tube.

7.21. With a 1000 µL micropipette, add 500 µL of Solution C5 into each MB Spin Column.

7.22. Centrifuge at 15000 x g for 1 min.

7.23. Discard the flow-through and place the MB Spin Column into a new 2 mL Collection Tube.

7.24. Centrifuge at up to 16000 x g for 2 min.

7.25. Carefully transfer the MB Spin Column into a new 1.5 mL Elution Tube. Avoid splashing any of Solution C5 onto the Spin Filter.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

NOTE: *Solution C5 is an ethanol-based solution. Any trace of this solution can interfere with qPCR.*

7.26. With a 200 µL micropipette, add 100 µL of Solution C6 onto the center of each white filter membrane.

7.27. Centrifuge at 15000 x g for 1 min.

7.28. Discard the MB Spin Column. 7.29. Store purified DNA in the freezer. Use -20°C for short term storage and -80° C for long term storage.

8. QUALITY CONTROL

8.1. Negative Controls

8.1.1. Analyze a Method Blank (MB) control for each set of samples filtered (maximum of 13 per batch). Method blanks are filtered per SOP PCR-4.0 and are stored with sample batches. Dip sterile forceps into 80% ethanol and shake to remove excess. Place method blank into a PowerBead Pro tube and extract along with samples in the batch.

8.1.2. Analyze an Extraction Blank (EB) control for each set of samples extracted (maximum of 13 per batch). Place the filter into a PowerBead Pro tube and extract along with samples in the batch.

8.2. Extraction Control (EC)

8.2.1. Include an EC for each assay to be analyzed for the batch (maximum of 13 samples per batch). A vacuum manifold and sterile funnel base can be used to fold a sterile filter as described in section 7.1.2. Place the sterile filter into a PowerBead Pro tube and spike with 5 µL of E^6 HFG2GD spike, E^7 BRD3 spike, or E^7 C2H2 spike. Make separate ECs for each assay needed. Extract along with samples in the batch.

8.3. Sample Processing Control (SPC)

8.3.1. Add 2 µL of 1E4 SKETA standard to each PowerBead Pro Tube after adding the sample/control filter in Section 6.2. The SKETA standard is added along with Solution CD1.

8.4. Refer to SOP PCR-5.1 for further details on uploading quality control data into LIMS, including calculations and control limits.

9. SAFETY/HAZARDOUS WASTE MANAGEMENT

9.1. Wear gloves and avoid all skin contact with kit reagents.

9.2. All media/bacteria/plates/tubes must be autoclaved for 45 minutes before disposal.

9.3. Refer to the Health and Safety Plan. Read MSDS information for reagent handling and disposal.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

10. REFERENCES

10.1. Qiagen DNeasy PowerSoil® Pro Kit Instruction Manual (Catalog#47014)

10.2. FDEP Laboratory Quality Assurance Manual

10.3. Health and Safety Plan

10.4. DEP SOP IZ-09, *Preparation of 80% Ethanol*

10.5. DEP SOP PCR-3.1, *Preparation of SKETA assay standards and reagents for qPCR analysis*

10.6. DEP SOP PCR-3.2, *Preparation of HF183 Assay Reagents and Standards for qPCR analysis*

10.7. DEP SOP PCR-3.3, *Preparation of Gull2 Assay Reagents and Standards for qPCR analysis*

10.8. DEP SOP PCR-3.4, *Preparation of GFD Assay Reagents and Standards for qPCR analysis*

10.9. DEP SOP PCR-3.7, *Preparation of CowM2 Assay Reagents and Standards for qPCR analysis*

10.10 DEP SOP PCR-3.9, *Preparation of Hof597 Assay Reagents and Standards for qPCR analysis*

10.11. DEP SOP PCR-4.1, *Crude Lysate Extraction Method for qPCR Analysis*

10.12. DEP SOP PCR-5.1, *LIMS QC Manager Data Upload, Verification, and Sample Authorization Instructions for qPCR samples*

Appendix of Changes

05/11/2016 Added optional method of purification that includes the use of vacuum manifold.

06/05/2017 Added extraction controls for additional markers. Replaced bead beater and added instructions specific to new bead beater.

06/05/2018 Updated the equipment, supplies and reagents list. Update the Health and Safety plan link.

05/28/2019 Added extraction controls for DG3, CowM2 and HumM2 markers.

02/11/2020 Updated the procedure with new specifications for Qiagen DNeasy PowerSoil® Pro Kit.

01/27/2022 Updated spikes to HFG2GD and BRD3. Removed BRD3GD and Gull2 spikes.

11/16/2023 Reformatted. Added table in Section 4.

02/08/2024 Added section 6.0 Preparation including diagram of tube setup.

01/16/2025 Changed tables font.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

04/10/2026	Updated preparation information to allow set up of all tubes the day prior to extraction. Specified number of labels to print. Removed preparation of extraction controls via filtration manifold. Added reference for SOP PCR-3.9

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.