Preparation of Artificial Seawater (ASW) Medium for Culturing *Dunaliella tertiolecta*

1. SCOPE AND APPLICATION

1.1. This method is based on Maloney, T.E., W.E. Miller, and D.T. Specht. 1974. The Marine Algal Assay Procedure: Bottle Test. National Environmental Research Center Office of Research and Development. U.S. EPA, Corvallis, Oregon. 43 p. It is used to prepare medium for culturing of Dunaliella teriolecta.

2. SUMMARY OF THE METHOD

2.1. Nutrients are added to Basal medium. Solution is filtered and then stored in the refrigerator until use. A record of preparation is maintained in the *AGP Record of Preparation* logbook.

3. SAFETY/HAZARDOUS WASTE MANAGEMENT

3.1. Refer to the Bureau of Laboratories Health and Safety Plan at https://www.fldepnet.org/content/dear/bureau-laboratories-safety-documents

3.2. Read MSDS information for reagent handling and disposal.

3.3. Use proper safety equipment including the lab hood or chemical ventilation system, a lab coat, safety glasses, a respirator and latex gloves (or a suitable substitute if allergic to latex) when working with chemicals and unprocessed samples.

4. SUPPLIES AND EQUIPMENT

4.1. Supplies

4.1.1. Prewashed 0.45 µm membrane filters

4.1.2. Autoclaved deionized water

4.1.3. Volumetric pipets, 10 ml, 4 each

4.1.4. *AGP Record of Preparation* logbook

Note: To correct raw data entries, place a single line through the incorrect or illegible entry (do not overwrite it), write the corrected entry near the error, initial and date the correction. If the correction requires an explanation, write a comment and reference the error.

4.2. Equipment

4.2.1. Erlenmeyer flask, 4 L

4.2.2. Filter flask, 4 L

4.2.3. Filtering apparatus (filter funnel, clamp, filter holder)

4.2.4. Vacuum source with attached tubing

4.2.5. Graduated cylinder, 1 L

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

4.2.6. Filter forceps

5. REAGENTS AND STANDARDS

5.1. Basal medium, 30 ppth salinity, prepared per SOP TA-09.10

5.2. PAAP nutrients, prepared in advance per SOP TA-09.09

5.2.1. Sodium nitrate

5.2.2. Trace


5.2.3. Potassium phosphate

5.2.4. EDTA

6. PREPARATION (FOUR LITERS)

6.1. Use a new 10-mL pipet for each PAAP nutrient to add 4 mL of each nutrient to four liters of 30 ppth basal medium.

6.2. Filter the solution (from 6.1) through a pre-washed 0.45 µm membrane filter. Place the filtered solution in a 4 L flask.

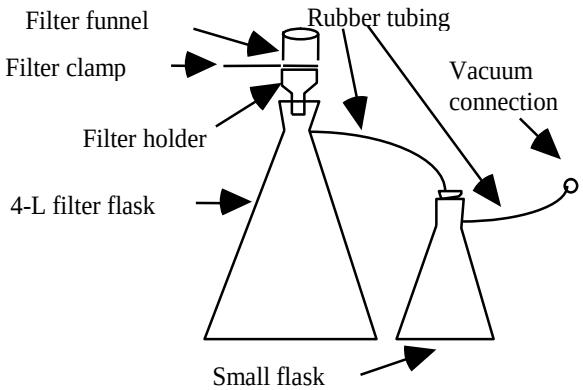

Note: Usually, the overflow flask and rubber tubing connections are already set up and it is only necessary to attach the filter flask, filter holder, and filter funnel. If not, set up the filter apparatus as indicated above. Flask size can vary.

6.3. Seal the flask (from 6.2) with para-film. Label the flask with *ASW medium, 30 ppth, date filtered, w/NEPT, batch #,* and *preparer's initials.* Add an expiration date of one (1) year from the date of preparation. Store the flask in the Biology Section's walk-in refrigerator.

7. REFERENCES

7.1. Maloney, T. E., W. E. Miller, and D. T. Specht. 1974. The Marine Algal Assay Procedure: Bottle Test. National Environmental Research Center Office of Research and Development. U. S. EPA, Corvallis, Oregon. 43 p

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

7.2. DEP SOP TA-06.16. Washing Filters for Algal Assay to Remove Trace Contaminants

7.3. DEP SOP TA-09.09. Preparation of PAAP-Medium Stock Solutions

7.4. DEP SOP TA-09.10. Artificial Seawater Basal Stock Solution for Culturing Marine Algae

7.5. Bureau Quality Manual, FDEP Bureau of Laboratories, January 2014

APPENDIX OF CHANGES:

07/25/2008	Added Scope and Application, Summary of the Method, References and Appendix of Changes
07/09/2009	Added "the" before "refrigerator" in Summary of the Method. Removed redundant Scope and Summary from Materials comment. Updated Quality manual reference
05/17/2010	Added definition of 'ppth' to Materials Item 1 comments. Removed () from 'nutrients in Methods Step 1. Added 'CAUTION Individuals allergic to latex should not use latex gloves. A suitable substitute should be found and used.' and deleted 'Vinyl gloves can be substituted for individuals allergic to latex.' to Safety Section comments. Formatted tabs under References. Updated reference to Biology Quality Manual. Capitalized 'Method' in Appendix of Changes 7/09/09.
04/26/2011	Updated reference to Biology Quality Manual
03/29/2012	**Updated reference to Bureau Quality Manual. Added "Do not overwrite entries" to Step 1 of Documentation and Corrective Action section**.
03/01/2013	Reformatted the SOP.
02/09/2015	Formatted for 508 compliance.
11/04/2015	Updated links.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.