DETERMINATION OF TRUE COLOR IN WATER USING THE AQ2 DISCRETE ANALYZER

TABLE OF CONTENTS

Page

1. SCOPE AND APPLICATION 1
2. SUMMARY OF THE METHOD 1
3. EQUIPMENT AND SUPPLIES 2
4. REAGENTS AND STANDARDS 2
5. SAMPLE PRESERVATION, PREPARATION AND HANDLING 4
6. SAMPLE ANALYSIS AND QUANTITATION 4
7. DATA ARCHIVAL 8
8. QUALITY CONTROL 8
9. SAFETY/HAZARDOUS WASTE MANAGEMENT 9
10. REFERENCES/DEFINITIONS/VALIDATION DATA 9

Appendix A: Appendix of Significant Changes 10
Appendix B: Example Tray Report/Bench Sheet 11
Appendix C: Example Color pH Sheet 12

1. SCOPE AND APPLICATION

1.1. The test covered by this SOP is W-COLOR. This method is applicable to waters and wastes that have color characteristics derived from naturally occurring materials.

1.2. The calibration range is from the method detection limit (MDL) of 2.5 to 150 Platinum Cobalt Units (PCU). Sample values above 150 PCU may be determined by quantitative dilution.

2. SUMMARY OF THE METHOD

2.1. True sample color is determined using Standard Methods 2120 B which allows for the use of a spectrophotometer.

2.2. Color is determined spectrophotometrically at a wavelength of 450 nm and calibrated with platinum cobalt standards. One unit of color is that produced by 1 mg/L platinum in the form of the chloroplatinate ion.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

2.3. Interferences: Turbidity and the presence of colloidal and suspended particles will interfere with this analysis. All samples analyzed are lab filtered before true color is measured. The results are pH dependent. Highly colored industrial wastes may have hues far different from those of the platinum-cobalt standards. Their color characteristics may not be derived from naturally occurring materials such as humic acids.

3. EQUIPMENT AND SUPPLIES

3.1. Spectrophotometer - SEAL AQ2 Discrete Analyzer (AQ2 3 (primary) and AQ2 2)

3.2. Sample cups - 1.2 mL for 100 position tray (SEAL part no. 5385)

3.3. Assorted pipettes

3.4. 5 mL Syringes

3.5. 0.45 µm and 0.20 µm pore diameter membrane syringe filters

3.6. Plastic pH strips:

3.6.1. 0.0 - 6.0 range (Fisher catalog no. 14-850-114)

3.6.2. 5 .0 - 9.0 range (Fisher catalog no. 14-850-113)

3.7. Reaction Segments (SEAL part no. 5000)

3.8. Reagent Containers (SEAL part no. 5010)

3.9. Class A volumetric flasks

3.10. 20 mL Disposable beakers

3.11. Fisher accumet XL150 pH meter

4. REAGENTS AND STANDARDS (SEE DEP SOP NU-042)

Note: Before preparing standards or running samples, please read the safety notes in the safety/hazardous waste section of this SOP and the MSDS for the chemicals used in this analysis.

Note: Standard preparations must be documented in the Laboratory Information Management System (LIMS) 'Standard Preparation Tracker' module. Unique ID numbers for each standard are generated in the module and documented in the Reagent and Standard Logsheet (NU-RGTSTD) located on the network at \\FLDEP1\DEAR\Labs\lab3\Chem_Data\Inorganic\Nutrients\REAGENT_PREP\Color\. This excel sheet needs to be updated and submitted as part of the analysis report folder. If available, a Certificate of Analysis for each purchased standard or reagent is also maintained in the REAGENT_PREP\Color folder. The original Certificate of Analysis for the standard reference material is kept in a binder in C307. All containers holding standards or reagents must be labeled with their unique ID and expiration date.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

Note: The volumes prepared for the reagents and standards may be altered by the analyst as necessary. The final concentrations of the constituents of a standard or reagent must not be changed.

4.1. Stock Standard Chloroplatinate Solution (500 American Dye Manufacturers Institute (ADMI) color units e.g., Hach catalog no. 141453). Check the bottle label for the manufacturer's expiration date. Store in the dark when not in use.

4.1.1. Calibration (Top) standard (150 PCU) - To a 250 mL volumetric flask containing approx. 100 mL of deionized (DI) water, add 75 mL of 500 PCU stock standard (4.1) and bring to volume. Store in the dark when not in use. Hold time: 30 days.

4.1.2. The remaining standards are prepared from the 150 PCU calibration standard (4.1.1.). Using a class A glass volumetric pipette, place the appropriate volume of top standard in a volumetric flask partially filled with DI water and dilute to the mark with DI water. Use the below table as a guide. Store in the dark when not in use. Hold time: 30 days.

Standard Concentration (PCU)	Final Volume (mL)	Volume of TOP Standard (4.1.1.) (mL)
120	50	40
90	50	30
60 **(CCV)**	200	80
45	50	15
30	50	10
15	50	5
10.5	50	3.5
6 **(PQL)**	100	4
0	0	0

4.2. A certified quality control sample, ALTQC (e.g. ERA catalog no. 070). Check the bottle label for the manufacturer's expiration date. Store in the dark when not in use.

4.3. Hydrochloric Acid, HCl, analytical grade (for pH adjustments)

4.4. Sodium Hydroxide, NaOH, analytical grade (for pH adjustments)

4.5. Cuvette cleaning solution (Seal part #5676)

4.6. Bleach Cleaning Solution - To a 500 mL volumetric flask, add 100 mL of Clorox (approximately 8.25 % sodium hypochlorite) and bring to volume.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

5. SAMPLE PRESERVATION, PREPARATION AND HANDLING

Since biological activity may change the color characteristics of a sample, the color determination should be made as soon as possible. Samples are stored refrigerated ≤ 6°C (not frozen) and must be analyzed within 48 hours of sampling.

6. SAMPLE ANALYSIS AND QUANTITATION (SEE DEP SOP NU-033)

Note: Observe routine maintenance as appropriate to keep the instrument in good working condition and enter any maintenance performed in the maintenance logbook.

6.1. Initial Preparation:

Remove standards from their designated storage area and check the volumes and expiration dates. Discard any expired standards and reprepare if necessary. Make sure there is sufficient volume of all standards to complete the analysis.

6.2. Calibration:

6.2.1. The AQ2 is manually calibrated at least every 30 days using the calibration standards in Section 4.1.

6.2.2. The standards in the calibration are as follows:

STD	Conc. (PCU)	Source
S1	0.0	DI water
S2	6.0 (PQL)	See 4.1.2.
S3	10.5	See 4.1.2.
S4	15.0	See 4.1.2.
S5	30.0	See 4.1.2.
S6	45.0	See 4.1.2.
S7	60.0 (CCV)	See 4.1.2.
S8	90.0	See 4.1.2.
S9	120.0	See 4.1.2.
S10	150.0	Top Standard (4.1.1.)
S0	0.0	DI water

6.3. Controls and Calibration Standard Loading:

6.3.1. Fill two separately labeled AQ2 reagent containers with the following: continuing calibration blank (CCB/DI water) (refill daily) and continuing calibration verification standard (CCV).

6.3.2. The reagent tray, with plastic wedge shaped containers, is located on the right side of the injector arm. It holds the CCV in position 14 and the CCB in position 15. When a calibration is done, the sample tray is loaded from S1 to S0 as in section 6.2.2. starting at cup 1. If calibration is not necessary, the run starts with a CCV and CCB (see Appendix B).

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

6.4. Analysis Setup/Sequence Setup:

6.4.1. Turn on the instrument computer.

6.4.2. Double click on the AQ2 software icon. Select '**Instrument 2 or 3**' as the Log-in name from the drop down list and enter the password to display the Main Operation Screen.

6.4.3. Turn on the main power switch (large switch) at the rear of the instrument. **Allow the instrument to warm-up for at least 15 minutes before performing the Daily Start-up procedure.**

6.4.4. AQ2 System Flush (performed before the Daily Start-up when necessary):

6.4.4.1. Remove AQ2 water inlet sensor and place in a beaker containing the cuvette cleaning solution solution (4.5.).

6.4.4.2. From the maintenance menu perform the following:

a. Prime the syringe - 10 primes (3 times).

b. Verify flow through the probe wash assembly.

c. Wash the cuvette - Auto Wash (at least 10 times).

d. Flush the beaker and the water inlet sensor with DI water.

e. Using fresh DI water, repeat steps a-c with fewer replicates.

f. Repeat d and e to ensure the cleaning solution is removed completely.

Note: A bleach cleaning solution (4.6.) in place of the cuvette cleaning solution may be used monthly or as needed.

6.4.5. Daily Start-up:

a. Empty and refill the large wash container located on the side of the instrument with DI water.

b. Empty the waste containers also located on the side of the instrument (if necessary).

c. Load the instrument with clean reaction well segments. Zero the reaction wells by selecting *Settings→Maintenance* and clicking on the picture of the reaction well.

d. Perform the Daily Start-up procedure by selecting *Settings→Daily Start-Up*.

Note: This procedure primes the syringe, washes the cuvette, and performs water baselines. The water baselines should read 'OK'. If not, contact a supervisor.

6.4.6. Analytical Setup - Select *Scheduling* from the Main Operations Screen. Select a *Tray Number* (**first FREE #**) and *Reagent Set 1* (use *Reagent Set 5* for AQ2 2). Enter the tray name which will be the run file name as follows: CTYYMMDDR# for Instrument AQ2 3 and CTYYMMDDQ#

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

for instrument AQ2 2, where # represents a number1, 2 etc. to indicate the 1st, 2nd etc. run of the day. Press *OK*. The sequence screen should appear. Starting at line 1, type the sample sequence in the 'Sample ID' column. Next, highlight the typed lines and select the **MAN_CLR** test. Click on the *Save* icon to save the schedule and the *Print* icon to print the tray report/bench sheet. Filter the samples using the 0.45 µm filter and syringe apparatus into **labeled** disposable beakers.

6.4.7. If filtration with a 0.45 µm filter does not remove all particulate matter from a sample, or if the Turbidity result for the sample is >15 NTU (see DEP SOP NU-072), it should be filtered using a 0.20 µm filter. If particulate matter is still present, then centrifuge an aliquot of the sample at the maximum speed, and re-filter the supernatant through a 0.20 µm filter. Samples prepared with the 0.20 µm filter are batched together. Record these steps on the tray report/bench sheet where applicable. Samples with results affected by particulate matter may have to be qualified.

Note: In some rare cases precipitation may occur over time in the already filtered aliquot of the sample. Samples should be analyzed as soon as possible after filtration. Always document visible changes of the sample appearance on the tray report/bench sheet.

6.4.8. Pour filtered samples from beakers into the 1.2 mL instrument sample cups according to the schedule. See Appendix B for an example run schedule. Note that each batch begins with an MBLK that is filtered. All samples and QC that appear on the run schedule are poured in individual sample cups. The method is set up to automatically schedule bracketing CCVs and CCBs every 10 samples and at the end of the run. These do not appear in the run schedule and the solutions are picked up from the reagent wedges (6.3.2.). The CCV and CCB may also be poured individually on the run instead of using the scheduling routine described above.

6.4.9. Start the run by double-clicking on the *Run* icon. A dialog box should appear. Select the tray (run file) to be analyzed and click *OK*. On the *Options* dialog box, check the auto-standardize box **only** if you intend to re-calibrate the instrument. Click *OK*.

6.4.10. After the analysis is complete, determine the pH of the samples including any diluted sample aliquots and record on the NU-COLORpH sheet (Appendix C).

First, check the pH of the samples using pH strips with range 5.0 – 9.0 (3.6.2). If the pH of a sample is observed as 5 on this strip range, verify the sample pH using an indicator strip with range 0.0 – 6.0 (3.6.1).

If the pH of a sample aliquot is < 4 or > 9 it must be verified with a pH meter (See DEP SOP NU-026). If the pH is verified, adjust a sufficient aliquot of the sample to pH 7 with NaOH or HCl. The color measurement must be taken from the pH adjusted aliquot. A test comment is entered in

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

the QC Manager indicating that the sample result was obtained from a sample aliquot adjusted to pH 7.

If the pH of a sample is < 2, perform a $BaCl_2$ test to check for improper preservation (see your supervisor if the sample has been preserved with acid).

6.5. Quantitation/Calculations:

6.5.1. The instrument will give the analysis results for Color in PCU.

6.5.2. The QC Manager calculates relative percent differences (RPDs) for duplicate samples as well as recoveries for QC standards. For general calculations including RPDs, see Section 5.9 of the FDEP Quality Manual.

6.5.3. The high calibration standard is prepared according to Section 4.1.1. The other standards used in the linear calibration are prepared from the high standard. The correlation coefficient for the calibration is calculated by the instrument software.

6.6. Note any equipment/field/trip blanks or field duplicates on the backlog. The backlog fields Matrix, Field ID, and Sampling Location are useful for identifying these sample types.

6.6.1. All equipment/field/trip blanks should have a color reading of < 2.5 and > - 2.5 PCU. If the measured color is not in this range, rerun the blank for confirmation on the same run (if possible) and add a comment for the sample in QC Manager.

6.6.2. The RPD for the field duplicates should be less than 20 % if the color values are above the practical quantitation limit (PQL). If the RPD is 20 % or greater, both field duplicates should be rerun for confirmation on the same run, if possible, and commented in QC Manager. Also note any visible differences between the field duplicates on the tray report and in the comments section of QC Manager.

6.7. Post-Run Procedure:

6.7.1. Click on the *Acceptance* menu located on the left side of the screen.

6.7.2. Make any necessary changes to the sample list (dilutions, additional samples, etc.). Click on the *Accept All* icon. This adds 'ACC' to the results. A pop-up window appears where the file name needs to be entered. Follow the naming convention described in 6.4.6. Accepting the run in AQ2 creates the excel transfer file used for uploading the data to the Laboratory Information Management System (LIMS). This file should be saved in the current year's subfolder in \\FLDEP1\DEAR\Labs\lab3\Chem_Data\Inorganic\Nutrients\AQ2\COLOR\.

6.7.3. Print the AQ2 Test Parameters Report from the 'Settings' menu on the top of the screen.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

6.7.4. Click on the *Reporting* icon on the left side of the screen. Check the reporting format 'Color' and *Generate Report.* Print the run report. It contains the calibration and all the data. Select the chronological order report function to print your report.

6.7.5. Convert the excel file into QC Manager raw data. Use instrument 'GenChem: Seal AQ2 3' or 'Nutrients: Seal AQ2 2'.

7. DATA ARCHIVAL

7.1. The raw data is archived in paper format. A report folder is created for each analytical run. Occasionally more than one run may be in one folder. The folder is labeled with the test (w-color), filename/s (CTYYMMDDR# (Instrument AQ2 3), or CTYYMMDDQ# (Instrument AQ2 2)), and initials of the analyst.

7.2. Each report folder should contain the following items - AQ2 run report, AQ2 test parameters, AQ2 tray report/bench sheet, NU-RGTSTD_W-Color sheet, NU-COLORpH sheet, backlog, and a copy of the data and QC assessments uploaded from QC Manager to LIMS at the time of test authorization.

7.3. The report folders are on file with the nutrients and general chemistry supervisory personnel. Periodically, the files are moved to the State Records Center for long time storage according to DEP SOP LB-006.

7.4. Color instrument data excel transfer files for upload to QC manager are stored in electronic format in the current year's subfolder located in: \\FLDEP1\DEAR\Labs\lab3\Chem_Data\Inorganic\Nutrients\AQ2\COLOR\.

8. QUALITY CONTROL

8.1. For method performance and acceptance criteria for QC measures and corrective action for out-of-control data and contingencies for handling out-of-control or unacceptable data, refer to the following DEP SOPs, NU-038 and NU-043.

8.2. A CCV and a CCB are run as check standards at the beginning of the run, every 10 samples, and at the end of each run. After the initial CCV and CCB a PQL level standard and an ALTQC (4.2.) standard are analyzed. Acceptance limits for recoveries of the standards are 70 - 130 % for the PQL and 85 - 115 % for the CCV. The ALTQC recovery should be within the manufacturer's certified limits. Every MBLK (filtered) marks the beginning of a new batch of ten (or less) samples. Acceptance limits for the CCB and MBLK (filtered) are < 2.5 PCU and > -2.5 PCU.

8.3. Analyze a duplicate sample for every batch. The RPD for the duplicates should be less than 20 % provided the duplicate results are above the PQL (6 PCU).

8.4. Calibration Range - See QA Manual Appendix B Section 7.1.1, 'Protocol for Determining the Test Method Range of Applicability'. Acceptance limits for recoveries of calibration standards above the PQL are 90 - 110 % and 70 - 130 %

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

at the PQL level. The correlation coefficient of the calibration must be equal to or greater than 0.995 for the calibration to be acceptable.

8.5. The pH of each aliquot analyzed for color is recorded on the NU-COLORpH sheet (Appendix C). If the pH of a sample is verified to be below 4 or above 10 units, it needs to be adjusted to pH 7. The adjusted aliquot is analyzed and a test comment is added in QC Manager indicating that the aliquot analyzed was adjusted (see 6.4.10.).

9. SAFETY/HAZARDOUS WASTE MANAGEMENT

9.1. Review the Laboratory Health and Safety Plan (accessible on the intranet at https://www.fldepnet.org/content/dear/bureau-laboratories-health-and-safety-plan).

9.2. Safety glasses, gloves and lab coats should be worn at all times whilst handling samples and standards.

9.3. Before performing the color analysis for the first time, read the MSDS sheet on the platinum-cobalt standard used in this analysis.

9.4. All purchased color standards and wastes thereof may be corrosive. Wastes exhibit the characteristic of corrosivity if they are aqueous and have a pH ≤ 2 or $\geq$ 12.5. Test the pH of any instrument waste before discarding it down the sink. If it is determined to be corrosive it must be discarded as hazardous waste (acidic waste drum in the waste room if the pH is below 2). The standard described in section 4.1 must be discarded as hazardous waste (pH < 2).

10. REFERENCES/DEFINITIONS/VALIDATION DATA

For definitions see section 3 of the FDEP Quality Manual. All abbreviations used throughout the document are defined at least once in the document usually at the first time of use. Method validation data is on file with the supervisor.

10.1. Standard Methods Online 2120 B-01.

10.2. DEP SOP NU-026, Manual Determination of pH.

10.3. DEP SOP NU-033, Summary of Sample Analysis Procedures for the General Chemistry and Nutrients Workgroups.

10.4. DEP SOP NU-038, Nutrients and General Chemistry: Quality Assurance and Quality Control.

10.5 DEP SOP NU-042, Standard Preparation Procedure

10.6. DEP SOP NU-043, Nutrients and General Chemistry Quality Assurance and Data Authorization Guidelines.

10.7. Quality Manual, Bureau of Laboratories, FDEP.

10.8. Laboratory Health and Safety Plan.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

10.9. DEP SOP LB-006, Standard operating procedure for records maintenance and storage.

Appendix A


Appendix of Significant Changes

11/17/2011	Section 6.4.5: Added Bleach Cleaning Procedure. Section 6.4.10. and 8.4.: Procedure for recording and adjusting the pH.
01/13/2012	MDL was adjusted to 2.5 PCU. Clarified waste disposal procedure. A filtered MBLK is run at the beginning of every batch (6.4.7.).
02/01/2012	Standard preparation was corrected in Section 4. Naming convention corrected in Section 6.4.6.
05/24/2012	Calibration preparation procedure changed, see Section 4. Calibration and analysis procedure edited in Section 6.
01/08/2015	File naming convention was modified for consistency and to specify the instrumentation used (see 6.4.6. and 7.1.). Data archival information was clarified in section 7.
01/07/2015	Filtering procedure in Section 6.4.7. and pH determination procedure in Section 6.4.10. modified. Appendix C example NU-COLORpH sheet updated.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

Appendix B

Example Tray Report/Bench Sheet:

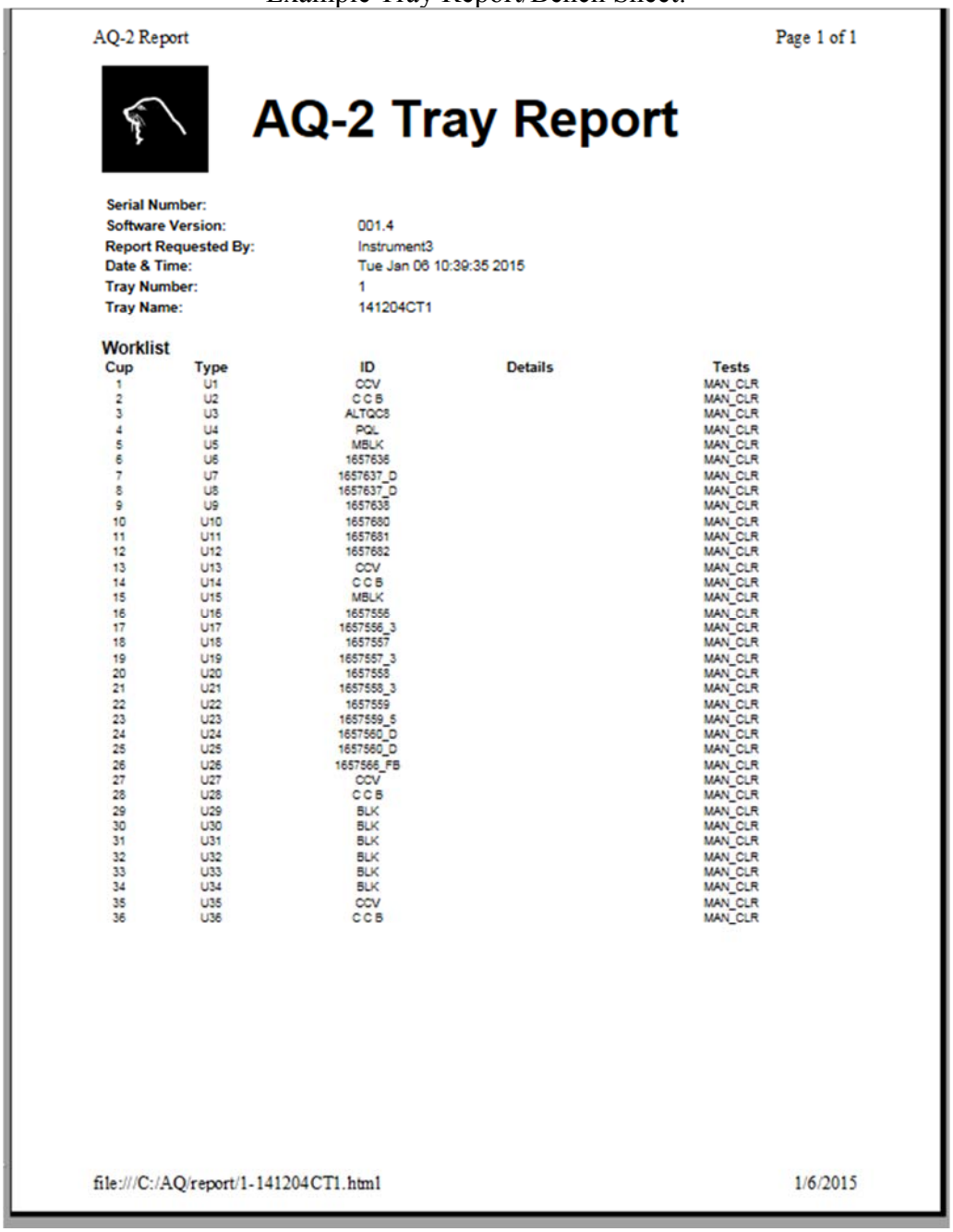

AQ-2 Report Page 1 of 1

AQ-2 Tray Report

Serial Number:
Software Version: 001.4
Report Requested By: Instrument3
Date & Time: Tue Jan 06 10:39:35 2015
Tray Number: 1
Tray Name: 141204CT1

Worklist

Cup	Type	ID	Details	Tests
1	U1	CCV		MAN_CLR
2	U2	C C B		MAN_CLR
3	U3	ALTQC8		MAN_CLR
4	U4	PQL		MAN_CLR
5	U5	MBLK		MAN_CLR
6	U6	1657636		MAN_CLR
7	U7	1657637_D		MAN_CLR
8	U8	1657637_D		MAN_CLR
9	U9	1657638		MAN_CLR
10	U10	1657680		MAN_CLR
11	U11	1657681		MAN_CLR
12	U12	1657682		MAN_CLR
13	U13	CCV		MAN_CLR
14	U14	C C B		MAN_CLR
15	U15	MBLK		MAN_CLR
16	U16	1657556		MAN_CLR
17	U17	1657556_3		MAN_CLR
18	U18	1657557		MAN_CLR
19	U19	1657557_3		MAN_CLR
20	U20	1657558		MAN_CLR
21	U21	1657558_3		MAN_CLR
22	U22	1657559		MAN_CLR
23	U23	1657559_5		MAN_CLR
24	U24	1657560_D		MAN_CLR
25	U25	1657560_D		MAN_CLR
26	U26	1657566_FB		MAN_CLR
27	U27	CCV		MAN_CLR
28	U28	C C B		MAN_CLR
29	U29	BLK		MAN_CLR
30	U30	BLK		MAN_CLR
31	U31	BLK		MAN_CLR
32	U32	BLK		MAN_CLR
33	U33	BLK		MAN_CLR
34	U34	BLK		MAN_CLR
35	U35	CCV		MAN_CLR
36	U36	C C B		MAN_CLR

file:///C:/AQ/report/1-141204CT1.html 1/6/2015

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

Appendix C

Example NU-COLORpH sheet:

Analyst: ____________ **Date:** ____________ **File:** ____________

Sample ID #	Filtered in Lab	pH reading with pH strips	Reading with pH meter		Correction - drops acid or base needed		Adjusted pH reading with pH meter	Notes
			Sample less 4	sample greater than 10	-OH	+H		

	Unique ID	Exp. date
1M HCL:		
1N NaOH:		

pH Range:		pH Range:	
Lot #:		**Lot #:**	
Exp date:		**Exp date:**	

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.