Effective Date: 02/14/2020
Review/Revision Date: 02/14/2020
Revision Author: J.Sherman
MB-01.1-7.15

Standard Operating Procedure for:
Membrane Filter Method for Fecal Coliforms

DEPARTMENT OF ENVIRONMENTAL PROTECTION
BUREAU OF LABORATORIES
BIOLOGY SECTION
TALLAHASSEE FLORIDA

TABLE OF CONTENTS

1. SCOPE AND APPLICATION
2. SUMMARY OF METHOD
    2.1 Interferences
    2.2 Definitions
3. EQUIPMENT AND SUPPLIES
4. REAGENTS AND STANDARDS
5. CALIBRATION AND STANDARDIZATION
6. SAMPLE COLLECTION, PRESERVATION, AND HANDLING
7. SAMPLE PREPARATION
8. SAMPLE ANALYSIS
9. VERIFICATION
10. DATA ARCHIVAL
11. QUALITY CONTROL
12. SAFETY/HAZARDOUS WASTE MANAGEMENT
13. REFERENCES
14. TABLES, DIAGRAMS

Appendix of Changes

1. SCOPE AND APPLICATION

1.1. This method is used to detect and count fecal coliform bacteria in non-potable water and is based on Standard Methods SM9222D. Since fecal coliforms are commonly found in the feces of humans and warm-blooded animals, their presence can indicate the potential presence of bacterial or viral pathogens.

1.2. This procedure is used for samples with Laboratory Information Management System (LIMS) test IDs: FCOLI-MF, FC-MF-CL.

2. SUMMARY OF METHOD

2.1. A water sample is filtered through a membrane that has a 0.45 µm pore size to capture bacteria. The membrane filter is placed on an mFC dish, which is a selective medium for fecal coliforms. The dish is then incubated for 24 ± 2 hours at 44.5 ± 0.2 °C . Positive colonies, which have a blue color, are counted and recorded.

2.2. Interferences

2.2.1. Suspended particulate materials in sample water can clog the membrane filter and cause spreading of bacterial colonies. Also, the incubation temperature is critical to the specificity of the test.

2.3. Definitions

2.3.1. **Fecal Coliforms** - Gram-negative, non-spore forming, rod-shaped bacteria that ferment lactose in 24 ± 2 hours at 44.5 ± 0.2 °C and produce blue colonies on mFC media. Fecal coliforms are a subgroup of total coliforms that includes *Escherichia* and *Klebsiella* from the Enterobacteriaceae family. The predominant species is *E. coli*.

2.3.2. **Target Organism** – The target organism is the specific bacteria for which the method is testing. In this method the target organism is fecal coliforms.

2.3.3. **Positive Growth** – Bacterial growth that resembles a typical positive colony which is a blue colored colony.

2.3.4. **Negative Growth** – Bacterial growth that does not have the characteristics of a typical positive colony and is considered a non-coliform.

3. EQUIPMENT AND SUPPLIES

Record any equipment repairs or maintenance in LIMS Laboratory Logbook *Maintenance_Log*

3.1. Sterile gloves, latex or nitrile

3.2. Petri dishes, 47 mm, sterilized, with tight fitting lids

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

3.3. Flask, Erlenmeyer, 2 L

3.4. Weighing Dish, Pour Boat Style, medium size

3.5. Weighing paper, medium size

3.6. Graduated cylinder, 1 L

3.7. Magnetic stir bar, medium size

3.8. Balance, Mettler Toledo MS303S (S/N B247544008) - Calibrate per SOP PCR 2.1

3.9. Hot plate with stirrer

3.10. Autoclave, Steris Amsco 250 LS Model 20VS (S/N 031761851)

3.11. pH meter, Accumet Excel XL 25 (S/N XL94005930) - Calibrate per SOP MB-03.4

3.12. Repeater pipette with 50 mL tips - Calibrate per SOP MB-03.3

3.13. Isotemp Refrigerator (S/N 0162177801170112), maintained at 1-5 °C

3.14. Chlorine test strips

3.15. Vacuum Filtration system (See Figure 14.1)

3.16. Membrane filtration unit, sterile wrapped in foil

3.17. UV Sterilizer unit

3.18. Forceps, sterile, straight or curved with smooth tips

3.19. Incinerator

3.20. Beaker, 30 mL

3.21. Serological disposable pipettes (1, 2, 5, 10, 50 mL as needed), sterile and individually wrapped, with pipette pump. The volume of each new lot is checked per SOP MB-03.3

3.22. Membrane filters, sterile, white gridded, 47 mm diameter, 0.45±0.02 µm pore size

3.23. Inoculating loops, sterile

3.24. Test tubes, 16 X 125 mm, with screw caps

3.25. Stereo microscope, Fisher Scientific

3.26. Heratherm IMH incubator (S/N 41784672), maintained at 35.0 ± 0.5 °C

3.27. Thermo Precision Water Bath (S/N 300263160), maintained at 44.5 ± 0.2 °C

4. REAGENTS AND STANDARDS

4.1. Deionized Water

4.2. Buffered Dilution Water - Prepare per SOP MB-02.24.

4.3. Media - Prepare per SOP MB-02.33.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

4.3.1. mFC broth (mFC) or mFC agar (mFCA)

4.3.2. Lauryl Tryptose Broth (LTB)

4.3.3. EC medium (EC)

4.4. Control cultures (See SOP MB-03.1.)

4.4.1. *Escherichia coli* (*E. coli*) ATCC #11775

4.4.2. *Pseudomonas aeruginosa* (*P. aeruginosa*) ATCC #10145

5. CALIBRATION AND STANDARDIZATION

5.1. See the Laboratory Quality Manual.

6. SAMPLE COLLECTION, PRESERVATION, AND HANDLING

6.1. See DEP Field SOP FS 1000 for sample collection and preservation.

7. SAMPLE PREPARATION

7.1. See DEP SOP BB-30 for sample custody and labels.

7.2. Turn on UV sterilizers and let warm up for at least 10 minutes.

7.3. At the beginning of each day, gather 6 funnels (3 with bases), unwrap and place 3 funnels in the UV sterilizer. Place bases in vacuum manifold with the 3 remaining funnels. Use a new set of funnels daily.

7.4. Remove any bottles of Buffered Dilution Water from the refrigerator and allow them to reach room temperature before using.

> **Note:** Any opened bottles of Buffered Dilution Water must be refrigerated after sample filtration has been completed for the day.

7.5 Place forceps in a small beaker of 95% ethanol and turn on the incinerator. Fill a new, clean, dry rinse bottle with buffered dilution water.

7.6 Select sample volumes that will produce 20-60 colonies based on expected pollution level. A minimum of three dilutions should be analyzed for each sample, with the exception of field or equipment blanks, to ensure that a countable dish is obtained. For field or equipment blanks, one 100 mL aliquot should be filtered. If sample is turbid, filter 0.1, 1, and 10 mL. For most samples, 1, 10, and 50 mL dilutions should be sufficient. 50 mL of buffered dilution water should be filtered for method blanks. Also, consider the sample's previous results when choosing volumes or refer to Standard Methods Table 9222:III.

7.7 Place sample label on the middle section of each mFC dish lid. Write the volume to be filtered on the lower section of the lid (below the label) with a black marker. If using mFC broth, add 2 mL of the broth to the absorbent pad in each dish.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

8. SAMPLE ANALYSIS

8.1. Use sterile forceps to grab and place a sterile membrane filter, grid-side up, on the filter base and attach the funnel to the base. If analyzing sample volume <20 mL, add 20-30 mL of buffered dilution water to the funnel **prior** to adding sample.

Note: Keep forceps in a small beaker of 95% ethanol during filtration process. Dip in the ethanol and flame-sterilize it by passing it through the flame of the incinerator before grasping each filter. Use caution. If forceps are accidentally dropped or placed on counter, replace with a new pair of sterile forceps from the UV sterilizer.

8.2. Shake the sample container vigorously, approximately 25 times to uniformly distribute bacteria. Carefully open the container. Avoid touching the top rim of the container to prevent contamination of the sample.

8.2.1 If the test ID includes –CL, check the sample for chlorine with a chlorine test strip. Pour a small volume of sample into a small container (such as a Durham tube), just enough to submerge the pads at one end of the strip. Swirl the test strip in the sample for 30 seconds. If the pads remain colorless, then the sample has already been dechlorinated and should be analyzed. If the pads turn purple, chlorine is present. Do not run the sample. Report as "No Result" with a "Y" qualifier for improper preservation. Record the lot number and result of the chlorine test strip in the **Chlorine Residual** column in Prep Worksheet Manager. See SOP BB-30 for instructions on using Prep Worksheet Manager.

8.3. Measure appropriate volume of sample with a sterile pipet into the funnel and filter. Rinse the walls of the funnel with 20-30 mL of buffered dilution water (rinse all the way around 3 times).

8.4. Turn vacuum off and aseptically remove filter with sterile forceps.

8.5. Place filter, grid side up, on the petri dish by rolling it onto the mFC agar or broth soaked pad. Avoid creating bubbles between the membrane and the medium. If bubbles do occur, lift and reset the filter. Also, run forceps along outside edge of filter to help it bind to the surface of the medium. If using mFC broth, make sure to carefully decant excess liquid from the plate before placing the membrane on the pad. Close dish, invert and set to the side.

8.6. Place the three used funnels into the UV sterilizer and place the other three that have been sterilized onto the funnel bases.

8.7. Repeat Steps 8.1.-8.6. until all dishes for all samples are filtered. Make sure funnels have been UV sterilized for at least 2 minutes before starting each new sample.

8.8. Incubate inverted dishes at 44.5 ± 0.2 °C for a period of 24 ± 2 hours. Dishes should be incubated no later than 30 minutes after they are filtered. Do not wait for all the samples to be filtered before incubation.

8.8.1. Place inverted dishes into a waterproof plastic bag and seal. Write on the outside of the bag with permanent marker the date and time of

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

incubation. Submerge bag in the Thermo Scientific Water Bath 1 (TSWB 1) and add a weight to secure them below the water line to maintain temperature requirements.

8.9. After incubation, count and record on dish above sample label all colonies that are various shades of blue. Negative growth may produce grey to cream-colored colonies. Verify any colonies that are questionable. Use a fluorescent lamp and dissecting microscope for better visibility of colonies.

9. VERIFICATION

9.1. Before verification, remove LTB broth and EC broth tubes from the refrigerator and incubate for 24 hours at 35.0 ± 0.5 °C and 44.5 ± 0.2 °C respectfully. Wrap the plate that will be verified with aluminum foil and place in the refrigerator. Let the plate reach room temperature before inoculating. After the broth tubes have been incubated, use a sterile inoculating loop to transfer a loopful of growth from a well-isolated colony to an LTB broth tube and into an EC broth tube (discard each inoculating loop after it has been transferred from both LTB to EC broth tube before moving onto the next set of tubes). Document the type of isolated colonies that will be transferred into the broth tubes in the *Monthly Verification Tests Logbook* (for example: colonies #1-3 are smooth, round, large colonies; colonies #4-5 are elevated, irregular, small colonies). Incubate LTB tubes for 48 ± 3 hours at 35.0 ± 0.5 °C. Incubate EC tubes at 44.5 ± 0.2 °C for 24 ± 2 hours. Observe for growth (turbidity) and gas formation in the Durham tube.

9.2. Colonies that produce gas in both LTB within 48 ± 3 hours and EC within 24 ± 2 hours of incubation are verified as positive for fecal coliforms.

9.3. Record results in *Monthly Verification Tests Logbook.*

10. DATA ARCHIVAL

10.1. Record filter, incubation, read dates/times, analyst, buffered water and media used, QC record numbers for consumables used, and blanks and dups with their respective Test IDs for each sample in Prep Worksheet Manager per SOP BB-30.

10.2. When dishes are read, open the worksheet and record the read dates/times and analyst. Review the prep date/time, incubation dates/times, read dates/times, prep and read analysts, and media and reagents for accuracy and ensure the correct number of blanks and dups are on the worksheet and are labelled correctly. Comment as necessary. Enter the Date/Time Reviewed and select the Reviewed By analyst. Click the **Complete** button and close the worksheet.

10.3. Enter counts into the *MICRO_DATA_LOG* in the LIMS Laboratory Logbook. Upload results to the QC Manager per SOP BB-31.

10.4. After results are uploaded to QCM, refer to SOP BB-32 for data authorization.

10.5. Document all verification of colonies in the *Monthly Verification Tests Logbook.*

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

10.6. Calculations (based on SM 9222B.6):

10.6.1. If at least one dilution is within the ideal colony count range of 20-60 colonies, use only the dilutions that fall within this range to calculate the final result in the equation below.

10.6.2. If no dilutions fall within the ideal colony count range, use all the dilutions to calculate the final result in the equation below and add a "B" qualifier code. If a 100 mL dilution was filtered, and it has a count less than 20 colonies use only this dilution and do not add a "B" qualifier.

10.6.3. If all dilutions have a colony count of <1, use 1 as the total number of colonies counted and use all the dilutions for the total volume filtered to calculate the final result in the equation below and add a "U" qualifier code.

10.6.4. If there are >200 target colonies in all dilutions, use the upper limit of the ideal range (60) and the smallest filtration volume in mL to calculate an estimated final result. Report with a "Z" qualifier.

10.6.5. If there are >200 non-target colonies, or if the colonies are not distinct enough to count (confluent growth) in all dilutions, report as "No Result" with a "Z" qualifier. A comment can be added, such as "Smear" or "Confluent", to provide a reason for the "No Result".

$$CFU/100mL = \frac{C \times 100}{V}$$

where,
CFU/100 mL = Fecal coliform colony forming units per 100 mL of sample.
C = Total number of positive colonies counted on all acceptable dilutions.
V = Total volume (sum of all acceptable dilutions)

11. QUALITY CONTROL

11.1. Sample Replicates

11.1.1. Analyze a laboratory replicate for 10% of samples with at least one replicate per filtration batch.

11.1.2. Maintain a control chart of duplicate results. Convert results to logarithms and calculate the range (difference of logarithms of a replicate set). The range should be < 3.27 times the mean range of the last 15 replicate sets.

11.1.3. Corrective Actions: Report laboratory replicate results that do not meet the acceptance criteria with a "J" qualifier and the comment "Result is J-qualified due to a lab replicate precision failure".

11.2. Method Blanks

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

11.2.1. Filter a 50 mL volume of buffered dilution water for each funnel used at the beginning and end of a sample filtration batch. Place the filter on an mFC dish and label the lid with "CB" (Control Before) or "CA" (Control After), funnel number, and date/time. If more than 30 minutes elapses between filtration batches, run a new set of controls. Incubate blanks with samples at 44.5 ± 0.2 °C for 24 ± 2 hours.

11.2.2. Absence of growth in method blanks demonstrates that media, filtration equipment, filters, and dilution water are free of contamination of the target organism from improper handling, inadequate sterilization or environmental exposure.

11.2.3. Corrective Actions: If method blanks have growth, determine the cause of contamination. The rinse water batch, media, and equipment can be checked for sterility. If the cause of contamination has also affected samples, report samples in the filtration batch that are greater than MDL with a "V" qualifier and the comment "Result is V-qualified due to control blank failure."

11.3. Negative Controls

11.3.1. For each batch of prepared media, transfer a small amount of *P. aeruginosa* to a 99 mL phosphate buffered water bottle and shake vigorously. Filter a 1 mL volume of this suspension and place filter on an mFC dish. Place a media label on the dish and mark it with a "-" sign for negative control. Incubate at 44.5 ± 0.2 °C for 24 ± 2 hours.

11.3.2. To check viability of negative control, also add a small amount of the suspension to a BHI broth tube and incubate at 35.0 ± 0.5 °C for 24 ± 2 hours. If negative control does not grow in BHI broth, then the *P. aeruginosa* may not be viable and the negative control test on the media should be repeated with a new culture.

11.3.3. Absence of growth in negative control indicates that the medium does not support growth of non-target organisms. Any negative growth indicates that non-target organisms will not resemble typical positive colonies.

11.3.4. Corrective Actions: If there is positive growth on a negative control, determine cause of contamination. If the cause of contamination has affected the media, discard media batch and prepare a new batch.

11.4. Positive Controls

11.4.1. For each batch of prepared media, transfer a small amount of *E. coli* to a 99 mL phosphate buffered water bottle and shake vigorously. Filter a 1 mL volume of this suspension and place filter on an mFC dish. Place a media label on the dish and mark it with a "+" sign for positive control. Incubate at 44.5 ± 0.2 °C for 24 ± 2 hours.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

11.4.2. Positive growth indicates that the medium does support growth of target organisms and produces the expected reaction (growth resembles typical positive colonies).

11.4.3. Corrective Actions: If no growth occurs, re-run the positive control with a new culture of *E. coli*. Also, run a pour plate with HPC (See SOP MB-01.5.) on the new culture to check for viability. If there is growth on the HPC medium but no growth on the mFC medium, then the media batch does not support growth of target organisms. Discard the batch and prepare a new media batch.

11.5. Sterility Checks

11.5.1. Check sterility on each batch of media by incubating an mFC dish along with positive and negative control dishes. Place a media label on the dish and mark it with and "S" for sterility control. If growth occurs, discard media batch and prepare a new batch.

11.5.2. Check sterility on each lot of membrane filters, each lot of sample containers, and each sterilized batch of glassware by incubating with a small amount of TSB 2X (see SOP MB-02.33) that has been swirled around to contact all surfaces at 35.0 ± 0.5 °C for 24 ± 2 hours. If there is growth, recheck lot or batch. If recheck also fails (shows growth), do not use the lot or rewash batch of glassware and contact a supervisor. Record results of sterility checks in the appropriate LIMS **Standard Preparation Tracker** module or **Laboratory Logbook** log. **Sterility of sample containers should be checked before they are packed into kits and shipped to customers.**

11.6. To enter media QC results, go to LIMS Laboratory/Technical/Standard Preparation Tracker. Choose the Group "BIO-MICRO". Click Media → Media QC Pass/Fail. Highlight the media preparation, enter the Date QC Test Read, QC Test Read By, Sterility Results, and click Pass or Fail.

11.7. Sodium Thiosulfate Dechlorination Check

11.7.1. Check each lot of sample containers that have a sodium thiosulfate tablet for adequate dechlorination.

11.7.1.1. Dilute bleach to 0.1%. Check the concentration with the Hach Colorimeter II per SOP TA-06.22 and record in Standard Preparation Tracker.

11.7.1.2. Dilute the 0.1% solution to 16 mg/L Cl_2 (also check the concentration with the Colorimeter II) and record the new solution in Standard Preparation Tracker.

11.7.1.3. Add 100 mL of the 16 mg/L Cl_2 solution to the sample container with the tablet that is being tested. Shake until the tablet is dissolved. Use a chlorine test strip to determine if the solution has been dechlorinated. Record the date, analyst, SN of the 16 mg/L Cl_2

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

solution and its concentration reading from the Colorimeter II, the final concentration of the solution after dechlorinating, and the lot number of the chlorine test strips in the ENTRY_COMMENTS of the corresponding sterility check record for the sample container lot in LIMS **Laboratory/Technical/Laboratory Logbook** *Lab_Supplies_QC*. If the chlorine test strip shows that the solution still has chlorine, repeat the test. If the repeat test fails (does not dechlorinate), do not use the lot and contact a supervisor.

11.8. Colony Verification

11.8.1. Verify at least 10 typical colonies per month. Verify atypical colonies to check for false negatives. See Section 9 for verification procedure.

11.8.2. Adjust final colony counts for verified samples.

11.9. Colony Counting Variability

11.9.1 Each trained analyst must compare colony counts against another trained analyst from one positive sample monthly. Counts should be within 10% of each other. Comparisons must be performed by all trained analysts performing this test within one hour after the dish has been removed from incubation. When comparing counts, each trained analyst will write their count on the back of the plate along with their initials. The counts will then be recorded into the MB-PARCOUNTLOG.

12. SAFETY/HAZARDOUS WASTE MANAGEMENT

12.1. All media/bacteria must be autoclaved for 45 minutes before disposal.

12.2. Refer to the Laboratory Health and Safety Plan.

12.3. Read MSDS information for reagent handling and disposal.

12.4 The UV sterilizer unit has an automatic shutoff of the lamp when opened. In the event this shutoff malfunctions or needs to be overridden to perform maintenance or repairs on the unit, wear UV blocking safety glasses. Do not look at the UV bulb without eye protection as serious eye damage can result.

13. REFERENCES

13.1. Standard Methods 9222 D and 9222 B -2006

13.2. Standard Methods 9020 -2017

13.3. NELAC Quality Systems Manual

13.4. Difco/BBL Manual

13.5. FDEP Laboratory Quality Manual e
http://www.dep.state.fl.us/labs/library/progplan.htm

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

13.6. FDEP Laboratory Health and Safety Plan
http://www.fldepnet.org/content/dear/bureau-laboratories-safety-documents

13.7. DEP SOP-001/01 FS 1000, *General Sampling Procedures*

13.8. DEP SOP PCR-2.1, *Calibration, Use and Maintenance of Mettler Toledo MS303S Balance*

13.9. DEP SOP MB-01.5, *Pour Plate Method for Heterotrophic Plate Count*

13.10. DEP SOP MB-02.24, *Preparation of Buffered Dilution Water and Stock Components*

13.11. DEP SOP MB-02.33, *Preparation of Dehydrated Culture Media*

13.12. DEP SOP MB-03.1, *Maintenance of Bacterial Cultures*

13.13. DEP SOP MB-06.18, *AMSCO Renaissance Series 3021 Autoclave Calibration and Use*

13.14. DEP SOP MB-03.3, *Calibration and Use of Mechanical and Disposable Pipettes*

13.15. DEP SOP MB-03.4, *Accumet Excel XL 25 pH Meter Calibration and Use*

13.16. DEP SOP BB-30, *Sample Custody, Preparation Labels and Worksheet Instructions for Bench Biology Samples*

13.17. DEP SOP BB-31, *LIMS QC Manager Data Uploading Instructions for Bench Biology and Microbiology Samples*


13.18. DEP SOP BB-32, *Test Level Authorization of Uploaded Results for Bench Biology Samples*

13.19. DEP SOP TA-06.22, *Measuring Total Residual Chlorine using the Hach Pocket Colorimeter and Colorimeter II Test Kits*

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

14. TABLES, DIAGRAMS

Figure 14.1 Set up for Vacuum filtration system

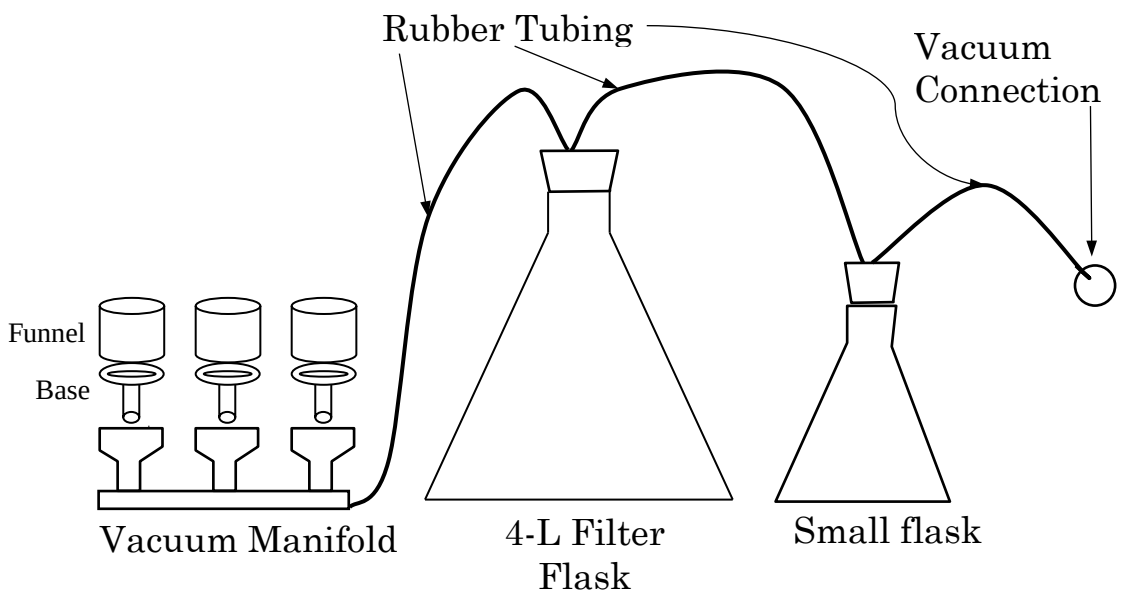

APPENDIX OF CHANGES

08/21/2008	Converted to new format and updated to Standard Methods Online Edition.
08/21/2009	Added new LIMS test ID and prep worksheet documentation.
08/19/2010	SOP reviewed and edited to update according to current procedures. Removed reference to older version of safety plan and added new version.
02/02/2011	Updated references to current version of quality manual
01/31/2011	Updated incubator and microscope information. Added corrective actions and removed reference to QAM for corrective actions.
01/28/2013	Added decimal place to some temperatures. Added items to supplies list.
01/15/2014	Removed old balance and added new balance to Equipment and Supplies; replaced “whirl-pak” with “sample container;” added a separate step for changing out dirty and sterilized funnels in Sample Preparation; added comment about checking sample container sterility in Sterility Checks.
12/16/2014	Corrected the Sterilmatic autoclave Model ID. Updated DATA ARCHIVAL to reflect the changes made in LIMS Prep Worksheet Manager. Clarified “V” qualifier in Corrective Actions for Method Blanks failures.
12/23/2014	Formatted for 508 compliance.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.

03/10/2015 Added procedure for checking volume on funnels and doing QC on sodium thiosulfate tablets.

03/14/2017 Updated equipment information, removed funnel calibration, added prep worksheet review, removed "> 10 x blank" instructions for blank failures, updated parallel count instructions, and corrected hyperlinks.

03/12/2018 Specified the range of incubation time and temperature on the Summary of Method section. Paragraph revision on section 8.5 to include decanting of excess media if using mFC broth on petri dishes. Added paragraph 8.8.1 to mention the incubators available for Fecal Coliform samples. Instructions of incubation in the water bath also specified. Updated the Data Archival paragraph and included the addition of the "incubator ID" information to the work sheet. Updated equipment information in sections 3 and 13.

09/12/2018 Updated Standards Methods version numbers in section 13.

09/04/2019 Updated autoclave and water bath information. Removed equipment that is no longer in service. Updated section 7. Sample Preparation and 11.8 Colony Verification.

01/10/2020 Added procedure for open bottles of Buffer water in section 7.3. Added information for broth tubes used for verification before inoculation; added procedure of discarding inoculating loops during each transfer; added documenting type of isolated colonies used for verification in section 9. In section 11.9, removed "and signed off on in Microbiology QC checklist;" corrected on where analyst will write parallel count; added which analyst count are excluded into the MB-PARCOUNTLOG.

The OFFICIAL COPY is the on-line version. All other copies are considered unofficial and uncontrolled.